

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

> [CAMRY](#) /

> [CAMRY Smart Bluetooth Body Fat Scale iF1060D-B User Manual](#)



CAMRY iF1060D-B

CAMRY Smart Bluetooth Body Fat Scale

Model: iF1060D-B

1. INTRODUCTION

Thank you for choosing the CAMRY Smart Bluetooth Body Fat Scale. This device is designed to provide comprehensive body composition analysis, helping you monitor your health and fitness progress. It measures various metrics including body weight, body fat, muscle mass, and more, and synchronizes data with a dedicated app on your smartphone.

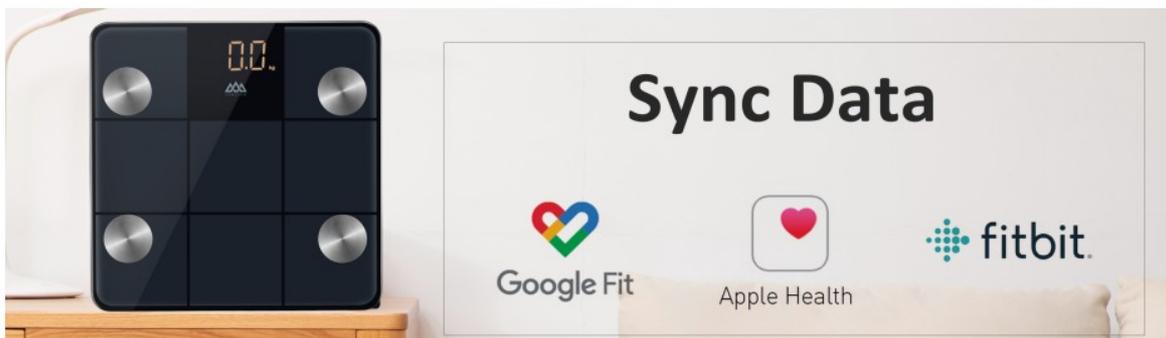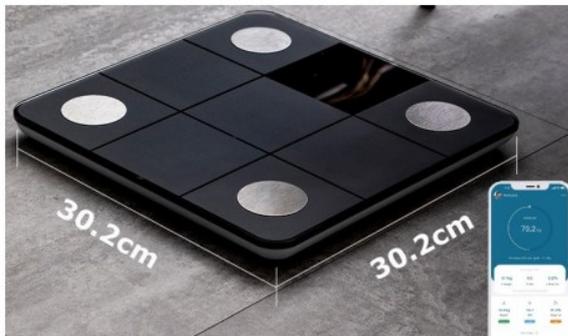

The CAMRY Smart Bluetooth Body Fat Scale, designed for comprehensive body analysis and seamless app integration.

2. SAFETY INFORMATION

- Do not use if you have a pacemaker or other internal medical device.
- Do not use if you are pregnant.
- Keep the scale dry. Do not immerse in water or use in extremely humid environments.
- Place the scale on a hard, flat surface for accurate measurements. Avoid carpets or uneven floors.

- The scale surface can be slippery when wet. Step on and off carefully.
- This scale is for household use only and should not be used for medical diagnosis or treatment.

3. PRODUCT OVERVIEW

The CAMRY Smart Bluetooth Body Fat Scale features advanced technology for precise body composition analysis. It includes high-sensitive electrodes and sensors for accurate readings.

Key Features:

- **13 Body Composition Functions:** Measures body weight, body fat, muscle, protein, bone, BMI, BMR, hydration, fat-free body mass, AMR, subcutaneous fat, skeletal muscle, and body shape.
- **High Precision:** Equipped with 4 highly sensitive electrodes and sensors for accurate measurements up to 180 kg (0.1 kg increments).
- **LED Display:** Clear LED monitor display for instant weight readings.
- **Smart App Synchronization:** Connects via Bluetooth 4.0 to the SenssunLife app, compatible with Android 4.3+ and iOS 7.0+.
- **Fitness App Integration:** Syncs data with Apple Health, Google Fit, and Fitbit.
- **Automatic Calibration & Shut-off:** Ensures ease of use and battery efficiency.

Components:

- Tempered Glass Platform
- High Sensitive Electronic Pads (Electrodes)
- LED Monitor Display
- High Precision Sensors
- Anti-slip Stable Pads



Internal structure of the scale, highlighting high-sensitive electronic pads, LED display, high precision sensors, and anti-slip stable pads.



The scale measures 15 key body composition values, including weight, BMI, body fat, and more.


4. SETUP

4.1 Battery Installation

The scale requires 3 AAA batteries (included). Open the battery compartment on the underside of the scale and insert the batteries, ensuring correct polarity (+/-).

4.2 App Download and Pairing

1. Download the **SenssunLife** app from the App Store (for iOS devices) or Google Play (for Android devices).
2. Enable Bluetooth on your smartphone or smart device.
3. Open the SenssunLife app and follow the on-screen instructions to create a user profile.
4. Step onto the scale to activate it and pair it with your device. Your body composition analysis will be measured and synchronized automatically.

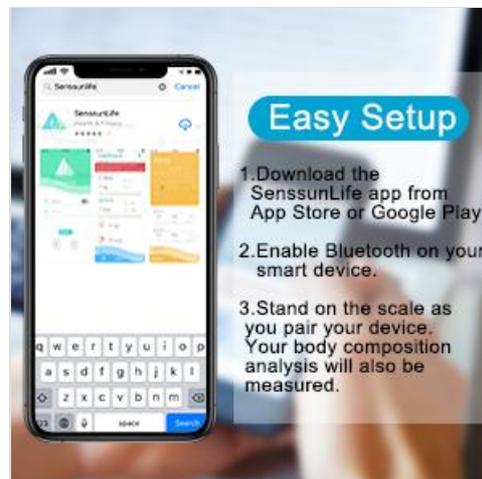

Easy setup process: Download the app, enable Bluetooth, and step on the scale to pair.

5. OPERATING INSTRUCTIONS

5.1 Taking a Measurement

1. Place the scale on a hard, flat surface.
2. Step onto the scale barefoot. Ensure your feet are in contact with the metal electrodes for accurate body composition analysis.
3. Stand still until your weight is displayed on the LED screen.
4. If the app is open and Bluetooth is enabled on your device, the data will automatically synchronize. If not, the scale will store the weight, and it will sync once connected to the app again.



The large backlit LED display provides clear weight readings, even in low light conditions.

6. APP INTEGRATION AND DATA TRACKING


The SenssunLife app is your central hub for tracking and analyzing your health data. It provides detailed insights into your body composition and allows you to monitor trends over time.

Key App Features:

- **Data Storage:** All measurement data is stored within the app.
- **Trend Charts:** Visualize your progress with body composition trend charts.
- **Goal Setting:** Set personal fitness goals and track your achievements.
- **User Profiles:** Supports multiple user profiles, allowing the whole family to track their data separately.
- **Third-Party App Sync:** Seamlessly integrates with Apple Health, Google Fit, and Fitbit.

User-Friendly App

SenssunLife APP allows you to set a goal, monitor your fitness trends with an automatic chart.

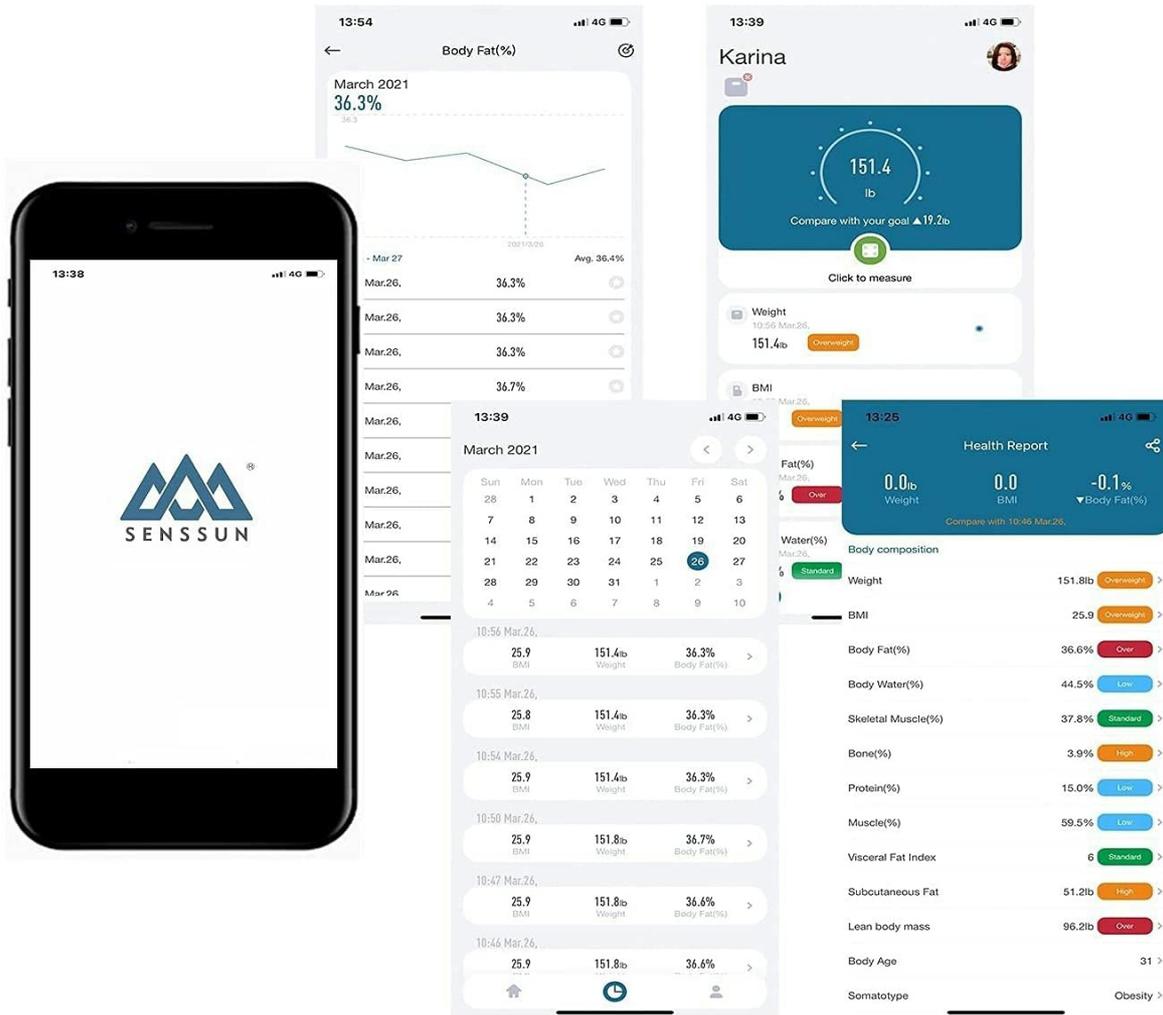

The SenssunLife app provides a user-friendly interface to monitor body fat, weight, and historical data.



Track your progress with detailed trend charts for body weight, BMI, and body fat within the app.

7. MAINTENANCE

- Clean the scale with a damp cloth and mild detergent. Do not use abrasive cleaners or immerse the scale in water.
- Store the scale in a cool, dry place.
- Remove batteries if the scale will not be used for an extended period to prevent leakage.

8. TROUBLESHOOTING

Inaccurate Readings:

- Ensure the scale is on a hard, flat surface.
- Make sure your feet are clean, dry, and in full contact with the electrodes.
- Recalibrate the scale by stepping on it briefly, then stepping off and waiting for it to display '0.0' before taking a measurement.
- Check battery levels. Low batteries can affect accuracy.

App Connection Issues:


- Ensure Bluetooth is enabled on your smartphone.
- Make sure the SenssunLife app is open and running in the foreground or background.
- Restart your smartphone and the scale (by removing and reinserting batteries).
- Ensure your smartphone meets the minimum OS requirements (Android 4.3+ or iOS 7.0+).

Scale Not Turning On:

- Check if batteries are inserted correctly.
- Replace batteries with new ones.

9. SPECIFICATIONS

Feature	Detail
Model Number	iF1060D-B
Dimensions	30 x 30 x 0.6 cm
Weight	1.8 kg
Weight Limit	180 Kilograms
Measurement Increments	0.1 kg
Display Type	LED
Connectivity	Bluetooth 4.0
Power Source	3 x AAA Batteries (included)
Material	ABS
Certifications	CE, RoHS approved

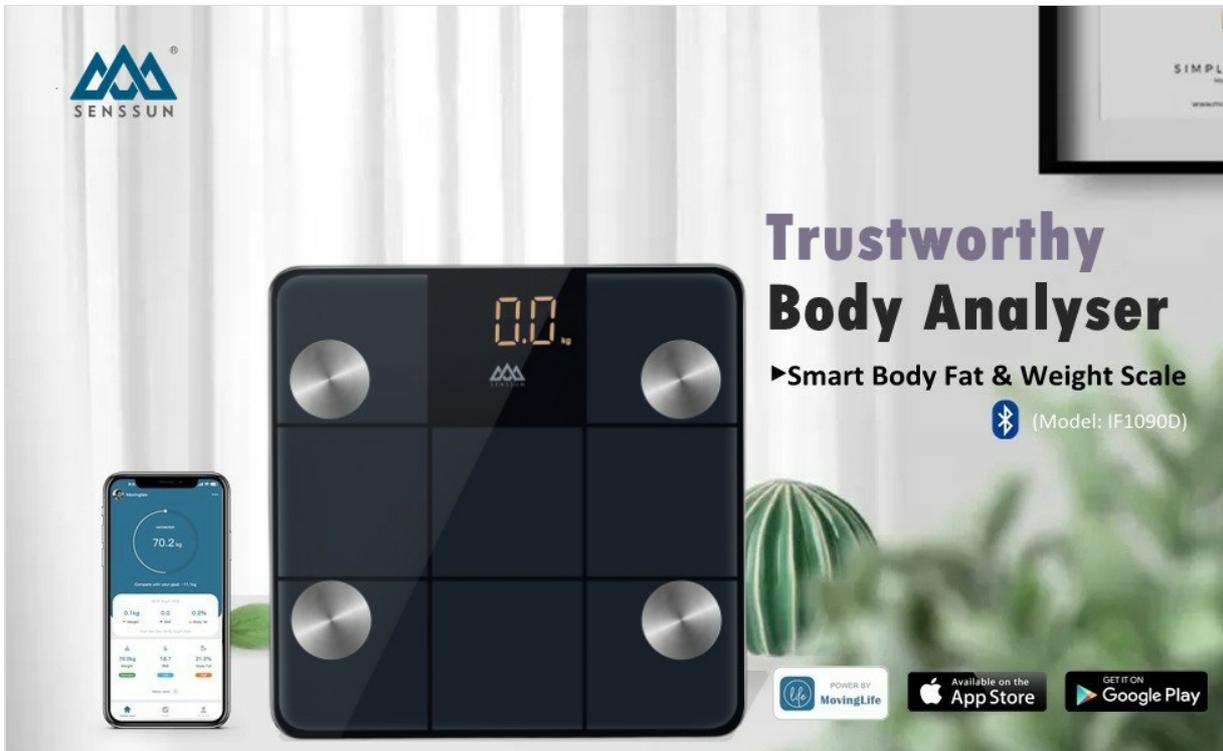

Detailed dimensions of the scale: 30.2 cm x 30.2 cm x 2 cm (height).

10. WARRANTY AND SUPPORT

CAMRY provides the following for your Smart Bluetooth Body Fat Scale:

- **30-Day Money-Back Guarantee:** If you are not satisfied with your purchase.
- **1-Year Warranty:** Covers manufacturing defects and malfunctions.
- **Lifetime Support:** Provided by SENSSUN.
- **Customer Service:** Expect a mail reply service within 24 hours for any inquiries or support needs.

For support, please refer to the contact information provided with your purchase or visit the official SENSSUN website.