

Manuals.plus /

- › Philips Sonicare /
- › Philips Sonicare ExpertClean 7500 Electric Toothbrush User Manual HX9690/06

## Philips Sonicare HX9690/06

# Philips Sonicare ExpertClean 7500 Electric Toothbrush User Manual

MODEL: HX9690/06  
Brand: Philips Sonicare

## INTRODUCTION

---

The Philips Sonicare ExpertClean 7500 electric toothbrush is designed to enhance your oral health routine. It offers advanced features to guide you towards improved brushing habits and achieve a superior clean. This manual provides essential information for setting up, operating, and maintaining your device to ensure optimal performance and longevity.

## WHAT'S INCLUDED

---

Your Philips Sonicare ExpertClean 7500 package includes the following components:

- 1 Philips Sonicare handle
- 1 charging travel case
- 1 charger
- 1 Premium Plaque Control brush head
- 1 Premium Gum Health brush head

## SETUP

---

### 1. Charging the Toothbrush

Before first use, fully charge your Philips Sonicare ExpertClean 7500. Place the handle on the charging base. The battery indicator lights will illuminate to show charging progress. A full charge typically takes 24 hours and provides up to 2 weeks of regular use.

For travel, utilize the included charging travel case. Connect the travel case to a power source using the provided USB cable to charge the toothbrush while on the go.

### 2. Attaching the Brush Head

Align the brush head with the metal shaft on the handle and push it down firmly until it clicks into place. Ensure a secure fit before use.

# Fits all Philips Sonicare click-on handles



Image: The Philips Sonicare handle is compatible with all Philips Sonicare click-on brush heads, which are easily interchangeable.

## OPERATING INSTRUCTIONS

### 1. Powering On/Off and Brushing

Wet the bristles and apply a small amount of toothpaste. Place the brush head against your teeth at a slight angle towards the gumline. Press the power button to start brushing. Gently move the brush head slowly across your teeth, allowing the bristles to do the work. The built-in timer will guide you through the recommended 2-minute brushing cycle, pausing briefly every 30 seconds to indicate when to move to the next quadrant of your mouth.

### 2. Pressure Sensor

The toothbrush features a pressure sensor that alerts you if you are brushing too hard. If excessive pressure is applied, the handle will vibrate differently, and a light ring may illuminate (depending on the model). Reduce pressure to protect your gums and teeth.

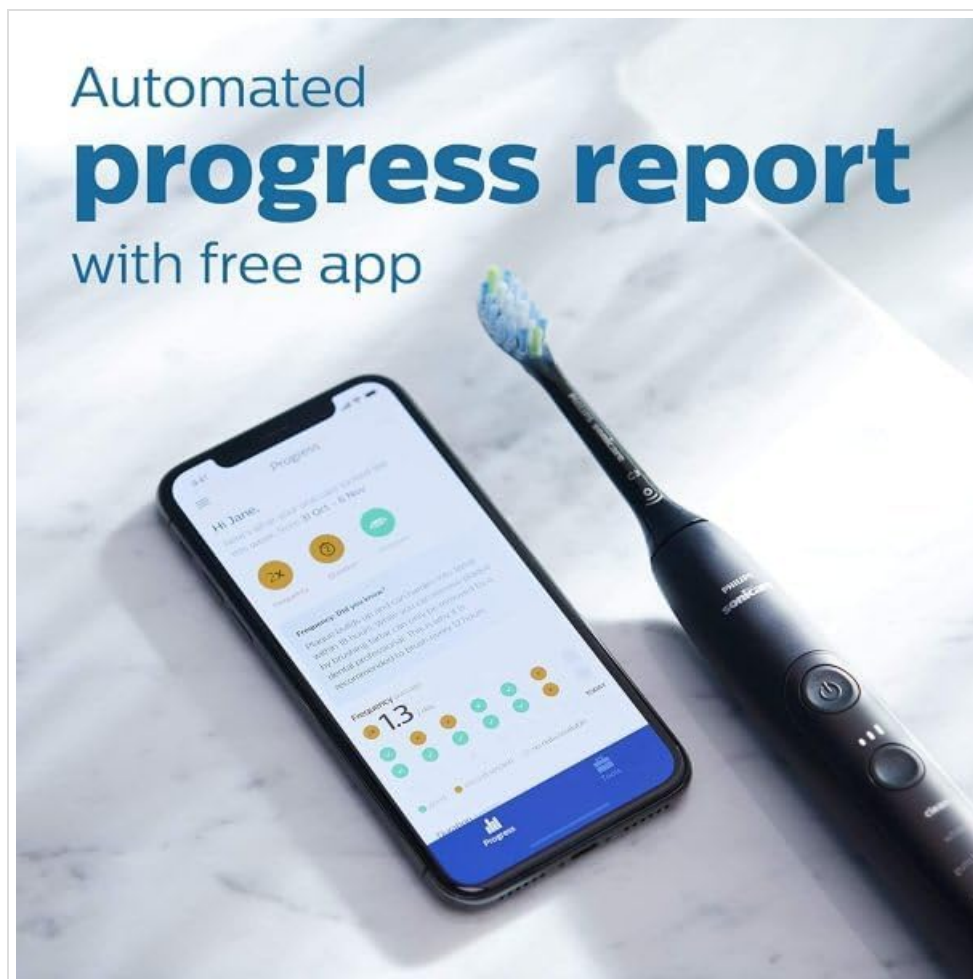


Image: The pressure sensor alerts you when you brush too hard, helping to protect your gums and teeth.

### 3. Modes and Intensities

The ExpertClean 7500 offers 4 brushing modes and 3 intensity settings to personalize your brushing experience:

- **Clean:** For an exceptional everyday clean.
- **White+:** To remove surface stains and whiten teeth.
- **Gum Health:** Provides a gentle yet effective clean for gums.
- **DeepClean+:** For an invigorating deep clean.

You can adjust the intensity (low, medium, high) by pressing the intensity button on the handle. The selected mode and intensity will be indicated by lights on the handle.

Pressure sensor alerts you  
when you brush too hard

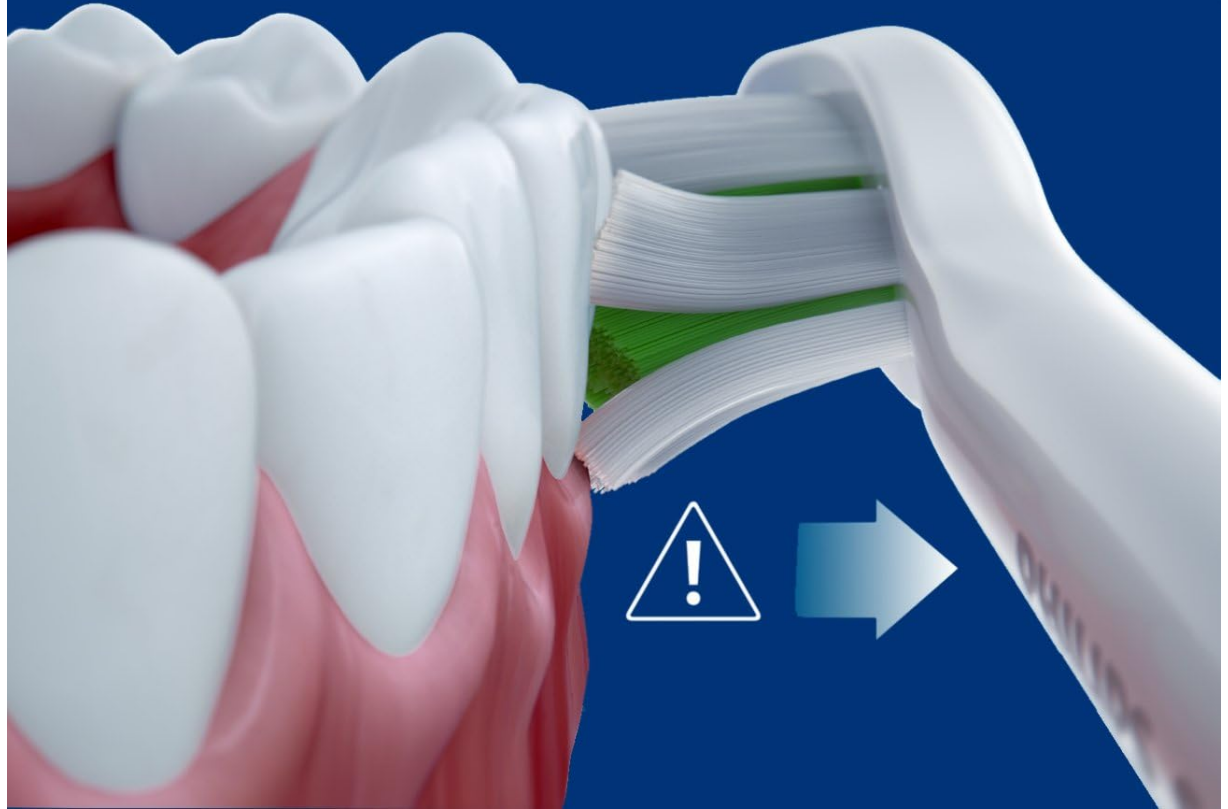


Image: Personalize your brushing with 3 intensities and 4 modes for a tailored oral care experience.

## MAINTENANCE

---

### 1. Cleaning the Toothbrush

After each use, rinse the brush head and bristles thoroughly under running water. Wipe the handle with a damp cloth. Do not immerse the handle in water or clean it in a dishwasher.

### 2. Brush Head Replacement

The Philips Sonicare ExpertClean 7500 features BrushSync Technology, which tracks your brush head usage and reminds you when it's time for a replacement. A light indicator on the handle will alert you when a new brush head is needed, typically every 3 months for optimal performance.

# BrushSync replacement reminder lets you know when to replace brush head

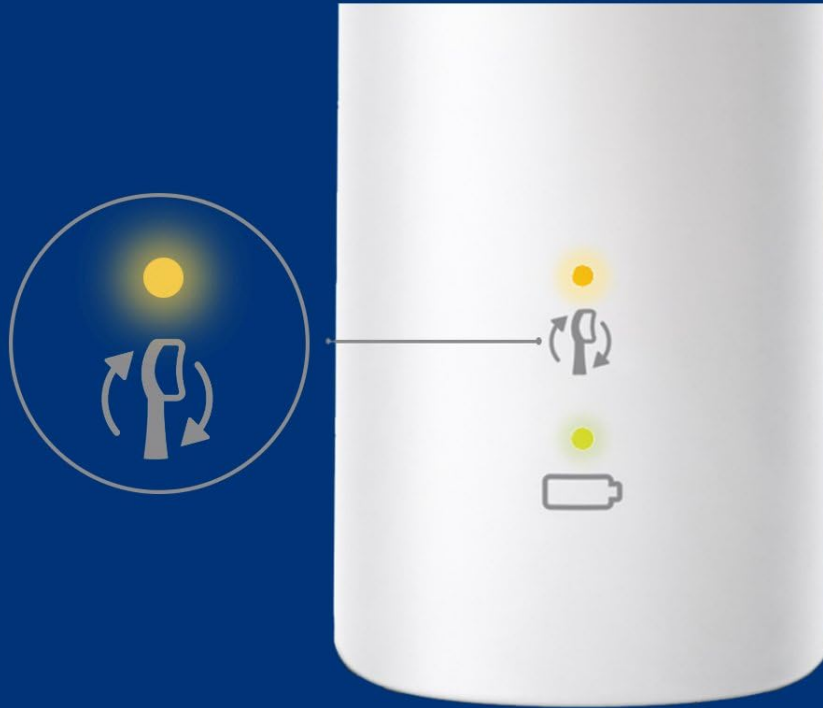


Image: The BrushSync replacement reminder light indicates when to replace your brush head for continued effective cleaning.

## TROUBLESHOOTING

---

If you encounter any issues with your Philips Sonicare ExpertClean 7500, please refer to the following common solutions:

- **Toothbrush not turning on:** Ensure the battery is charged. Place the handle on the charger for at least 24 hours.
- **Reduced brushing power:** The brush head may be worn out. Replace the brush head.
- **Charging issues:** Check that the charger is properly plugged into a working outlet. Ensure the toothbrush handle is correctly seated on the charging base.
- **Excessive vibration/noise:** Ensure the brush head is securely attached to the handle.

If problems persist, please contact Philips Sonicare customer support for further assistance.

## SPECIFICATIONS

---

Feature	Detail
Product Dimensions	3.07 x 6.81 x 9.41 inches
Item Weight	1.08 Pounds
Model Number	HX9690/06
Power Source	Corded Electric
Battery Type	1 A battery (included)
Number of Settings	4 Modes, 3 Intensities
Special Features	Brushing Timer, Customizable, Pressure Sensor
Handle Material	Plastic
Bristle Material	Nylon
Color	White

## WARRANTY AND SUPPORT

---

For warranty information, product registration, or technical support, please visit the official Philips Sonicare website or contact their customer service. Keep your proof of purchase for warranty claims.

**Important Safety Information:** Choking hazard: Keep out of reach of children under 3 years. Do not use in shower or near water to prevent electrical shock.