

[Manuals.plus](#) /

> [FIGERM](#) /

> FIGERM Foot Warmer Electric Heating Pad (Model B07TVGMM32) Instruction Manual

## FIGERM B07TVGMM32

# FIGERM Foot Warmer Electric Heating Pad (Model B07TVGMM32) Instruction Manual

## 1. PRODUCT OVERVIEW

---

This instruction manual provides essential information for the safe and effective use of your FIGERM Foot Warmer Electric Heating Pad, Model B07TVGMM32. This king-size, ultra-soft flannel heating pad is designed for versatile use on feet, legs, back, shoulders, and abdomen to provide warmth and comfort.



The FIGERM Foot Warmer Electric Heating Pad in blue, showcasing its soft flannel material and compact design.

## 2. IMPORTANT SAFETY INSTRUCTIONS

---

Please read all instructions carefully before using this heating pad. Failure to follow these instructions may result in fire, electric shock, or injury.

- Do not immerse the heating pad or its controller in water.
- Keep out of reach of children.
- Avoid direct contact with heating elements for prolonged periods to prevent burns.
- Do not use while sleeping or on an unconscious person.
- Do not use with ointments, liniments, or salves, as skin burns may result.
- Always unplug the heating pad when not in use.
- Do not use if the cord or plug is damaged.
- Do not use if the pad is wet.
- Do not puncture the pad with pins or other sharp objects.
- The heating pad features an automatic 2-hour shut-off to prevent overheating and conserve energy.

## 3. PRODUCT FEATURES

---

- **Full Body Application:** Suitable for feet, legs, back, shoulders, and abdomen.
- **Fast Heating Technology:** Provides rapid warmth.

- **Adjustable Temperature Settings:** 3 levels (Low, Medium, High).
- **Automatic Shut-Off:** Powers off after 2 hours for safety and energy saving.
- **Ultra-Soft Flannel Material:** Ensures comfort and warmth.
- **Machine Washable:** Detachable controller allows for easy cleaning.
- **Extra-Long Cord:** 10-foot power cord for convenience.
- **ETL Certified:** Ensures product safety standards are met.
- **Dimensions:** Approximately 22 inches x 20 inches.



Illustrates the versatility of the heating pad, showing it applied to feet, abdomen, and back for targeted warmth.

#### 4. SETUP

1. Unpack the heating pad and controller from its packaging.
2. Ensure the controller cable is securely plugged into the heating pad's connector.
3. Place the heating pad on a flat, stable surface or the desired body area.

# Temperatures Control Heating pad



**ETL Intertek  
Certified Safe**




**Fast Heat  
in Seconds**



**Machine  
Washable**



**10ft  
Extra-  
Long Cord**

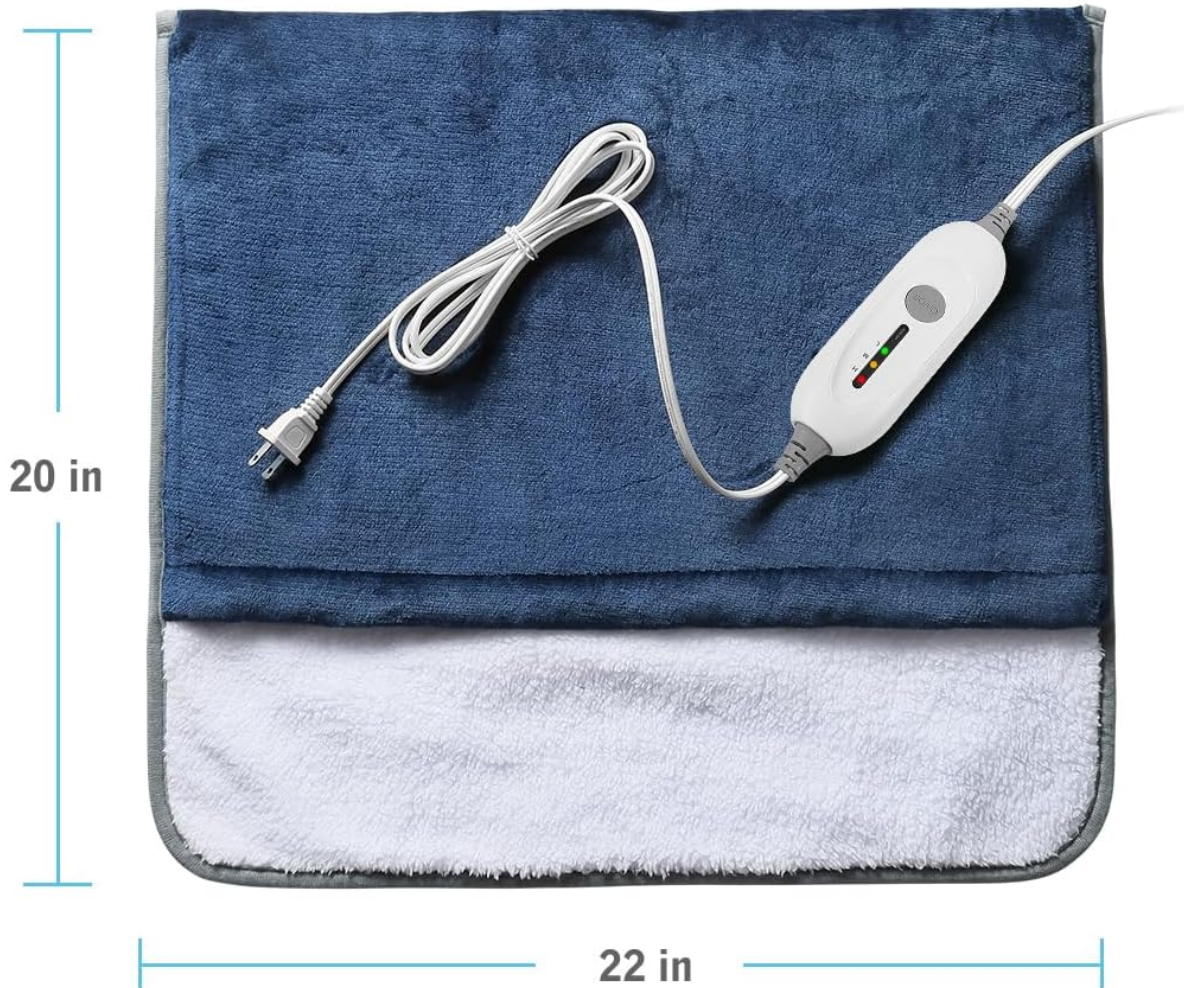

A view of the heating pad laid flat, showing the detachable controller and the 10-foot power cord. Dimensions of 22x20 inches are indicated.

## 5. OPERATING INSTRUCTIONS

1. Plug the power cord into a standard electrical outlet.
2. Press the "On/Off" button on the controller to turn on the heating pad.
3. Select your desired heat level by pressing the button:
  - Low (Green Light)
  - Medium (Orange Light)
  - High (Red Light)
4. Place your feet inside the pocket or apply the pad to the desired body area.
5. The heating pad will automatically turn off after 2 hours. To restart, press the "On/Off" button again.

# 2 Hours Auto-off

## Protect you from excessive heating



A person comfortably using the FIGERM Foot Warmer to warm their feet while seated on a couch.

## 6. CARE AND MAINTENANCE

1. Before cleaning, always unplug the heating pad from the electrical outlet and disconnect the controller from the pad.
2. The fabric portion of the heating pad is machine washable. Wash on a gentle cycle with cold water.
3. Tumble dry on low heat or air dry. Do not iron.
4. Ensure the heating pad is completely dry before reconnecting the controller and plugging it back in.
5. Store the heating pad in a cool, dry place when not in use.



An image demonstrating the machine washable feature of the heating pad, with the pad being placed into a washing machine.

## 7. TROUBLESHOOTING

Problem	Solution
Heating pad does not turn on or heat up.	Ensure the power cord is securely plugged into a working electrical outlet. Check that the controller is firmly connected to the heating pad.
Heating pad is not warm enough.	Adjust the heat setting to Medium or High using the controller. Allow a few minutes for the pad to reach the desired temperature.
Heating pad turns off unexpectedly.	This is the automatic 2-hour shut-off safety feature. To continue use, press the "On/Off" button on the controller to restart.

## 8. SPECIFICATIONS

Feature	Detail
Brand	FIGERM
Model	B07TVGMM32
Dimensions	22 inches x 20 inches
Material	Ultra Soft Flannel
Heat Settings	3 (Low, Medium, High)
Auto Shut-Off	2 hours
Cord Length	10 feet
Certification	ETL Certified
Item Weight	0.84 Kilograms (1.85 Pounds)

## 9. WARRANTY AND SUPPORT

---

The FIGERM Foot Warmer Electric Heating Pad comes with **lifetime warranty**. For any product-related issues or support, please contact customer service through the retailer or manufacturer's official channels. Please retain your proof of purchase for warranty claims.