

TUREWELL FC165, FC168, FC169, FC166, FC163, FC162, FC188, FC288

TUREWELL Water Flosser Replacement Nozzles Instruction Manual

For Models: FC165, FC168, FC169, FC166, FC163, FC162, FC188, FC288

INTRODUCTION

This manual provides comprehensive instructions for the TUREWELL 7-piece replacement jet tip set, designed for use with TUREWELL oral irrigators including models FC165, FC168, FC169, FC166, FC163, FC162, FC188, and FC288. This set includes various tips to support a complete oral hygiene routine.

IMPORTANT SAFETY INFORMATION

WARNING:

- Do not point the nozzle towards the eyes while in use.
- Children must use this product under the guidance of parents.
- Individuals with serious dental conditions should consult a dentist before use.
- These nozzles are not suitable for cordless water flossers.

PRODUCT OVERVIEW: INCLUDED NOZZLES

Your TUREWELL replacement set includes 7 specialized jet tips for various oral care needs:



Image: A visual representation of the seven distinct TUREWELL water flosser replacement nozzles included in the set. From left to right, a blue tongue cleaner, followed by two standard jet tips with blue bases, one standard jet tip with a pink base, one standard jet tip with a teal base, an orthodontic tip with a blue base, a plaque brush tip with a blue base, and a periodontal tip with a blue base.

7 NOZZLES With Different Functions



Image: A detailed diagram showcasing the seven different TUREWELL nozzles and their specific functions. It labels the Classic Jet Tip for general cleaning, the Tongue Cleaner for removing bacteria buildup, the Orthodontic Tip for removing food debris in braces, the Dental Plaque Tip for removing plaque and food debris in hard-to-reach areas, and the Periodontal Tip for rinsing and massaging deep below the gumline.

- **3 x Standard Jet Tips:** Suitable for daily dental cleaning, offering effective plaque removal and gum massage to promote blood circulation and strengthen teeth.
- **1 x Orthodontic Tip:** Designed for daily dental care for individuals with braces, orthodontic appliances, or implants. It helps clean around brackets and wires, massaging gums and promoting gum health.
- **1 x Plaque Brush Tip:** Ideal for cleaning implants, dentures, bridges, and crowns. Its novel design effectively removes plaque from hard-to-reach areas.
- **1 x Periodontal Pocket Tip:** Specifically for cleaning periodontal pockets. This fine nozzle delivers a pulsed water column 6-7mm deep into the gumline, cleaning harmful substances with low pressure and massaging gums comfortably.
- **1 x Tongue Cleaner:** Used for cleaning the tongue, removing bad breath, and clearing tongue waste. Its circular arc design conforms to the tongue's shape, and pulsed water helps clean food residue for fresh breath.

SETUP

1. **Ensure Safety:** Before attaching or detaching any nozzle, ensure your TUREWELL oral irrigator is turned off and unplugged from the power source.
2. **Select Nozzle:** Choose the appropriate nozzle for your specific oral care need.
3. **Attach Nozzle:** Firmly insert the desired nozzle into the opening on the irrigator handle. Push down until you hear a click, indicating it is securely locked in place.
4. **Remove Nozzle:** To remove a nozzle, press the nozzle release button (if applicable on your irrigator model) and gently pull the nozzle straight out from the handle.

OPERATING INSTRUCTIONS

1. **Fill Reservoir:** Fill the water flosser reservoir with lukewarm water. You may add a small amount of mouthwash if desired.
2. **Insert Nozzle:** Attach the chosen nozzle as described in the Setup section.
3. **Position Nozzle:** Lean over a sink. Place the nozzle tip into your mouth, aiming it at your gumline at a 90-degree angle. Keep your mouth slightly open to allow water to flow into the sink.
4. **Start Flossing:** Turn on the water flosser. Start with the lowest pressure setting and gradually increase it to a comfortable level.
5. **Guide Water Jet:** Guide the water jet along your gumline, pausing briefly between teeth. Ensure you clean both the front and back surfaces of your teeth.
6. **Tongue Cleaning (if applicable):** If using the tongue cleaner, set the pressure to low. Place the cleaner on the back of your tongue and move it forward, cleaning the entire surface.
7. **Finish:** Turn off the water flosser when finished. Remove the nozzle and rinse it.

USING PULSING WATER FLOWS TO CLEAN THE TEETH.

Use 1700 times per minutes frequency conversion motor, and the maximum water pressure can reach 100psi.

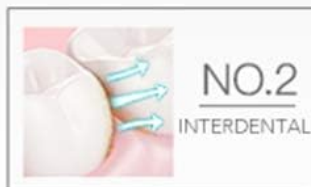
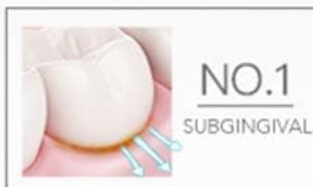


Image: An illustration demonstrating how pulsing water flows clean teeth. It highlights three key areas of cleaning: subgingival (below the gumline), interdental (between teeth), and cariosity (cavity-prone areas), showing the effectiveness of the water flosser in reaching these regions.

MAINTENANCE

- **Rinsing:** After each use, remove the nozzle and rinse it thoroughly under running water to remove any debris or residue.
- **Cleaning:** Periodically, clean the nozzles with a mild soap solution and rinse them thoroughly with clean water.
- **Replacement Schedule:**
 - Standard Jet Tips and Tongue Cleaner: Replace every six months.
 - Orthodontic Tip, Plaque Brush Tip, and Periodontal Pocket Tip: Replace every three months.
- **Storage:** Store nozzles in a clean, dry place when not in use.

TROUBLESHOOTING

- **Nozzle Does Not Fit:** Ensure you are using TUREWELL replacement nozzles specifically designed for your irrigator model (FC165, FC168, FC169, FC166, FC163, FC162, FC188, FC288). These nozzles are not compatible with cordless water flossers.
- **Low Water Pressure:** Check the water level in your irrigator's reservoir. Ensure the nozzle is securely attached to the handle and that there are no blockages in the nozzle tip.
- **Water Leakage:** Verify that the nozzle is fully inserted and clicked into place. Check for any visible damage to the nozzle or the irrigator handle's connection point.

SPECIFICATIONS

Feature	Detail
Product Dimensions	3.15 x 3.15 x 0.39 inches
Item Weight	1.41 ounces
Unit Count	7 Count
Brand	TUREWELL
Compatible Oral Irrigator Models	FC165, FC168, FC169, FC166, FC163, FC162, FC188, FC288

WARRANTY AND SUPPORT

For warranty information regarding your TUREWELL oral irrigator or replacement parts, please refer to the documentation included with your original irrigator purchase or contact TUREWELL customer support directly. For further assistance or product inquiries, please visit the official TUREWELL Store on Amazon:[TUREWELL Amazon Store](#).