



Manuals.plus /

> **Burst** /

> Burst Original Sonic Electric Toothbrush Instruction Manual (Model SG-987)

Burst SG-987

Burst Original Sonic Electric Toothbrush Instruction Manual

Model: SG-987

Brand: Burst

[Overview](#) [Setup](#) [Operation](#) [Maintenance](#) [Safety Information](#) [Product Troubleshooting](#) [Specifications](#) [Warranty](#) [Support](#)

IMPORTANT SAFETY INFORMATION

For adult use only. Always follow usage instructions for safe operation. Keep out of reach of children. Do not immerse the charging base in water. Consult a dental professional if you experience persistent discomfort or bleeding gums.

PRODUCT OVERVIEW

The Burst Original Sonic Electric Toothbrush is designed to provide a deep clean with 33,000 sonic vibrations per minute. It features ultra-soft, multi-length, interdental-cut PBT nylon bristles that are gentle on gums and enamel while effectively removing plaque and stains. The toothbrush includes a 2-minute smart timer with 30-second quadrant reminders to ensure thorough brushing.



Image: Burst Original Sonic Electric Toothbrush in Black, showcasing its sleek design and brush head.

burst[®] Original Sonic Toothbrush

Clinically proven* results

- Removes 10x more plaque than a manual toothbrush
- Reduce gum bleeding up to 3x more than a manual toothbrush
- Proven healthier gums with twice daily brushing

*Compared to a manual toothbrush, brushing twice daily in whitening mode for two minutes for a period of 30 days. Data on file at <https://www.burstoralcare.com/clinical>.

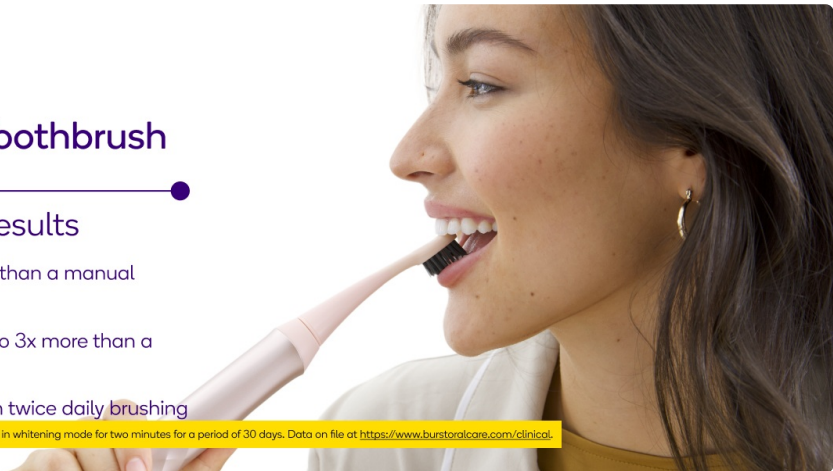


Image: Infographic highlighting clinically proven results, including 10x more plaque removal and reduced gum bleeding.

SETUP

1. Initial Charging

Before first use, fully charge your Burst Sonic Electric Toothbrush. The toothbrush features a long-

lasting battery, providing up to one month of use on a single overnight charge. Use the included USB charger and wireless charging base.

Power That Lasts

- Up to 4 weeks on a single overnight charge using the included wireless charging base.
- (USB cord included; wall adapter not included.)

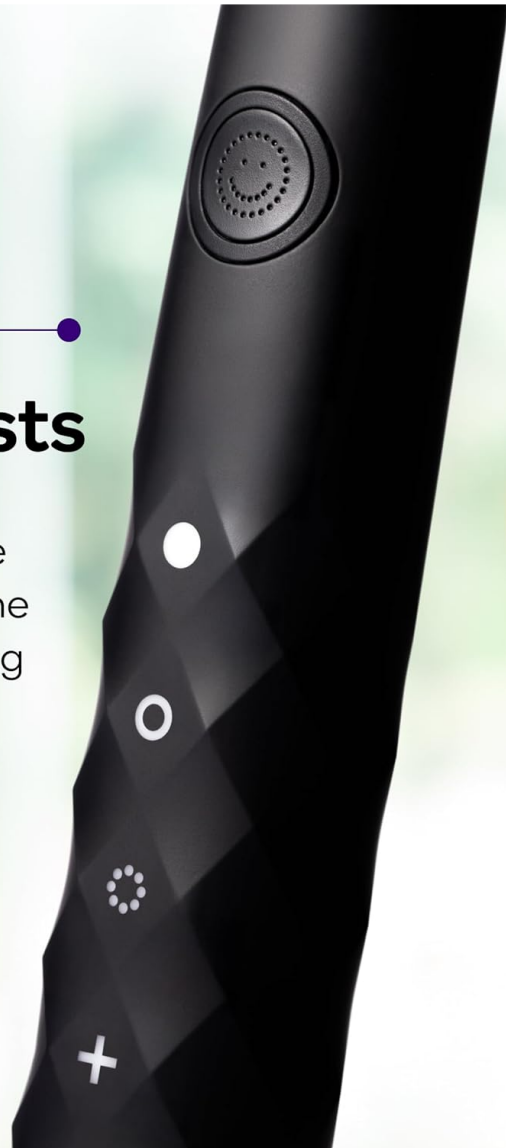


Image: The Burst Sonic Electric Toothbrush handle placed on its wireless charging base, connected via USB.

2. Attaching the Brush Head

Align the brush head with the metal shaft on the toothbrush handle. Gently push the brush head down until it clicks securely into place. Ensure it is firmly attached before use.

OPERATING INSTRUCTIONS

1. Selecting a Brushing Mode

The Burst toothbrush offers three brushing modes: Whitening, Sensitive, and Massage. Press the power button once to turn on the toothbrush. Press it again within a few seconds to cycle through the modes until you reach your desired setting. The indicator lights on the handle will show the active mode.

2. Brushing Technique

Apply a small amount of toothpaste to the brush head. Place the brush at a 45-degree angle to your gumline. Gently move the brush slowly from tooth to tooth, allowing the sonic vibrations to do the

cleaning. Do not scrub vigorously.

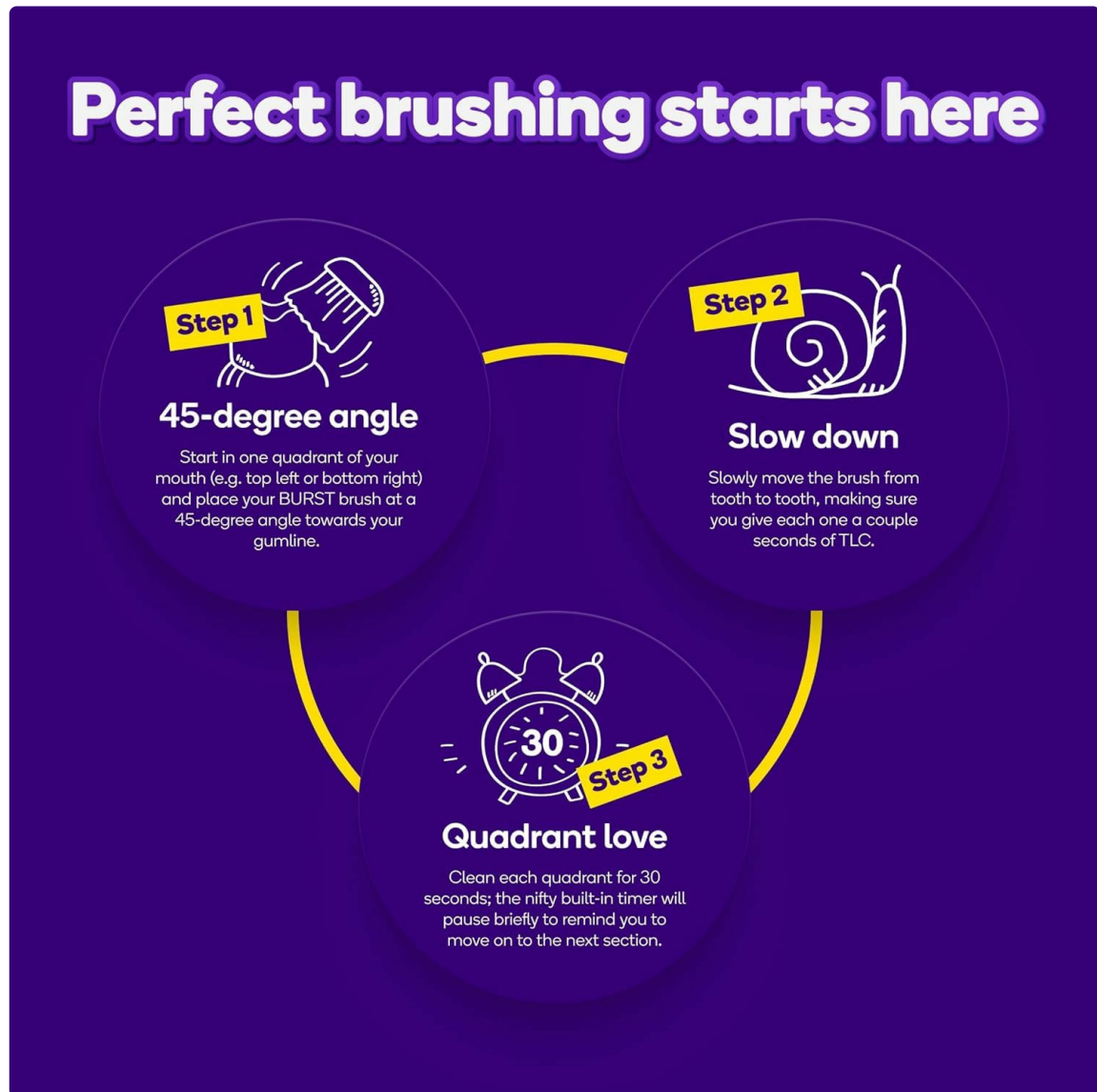


Image: Visual guide demonstrating the recommended 45-degree brushing angle, slow movement, and 30-second quadrant timing.

3. Smart Timer and Quad Pacer

The toothbrush features a 2-minute smart timer, which is the dentist-recommended brushing duration. It will pulse every 30 seconds to remind you to move to a different quadrant of your mouth (upper right, upper left, lower left, lower right), ensuring an even and thorough clean.

Your browser does not support the video tag.

Video: Official product video demonstrating the Burst Sonic Electric Toothbrush in use, highlighting its features and brushing action.

MAINTENANCE

1. Replacing Brush Heads

Dental professionals recommend replacing your brush head every 3 months, or sooner if the bristles appear frayed. To replace, firmly pull the old brush head straight off the handle. Align the new brush head and push it down until it clicks into place.

BURST Signature Soft Bristles

Loved by dentists and dental hygienists, these multi-length, tapered bristles reach deep & whiten as you brush



Image: A close-up view of the Burst Signature Soft Bristles, emphasizing their multi-length and tapered design.

2. Cleaning the Toothbrush

After each use, rinse the brush head thoroughly under running water. Wipe the handle with a damp cloth. Do not use abrasive cleaners or immerse the handle in water for extended periods.

TROUBLESHOOTING

Toothbrush Not Turning On or Low Power

- Ensure the toothbrush is fully charged. Place it on the charging base and confirm the charging indicator light is active.
- Check the charging cable and adapter for any damage.
- If the issue persists after charging, try a different power outlet.

Brush Head Feels Loose

- Remove the brush head and reattach it, ensuring it clicks firmly into place.
- If the brush head remains loose, it may be worn and require replacement.

For any other issues not covered here, please refer to the Customer Support section.

SPECIFICATIONS

- **Product Dimensions:** 12 x 3 x 10 inches; 8.11 ounces
- **Item Model Number:** SG-987
- **Batteries:** 1 Lithium Ion batteries required (included)
- **UPC:** 860203000708
- **Manufacturer:** Burst Oral Care
- **Brand:** Burst
- **Age Range (Description):** Adult
- **Special Feature:** Brushing Timer
- **Power Source:** Battery Powered
- **Item Firmness Description:** Soft, Ultra Soft

WARRANTY

BURST offers a one-year warranty on any BURST device. Additionally, a limited lifetime warranty is provided to customers enrolled in the device's replacement subscription program. For detailed terms and conditions, please refer to the official Burst Oral Care website or contact customer support.

CUSTOMER SUPPORT

For further assistance, product inquiries, or warranty claims, please visit the official **Burst Store on Amazon** or contact Burst Oral Care directly through their official channels.