

Medisana 60258

Medisana FWS Foot Warmer Model 60258 Instruction Manual

Your guide to safe and effective use

1. INTRODUCTION

Thank you for choosing the Medisana FWS Foot Warmer, Model 60258. This instruction manual provides important information regarding the safe operation, setup, maintenance, and troubleshooting of your device. Please read this manual thoroughly before first use and keep it for future reference.

2. SAFETY INSTRUCTIONS

Please read and understand all safety warnings before use.

- **Electrical Safety:** Ensure the voltage indicated on the appliance corresponds to your local mains voltage before connecting. Do not use if the power cord or plug is damaged.
- **Water and Moisture:** Do not use the foot warmer in wet conditions or immerse it in water. Keep the control unit dry.
- **Children and Pets:** This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Supervise children to ensure they do not play with the appliance.
- **Proper Use:** Do not fold, crease, or bunch the foot warmer during operation. Do not use while sleeping. Do not insert pins or sharp objects into the foot warmer.
- **Overheating Protection:** The device is equipped with an automatic shut-off function after 90 minutes and overheat protection. If the device overheats, it will automatically switch off.
- **Maintenance:** Always disconnect the power plug from the socket before cleaning or when the device is not in use.
- **Damage Check:** Regularly inspect the foot warmer for any signs of wear or damage. If damage is detected, discontinue use and contact customer support.

3. PRODUCT OVERVIEW

Components:

- Medisana FWS Foot Warmer Unit
- Detachable Control Unit with Power Cord

Key Features:

- **Adjustable Heat:** 6 temperature levels for personalized comfort.
- **Timer Function:** 6 timer settings for automatic shut-off.
- **Fast Heating:** 100 watts of power for quick warmth.
- **Washable Liner:** Removable, fluffy inner lining is machine washable at 30°C.
- **Safety Features:** Automatic shut-off after 90 minutes and integrated overheat protection.
- **Spacious Design:** Accommodates feet up to size 46.



The Medisana FWS Foot Warmer, shown in gray, features a soft, plush exterior and a detachable control unit for easy operation.

Ergonomically designed controller with LED control display

Automatic switch-off after selected timer

Safety system with overheating protection



This image details key features: 100 watts for quick heating, a spacious interior for feet, and the control unit with options for 6 temperature and 6 timer settings.

4. SETUP

1. **Unpack:** Carefully remove the foot warmer and all components from the packaging. Retain packaging for storage or disposal.
2. **Inspect:** Check the foot warmer and control unit for any visible damage. Do not use if damaged.
3. **Connect Control Unit:** Ensure the control unit is securely plugged into the foot warmer unit. The connector should click into place.
4. **Position:** Place the foot warmer on a flat, stable surface, such as the floor, where you intend to use it. Ensure there is adequate space around the unit for air circulation.
5. **Power Connection:** Plug the power cord of the control unit into a suitable electrical outlet.

23 cm
9.1 in



30 cm
11.8 in



30 cm
11.8 in



This image illustrates the dimensions of the foot warmer, measuring approximately 23 cm (9.1 inches) in height and 30 cm (11.8 inches) in width, providing ample space for feet.


5. OPERATION

The Medisana FWS Foot Warmer is designed for ease of use with its intuitive control unit.

1. **Power On:** Press the power button () on the control unit to turn on the foot warmer. The LED display

will illuminate.

2. **Adjust Temperature:** Press the temperature button (⊛) repeatedly to cycle through the 6 available heat settings. The selected level will be indicated on the LED display. Start with a lower setting and increase as desired.
3. **Set Timer:** Press the timer button (🕒) to select one of the 6 available timer settings. The foot warmer will automatically switch off after the set time has elapsed. If no timer is set, the device will automatically shut off after 90 minutes.
4. **Enjoy:** Place your feet inside the foot warmer.
5. **Power Off:** To turn off the device manually, press the power button (🔌) again. Always unplug the device from the mains when not in use.

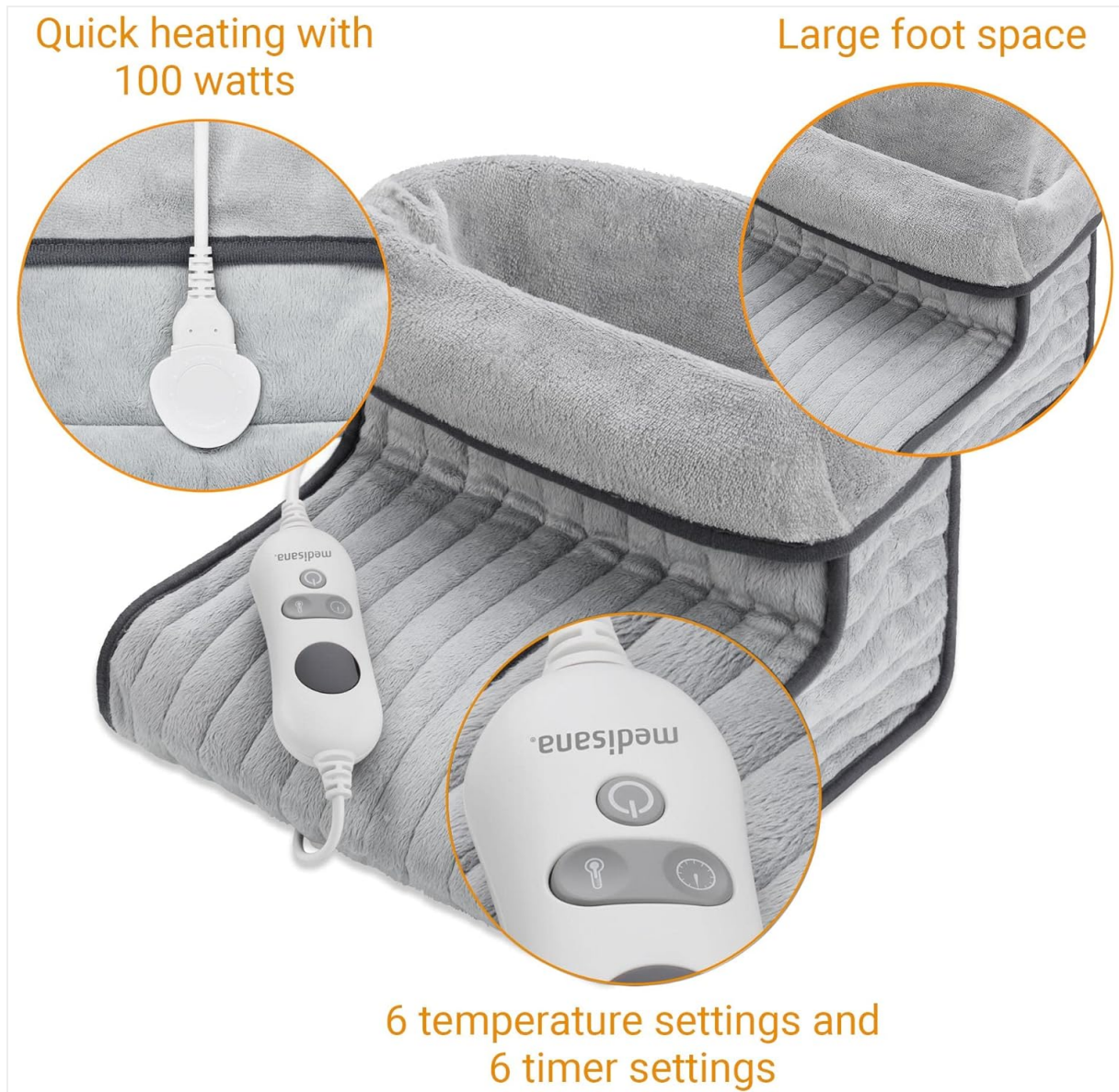

The ergonomic control unit features an LED display, a power button, a timer button, and a temperature adjustment button. It also includes a safety system with overheating protection and automatic switch-off.



A user is shown with their feet inside the Medisana FWS Foot Warmer, demonstrating its comfortable and practical use in a home setting.


6. CLEANING AND MAINTENANCE

Regular cleaning ensures hygiene and extends the lifespan of your foot warmer.

1. **Disconnect Power:** Always unplug the foot warmer from the electrical outlet before cleaning.
2. **Remove Inner Lining:** The soft, fluffy inner lining can be easily removed from the main unit.
3. **Wash Inner Lining:** The inner lining is machine washable at 30°C (86°F) on a gentle cycle. Do not bleach, iron, or dry clean. Allow it to air dry completely before reinserting.
4. **Clean Outer Unit:** Wipe the outer surface of the foot warmer and the control unit with a soft, damp cloth. Do not use abrasive cleaners or immerse the control unit in water.
5. **Reassemble:** Once all parts are completely dry, reinsert the inner lining into the foot warmer unit.

Removable fluffy soft lining

Inner lining machine washable at 30°C

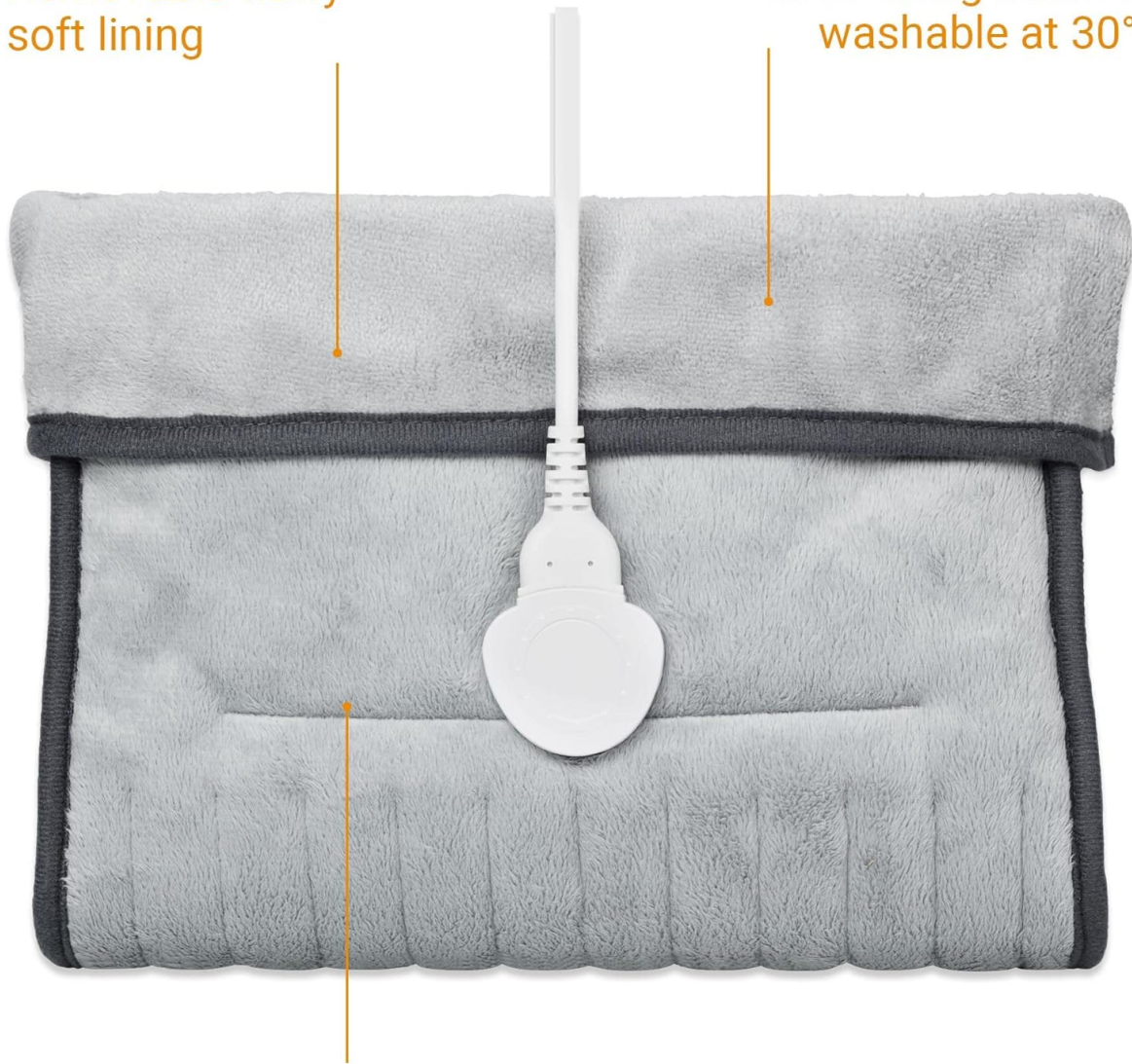

Extra-soft breathable core

The foot warmer's design includes a removable, fluffy soft inner lining that is machine washable at 30°C, and an extra-soft breathable core for comfort.

7. STORAGE

When not in use, store the Medisana FWS Foot Warmer in a clean, dry place. Ensure the device is completely cool and dry before storing. Avoid placing heavy objects on top of the foot warmer to prevent damage.

8. TROUBLESHOOTING

If you encounter any issues with your Medisana FWS Foot Warmer, please refer to the following common problems and solutions:

Problem	Possible Cause	Solution
Foot warmer does not heat up.	Not plugged in correctly, power outage, control unit not connected, device damaged.	Check power plug and outlet. Ensure control unit is securely connected. Check household fuse. If problem persists, contact customer support.

Problem	Possible Cause	Solution
Foot warmer switches off unexpectedly.	Automatic shut-off activated (90 minutes or timer setting), overheat protection activated.	This is normal operation. Wait a few minutes for the device to cool down before attempting to restart. If it continues to shut off prematurely, contact customer support.
LED display not working.	Power issue, control unit malfunction.	Check power connection. If the device is heating but the display is off, or if neither works, contact customer support.

If the problem cannot be resolved using the troubleshooting guide, do not attempt to repair the device yourself. Contact Medisana customer support for assistance.

9. SPECIFICATIONS

- **Model:** FWS (60258)
- **Brand:** Medisana
- **Power:** 100 W
- **Temperature Levels:** 6
- **Timer Settings:** 6
- **Automatic Shut-off:** After 90 minutes (default)
- **Upper Temperature Rating:** 45°C
- **Dimensions (Product):** Approximately 12.6 x 33.1 x 12.5 cm
- **Weight:** Approximately 0.68 kg
- **Inner Lining Washability:** Machine washable at 30°C
- **Maximum Foot Size:** Up to EU size 46

10. WARRANTY AND SUPPORT

For warranty information, please refer to the warranty card included with your product or visit the official Medisana website. If you require technical support or have questions not covered in this manual, please contact Medisana customer service through their official channels.