

## Oral-B Vitality FlossAction

# Oral-B Vitality FlossAction Electric Toothbrush Instruction Manual

## INTRODUCTION

---

The Oral-B Vitality FlossAction Electric Toothbrush is designed to provide superior plaque removal compared to a manual toothbrush. Its unique FlossAction brush head features micropulse bristles that clean deep between teeth. This manual provides essential information for the proper setup, operation, and maintenance of your new electric toothbrush.



Oral-B Vitality FlossAction Electric Toothbrush in its retail packaging.

## 1. INITIAL SETUP

### Charging the Toothbrush

Before using your Oral-B Vitality electric toothbrush for the first time, it is essential to charge it fully. This ensures optimal performance and battery life.

1. Place the toothbrush handle onto the charging base.
2. Plug the charging base into a suitable electrical outlet.
3. Allow the toothbrush to charge for the recommended duration as indicated in the full product manual (typically 16-24 hours for the first charge).



The Oral-B Vitality FlossAction Electric Toothbrush positioned on its charging base.

This video demonstrates the Oral-B Vitality Electric Toothbrush, showing its operation, the rotating brush head under water, and how it sits on its charging base.

A demonstration of the Oral-B Vitality Electric Toothbrush, highlighting its simple on/off operation and placement on the charging stand.

## 2. OPERATING INSTRUCTIONS

---

### Attaching the Brush Head

Ensure the brush head is securely attached to the handle before use.

1. Align the brush head with the metal shaft on the toothbrush handle.
2. Push the brush head down firmly until it clicks into place.

### Brushing Your Teeth

Follow these steps for an effective brushing experience:

1. Wet the brush head and apply a small amount of toothpaste.


2. Place the brush head against your teeth, holding it at a slight angle towards the gum line.
3. Press the power button to turn on the toothbrush.
4. Gently guide the brush head from tooth to tooth, allowing the rotating bristles to clean each surface. Do not scrub aggressively.
5. The built-in 2-minute timer will signal with a brief pause or vibration when you have brushed for the dentist-recommended duration.
6. After brushing, rinse your mouth and the brush head thoroughly.



Close-up view of the Oral-B FlossAction brush head, showcasing its micropulse bristles for effective cleaning.



# ROUND BRUSH HEAD CLEANS BETTER\*

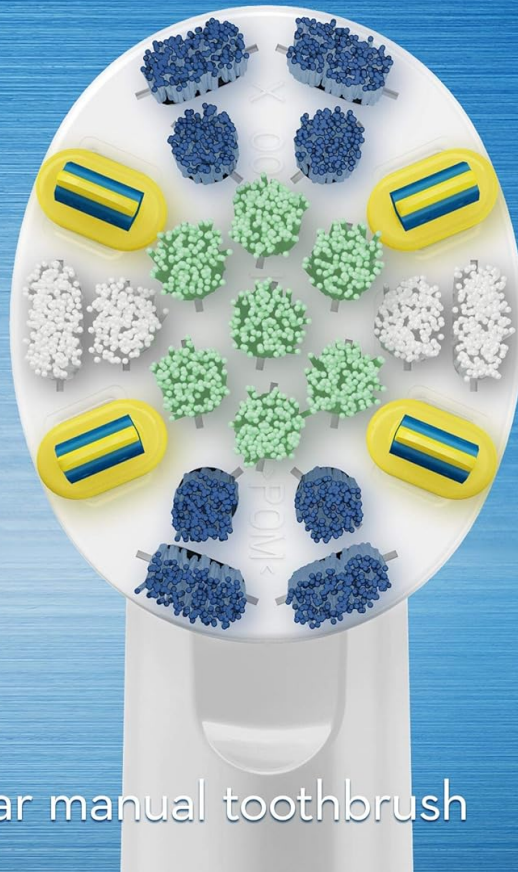

\*vs a regular manual toothbrush

Visual representation of the built-in 2-minute timer, indicating the recommended brushing duration.

## 3. MAINTENANCE

### Replacing the Brush Head

For optimal hygiene and cleaning effectiveness, replace your brush head regularly.

- Dentists recommend replacing your toothbrush head every 3 months, or sooner if the bristles become frayed or worn.
- To remove, simply pull the brush head straight off the handle.
- To attach a new brush head, align it with the shaft and push down until it clicks.

# COMPATIBLE WITH A VARIETY OF ORAL-B BRUSH HEADS\*



\*all brush heads sold separately

The Oral-B Vitality handle is compatible with a variety of Oral-B replacement brush heads, allowing for personalized oral care.

## Cleaning the Toothbrush

Regular cleaning helps maintain hygiene and extends the life of your toothbrush.

- After each use, rinse the brush head thoroughly under running water while the handle is still on.
- Remove the brush head and rinse both the brush head and the handle separately under running water.
- Wipe the handle dry with a clean cloth.
- Store the toothbrush in an upright position to allow it to air dry.

## 4. TROUBLESHOOTING

If you encounter any issues with your Oral-B Vitality electric toothbrush, refer to the following common troubleshooting tips:

- **Toothbrush not turning on:** Ensure the toothbrush is fully charged. Place it on the charger for at least 30 minutes and try again. Verify the power outlet is functional.
- **Reduced brushing power:** The battery may be low; place the toothbrush on the charger. If the issue

persists after charging, the brush head might be worn and require replacement.

- **Brush head not rotating effectively:** Check if the brush head is securely attached to the handle. If it is, and the power is on, consider replacing the brush head.
- **Unusual noise or vibration:** Ensure the brush head is correctly attached. If the noise persists, the internal mechanism might be affected. Contact customer support if necessary.

## 5. PRODUCT SPECIFICATIONS

Feature	Detail
Brand	Oral-B
Model	Vitality FlossAction
Power Source	Battery Powered (Rechargeable)
Brush Head Type	FlossAction (Micropulse bristles)
Timer	Built-in 2-minute timer
Compatibility	Compatible with CrossAction, 3D White, Sensitive Clean, Precision Clean, FlossAction, Deep Sweep, Ortho, Power Tip, and Dual Clean Oral-B replacement brush heads.
Color	Black
Product Dimensions	4.53 x 1.85 x 9.17 inches; 8.8 ounces
UPC	069055128275
ASIN	B07Y45NJ83

## 6. WARRANTY AND SUPPORT

### Warranty Information

The Oral-B Vitality FlossAction Electric Toothbrush typically comes with a limited warranty. For specific details regarding your product's warranty coverage, including duration and terms, please refer to the documentation included in your product packaging or visit the official Oral-B website.



# BUY WITH CONFIDENCE



Includes replacement or repair. Proof of date of purchase is required. US consumer is responsible for product return costs. (Full details in pack.)



A 2-year limited warranty is typically offered for the Oral-B Vitality Electric Toothbrush. Proof of purchase is required for replacement or repair.

## Customer Support

For further assistance, product registration, or to inquire about replacement parts, please contact Oral-B customer support through their official website or the contact information provided in your product manual.