

COMFIER CF-2307A-US

COMFIER Shiatsu Neck & Back Massager Instruction Manual

Model: CF-2307A-US | Brand: COMFIER

INTRODUCTION & OVERVIEW

The COMFIER Massager Chair with Heat is a portable shiatsu neck and back massager pad designed to provide a comprehensive massage experience. It combines multiple massage techniques and features to help you relax and alleviate muscle tension.

- **Full Body Massage:** Combines Shiatsu, Kneading, Rolling, Compression, Vibration, and Heat functions for complete fatigue and stress elimination.
- **Shiatsu Neck & Full Back Massager:** Features 4 Shiatsu nodes for neck and shoulder massage with adjustable positions, and 4 Shiatsu massage nodes for deep tissue kneading across the entire back.
- **Optional Heat, Rolling & Spot Massage:** Provides soothing warmth to loosen stiff muscles. Rolling massage gently glides along the spine. Spot massage allows targeting specific areas.
- **Adjustable Air Compression:** Offers 3 adjustable intensity levels on the waist and hips to relieve tension.
- **Ultimate Comfort & Versatility:** Designed to be strapped to various chairs, sofas, recliners, or office chairs for use at home or in the office.



Figure 1: The COMFIER Shiatsu Neck & Back Massager with its remote control and packaging.

SETUP GUIDE

Setting up your COMFIER massager is quick and easy. Follow these steps to begin your relaxation experience:

1. **Unpack:** Carefully remove the massager and all accessories from the packaging.
2. **Placement:** Place the massager on your desired chair, sofa, or recliner. Ensure it is securely positioned. The massager includes straps on the back to help secure it to a chair.
3. **Connect Power:** Connect the power adapter to the massager's input port, then plug the adapter into a standard electrical outlet.
4. **Remote Control:** The remote control is tethered to the massager and ready for use.



Figure 2: Key components including the removable neck cover, detachable back flap, and UL approved home adapter.

Unboxing Video

Your browser does not support the video tag.

Video 1: Detailed unboxing and initial setup of the COMFIER Shiatsu Neck and Back Massager.

OPERATING INSTRUCTIONS

Your massager is controlled via a user-friendly remote. Familiarize yourself with its functions to customize your massage experience.

Remote Control Overview

CONTROLLER

Vibration Button:

Press to turn on the seat Vibration, you may cycle through High, Medium, Low or off settings

Back Shiatsu Button:

press the direction button to switch the rotation

Spot Button:

Press to keep shiatsu function in one area.
Press **arrow** buttons to adjust the shiatsu nodes position up or down.

Compress Intensity Button:

Press to set both Compress massage intensity at Low, Medium or High level

Compress Back Button:

Press to turn on or off Compress massage on the back

Shiatsu Area Button:

Select Full, Upper or Lower back to designate the area covered by the shiatsu massage.

Auto Button:

Quick demonstration of the different functions. The L.E.D. lights will blink to indicate the active functions. Once complete, the device will automatically turn off.

Neck Shiatsu Button:

Press to turn on Shiatsu neck massage, Press again to shift rotation direction, Press once again to turn off
Press **arrow** buttons to adjust the shiatsu nodes position up or down

Rolling Button:


press to activate Rolling function, the massage nodes travel up & down across the muscle gently. Continue pressing the button to adjust the width of the massage nodes to custom fit the body.

Compress Seat Button:

Press to turn on or off Compress massage on the seat

Back Heat Button:

Press to turn on or off Heat
The heat can't be activated if back shiatsu function is not turned on

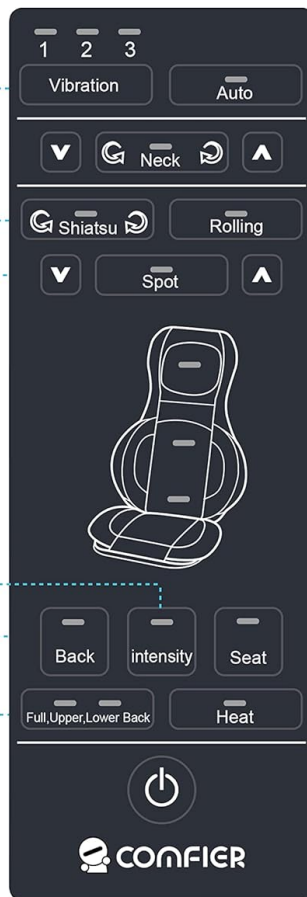

Press to turn on or off the device

Figure 3: Detailed diagram of the remote control and its various functions.

- **Power Button:** Press to turn the device on or off.
- **Vibration Button (1, 2, 3):** Cycles through High, Medium, Low, or Off settings for seat vibration.
- **Auto Button:** Initiates a quick demonstration of different functions. LED lights indicate active functions. The device automatically turns off upon completion.
- **Neck Shiatsu Button:** Press to turn on Shiatsu neck massage. Press again to switch rotation direction. Press a third time to turn off.
- **Neck Adjustment Arrows (Up/Down):** Adjusts the position of the Shiatsu nodes for neck massage.
- **Back Shiatsu Button:** Press to turn on Shiatsu back massage. Press again to switch rotation direction.
- **Rolling Button:** Activates the Rolling function, where massage nodes travel up and down the back. Pressing the button adjusts the width of the massage nodes.
- **Spot Button:** Keeps the Shiatsu function focused in one area. Use arrow buttons to adjust the Shiatsu nodes' position up or down.
- **Air Back Button:** Turns on or off compression massage on the back.
- **Air Intensity Button:** Sets compression massage intensity at Low, Medium, or High level.
- **Air Seat Button:** Turns on or off compression massage on the seat.
- **Full/Upper/Lower Back Button:** Selects Full, Upper, or Lower back to designate the area covered by the Shiatsu massage.

- **Heat Button:** Turns on or off the heat function. Heat cannot be activated if Shiatsu function is not turned on.

Massage Modes & Features

The massager offers a variety of modes to target different areas and provide diverse massage sensations.

- **Deep Kneading Shiatsu:** Targets neck, shoulders, back, waist, and thighs to relax tense muscles.



Figure 4: Deep kneading Shiatsu massage targets various body parts for muscle relaxation.

- **Pressure Massage:** Adjustable air compression helps relax waist and thigh muscles.



PRESSURE MASSAGE

Relax your waist and thighs muscles

Figure 5: Pressure massage provides relief to the waist and thigh areas.

- **Soothing Heat:** Provides warmth to the full back with a heating pad, aiding in muscle relaxation.



Figure 6: The soothing heat function helps to relax muscles in the back.

Product Demonstration Video

Your browser does not support the video tag.

Video 2: Official demonstration of the COMFIER Shiatsu Neck Back Massager's features and functions.

PRODUCT FEATURES

- **Adjustable Neck Rollers:** Suitable for various heights. For users under 5'1", a seat cushion or pillow can be used on the sitting part to ensure the neck rollers reach the desired area.

ADJUSTABLE NECK MASSAGE

Is suitable for Women/ Men of different Heights



Figure 7: The neck massage nodes are adjustable to accommodate users of different heights.

- **Versatile Use:** Can be easily strapped to various types of seating, including armchairs, office chairs, sofas, and recliners, providing comfort anywhere.



Figure 8: The massager is versatile and can be used on various types of chairs for ultimate comfort.

MAINTENANCE

Proper maintenance ensures the longevity and hygiene of your massager.

- **Cleaning:** The neck cover is removable and washable for easy cleaning.
- **Detachable Flap:** The massager includes a detachable flap that can be removed for a more intense massage or kept on for a softer experience.



Figure 9: The massager features a removable and washable neck cover and a detachable back flap for customized massage intensity.

TROUBLESHOOTING

If you encounter any issues with your massager, please refer to the following common solutions:

- **Device Not Turning On:** Ensure the power adapter is securely plugged into both the massager and a working electrical outlet. Check if the power button on the remote is pressed.
- **Massage Intensity Too Strong/Weak:** Adjust the intensity settings using the remote control. For Shiatsu, consider using or removing the detachable back flap.
- **Neck Rollers Not Reaching:** If the neck rollers do not align with your neck, try placing a seat cushion or pillow on the sitting part of the massager to adjust your height.
- **Overuse Protection:** For safe use, avoid continuous massage exceeding 15 minutes per session and no more than 3 times per day to prevent muscle strain. The device may automatically shut off after 15 minutes.

SPECIFICATIONS

Feature	Detail
Product Dimensions	30 x 17 x 18 inches
Item Weight	16 pounds
Manufacturer	COMFIER
Model Number	CF-2307A-US
Material	Leather
Power Source	Corded Electric
Use For	Shoulder, Back, Neck, Body
Date First Available	November 15, 2019



Figure 10: Key dimensions of the massager for proper fit and placement.

WARRANTY & SUPPORT

COMFIER is committed to providing high-quality products and customer satisfaction.

- **Return Policy:** If this chair massage pad does not meet your expectations for any reason, you may return it for a refund within 30 days of purchase.
- **Shipping for Returns:** COMFIER covers shipping only for quality-related returns.
- **Contact Support:** For further assistance, troubleshooting, or warranty claims, please visit the official COMFIER website or contact their customer service directly.

