

COMFIER B082CXMSVR

COMFIER Cordless Foot and Leg Massager User Manual

Model: B082CXMSVR

IMPORTANT SAFETY INFORMATION

Please read all instructions carefully before using the COMFIER Cordless Foot and Leg Massager. Retain this manual for future reference.

- Always follow the instructions provided with the massager.
- Keep the massager away from water and heat sources.
- Do not use the massager while sleeping or taking a bath.
- Always unplug the massager from the power source when not in use or during charging.
- Do not use the massager if you have any medical conditions that could be affected by heat or massage. Consult a physician before use if you have pre-existing conditions.
- Do not attempt to disassemble or repair the unit yourself. Contact qualified service personnel.
- Keep out of reach of children.

PRODUCT OVERVIEW

The COMFIER Cordless Foot and Leg Massager with Heat is designed to provide a relaxing and therapeutic massage experience for your feet, ankles, and calves. Its portable and rechargeable design allows for use anywhere, anytime.

Key Features:

- **Wireless & Portable:** Powered by a rechargeable lithium-ion battery, offering 80-90 minutes of use per full charge.
- **Customized Massage Settings:** Features 3 pre-programmed massage modes and 3 intensity levels to thoroughly massage from calf to ankles and toes.
- **Soothing Heat Function:** Integrated heat on the sole for foot toe and heel, with 2 heating settings and an over-heat protection system.
- **Flexible & Adjustable Size:** Ergonomic design with adjustable Velcro straps to fit various leg and calf sizes, up to a maximum calf size of 28 inches.

Components:



The main unit of the COMFIER Cordless Foot and Leg Massager, showing its full length and integrated control panel.



Contents of the package: the massager unit, USB charging cable, user manual, and packing box.



Diagram illustrating the various massage areas and benefits, including calf massage, instep relaxation, muscle tension relief, ankle shiatsu, and heat on the sole of the foot.

SETUP

1. **Unpack:** Carefully remove all components from the packaging. Ensure the massager unit, USB charging cable, and user manual are present.
2. **Initial Charge:** Before first use, fully charge the massager. Connect the USB charging cable to the massager's charging port and plug the other end into a compatible USB power source (e.g., laptop, power bank, or phone adapter). The charging indicator light will show the charging status.



Illustration of various USB charging options for the massager, including laptop, power bank, and phone adapter.

3. **Adjust Wraps:** The massager features adjustable Velcro straps. Before use, adjust the wraps to comfortably fit your feet and legs. Ensure the fit is snug but not overly tight to allow for proper compression and circulation. The maximum calf size accommodated is 28 inches.

ADJUSTABLE WRAPS TO FIT ANY SIZE

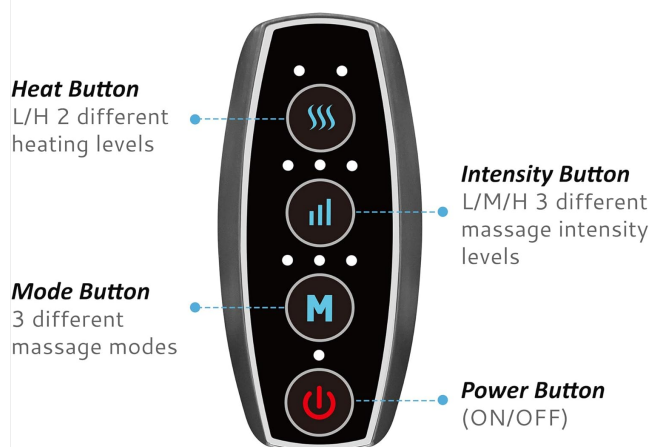


Diagram showing the adjustable wraps with Velcro and the overall dimensions of the massager, indicating a maximum calf circumference of 28 inches.

OPERATING INSTRUCTIONS


The control panel is located on the side of the massager. Familiarize yourself with the buttons before use.

CONTROLLER



Detailed view of the massager's control panel, highlighting the Heat Button, Intensity Button, Mode Button, and Power Button.

1. **Power On/Off:** Press the **Power Button** (red circle with power symbol) to turn the massager on or off.
2. **Select Massage Mode:** Press the **Mode Button** (blue 'M' circle) to cycle through the 3 different pre-programmed massage modes. Each press will change the mode, indicated by corresponding lights.
3. **Adjust Massage Intensity:** Press the **Intensity Button** (three vertical bars) to select from 3 different massage intensity levels (Low, Medium, High). The number of illuminated bars will indicate the current intensity.

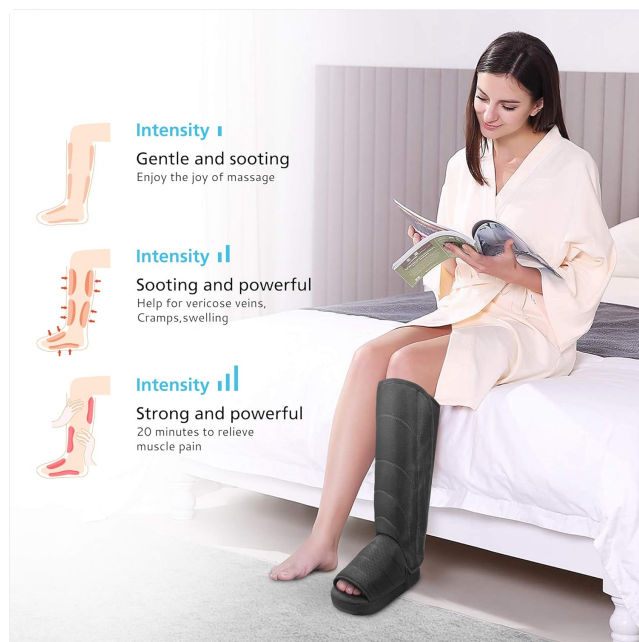

Visual representation of the three massage intensity levels: Gentle and soothing, Soothing and powerful, and Strong and powerful.

4. **Activate Heat Function:** Press the **Heat Button** (three wavy lines) to turn on the heat function. Press again to cycle between 2 heating levels (Low, High) or to turn off the heat.
5. **Enjoy Your Massage:** Relax and let the massager work. The device is designed for comfortable use while sitting.

MAINTENANCE AND CARE

- **Cleaning:** Always unplug the massager before cleaning. Wipe the surface of the massager with a soft, damp cloth. Do not use abrasive cleaners, solvents, or harsh chemicals. Do not immerse the unit in water.
- **Storage:** Store the massager in a cool, dry place away from direct sunlight and extreme temperatures. Avoid placing heavy objects on top of the massager.
- **Battery Care:** To prolong battery life, avoid fully discharging the battery frequently. Charge the unit regularly, even if not in frequent use.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Massager does not turn on.	Battery is low or depleted.	Charge the massager fully using the provided USB cable.
No heat or weak heat.	Heat function not activated or set to low.	Press the Heat Button to activate and cycle through heat levels. Allow a few minutes for the heat to build up.
Massage feels too weak or too strong.	Intensity level is not suitable.	Adjust the massage intensity using the Intensity Button. Ensure the wraps are adjusted correctly.
Massager stops working unexpectedly.	Battery depleted or automatic shut-off engaged.	Recharge the massager. The device may have an automatic shut-off timer for safety; simply restart if desired.

SPECIFICATIONS

Model: B082CXMSVR
Brand: COMFIER
Product Dimensions: 11 x 4.3 x 17.7 inches
Weight: 5.07 Pounds
Power Source: Battery Powered (Rechargeable Lithium-ion)
Material: Acrylonitrile Butadiene Styrene, Polyethylene (PE)
Color: Gray
Intended Use: Feet, Legs

WARRANTY AND SUPPORT

For warranty information or technical support, please contact COMFIER customer service. Refer to the contact details provided on the product packaging or the official COMFIER website.



Website: [Visit the COMFIER Store on Amazon](#)
Please have your model number (B082CXMSVR) and purchase date ready when contacting support.





OFFICIAL PRODUCT VIDEOS

No official product videos from the seller are available for this model at this time.

© 2023 COMFIER. All rights reserved.

Related Documents - B082CXMSVR

	<p>COMFIER FE-7204B Compression Leg Massager User Manual and Safety Guide</p> <p>Comprehensive user manual and safety guide for the COMFIER FE-7204B Compression Leg Massager, covering features, operation, maintenance, troubleshooting, and warranty information.</p>
	<p>Comfier CF-5913 Shiatsu Foot Massager User Manual</p> <p>User manual for the Comfier CF-5913 Shiatsu Foot Massager Machine. Learn how to use this kneading foot massager with heat for pain relief, plantar fasciitis, and general relaxation. Includes operation, safety, maintenance, and warranty details.</p>

<div><div><div>COMFIER</div><div>Leg and Foot Massager</div></div><div><div>User Manual CF-F25H</div></div></div>	<p>Comfier CF-F25H Leg and Foot Massager User Manual</p> <p>Comprehensive user manual for the Comfier CF-F25H Leg and Foot Massager, detailing safety instructions, product features, setup and operation, maintenance, troubleshooting, warranty information, and FCC compliance.</p>
<div><div><div>COMFIER</div><div>Air Compression Leg Massager</div></div><div><div>CE User Manual FE-7204B</div></div></div>	<p>COMFIER FE-7204B Air Compression Leg Massager User Manual</p> <p>Comprehensive user manual for the COMFIER FE-7204B Air Compression Leg Massager, detailing safety instructions, operating procedures, features, technical specifications, and troubleshooting for optimal use and care.</p>
<div><div><div>COMFIER</div><div>Heated Massage Footrest Reposero piedi de massage chauffant Ripoggiapedi massaggianti riscaldato Reposapies de massage con calefacción Vierușoare de masaj cu încălzire</div></div><div><div>CE User Manual CF-5420</div></div></div>	<p>Comfier CF-5420 Heated Massage Footrest User Manual</p> <p>User manual for the Comfier CF-5420 Heated Massage Footrest, providing detailed instructions on safety, operation, features, maintenance, and warranty information.</p>
<div><div><div>COMFIER</div><div>Cordless Neck Massager</div></div><div><div>User Manual CF-6026</div></div></div>	<p>Comfier Cordless Neck Massager CF-6026 User Manual</p> <p>Comprehensive user manual for the Comfier Cordless Neck Massager (Model CF-6026), detailing safety instructions, features, technical specifications, usage guidelines, troubleshooting tips, care and maintenance procedures, and warranty information.</p>