



Manuals.plus /

› Medisana /

› Medisana Electric Acupuncture Plus Instruction Manual

Medisana Electric Acupuncture Plus

Medisana Electric Acupuncture Plus Instruction Manual

Your guide to safe and effective needle-free acupuncture therapy.

1. INTRODUCTION

The Medisana Electric Acupuncture Plus is a compact, light, and user-friendly device designed for needle-free acupuncture therapy. It helps to locate and stimulate acupuncture points on the body, offering a natural and comfortable treatment experience. This manual provides essential information for the safe and effective use of your device.



Image 1: The Medisana Electric Acupuncture Plus device shown with its retail packaging, illustrating its compact design and how it is held for use.

2. SAFETY INFORMATION

Please read all safety instructions carefully before using the device. Failure to follow these instructions may result in injury or damage to the device.

- Do not use this device if you have a pacemaker, defibrillator, or any other implanted electronic device.
- Consult a physician before use if you are pregnant, have epilepsy, heart disease, or any other serious medical condition.
- Do not apply the device to open wounds, irritated skin, or sensitive areas such as the eyes or throat.
- Keep out of reach of children.
- This device is for personal use only and should not be shared without proper cleaning.
- Discontinue use immediately if you experience any discomfort or adverse reactions.

3. PACKAGE CONTENTS

Upon opening the package, please verify that all components are present and undamaged:

- Medisana Electric Acupuncture Plus device
- Instruction Manual
- Therapy Manual
- Button Batteries (typically 2 units)

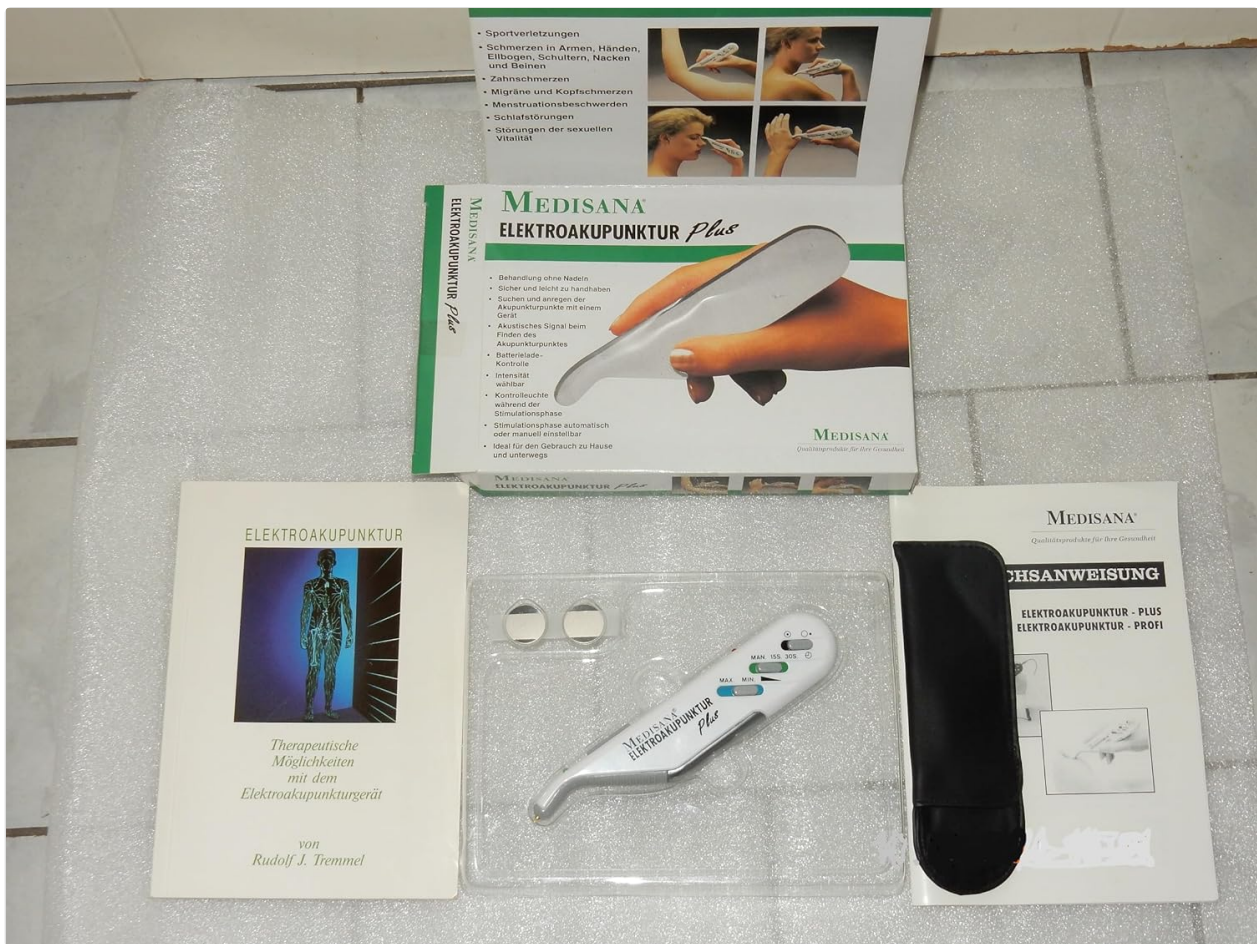


Image 2: An overhead view of the Medisana Electric Acupuncture Plus package contents, showing the device, the instruction manual, the therapy manual, and two button batteries.

4. DEVICE COMPONENTS

Familiarize yourself with the different parts of your Medisana Electric Acupuncture Plus device:

- **Stimulation Tip:** The contact point for applying stimulation to acupuncture points.
- **Grounding Plate:** Located on the side of the device, this must be in contact with your skin during use.
- **MAX/MIN Buttons:** Adjust the intensity of the electrical stimulation.
- **15S/30S Buttons:** Select the duration of the stimulation (15 or 30 seconds).
- **Indicator Light:** Illuminates to indicate device operation and battery status.
- **Battery Compartment:** Located at the rear of the device for battery insertion.



Image 3: A detailed view of the Medisana Electric Acupuncture Plus device, highlighting the MAX/MIN intensity buttons, 15S/30S duration buttons, and the stimulation tip.

5. SETUP

1. **Insert Batteries:** Open the battery compartment cover. Insert the two button batteries, ensuring correct polarity (+/-). Close the compartment securely.
2. **Prepare Skin:** Ensure the skin area where you intend to apply the device is clean and dry.

6. OPERATING INSTRUCTIONS

Follow these steps for effective use of your Electric Acupuncture Plus device:

1. **Hold the Device:** Grip the device firmly, ensuring your thumb or finger is in full contact with the grounding plate on the side.
2. **Locate Acupuncture Point:** Gently move the stimulation tip over the skin in the area where you wish to apply therapy. The device is designed to emit an acoustic signal when an acupuncture point is detected.
3. **Adjust Intensity:** Once an acupuncture point is located, press the 'MAX' or 'MIN' buttons to adjust the stimulation intensity to a comfortable level. Start with the lowest setting and gradually increase.
4. **Set Duration:** Select the desired stimulation duration by pressing the '15S' (15 seconds) or '30S' (30 seconds) button. The indicator light will illuminate during stimulation.
5. **Apply Stimulation:** Keep the stimulation tip on the acupuncture point for the selected duration. The device will automatically stop after the set time.
6. **Repeat:** Move to other acupuncture points as needed, following the guidance in the included Therapy Manual.

7. MAINTENANCE

- **Cleaning:** After each use, wipe the stimulation tip and the device body with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the device in a cool, dry place, away from direct sunlight and extreme temperatures.
- **Battery Replacement:** Replace batteries when the indicator light dims or the device fails to operate. Dispose of used batteries according to local regulations.

8. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Dead or incorrectly inserted batteries.	Check battery polarity or replace batteries.
No acoustic signal when locating points.	Poor skin contact or no acupuncture point detected.	Ensure firm contact with the grounding plate and stimulation tip. Try a different area.
Weak or no stimulation.	Low battery power or intensity set too low.	Replace batteries. Increase intensity using the 'MAX' button.

9. SPECIFICATIONS

- **Brand:** Medisana
- **Model:** Electric Acupuncture Plus
- **Product Dimensions:** 15 x 3 x 1 cm
- **Weight:** 1 kg (including packaging)
- **Power Source:** Button Batteries (type not specified, typically CR2032 or similar)
- **ASIN:** B083PVKS36
- **Manufacturer:** Medisana

10. WARRANTY AND SUPPORT

Specific warranty information for this product is not provided in the available documentation. For warranty claims, technical support, or further assistance, please contact Medisana customer service or refer to the official Medisana website for your region.