

## Pure Daily Care LEDMASK

# Pure Daily Care Luma Mask LED Skincare Device User Manual

Model: LEDMASK

[Information](#)   [Components](#)   [Setup](#)   [Operation](#)   [Maintenance & Support](#)   [Introduction](#)   [Safety](#)   [Troubleshooting](#)   [Specifications](#)   [Warranty](#)

## 1. INTRODUCTION

The Pure Daily Care Luma Mask is an advanced LED skincare device designed for at-home use. It utilizes natural light energy through 150 light-emitting diodes (LEDs) to deliver 7 distinct wavelengths of light, each targeting specific skin concerns. This manual provides comprehensive instructions for the safe and effective use of your Luma Mask, helping you achieve a younger, smoother, brighter, clearer, and more radiant complexion.

The Luma Mask offers 7 color modes and 5 light intensity levels, providing a customizable light therapy experience without UV exposure.

## 2. IMPORTANT SAFETY INFORMATION

- Read all instructions before using the device.
- This device is for external use only. Do not use on broken or irritated skin.
- The Luma Mask emits light that is UV-free. However, it is recommended to keep your eyes closed during treatment or use additional eye protection for enhanced comfort, especially if you experience discomfort from the light.
- Keep the device away from water and moisture. Do not use in wet environments.
- Do not attempt to disassemble or repair the device. Contact customer support for assistance.
- Keep out of reach of children.
- Consult a healthcare professional before use if you have any pre-existing medical conditions, are pregnant, or are taking photosensitizing medications.
- The device must be plugged into a power outlet during use. Ensure the USB power cable is securely connected to both the mask and the controller, and the controller is connected to a power source.

### 3. PRODUCT COMPONENTS

---

Your Pure Daily Care Luma Mask package includes the following items:

- 1 LED Face Mask
- 1 Adjustable Elasticized Strap
- 1 Light Controller
- 1 USB Power Cable
- User Guide (this document)

## LUMA MASK

### LED SKIN TREATMENT MASK

150 Light Emitting Diodes (LEDs)

7 Light Wavelength Colors

5 Light Intensity Levels

Non-Invasive & Pain-Free

UV-Free and Safety Tested



Image: The Luma Mask with a list of its key features: 150 Light Emitting Diodes (LEDs), 7 Light Wavelength Colors, 5 Light Intensity Levels, Non-Invasive & Pain-Free, UV-Free and Safety Tested.



Image: The retail packaging for the Pure Daily Care Luma Mask, showing the brand logo and product name.

## 4. SETUP GUIDE

---

1. **Unpack:** Carefully remove all components from the packaging.
2. **Connect Controller:** Plug the USB power cable into the light controller.
3. **Connect Mask:** Connect the other end of the USB power cable to the port on the Luma Mask. Ensure a secure connection to prevent interruptions during use.
4. **Power Source:** Plug the light controller into a standard electrical outlet. The device requires continuous power during operation.
5. **Adjust Strap:** Attach the adjustable elasticized strap to the mask and adjust it for a comfortable fit around your head. Alternatively, you may lie down and rest the mask on your face without the strap.

## 5. OPERATING INSTRUCTIONS

---

### 5.1 General Usage

1. **Preparation:** Cleanse your face thoroughly and ensure it is dry before placing the mask.
2. **Position Mask:** Place the Luma Mask comfortably on your face. Adjust the elastic strap if using, or lie down to support the mask's weight. If the mask feels heavy on the bridge of your nose, try adjusting its position or supporting it with your hands.
3. **Power On:** Press the power button on the light controller to turn on the device.
4. **Select Color Mode:** Use the "Color" button on the controller to cycle through the 7 available light colors (Red, Blue, Green, Cyan, Yellow, Purple, White) or the 8th setting which cycles through all colors.

5. **Adjust Intensity:** Use the "Intensity" buttons (+/-) to select one of the 5 light intensity levels. Start with a lower intensity and gradually increase as comfortable.
6. **Set Timer:** Use the "Timer" buttons (+/-) to set your desired session duration. The mask will automatically shut off when the treatment is complete.
7. **Eye Protection:** While the mask is designed to be UV-free, it is recommended to keep your eyes closed during treatment. For additional comfort, especially if you find the light too bright, consider using tanning goggles or similar eye protection.
8. **Relax:** During the session, relax and allow the light therapy to work.
9. **Power Off:** After the session, the mask will automatically turn off. You can also manually turn it off by pressing the power button on the controller.

## 5.2 Recommended Usage Frequency and Duration

For optimal results, it is generally recommended to use the Luma Mask 3-5 times per week, with each session lasting 15-30 minutes. Individual results may vary. Consistency is key for visible improvements.

You may choose to focus on one color per session or combine multiple colors. For example, you could use Red light for 15 minutes followed by Blue light for 15 minutes in a single session.

## 5.3 Understanding Light Therapy Modes

Each of the 7 light colors penetrates the skin to a different depth, offering unique benefits:

# 7 LIGHT WAVELENGTHS

to address an array of skin concerns

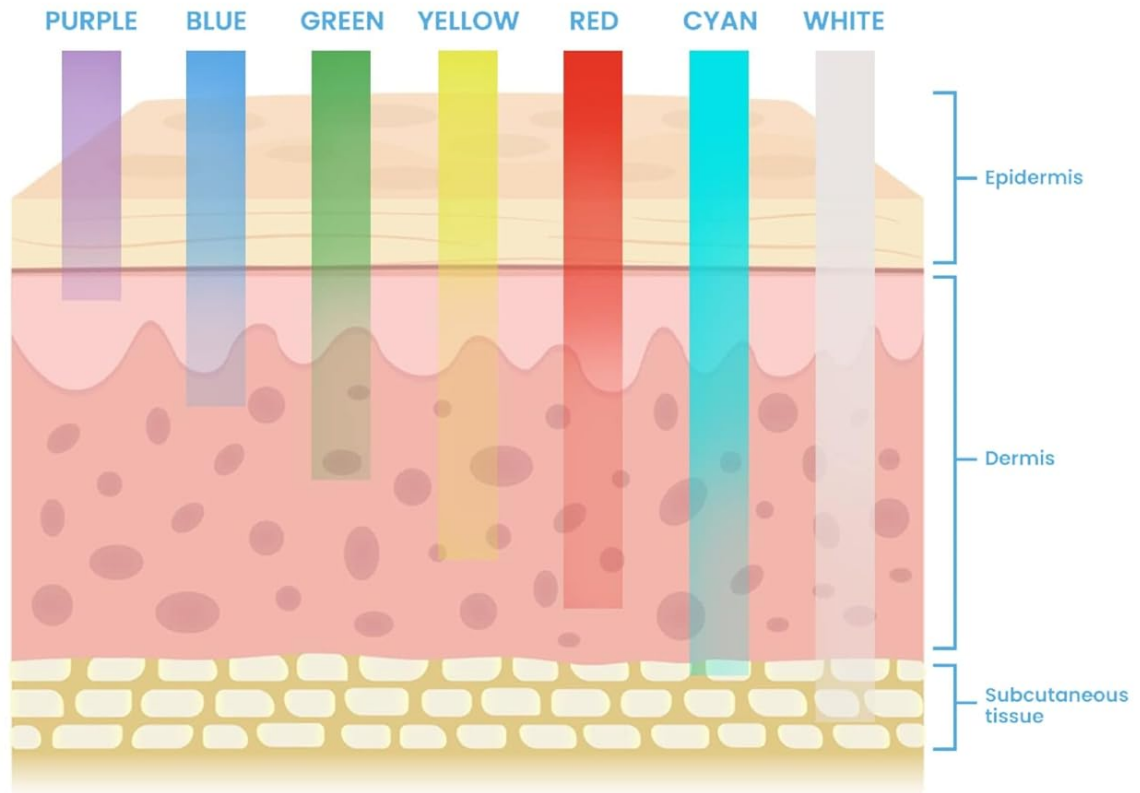


Image: A diagram illustrating how different light wavelengths (Purple, Blue, Green, Yellow, Red, Cyan, White) penetrate various layers of the skin (Epidermis, Dermis, Subcutaneous tissue).

- **Red Light (630nm):** Targets the epidermis and dermis. Promotes a more youthful appearance, reduces the look of fine lines and wrinkles, and supports skin revitalization.

# RED LIGHT

for younger-looking skin



Image: A woman holding the Luma Mask, which is emitting red light, indicating its use for a younger-looking complexion.

- **Blue Light (470nm):** Primarily targets the epidermis. Helps in balancing skin and clarifying the complexion.

# BLUE LIGHT

to soothe & clarify the complexion



Image: A woman holding the Luma Mask, which is emitting blue light, suggesting its use for soothing and clarifying the complexion.

- **Green Light (520nm):** Penetrates the epidermis. Contributes to a more radiant and even-toned complexion.

# GREEN LIGHT

for skin that looks brighter



Image: A woman holding the Luma Mask, which is emitting green light, indicating its benefit for brighter-looking skin.

- **Cyan Light (490nm):** A blend of blue and green light. Offers benefits for calming and soothing the skin.
- **Yellow Light (590nm):** Targets the epidermis. Helps improve skin tone and texture.
- **Purple Light (420nm):** A combination of red and blue light. Supports skin revitalization and clarity.
- **White Light:** Combines all wavelengths. Provides comprehensive benefits for overall skin health and appearance.

## 6. CARE AND MAINTENANCE

---

- **Cleaning:** After each use, gently wipe the mask with a soft, dry cloth. If necessary, use a slightly damp cloth with mild soap, then wipe dry immediately. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the Luma Mask in a cool, dry place, away from direct sunlight and extreme temperatures. Keep it in its original packaging or a protective case to prevent damage.
- **Cable Care:** Handle the USB power cable and controller with care. Avoid bending or twisting the cable excessively, as this can lead to internal damage and connection issues. Ensure the cable is securely plugged into the mask and controller during use.

## 7. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	<ul style="list-style-type: none"><li>Power cable not securely connected.</li><li>Power outlet not functioning.</li><li>Controller or mask malfunction.</li></ul>	<ul style="list-style-type: none"><li>Ensure the USB power cable is firmly connected to both the mask and the controller, and the controller is plugged into a working power outlet.</li><li>Try a different power outlet.</li><li>If the issue persists, contact customer support.</li></ul>
Lights flicker or turn off unexpectedly.	<ul style="list-style-type: none"><li>Loose connection of the USB power cable.</li><li>Intermittent power supply.</li></ul>	<ul style="list-style-type: none"><li>Check all cable connections to ensure they are secure.</li><li>Ensure the power source is stable.</li><li>Avoid excessive movement of the mask or cable during use.</li></ul>
Mask feels uncomfortable or heavy on the face.	<ul style="list-style-type: none"><li>Improper strap adjustment.</li><li>Pressure on nose/forehead.</li></ul>	<ul style="list-style-type: none"><li>Adjust the elastic strap for a looser fit.</li><li>Consider lying down during treatment to alleviate pressure.</li><li>Gently hold the mask to reduce pressure points if needed.</li></ul>
Only one color or limited colors are working.	<ul style="list-style-type: none"><li>Device malfunction.</li></ul>	<ul style="list-style-type: none"><li>Contact customer support for further assistance.</li></ul>

## 8. PRODUCT SPECIFICATIONS

Feature	Detail
<b>Model Number</b>	LEDMASK
<b>Brand</b>	Pure Daily Care
<b>Light Source</b>	150 LED Diodes
<b>Color Modes</b>	7 (Red, Blue, Green, Cyan, Yellow, Purple, White) + All Colors Cycle
<b>Intensity Levels</b>	5
<b>Power Source</b>	Wired (requires continuous power via USB cable and adapter)
<b>UV Emission</b>	None (UV-free)
<b>Package Dimensions</b>	10.63 x 8.9 x 5.55 inches
<b>Item Weight</b>	2.05 Pounds
<b>UPC</b>	810016711254
<b>Manufacturer</b>	Eyco Technology Co

## 9. WARRANTY AND SUPPORT

For information regarding product warranty, please refer to the warranty card included with your purchase or visit the official

Pure Daily Care website. Warranty terms and conditions may vary.

If you encounter any issues or have questions not covered in this manual, please contact Pure Daily Care customer support for assistance. You can typically find contact information on the product packaging or the brand's official website.

**Pure Daily Care Website:** [www.puredailycare.com](http://www.puredailycare.com)