YUNMAI PRO

YUNMAI PRO Smart Body Fat Scale User Manual

Model: PRO

1. INTRODUCTION

This manual provides essential information for the safe and effective use of your YUNMAI PRO Smart Body Fat Scale. Please read it thoroughly before first use and retain it for future reference.

The YUNMAI PRO Smart Body Fat Scale is designed to measure various body composition metrics, including weight, body fat percentage, muscle mass, and BMI, providing a comprehensive overview of your health and fitness progress. It connects via Bluetooth to the YUNMAI app for data tracking and analysis.

2. SAFETY INFORMATION

- Do not use this device if you have a pacemaker or other internal medical device.
- Do not use this device if you are pregnant.
- The scale surface can be slippery when wet. Ensure the surface is dry before stepping on.
- Keep the scale away from strong electromagnetic fields.
- Do not jump on the scale or drop heavy objects on it.
- This scale is designed for household use only and is not intended for medical diagnosis or treatment.
- Maximum weight limit: 180 kg (396 lbs).

3. PACKAGE CONTENTS

- YUNMAI PRO Smart Body Fat Scale
- User Manual (this document)
- USB Charging Cable

4. SETUP

4.1 Charging the Scale

Before first use, ensure the scale is fully charged. Connect the provided USB charging cable to the charging port on the

side or back of the scale and to a standard USB power adapter (not included). The LED display will indicate charging status.

4.2 App Installation

The YUNMAI PRO Smart Body Fat Scale requires the YUNMAI app for full functionality and data tracking. The app is compatible with iOS and Android devices.

1. Search for "YUNMAI" in the Apple App Store (for iOS) or Google Play Store (for Android).
2. Download and install the YUNMAI app.
3. Create a new user account or log in if you already have one.

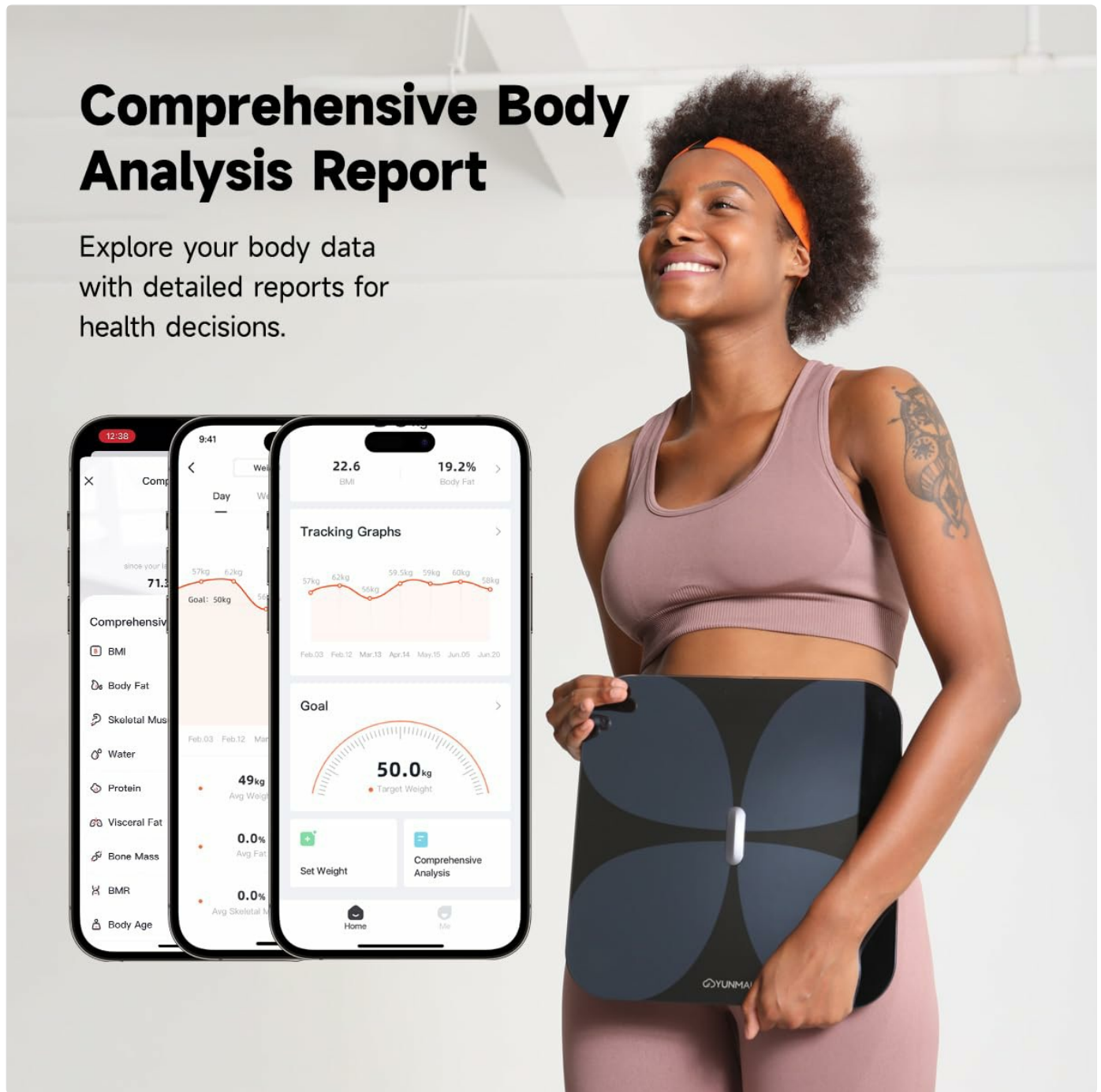
Image: The YUNMAI app integrates with popular fitness platforms like Apple Health, Google Fit, and Fitbit.

4.3 Pairing the Scale with the App

1. Ensure Bluetooth is enabled on your smartphone.
2. Open the YUNMAI app.

3. Step onto the scale briefly to activate it. The display will light up.
4. The app should automatically detect and connect to the scale. Follow any on-screen prompts within the app to complete the pairing process.
5. Once paired, your measurements will synchronize with the app after each weigh-in.

5. OPERATING INSTRUCTIONS

5.1 Taking a Measurement

1. Place the scale on a hard, flat surface. Avoid carpets or uneven floors, as this can affect accuracy.
2. Ensure your feet are clean and dry.
3. Step onto the scale barefoot, standing still with your feet evenly distributed over the ITO coating sensors.
4. The scale will display your weight first, then cycle through other body composition metrics.
5. For full body composition analysis, ensure the YUNMAI app is open and connected via Bluetooth. Your data will automatically sync.

Image: The YUNMAI PRO scale displaying weight, with the companion app showing detailed body analysis.

Track Your Progress Reach Your Goal

Monitor your health trends weekly, monthly and yearly, providing the insights to reach your fitness goals with YUNMAI App.

Image: The scale provides 13 key body metrics for comprehensive analysis.

5.2 Understanding Measurements

The YUNMAI app provides detailed explanations for each metric. Here are some common ones:

- **Weight:** Your total body mass.
- **BMI (Body Mass Index):** A measure of body fat based on height and weight.
- **Body Fat %:** The percentage of your total body weight that is fat.
- **Muscle Mass:** The total weight of muscle in your body.
- **Water %:** The percentage of your total body weight that is water.
- **Bone Mass:** The total weight of bone in your body.
- **Visceral Fat Index:** An indicator of fat stored around your organs.

Accurate Measurements

Detect slight changes in weight and body metrics,
accurate to 0.11lb.


4 High-precision Sensors

BIA Technology

0.11lb Increments

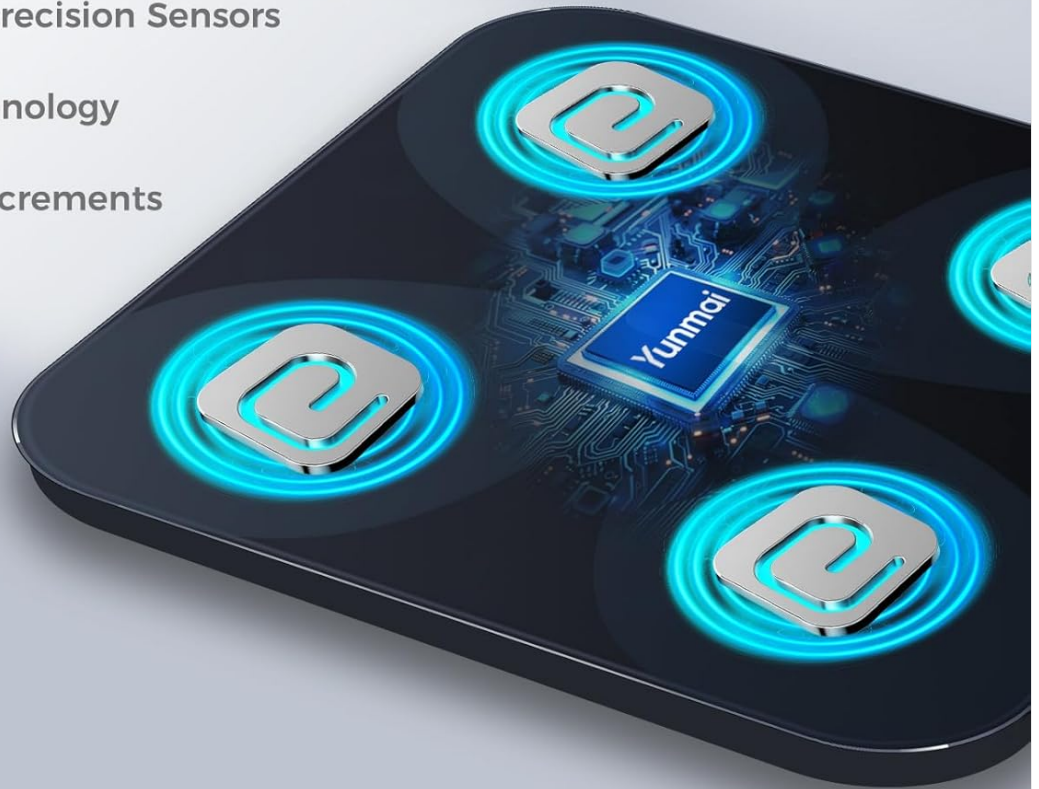*Designed for human use only.

Image: Track your progress over time with detailed graphs in the YUNMAI app.

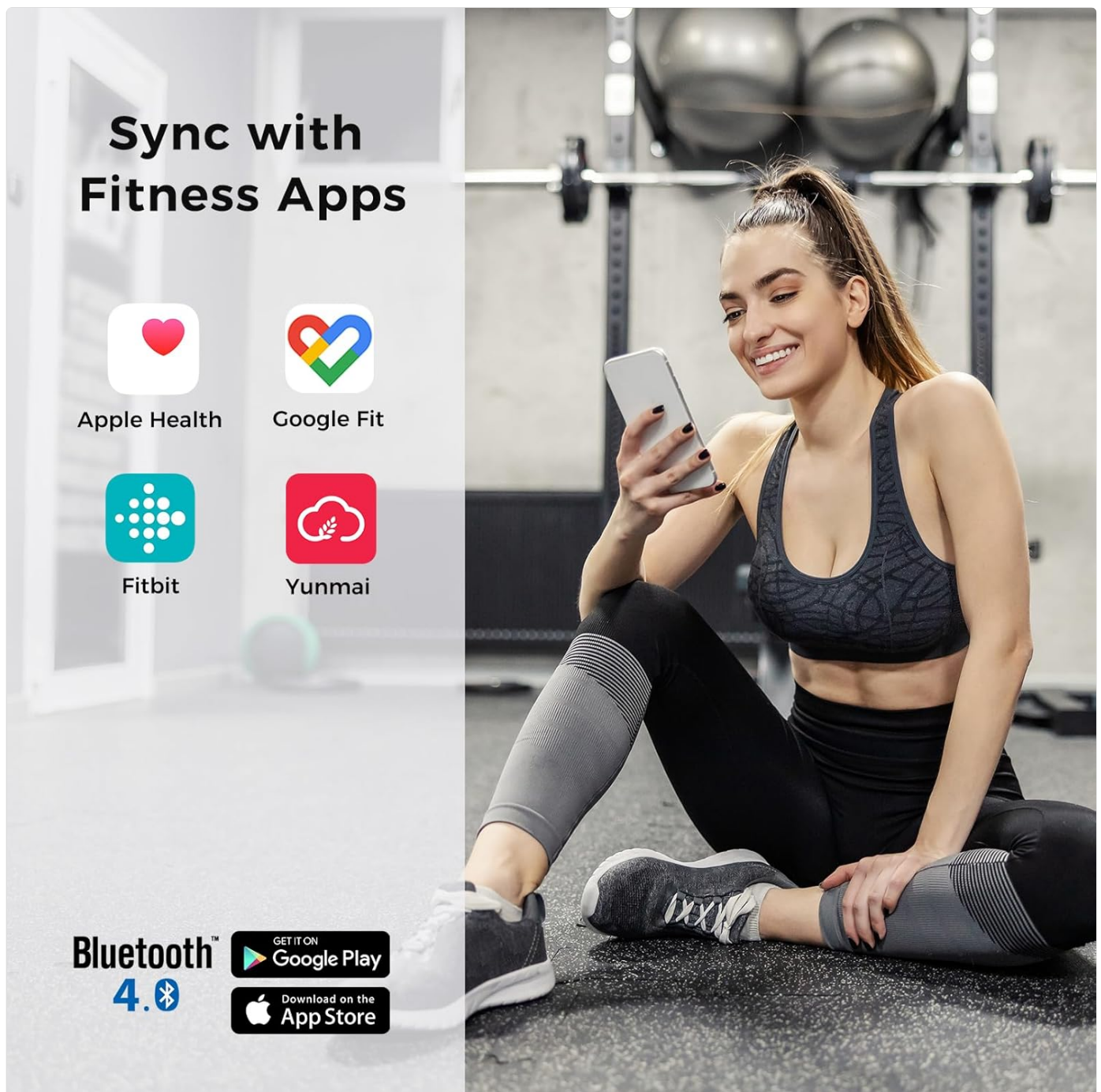
Image: The YUNMAI app provides a comprehensive body analysis report to help with health decisions.

5.3 Multi-User Management

The YUNMAI PRO scale supports up to 16 users, automatically recognizing individuals and syncing data to their respective profiles within the app. This feature is ideal for families or shared use.

1. Each user should create their own profile within the YUNMAI app.
2. Ensure each user enters their correct personal information (age, height, gender) for accurate body composition analysis.
3. When a user steps on the scale, the app will attempt to identify them based on their previous measurements. If multiple users have similar weights, the app may prompt you to select the correct profile.

Support 16 Users, Auto Recognition

Automatically syncs data to the correct profile for each user, perfect for families or shared use.

Image: The YUNMAI app supports up to 16 users with auto-recognition for family or shared use.

6. MAINTENANCE AND CARE

- Clean the scale surface with a damp cloth and mild detergent. Do not use abrasive cleaners or immerse the scale in water.
- Store the scale in a cool, dry place.
- Avoid exposing the scale to direct sunlight or extreme temperatures.
- The scale features a durable tempered glass surface. Handle with care to prevent damage.

Large Platform, Durable Design

Sleek Rounded Edges

Measurement Range:

11~ 400lb / 5 ~180kg

Image: The YUNMAI PRO scale features durable tempered glass and anti-slip pads for stability.

7. TROUBLESHOOTING

Problem	Possible Cause	Solution
Scale does not turn on.	Low battery or no power.	Charge the scale using the provided USB cable.
Inaccurate measurements.	Scale on uneven surface; feet not properly placed; scale not calibrated.	Place scale on a hard, flat surface. Ensure bare feet are centered on sensors. Step on and off quickly to calibrate to zero before weighing.
Cannot connect to app/Bluetooth.	Bluetooth off; app not open; scale not activated; interference.	Ensure Bluetooth is on. Open YUNMAI app. Step on scale to activate. Move closer to scale. Restart phone/app.

Problem	Possible Cause	Solution
Body composition data not showing.	Not standing barefoot; feet not making proper contact with sensors.	Ensure you are standing barefoot with feet centered on the ITO coating sensors. Ensure feet are clean and dry.

ITO Coating Technology

Larger sensor coverage, fit all size of feet.

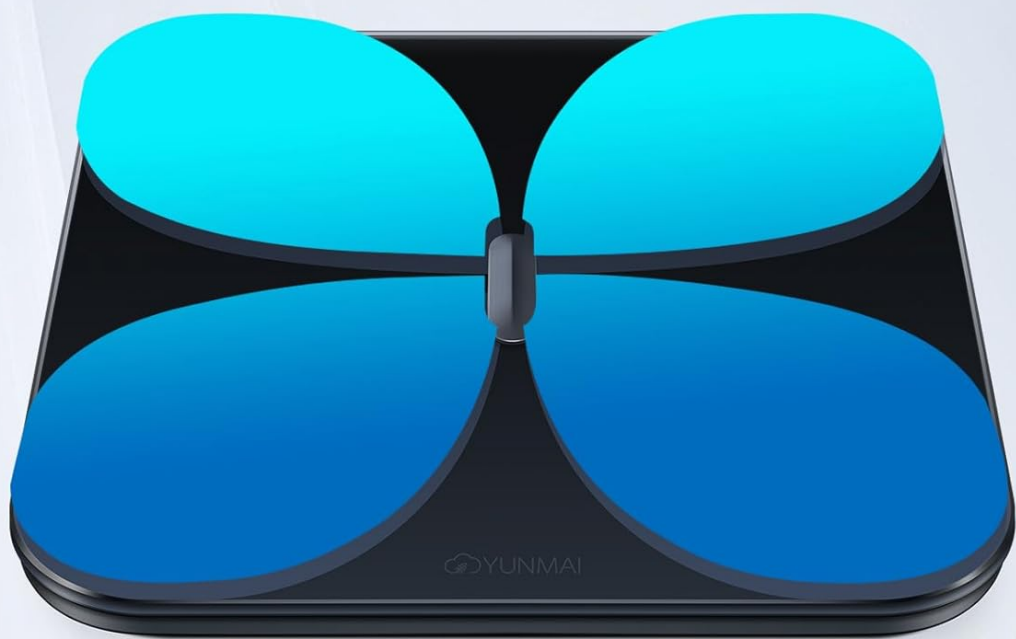
Image: The scale uses four high-precision sensors and BIA technology for accurate measurements.

13 Key Body Metrics Analysis

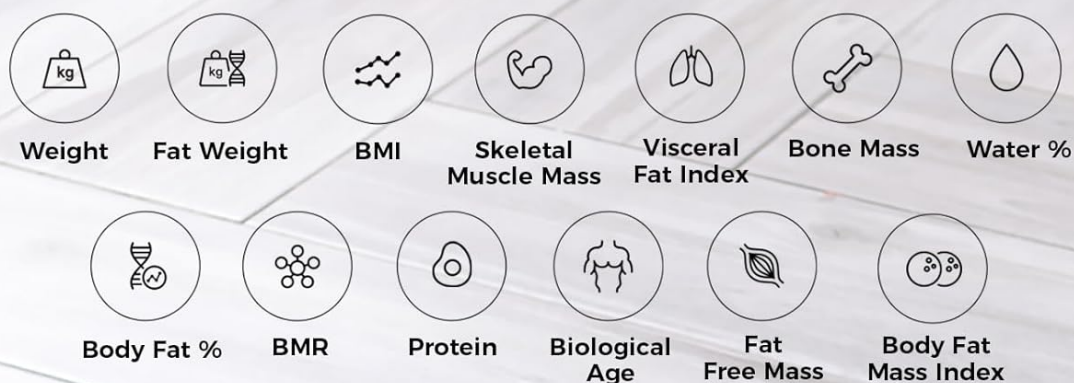
Image: ITO coating technology ensures larger sensor coverage for accurate readings across different foot sizes.

8. SPECIFICATIONS

- **Model:** YUNMAI PRO
- **Dimensions:** 30 x 30 x 2 cm (11.8 x 11.8 x 0.79 inches)
- **Weight:** 1.7 kg (3.75 lbs)
- **Weight Limit:** 5 - 180 kg (11 - 396 lbs)
- **Measurement Increments:** 0.05 kg (0.11 lbs)
- **Power Source:** 1 Lithium-polymer battery (included)
- **Connectivity:** Bluetooth 4.0
- **Material:** Tempered Glass, Acrylic
- **Color:** Black
- **Supported Users:** Up to 16


9. WARRANTY AND SUPPORT

YUNMAI products are designed and manufactured to the highest quality standards. For warranty information and technical support, please refer to the official YUNMAI website or contact their customer service directly. Keep your purchase receipt as proof of purchase.

For further assistance, visit the official YUNMAI support page or contact their customer service team.


© 2024 YUNMAI. All rights reserved.

Related Documents - PRO

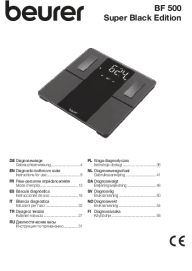	True Wireless Earphones User Manual - PRO 4 User manual for the PRO 4 True Wireless Earphones, including product specifications, operation instructions, pairing, calling, and troubleshooting.
	YUNMAI X Smart Scale YMBS-M268 User Manual - Setup and Operation Guide Comprehensive user manual for the YUNMAI X Smart Scale (Model YMBS-M268). Learn how to set up, connect via Bluetooth, use the YUNMAI App, and understand product features and specifications for a healthier lifestyle.
	Beurer BF 500 Super Black Edition Diagnostic Bathroom Scale - User Manual User manual for the Beurer BF 500 Super Black Edition diagnostic bathroom scale. Learn how to use, set up, and maintain your scale for accurate body composition analysis. Features Bluetooth connectivity and the HealthManager Pro app.
	YUNMAI Premium Smart Scale User Manual User manual for the YUNMAI Premium Smart Scale (M1301). Covers features, setup, usage, specifications, and warranty in English, Czech, Slovak, Hungarian, and German.

 YUNMAI S User Manual • Uživatelský manuál • Uživatelský manuál • Használati utasítás • Benutzerhandbuch	YUNMAI S Smart Scale User Manual Comprehensive user manual for the YUNMAI S Smart Scale, detailing setup, app integration, usage instructions, product specifications, and warranty information.
 YUNMAI Pro User Manual • Uživatelský manuál • Uživatelský manuál • Használati utasítás • Benutzerhandbuch	YUNMAI Pro Smart Scale User Manual and Guide Comprehensive user manual for the YUNMAI Pro Smart Scale, covering setup, usage, maintenance, specifications, and warranty information. Includes BMI tracking and app connectivity.


Documents - YUNMAI – PRO

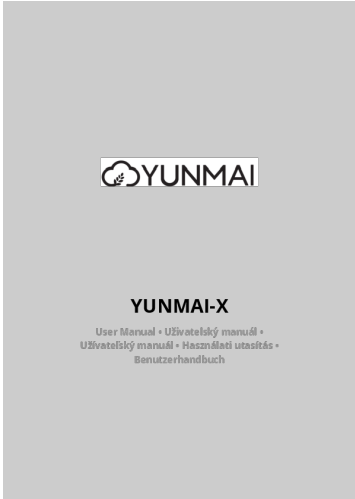[\[pdf\]](#) User Manual Instructions

Lukáš Chlebovský yunmai x Alza.cz server Please read the following instructions carefully before first use and keep this user manual for future reference Pay particular attention YUNMAI X dwn alza 100038 ||| YUNMAI-X User Manual Uživatelsk manu Uzvatesk manu Hasznlati utasts Benutzerhandbuch English 3 12 Cestina 13 22 Slovencina 23 32 Magyar 33 42 Deutsch 43 52 2 Dear customer, Thank you for purchasing our **product**. Please read the following instructions carefully before first use and kee...

lang:en score:40 filesize: 590.57 K page_count: 52 document date: 2021-12-07

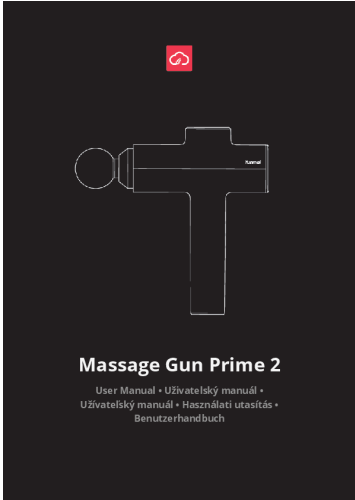[Yunmai Massage Gun Prime 2: Official User Manual & Operation Guide](#)

Official user manual for the Yunmai Massage Gun Prime 2 (YMFG-M418). Get detailed instructions on usage, safety, charging, maintenance, and troubleshooting for optimal muscle recovery.

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[YUNMAI S Smart Scale User Manual](#)

Comprehensive user manual for the YUNMAI S Smart Scale, detailing setup, app integration, usage instructions, product specifications, and warranty information.

lang:en **score:32** filesize: 686.39 K page_count: 47 document date: 2021-12-06

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
[YUNMAI Wrist Ball User Manual](#)

User manual for the YUNMAI Wrist Ball (Model YMGB-Z701/YMGB-Z702), detailing product features, usage instructions, safety precautions, and warranty information.

lang:en **score:31** filesize: 221.52 K page_count: 17 document date: 2021-12-07

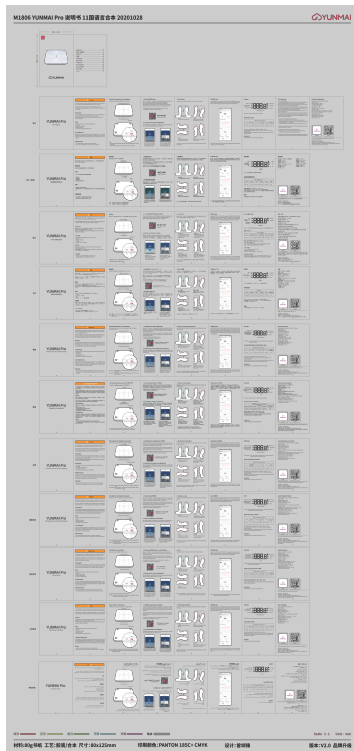
[\[pdf\]](#) User Manual

M1806 YUNMAI Pro 11 20201028 User manuals for Yunmai devices massagers smart scales jump ropes 11 Language drive google file d 181eylfvdiU3kLFkjc6kJF4qhm2z mmG5 view usp sharing u003eDownload |||

M1806 YUNMAI **Pro** 11 20201028 80mm0.5 80mm0.5 English 01 10 18 26 Deutsch 34 42 Fr ... eitungen fr die intelligente Waage 34 YUNMAI **Pro** Warning Do not disassemble, squeeze, puncture the **product** or put it into water, fire, or exposes it to temperature higher than 60 C 140 F . Please ...

lang:en **score:30** filesize: 12.16 M page_count: 1 document date: 2023-07-05

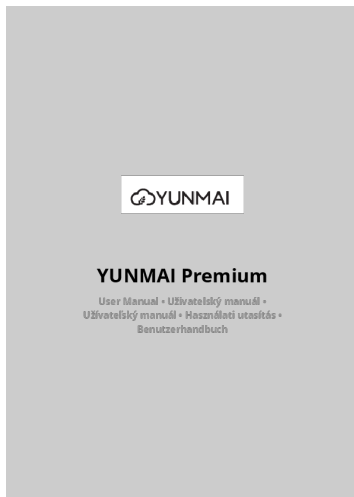[YUNMAI Premium Smart Scale User Manual](#)

User manual for the YUNMAI Premium Smart Scale (M1301). Covers features, setup, usage, specifications, and warranty in English, Czech, Slovak, Hungarian, and German.

lang:en **score:30** filesize: 383.05 K page_count: 52 document date: 2021-12-08

[\[pdf\]](#) User Manual


YUNMAI Jump Rope 20200921 Download User manual for Yunmai Fitness Pro manuals devices massagers smart scales jump ropes YMHR P701 drive google file d

1ekDP4Jl2gHz0vuq1gmXnpfKAWXv2FaBv view usp sharing |||

YUNMAI Jump Rope 20200921 125mm 80mm YUNMAI Jump Rope User Instruction **Product** Information Thank you for choosing the YUNMAI Jump Rope. Assembly 1. Remove the **protection** cap from the steel wire cable and insert the cable into the rope **protecting** tube Adjustment 1. Put your feet together and ...

lang:en **score:30** filesize: 260.21 K page_count: 1 document date: 2023-07-05

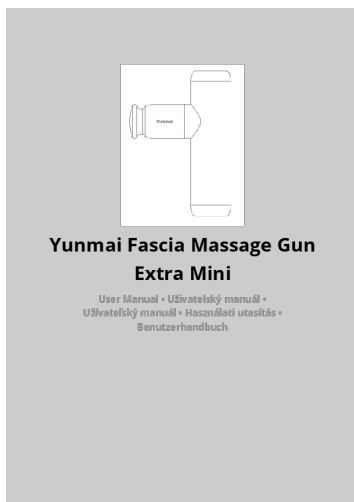
[\[pdf\]](#) User Manual Instructions Guide

Lukáš Chlebovský YUNMAI Extra Mini dwn alza cz manual 100030 |||

Yunmai Fascia Massage Gun Extra Mini User Manual Uživatelsk manu Uzvatesk manu Hasznlati utasts Benutzerhandbuch English 3 21 Cestina 22 39 Slovenscina 40 57 Magyar 58 75 Deutsch 76 - 93 2 The **product** user guide contains **product** features, how to use, and the operating **procedure**. Read th...

lang:en **score:29** filesize: 792.53 K page_count: 93 document date: 2021-12-06

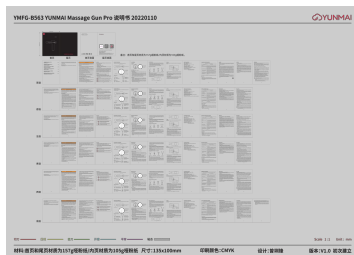
[\[pdf\]](#) User Manual

YMFGB563 YUNMAI Massage Gun Pro 20220106 User manuals for Yunmai devices massagers smart scales jump ropes 6 Language drive google file d 19pLPS 3xF HLTBv7D dPINBQIB1GbA02 view usp sharing u003eDownload |||

YMFGB563 YUNMAI Massage Gun **Pro** 20220110 100.00 mm YUNMAI Massage Gun **Pro** Made in China 135.00 ... 8 Chiwan 6th Road, Nanshan District, Shenzhen, China 157g/105g 01 Thank you for purchasing this **product**. Please read this manual carefully before use. For future reference, please ensure this manu...


lang:en **score:29** filesize: 1.1 M page_count: 1 document date: 2023-07-05

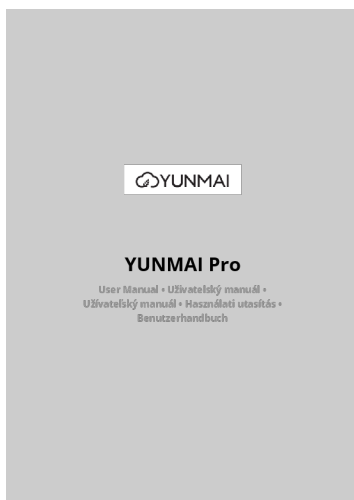[YUNMAI Pro Smart Scale User Manual and Guide](#)

Comprehensive user manual for the YUNMAI Pro Smart Scale, covering setup, usage, maintenance, specifications, and warranty information. Includes BMI tracking and app connectivity.


lang:en **score:28** filesize: 762.43 K page_count: 47 document date: 2021-12-06

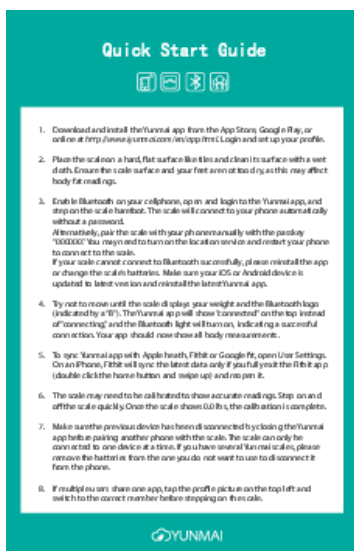[\[pdf\]](#) Quick Start Guide Installation Guide Frequently Asked Questions Guide

Yunmai Quick Start Guide0113 Amazon AWS Enable Bluetooth on your cellphone open and login to the app step scale barefoot The will connect phone automatically 1615868 Mini Smart Weight Scale Guide file yunmai color control bluetooth smart rrp 129 00 709080 02 v 637076190673913802 appinthe store s3 amazonaws 2016 globale FAQ data Home powerhousepc au media cdn mydeal 43780 98 APP 1. Download and install the Yunmai app from the App Store, Google Play, or online at <http://www.iyunmai.com/en/app.html>. Login and set up your **pro** le. 2. Place the scale on a hard, at surface like tiles and clean its surface with a wet cloth. Ensure the scale surface and your feet are not too...

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
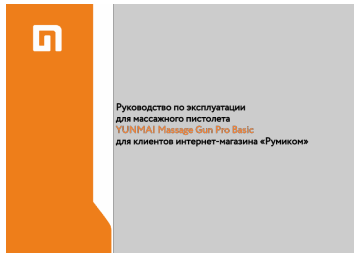
lang:so **score:27** filesize: 813.78 K page_count: 8 document date: 2021-01-21

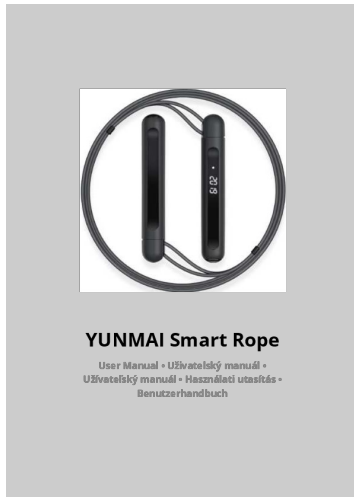[YUNMAI Smart Rope User Manual](#)

Comprehensive user manual for the YUNMAI Smart Rope, detailing product introduction, assembly, Bluetooth connection, startup, charging, and specifications. Includes warranty information and EU Declaration of Conformity.


lang:en **score:26** filesize: 612.31 K page_count: 47 document date: 2021-12-07

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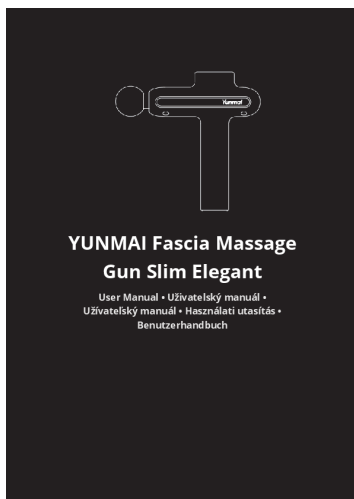
lang:hr **score:26** filesize: 946.09 K page_count: 19 document date: 2021-12-21

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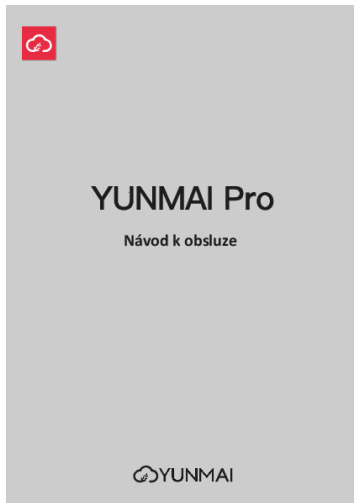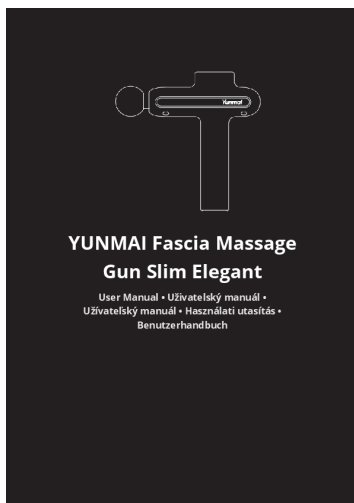[\[pdf\]](#) User Manual Guide

Lukáš Chlebovský YUNMAI Fascia Massage Gun Slim Elegant Dec 19 2021 — The red solid line in the figure is for alignment which bright black actual product When assembling or replacing massage head SLIM ELEGANT dwn alza cz manual 100035 ||| YUNMAI Fascia Massage Gun Slim Elegant User Manual Uzivatelsk manu Uzvatesk manu Hasznlati utasts Benutzerhandbuch English 3 20 Cestina 21 38 Slovencina 39 56 Magyar 57 74 Deutsch 75 92 2 The **product** user guide contains **product** features, how to use, and the operating **procedure**. Read t... lang:en score:24 filesize: 969.86 K page_count: 92 document date: 2021-12-07

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
lang:no **score:22** filesize: 1.14 M page_count: 12 document date: 2023-02-16

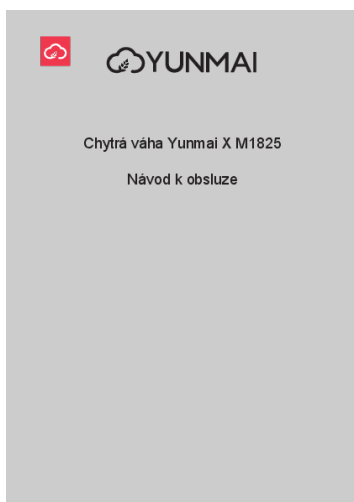[\[pdf\]](#)


tooooo docId 986ACFE4592FFC48FBA41CF464CDC17D Manuál 1 MB Chytrá váha Yunmai X M1825

bílá My Smart Home yunmai mysmarthome cz user related files |||

Chytr vha Yunmai X M1825 Nvod k obsluze Varovn Pstroj se nesm rozebrat, mackat, **propichovat**, umsovat do vody, ohn nebo vystavovat teplotm vyssm nez 60 C. Dlezit

Vha YUNMAI Smart Scale by nemla bt pouzvana k diagnostice nebo lcb jakhokoli zdravotnho stavu. Informace ped pouzitm - Pi prvnm pouzit vhy...

lang:hr **score:22** filesize: 1.14 M page_count: 10 document date: 2024-04-25

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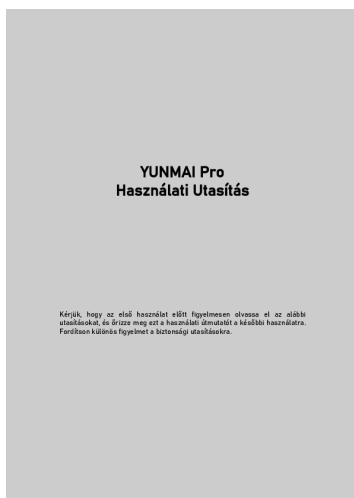
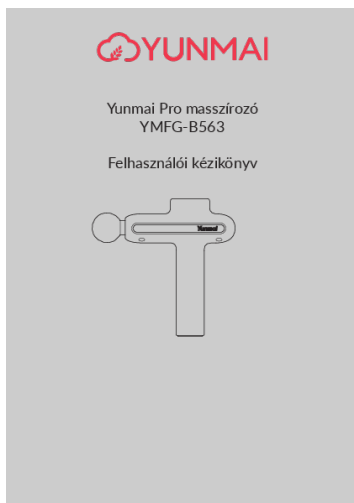
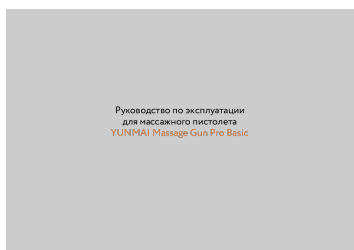
lang:pl score:21 filesize: 650.07 K page_count: 8 document date: 2022-09-14

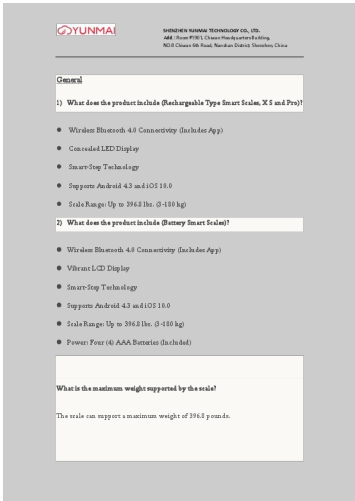
[pdf]

Microsoft Word YUNMAIÂ SMARTÂ SCALE Q A update 2019 10 25 My yunmai scale err1 Creation Date 2016 07 11 4 years 211 days left scales YUNMAI Official US Online Store Yunmai You will enjoy a more intelligent and interesting healthy life with App! Hurry only 1 left! original media s bol Rg9krm1gWxxO |||

SHENZHEN YUNMAI TECHNOLOGY CO., LTD. Add.: Room#1901, Chiwan Headquarters Building, NO.8 Chiwan 6th Road, Nanshan District, Shenzhen, China General 1 What does the **product** include Rechargeable Type Smart Scales, X S and **Pro** Wireless Bluetooth 4.0 Connectivity Includes App Concealed LED Displ...

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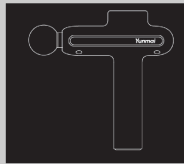
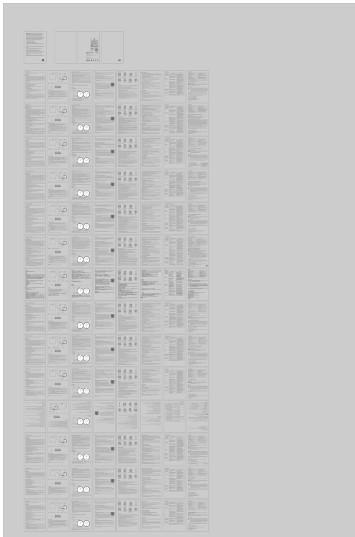
b2c i18n upload azams0 mi fds User Guide Upload Personal Care |||

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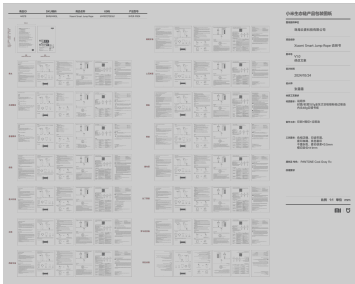Rope Blue UserGuide xiaomi b2c i18n upload azams0 mi fds User Guide Upload Personal Care |||

153mm 44578 105mm BHR6944GL Xiaomi Smart Jump Rope 6941812708361

44578 XMSR-P ... jumping, always wear shock-absorbing sneakers, such as thick-
soled and high-top sneakers, to be er **protect** your ankles and feet. Warm up before
jumping to help speed up the blood circulation of your ...

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
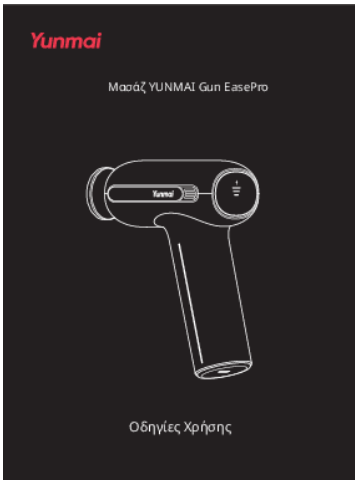
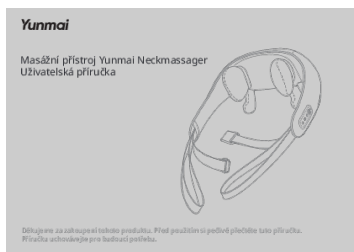
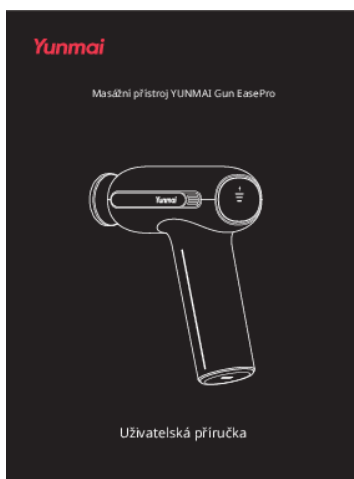
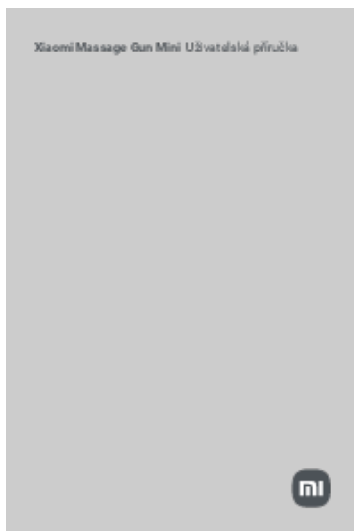
[\[pdf\]](#) User Manual Declaration of Conformity

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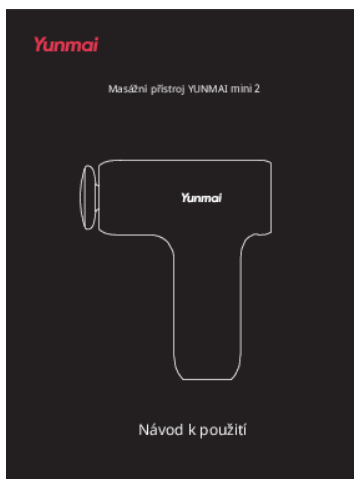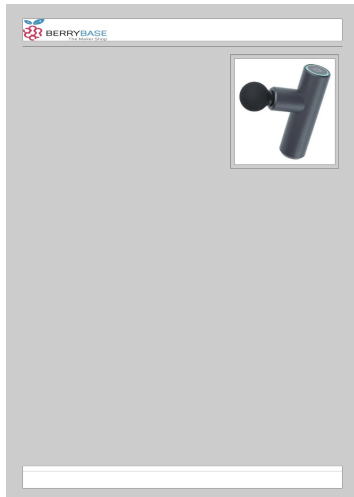
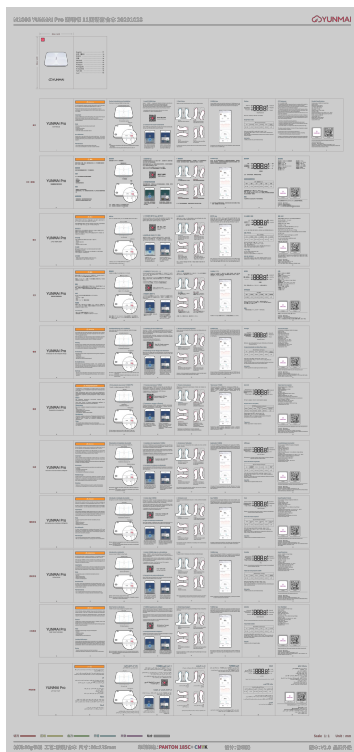
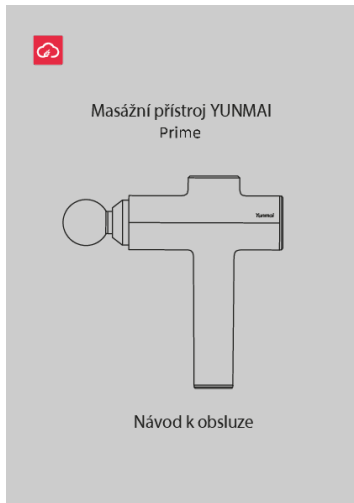
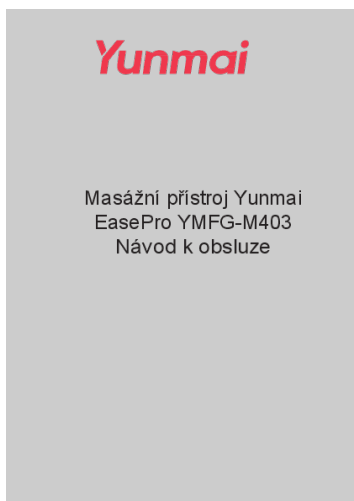
lang:hr **score:18** filesize: 1.12 M page_count: 13 document date: 2024-04-24

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