

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

› [COMFIER](#) /

› [COMFIER Shiatsu Full Back Massager User Manual](#)

COMFIER 2913-Light Gray

COMFIER Shiatsu Full Back Massager User Manual

MODEL: 2913-LIGHT GRAY

1. Product Overview

The COMFIER Shiatsu Full Back Massager is designed to provide deep tissue massage and soothing heat therapy for various parts of your back, including the neck, shoulders, and lower lumbar area. It features advanced 2D/3D shiatsu massage nodes and a vibrating seat cushion to enhance relaxation and comfort. This portable massage pad can be used on most chairs, making it ideal for home or office use.



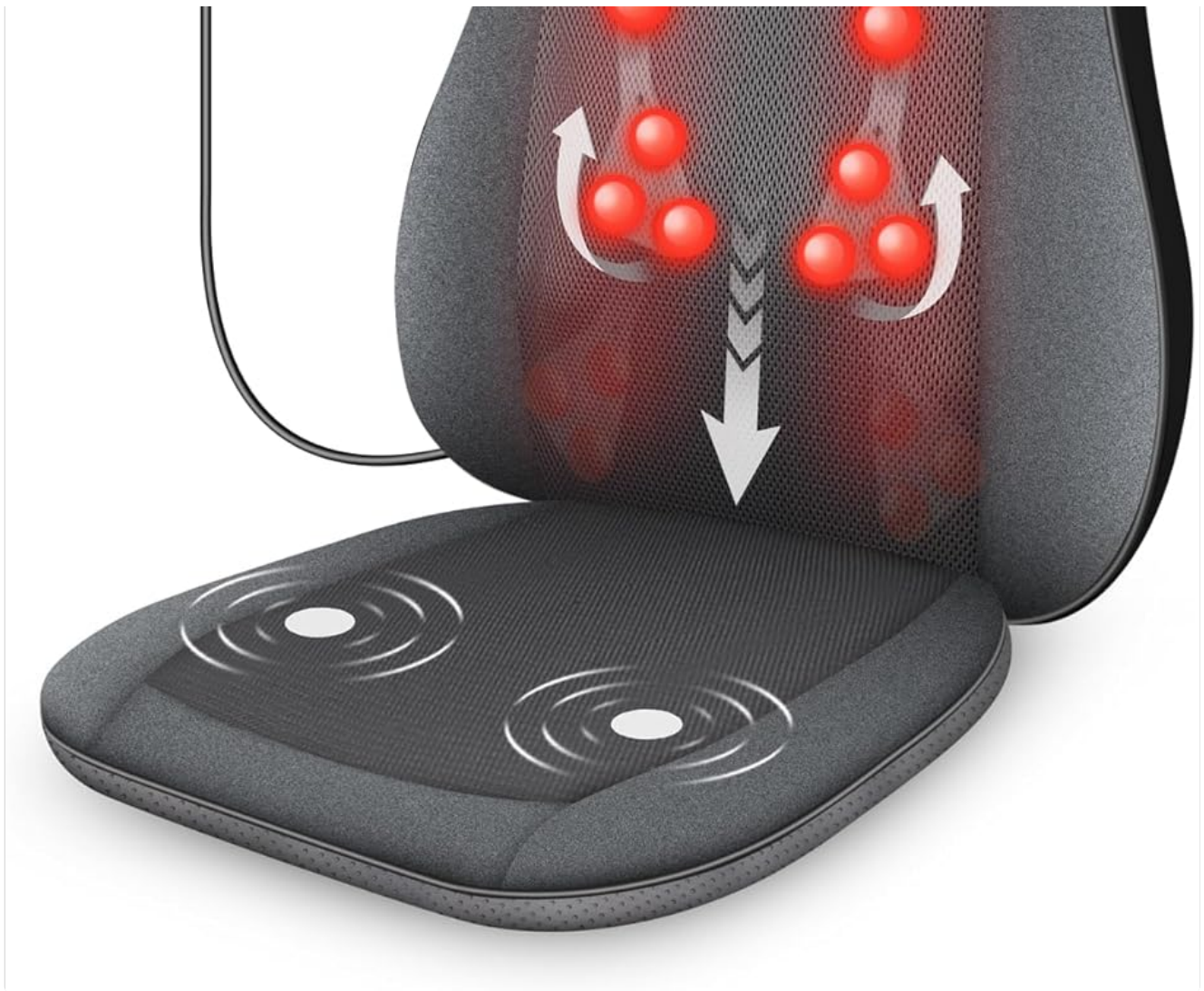


Image: The COMFIER Shiatsu Full Back Massager with its remote control.

2. Safety Information

Please read all instructions carefully before using this product. Failure to follow safety guidelines may result in injury or damage to the device.

- Avoid wearing the massager while sleeping or driving.
- Keep the massager away from water and do not use it in damp environments.
- Consult a doctor before use if you are pregnant, have a medical condition (e.g., heart disease, pacemakers, blood clots), or are recovering from surgery.
- Ensure the power adapter is connected securely and follow all electrical safety guidelines.
- Do not use the massager on swollen, inflamed, or injured areas of the skin.
- Do not use for more than 15 minutes continuously. Allow the device to cool down before another session.
- Keep out of reach of children.

3. Package Contents

Upon opening your COMFIER Shiatsu Full Back Massager package, please ensure all the following items are included:

- 1x COMFIER Shiatsu Full Back Massager with Heat
- 1x User Manual (this document)
- 1x UL Approval Home Adapter
- 1x Packing Box



Image: Package contents including the massager, user manual, and power adapter.

4. Key Features

The COMFIER Shiatsu Full Back Massager offers a range of features designed for optimal comfort and therapeutic benefits:

- **2-IN-1 Shiatsu Massager:** Provides unique 2D/3D shiatsu styles for deep tissue muscle massage with two intensity settings.
- **Full Coverage:** Features a shiatsu roller with 10 massage nodes and an extra-long massage track, covering the entire back from neck and shoulders to lower lumbar.
- **Adjustable Massage Zones:** Choose from Full back, Upper back, or Lower back massage zones.

- **Spot Massage Function:** Allows for pin-point relaxation on specific areas.
- **Rolling Back Massage:** Delivers gentle rolling massage along the spine, with adjustable width for customized comfort.
- **Soothing Heat Therapy:** Optional heating warms your back, lumbar, and shoulders, promoting blood circulation.
- **Vibration Seat Massage:** Three levels of intensity to relax your hips and thighs.

INNOVATIVE 3D/4D FINGER PRESSURE SHIATSU MASSAGE

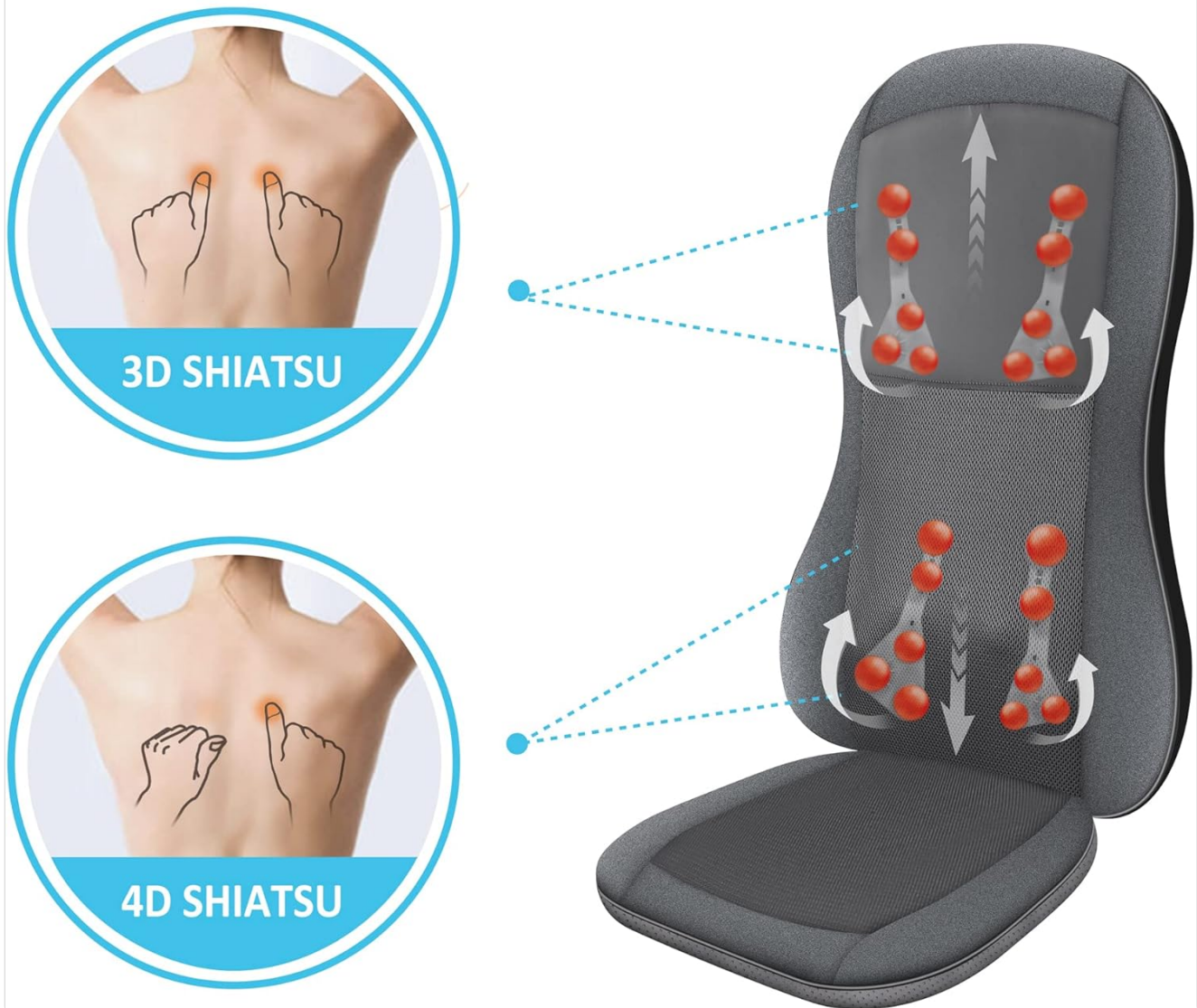


Image: Innovative 3D/4D Finger Pressure Shiatsu Massage.



The image shows a brown armchair with a grey massage pad on the backrest and seat. The pad features red nodes and white arrows indicating massage directions. Four circular callouts on the right show hands massaging specific areas: the shoulder, the upper back, the lower back/waist, and the thigh. A dark grey banner at the bottom contains a yellow circular icon with a white silhouette of a person, followed by the text 'DEEP KNEADING SHIATSU MASSAGE' and a descriptive sentence.

DEEP KNEADING SHIATSU MASSAGE
When Select Lower Back Area, You can Adjust Nodes to Reach Lower Lumbar Area as Expected

Image: Deep Kneading Shiatsu Massage targeting various body parts.



Image: Spot Massage function for targeted relaxation.



SOOTHING HEAT

It delivers gentle heat, not as hot as heating pad

Image: Soothing Heat Therapy for added comfort.

ALL THE BENEFITS

Comfier


2 Massage heads with **10** massage nodes



- ☒ 10 Massage nodes
- ☒ 2D/3D Shiatsu
- ☒ Flexible Massage heads
- ☒ Soothing heat therapy
- ☒ Seat vibration

Others

2 Massage heads with **4** massage nodes

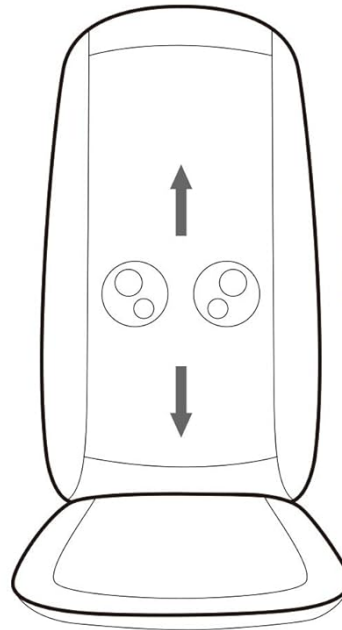

- ☒ 4 Massage nodes
- ☐ 2D/3D Shiatsu
- ☐ Flexible Massage heads
- ☐ Soothing heat therapy
- ☐ Seat vibration

VS

Image: Comparison of COMFIER's features versus other massagers.

5. Setup Instructions

Follow these steps to set up your COMFIER Shiatsu Full Back Massager:

1. **Unpack:** Carefully remove the massager and all accessories from the packaging.
2. **Position:** Place the massager on a suitable chair, sofa, or recliner. Ensure the chair has a backrest tall enough to support the massager. The massager includes straps to secure it to the chair.
3. **Connect Power:** Plug the provided UL approval home adapter into the massager's power input port. Then, plug the adapter into a standard electrical outlet.
4. **Remote Control:** The remote control is attached to the massager for easy access.

Video: Demonstrates setting up and using the COMFIER Back Massager for full back coverage.

6. Operating Instructions

Use the attached remote control to operate your massager. Familiarize yourself with the buttons:

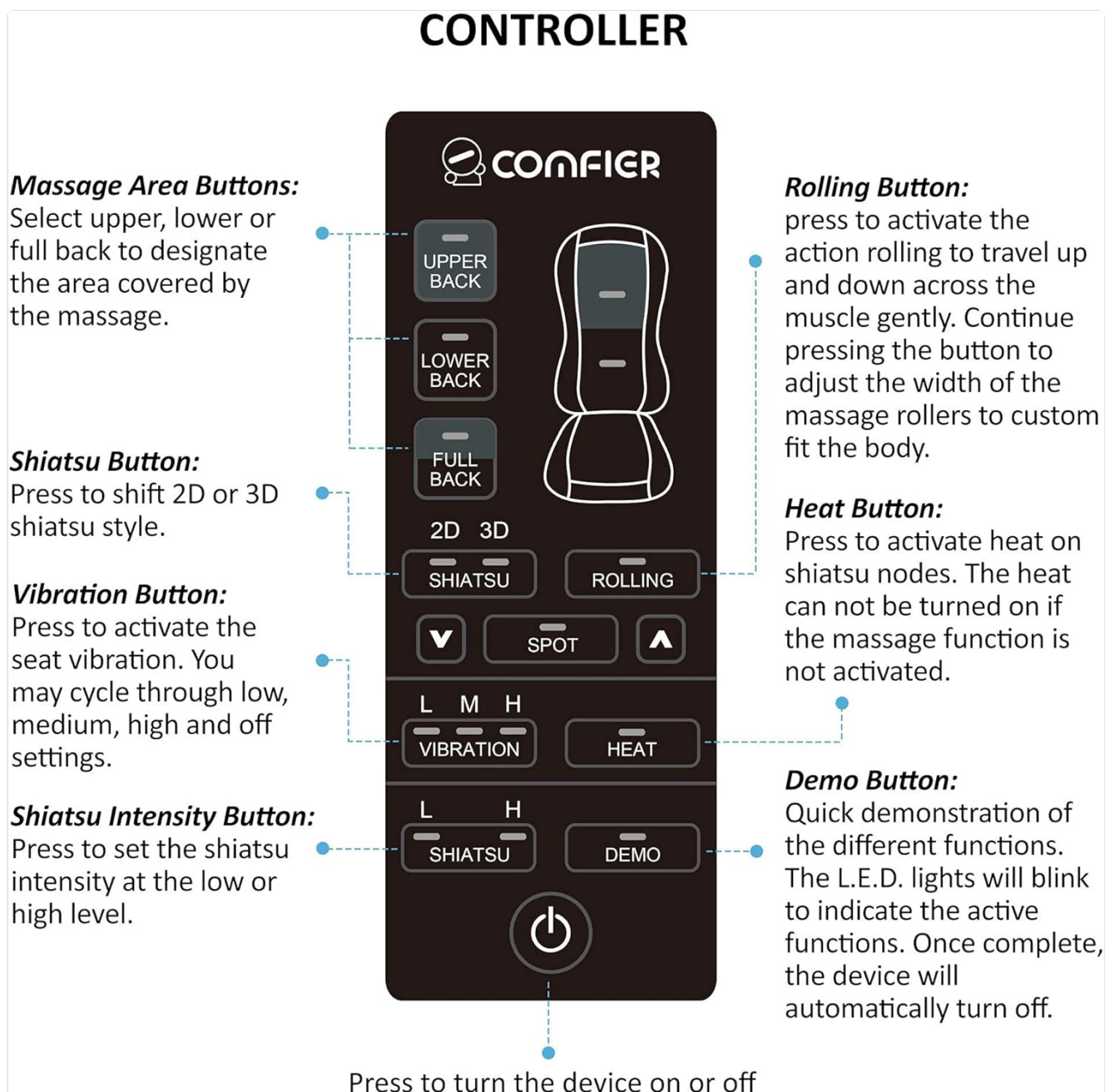


Image: Remote Controller Layout and Functions.

1. **Power On/Off:** Press the **Power Button** to turn the device on or off.
2. **Select Massage Area:** Use the **Massage Area Buttons** (UPPER BACK, LOWER BACK, FULL BACK) to designate the area covered by the massage.
3. **Choose Shiatsu Style:** Press the **Shiatsu Button** (2D, 3D) to shift between 2D or 3D shiatsu styles.
4. **Rolling Massage:** Press the **Rolling Button** to activate the action of rolling up and down the muscle gently. Continue pressing the button to adjust the width of the massage rollers to custom fit your body.
5. **Vibration Seat:** Press the **Vibration Button** (L, M, H, Off) to activate the seat vibration. You may cycle through low, medium, high, and off settings.
6. **Heat Function:** Press the **Heat Button** to activate heat on the shiatsu nodes. Heat cannot be turned on if the massage function is not activated.
7. **Shiatsu Intensity:** Press the **Shiatsu Intensity Button** (L, H) to set the shiatsu intensity at the low or high level.
8. **Demo Mode:** Press the **Demo Button** for a quick demonstration of the different functions. The L.E.D. lights will blink to indicate the active functions. Once complete, the device will automatically turn off.

9. **Automatic Shut-off:** The massager is equipped with an automatic shut-off feature after 15 minutes of continuous use for safety.

7. Maintenance

Proper maintenance ensures the longevity and optimal performance of your massager:

- **Cleaning:** Disconnect the massager from the power outlet before cleaning. Wipe the surface with a soft, slightly damp cloth. Do not use abrasive cleaners, brushes, or strong chemicals.
- **Storage:** Store the massager in a cool, dry place away from direct sunlight and extreme temperatures. Avoid placing heavy objects on the massager.
- **Care:** Do not attempt to repair the massager yourself. Refer to qualified service personnel for any repairs.

8. Troubleshooting

If you encounter any issues with your COMFIER Shiatsu Full Back Massager, please refer to the following common problems and solutions:

Problem	Possible Cause	Solution
Massager does not turn on.	Power adapter not connected properly or no power from outlet.	Check all power connections. Ensure the power outlet is functional by testing with another device.
Heat function is not working.	Heat function is only active when a massage mode is selected.	Activate a shiatsu or rolling massage mode before attempting to turn on the heat.
Massage feels too intense.	High intensity setting or direct contact with nodes.	Switch to a lower intensity setting. Place a towel or blanket between your back and the massager for a softer experience.
Massager stops unexpectedly.	Automatic shut-off activated.	The massager has a 15-minute auto shut-off. Wait a few minutes for the device to cool down, then restart if desired.

If the problem persists after trying these solutions, please contact COMFIER customer support.

9. Product Specifications

Specification	Detail
Brand	COMFIER
Model Number	2913
Product Dimensions	30.31 x 18.5 x 6.3 inches
Item Weight	9.92 Pounds (4500 Grams)

Specification	Detail
Material	Polyester
Power Source	Corded Electric
Massage Type	Shiatsu (2D/3D), Rolling, Vibration
Heat Function	Yes

10. Warranty and Support


COMFIER products are designed for durability and performance. This product comes with a standard manufacturer's warranty. For specific warranty details, please refer to the warranty card included in your product packaging or visit the official COMFIER website.






Customer Support:

- For technical assistance, troubleshooting, or warranty claims, please contact COMFIER customer support.
- You may also consider purchasing extended protection plans for additional coverage.

For the most up-to-date support information, please visit the official COMFIER website or contact their customer service directly.

Related Documents - 2913-Light Gray

<div><div><div>COMFIER</div><div>Shiatsu Neck and Back Massager with Heat</div></div><div></div><div><div>User Manual</div><div>CF-2303-KR</div></div></div>	<p>COMFIER CF-2303-KR Shiatsu Neck and Back Massager with Heat User Manual</p> <p>Comprehensive user manual for the COMFIER CF-2303-KR Shiatsu Neck and Back Massager with Heat, covering safety instructions, features, technical data, setup, operation, controller functions, care, cleaning, troubleshooting, and warranty information.</p>
<div><div><div>COMFIER</div><div>Shiatsu Massage Pillow with Heat</div></div><div></div><div><div>CE</div><div>User Manual</div><div>CF-1228</div></div></div>	<p>Comfier CF-1228 Shiatsu Massage Pillow with Heat User Manual</p> <p>User manual for the Comfier CF-1228 Shiatsu massage pillow with heat, detailing safety instructions, operating procedures, features, maintenance, disposal, technical specifications, and warranty information.</p>

<div><div><div>CONFIER</div><div>Shiatsu Neck and Back Massager</div></div><div><div>User Manual CF-2309A</div></div></div>	<p>Comfier Shiatsu Neck and Back Massager CF-2309A User Manual</p> <p>User manual for the Comfier Shiatsu Neck and Back Massager, model CF-2309A. Includes safety instructions, features, technical data, setup and operation, controller instructions, care and cleaning, troubleshooting, and warranty information.</p>
<div><div><div>CONFIER</div><div>Shiatsu massagecushion mit wärmefunktion DE Shiatsu massage cushion with heat GB Coussin de massage Shiatsu avec chaleur FR Cushion per massage Shiatsu riscaldato IT Cushy massager Shiatsu mit wärmefunktion ES</div></div><div><div>CE User Manual CF-2606MC</div></div></div>	<p>COMFIER CF-2606MC Shiatsu Massage Cushion with Heat User Manual</p> <p>Comprehensive user manual for the COMFIER CF-2606MC Shiatsu massage cushion with heat. Includes safety instructions, operating guide, features, maintenance, and warranty information.</p>
<div><div><div></div><div>User Manual</div></div><div><div>CF-6701 DEEP-KNEADING SHIATSU NECK & SHOULDER MASSAGER</div></div></div>	<p>Comfier CF-6701 Deep-Kneading Shiatsu Neck & Shoulder Massager User Manual</p> <p>User manual for the Comfier CF-6701 Deep-Kneading Shiatsu Neck & Shoulder Massager, covering product features, technical data, setup, operation, safety instructions, care, troubleshooting, and warranty information.</p>
<div><div><div>CONFIER</div><div>NECK & SHOULDER SHIATSU MASSAGER WITH HEAT DE NECK & SHOULDER SHIATSU MASSAGER WITH HEAT GB NECK & SHOULDER SHIATSU MASSAGEUR AVEC CHALEUR FR NECK & SHOULDER SHIATSU MASSAGER WITH HEAT IT NECK & SHOULDER SHIATSU MASSAJER S/ CALORIFIC ES</div></div><div><div>CE User Manual CF-6302</div></div></div>	<p>Comfier CF-6302 Neck & Shoulder Shiatsu Massager with Heat User Manual</p> <p>Comprehensive user manual for the Comfier CF-6302 Neck & Shoulder Shiatsu Massager with Heat. Learn about safety instructions, features, operation, maintenance, and warranty.</p>