

## Oral-B iO9 Power Handle

# Oral-B iO Series 9 Electric Toothbrush Instruction Manual

Model: iO9 Power Handle | Brand: Oral-B

## INTRODUCTION

The Oral-B iO Series Rechargeable Toothbrush is engineered to provide an advanced cleaning experience. It integrates powerful, yet gentle micro-vibrations with Oral-B's unique dentist-inspired round brush head design. This combination aims to deliver a professional clean feeling daily, promoting improved gum health and plaque removal.

Key features include:

- **Advanced Plaque Removal:** Designed to remove 100% more plaque and contribute to 2x healthier gums in just one week.
- **Smart Brushing Experience:** Features an interactive color display for selecting cleaning settings, tracking brushing time, monitoring brush head replacement, and checking battery level.
- **Oral-B App Integration:** Connects to the Oral-B app for real-time tracking of brushing coverage, ensuring no spot is missed.
- **Gum Protection:** Equipped with an iO pressure sensor that signals if you brush too hard, too soft, or just right, protecting gums and enamel.
- **Personalized Cleaning Settings:** Offers 7 distinct cleaning settings to cater to various oral care needs.
- **Convenient Brushing:** Includes a 3-hour fast charger and a charging travel electric toothbrush case for on-the-go use.

## SAFETY INFORMATION

Do not use for purposes other than oral hygiene.

## SETUP

1. Unbox the toothbrush and charging case.
2. Insert the brush head onto the handle.
3. Place the toothbrush on the charging base to fully charge before first use.



Image: Oral-B iO Series 9 Electric Toothbrush components including the handle, multiple brush heads, charging base, and travel case.

## OPERATING INSTRUCTIONS

---

To begin brushing:

1. Press the power button to turn on the toothbrush.
2. Use the mode selection button to cycle through the 7 available brushing modes: Daily Clean, Intense, Sensitive, Super Sensitive, Whiten, Gum Care, and Tongue Clean.
3. Apply the brush head to your teeth. The Smart Pressure Sensor will provide visual feedback: white for too soft, red for too hard, and green for just right.
4. The built-in 2-minute quadrant timer pulsates every 30 seconds to guide you through a complete clean of all areas of your mouth.

# Brush like your Dentist Recommends with iO™

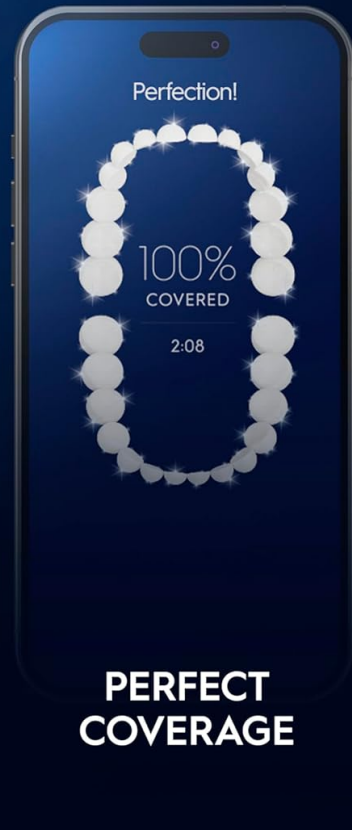


Image: The interactive display of the Oral-B iO Series 9, illustrating feedback for perfect brushing time, optimal pressure, and comprehensive coverage.



Image: Close-up of the interactive color display on the Oral-B iO Series 9 handle, showing battery percentage, current brushing mode (Smart Adapt), and a countdown for brush head replacement.

For enhanced guidance, connect your toothbrush to the Oral-B app via Bluetooth. The app provides 3D Teeth Tracking and AI Recognition to offer real-time feedback and ensure thorough cleaning across all areas of your mouth.

## MAINTENANCE

---

### Brush Head Replacement

For best results and optimal performance, replace your Oral-B iO Series Replacement Brush Heads every 3 months. The interactive display on your toothbrush will also track and remind you when it's time for a replacement.

### Cleaning the Toothbrush

After each use, rinse the brush head thoroughly under running water. Detach the brush head from the handle and rinse both parts separately. Wipe the handle dry with a clean cloth. Ensure the charging base is kept clean and dry.

## TROUBLESHOOTING

---

### Battery Not Charging

Ensure the charging base is properly plugged into a working electrical outlet. Place the toothbrush securely on the magnetic charging base. The light indicator on the toothbrush should illuminate to confirm charging. If the issue persists, try a different outlet.

### App Connectivity Issues

Ensure Bluetooth is enabled on your smartphone and the Oral-B app is updated to the latest version. Restart both your toothbrush and your smartphone, then attempt to reconnect through the app.

### Pressure Sensor Indicator

If the pressure sensor consistently shows red, you may be brushing too hard. Adjust your brushing technique to apply gentler pressure until the indicator turns green. If it consistently shows white, apply slightly more pressure.

### SPECIFICATIONS

Product Dimensions	7.01 x 3.94 x 9.96 inches
Item Weight	1.39 Pounds
Item Model Number	iO M9.4A1.1A
Batteries	1 Lithium Ion batteries required. (included)
UPC	069055128893
Manufacturer	Procter & Gamble
Brand	Oral-B
Age Range (Description)	Adult
Power Source	Corded Electric
Item Firmness Description	Soft
Color	Whitealabaster

### WARRANTY AND SUPPORT

For warranty information, technical support, or to purchase replacement parts, please refer to the official Oral-B website or contact their customer service directly. Keep your proof of purchase for warranty claims.