

Oral-B iO M7.1A1.1DH

# Oral-B iO Series 7G Electric Toothbrush User Manual

Model: iO M7.1A1.1DH | Brand: Oral-B

## INTRODUCTION

The Oral-B iO Series 7G Electric Toothbrush combines advanced technology with a dentist-inspired round brush head for an effective and gentle cleaning experience. This manual provides essential information for the proper setup, operation, and maintenance of your new electric toothbrush.



Image: The Oral-B iO Series 7G Electric Toothbrush, including the handle, brush head, magnetic charger, and travel case, presented in its packaging.

## IMPORTANT SAFETY INFORMATION

---

Please read all instructions carefully before using this product. Keep this manual for future reference.

- Do not use this appliance for purposes other than oral hygiene.
- Supervise young children to avoid choking hazards from small parts or brush heads.
- If the product is dropped, the brush head should be replaced before the next use, even if no visible damage is apparent.
- Do not modify or repair the product. This may result in fire, electric shock, or injury. Consult your dealer for repairs.
- Do not place the charger in water or liquid.
- This product contains batteries that are not replaceable.

## WHAT'S INCLUDED

---

Your Oral-B iO Series 7G package typically includes the following components:

- Oral-B iO Series 7G Toothbrush Handle
- Oral-B iO Ultimate Clean Replacement Brush Head
- Magnetic Charging Stand
- Travel Case



Image: A top-down view of the Oral-B iO Series 7G Electric Toothbrush components: the handle, a brush head, the magnetic charging base, and a soft travel pouch.

## SETUP

---

1. **Initial Charge:** Place the toothbrush handle on the magnetic charging stand. The display will show the battery level. A full charge typically takes approximately 3 hours.
2. **Attach Brush Head:** Align the brush head with the shaft on the toothbrush handle and push it down until it clicks into place.
3. **Download Oral-B App (Optional but Recommended):** For enhanced features like AI recognition and real-time tracking, download the Oral-B app from your device's app store and follow the in-app instructions to pair your toothbrush via Bluetooth.

# Experience Oral-B's smart brushing



Image: The Oral-B iO Series 7G Electric Toothbrush handle on its magnetic charging base, alongside an open travel case. The image highlights the fast charging capability.

## OPERATING INSTRUCTIONS

### Powering On/Off and Mode Selection

1. **Power On:** Press the power button ( ) to turn on the toothbrush.
2. **Select Mode:** Use the mode button (●) to cycle through the available cleaning modes. The iO Series 7G offers 5 modes:
  - Daily Clean
  - Sensitive
  - Intense
  - Whitening
  - Gum Care
3. **Start Brushing:** Once your desired mode is selected, press the power button again to begin brushing.
4. **Power Off:** Press and hold the power button to turn off the toothbrush.



# Personalize your clean

with 5 cleaning modes



Image: The interactive display of the Oral-B iO Series 7G toothbrush, illustrating the five available cleaning modes: Daily Clean, Intense, Sensitive, Gum Care, and Whiten.

## Smart Pressure Sensor

The toothbrush features an intelligent pressure sensor that provides real-time feedback to protect your gums and enamel:

- **Green Light:** Indicates optimal brushing pressure.
- **Red Light:** Signals that you are brushing too hard. Reduce pressure immediately.
- **Other Colors (e.g., Orange/White):** Suggests insufficient pressure. Apply slightly more pressure for an effective clean.

# Protect your gums



Image: A person brushing their teeth with the Oral-B iO Series 7G, showing the green light ring on the toothbrush handle, indicating "Perfect Pressure" for effective and safe cleaning.

## 2-Minute Timer and Quadrant Pacer

The iO Series 7G includes a built-in 2-minute timer to ensure you brush for the dentist-recommended duration. It also features a quadrant pacer that pulsates every 30 seconds, prompting you to move to the next section of your mouth for a thorough clean.

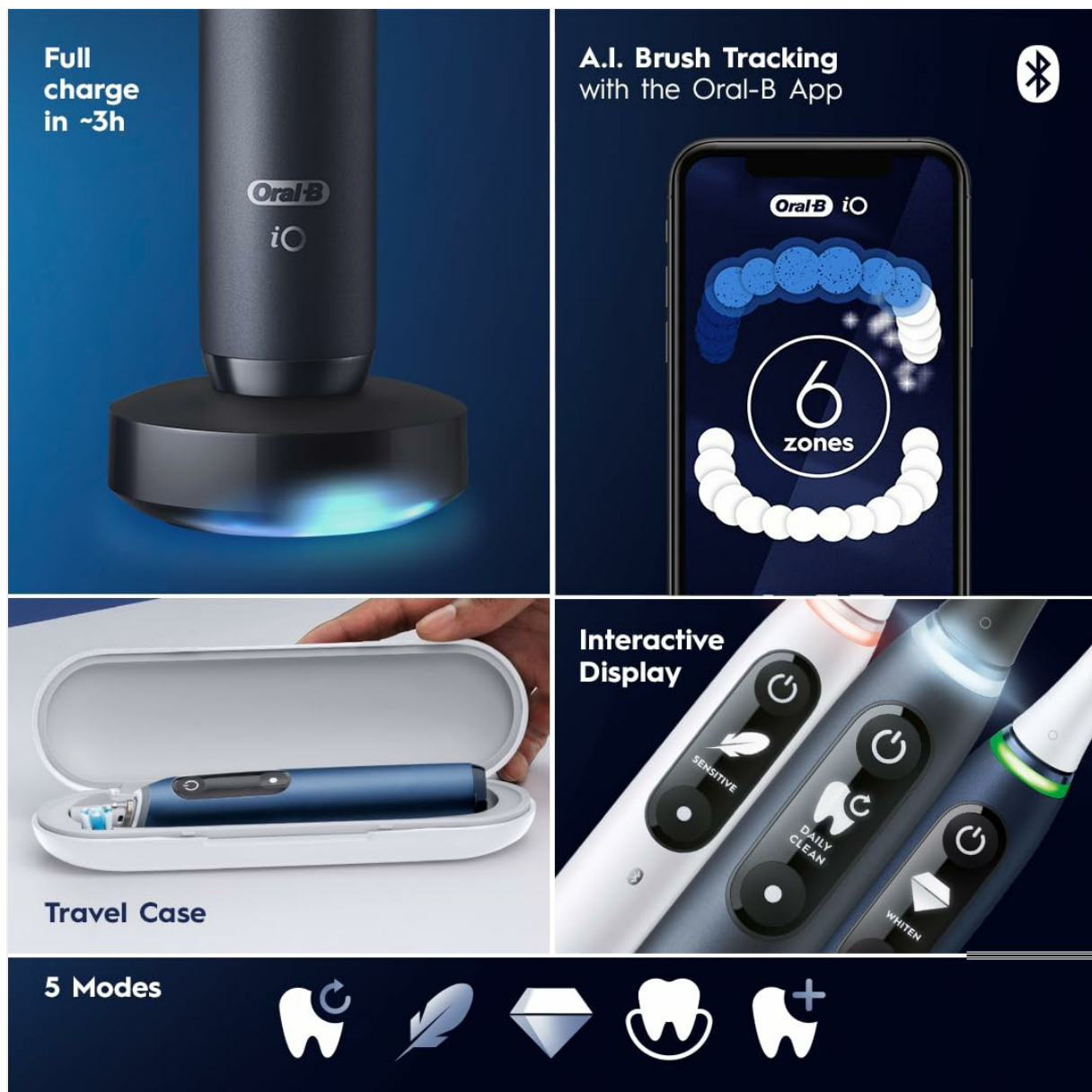


Image: A close-up of the Oral-B iO Series 7G display, highlighting smart brushing features such as 5 cleaning modes, a 2-minute timer, a brush head change alert, and a battery level indicator.

## Oral-B App Integration (AI Tracking)

Connect your iO Series 7G to the Oral-B app via Bluetooth for an enhanced brushing experience. The app provides:

- **Real-time AI Recognition:** Tracks how and where you brush in real-time, guiding you to cover all areas of your mouth.
- **Session Score:** Provides a score for each brushing session based on coverage, brush time, and pressure.
- **Personalized Guidance:** Offers insights and recommendations to improve your brushing habits.
- **Brush Head Replacement Reminders:** Alerts you when it's time to replace your brush head for optimal performance.

# Become a better brusher

With iO Artificial Intelligence



Image: The Oral-B iO Series 7G toothbrush next to a smartphone displaying the Oral-B app, which shows real-time brushing feedback and AI-powered tracking across 6 zones of the mouth.

*Note: No official product videos from the seller were available for embedding based on the provided data.*

## MAINTENANCE

### Cleaning Your Toothbrush

- After each use, rinse the brush head thoroughly under running water with the handle switched on.
- Remove the brush head and clean both the handle and the brush head separately under running water.
- Wipe the handle dry with a soft cloth.
- Clean the charging stand regularly with a damp cloth. Ensure it is unplugged before cleaning.

### Replacing Brush Heads

Dentists recommend replacing your brush head every 3 months, or sooner if the bristles are splayed or worn. The Oral-B iO Series 7G uses specific Oral-B iO Ultimate Clean replacement brush heads. Other Oral-B brush heads are not compatible.





Image: A single Oral-B iO Series replacement brush head, designed for optimal cleaning performance.

### Battery Care

To maintain battery life, avoid leaving the toothbrush on the charger continuously when not in use. Charge it when the battery indicator shows a low level.

### TROUBLESHOOTING

Problem	Possible Cause / Solution
Toothbrush does not turn on.	Battery may be depleted. Place on charger for at least 30 minutes. Ensure charger is plugged in correctly.
Battery does not hold charge for long.	Battery performance can decrease over time. Ensure full charges are completed. Avoid continuous charging when not in use.
Pressure sensor light is always red.	You may be applying too much pressure. Adjust your brushing technique to maintain a green light.

Problem	Possible Cause / Solution
Oral-B App not connecting or tracking correctly.	Ensure Bluetooth is enabled on your device and the app is open during brushing. Restart the app and/or toothbrush. Check for app updates.
Brush head feels loose.	Ensure the brush head is fully pushed down onto the shaft until it clicks. If it remains loose, the brush head may be worn or damaged and needs replacement.

## SPECIFICATIONS


- **Model Number:** iO M7.1A1.1DH
- **Product Dimensions:** 7.01 x 3.94 x 9.96 inches
- **Item Weight:** 1.08 Pounds
- **Power Source:** Corded Electric (via magnetic charging stand)
- **Battery Type:** 1 Lithium Ion battery (included)
- **Brush Head Firmness:** Soft
- **Color:** White Alabaster
- **UPC:** 069055128831

## WARRANTY AND SUPPORT

For warranty information, product registration, or technical support, please refer to the official Oral-B website or contact Oral-B customer service. Keep your proof of purchase for warranty claims.

You can visit the official Oral-B store for more information:[Oral-B Store](#)

### Related Documents - iO M7.1A1.1DH

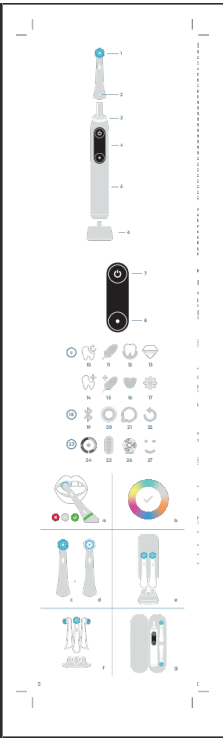

[Oral-B iO Electric Toothbrush User Manual](#)

Comprehensive user manual for the Oral-B iO electric toothbrush, covering setup, usage, features, maintenance, troubleshooting, and warranty information.



[Oral-B Vitality Electric Toothbrush User Manual](#)

Comprehensive user manual for the Oral-B Vitality electric toothbrush, covering features, usage instructions, safety guidelines, brush head information, cleaning, and warranty. Includes model information for Charger Type 3757 and Handle Type 3710.



**Customer Service**

**DE** 00 800 27 28 64 63  
**DE** 05 44 - 88 40 10  
**DE** 0 800 731 1792  
**DE** 1 800 003 448  
**DE** 0 800 844 803  
(service at Apple  
stores)  
**DE** 0 800 14 592  
**IT** 800 21172  
**ES** 900 814 208  
**PT** 808 20 00 33  
**IT** 800 124 600  
**NL** 0800 445 53 88  
**DK** 70 15 00 13  
**SE** 22 63 00 93  
**SE** 020 - 21 33 21  
**FI** 020 377 877

**Internet:**  
www.orb.com  
www.braun.com  
www.service.braun.com

Braun GmbH  
Friedrichstraße 145  
63476 Kronberg, Germany

DELUXEFLY FÜR KINDER  
1070046-01-0

[Oral-B Kids Electric Toothbrush - User Manual](#)

Comprehensive user manual for the Oral-B Kids electric toothbrush, powered by Braun. Includes safety instructions, charging, usage, modes, cleaning, and warranty information for children's oral hygiene.





© 2015 P&G



1-800-775-1752  
1-800-900-6466  
1-800-200-20-00

© 2015 P&G  
All rights reserved.  
P&G, the P&G logo, and the Oral-B logo are trademarks of P&G.

Oral-B  
P&G Consumer Goods  
23715 Northpark Drive  
Beverly Hills, CA 90210

Charger Type 3757  
Handle Type 3767

© 2015 P&G

### [Oral-B SMART/JUNIOR Electric Toothbrush User Manual and Guide](#)

Comprehensive user manual for the Oral-B SMART/JUNIOR electric toothbrush, covering setup, operation, features, safety, troubleshooting, and warranty information. Learn how to use your smart toothbrush effectively.



© 2015 P&G



1-800-775-1752  
1-800-900-6466  
1-800-200-20-00  
1-800-900-6466  
1-800-900-6466  
1-800-900-6466

© 2015 P&G  
All rights reserved.  
P&G, the P&G logo, and the Oral-B logo are trademarks of P&G.

Oral-B  
P&G Consumer Goods  
23715 Northpark Drive  
Beverly Hills, CA 90210

Charger Type 3757  
Handle Type 3767

© 2015 P&G

### [Oral-B SMART/Teen Electric Toothbrush User Manual and Guide](#)

Comprehensive user manual for the Oral-B SMART/Teen electric toothbrush (Charger Type 3757, Handle Type 3767). Learn about features, usage, charging, cleaning, troubleshooting, and warranty.



© 2015 P&G

### [Oral-B GENIUS 8000-9000 Electric Toothbrush User Manual](#)

User manual for the Oral-B GENIUS 8000-9000 electric toothbrush, covering operation, features, app connectivity, maintenance, and troubleshooting.