

mornwell F5106-F1

Mornwell Water Flosser User Manual

Model: F5106-F1

Brand: mornwell

PRODUCT OVERVIEW

The Mornwell Water Flosser is a portable dental oral irrigator designed to enhance your oral hygiene routine. It features a 330ml water tank, three cleaning modes, and a 360° rotatable nozzle to effectively clean teeth, braces, and bridges, removing food debris and plaque.



Image: The Mornwell Water Flosser, its packaging, various nozzles, and a USB charging cable.

IMPORTANT SAFETY INFORMATION

- Do not immerse the main unit in water. The device is IPX7 waterproof, but the charging port should remain dry during charging.
- Ensure the base cover is tightly attached to maintain waterproof integrity during use.
- Keep out of reach of children.
- Do not use if the device or charging cable is damaged.
- Use only clean water or dentist-recommended oral rinse solutions.
- Consult your dentist before use if you have severe periodontal disease or have undergone oral surgery recently.
- It is recommended to use the Soft mode for the first use to adapt to the water pressure.

IPX7 Waterproof



Image: The water flosser demonstrating its IPX7 waterproof capability, with a note to ensure the base cover is tightly attached.

PACKAGE CONTENTS

Upon opening your Mornwell Water Flosser package, you should find the following items:

- 1 x Cordless Water Flosser Handle
- 1 x Standard Nozzle
- 1 x Orthodontic Nozzle
- 1 x Pocket Nozzle
- 1 x Tongue Scraper
- 1 x USB Charging Cable
- 1 x Storage Bag
- 1 x User Manual (this document)

You May Receive

NEW Nozzle configuration

Standard Nozzle*3
Orthodontic nozzle*1



OLD Nozzle configuration

Standard Nozzle*3
Tongue Cleaner Nozzle*1



Image: A visual representation of the nozzles included with the water flosser, highlighting the new nozzle configuration.

SETUP

Charging the Device

Before first use, fully charge the water flosser. The built-in powerful lithium battery can last up to 14 days on a full charge. Use the provided USB cable to connect the device to any power source with a USB port, such as a 5V power adapter, computer notebook, car charger, or portable battery.

Powerful Lithium Battery & Universal USB 5V Charging



Image: Illustration of the versatile USB charging capabilities of the water flosser.

Filling the Water Tank

The water flosser features a 300ml large capacity detachable water tank. To fill, detach the tank from the main unit by twisting or pulling down (refer to specific model instructions if different). Fill the tank with lukewarm water or a dentist-recommended oral rinse. Reattach the tank securely to prevent leaks.



Image: The water flosser with its large, detachable 300ml water tank.

Attaching Nozzles

Select the desired nozzle (Standard, Orthodontic, Pocket, or Tongue Scraper) and insert it firmly into the opening at the top of the flosser handle until it clicks into place. The nozzle is 360° rotatable for easy access to all areas of your mouth.

OPERATING INSTRUCTIONS

Selecting a Mode

The Mornwell Water Flosser offers three cleaning modes to suit various oral care needs. Press the 'Mode' button to cycle through them:

- **Normal Mode:** Designed for experienced users, providing powerful cleaning.
- **Soft Mode:** Recommended for first-time users and individuals with sensitive teeth or gums.
- **Pulse Mode:** Ideal for teeth and gum massage, promoting gum health.

The device features a memory function that records the last used mode, so you don't need to adjust it each time.



Image: Explanation of the three distinct cleaning modes available on the water flosser.

Using the Water Flosser

Once the tank is filled and the nozzle is attached, place the nozzle tip into your mouth, aiming it at the gumline at a 90-degree angle. Lean over a sink to avoid splashing. Press the 'On/Off' button to start the water flow. Guide the water jet along your gumline and between your teeth, pausing briefly between each tooth. The 360° rotatable nozzle allows you to reach all areas effectively. Continue until you have cleaned all areas of your mouth. The device is highly effective at removing 99.99% of food debris and plaque.

Effectively Remove up to **99.99%** Plaque




Mornwell

0.5mm
Water column

1400-1800
Pulsations
per minute

30-120 psi

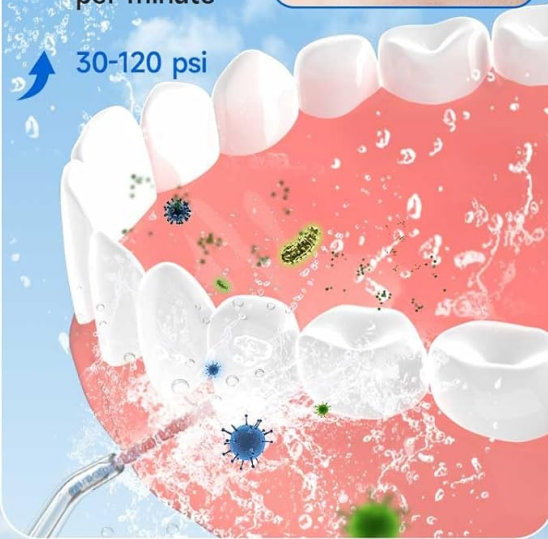

others

0.66mm
Water column

1200
Pulsations
per minute

20-105 psi

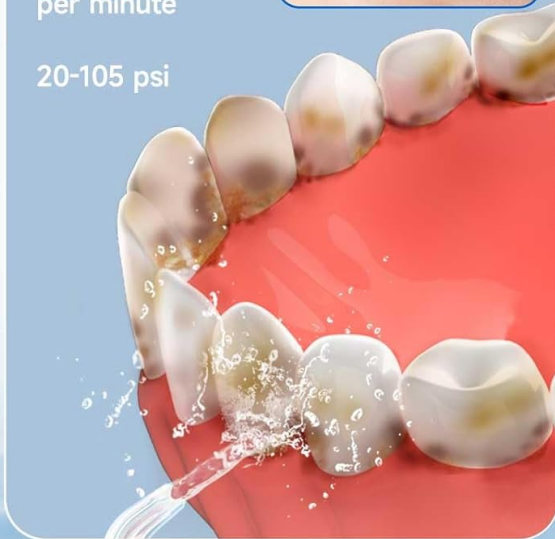


Image: A comparison demonstrating the Mornwell water flosser's effectiveness in plaque removal.

Solve different oral problems



Bad breath



Decay



Bleeding gums



Tooth plant



Dental plaque



Teeth stain



Image: Common oral problems that can be addressed by using the water flosser.

MAINTENANCE AND CLEANING

Cleaning the Water Tank

After each use, empty any remaining water from the tank. Detach the tank and rinse it thoroughly with clean water. For a deeper clean, you can use a mild soap solution and then rinse completely. Allow the tank to air dry before reattaching or storing.

Cleaning Nozzles

Remove the nozzle after each use and rinse it under running water. Periodically, you can soak the nozzles in a solution of 2 parts water to 1 part white vinegar for 5-10 minutes to remove mineral deposits, then rinse thoroughly.

Cleaning the Main Unit

Wipe the exterior of the main unit with a soft, damp cloth. Do not use abrasive cleaners or immerse the main unit in water.

Storage

Store the water flosser in a cool, dry place away from direct sunlight. Ensure the water tank is empty and dry before storing for extended periods. Use the provided storage bag for convenient and hygienic travel or storage.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery; Device not charged.	Charge the device fully using the provided USB cable.
Weak water pressure.	Water tank is empty or low; Nozzle is clogged; Incorrect mode selected.	Refill the water tank; Clean the nozzle; Switch to Normal or Pulse mode.
Water leakage.	Water tank not properly attached; Base cover not tightly closed.	Ensure the water tank is securely clicked into place; Verify the base cover is tightly sealed.
Device stops working during use.	Battery depleted; Overheating.	Recharge the device; Allow the device to cool down for a few minutes before resuming use.

SPECIFICATIONS

Feature	Detail
Model Number	F5106-F1
Product Dimensions	4.21 x 2.64 x 9.45 inches
Item Weight	1.08 Pounds
Water Tank Capacity	330ml
Cleaning Modes	Normal, Soft, Pulse
Nozzle Rotation	360° Rotatable
Battery Type	Lithium Ion (1 included)
Charging Method	USB Rechargeable
Waterproof Rating	IPX7
UPC	639510833923

WARRANTY AND SUPPORT

Mornwell electric water flossers come with a **30-day money-back guarantee** and a **12-month warranty**. If you have any questions or encounter issues with your device, please do not hesitate to contact Mornwell customer support. We aim to respond to all inquiries within 24 hours.

For support, please visit the official Mornwell website or refer to the contact information provided in your

product packaging.