

## realink B08B39K9LS

# realink Cordless Smart Neck Massager with Heat: User Manual

Model: B08B39K9LS

## 1. INTRODUCTION

---

The realink Cordless Smart Neck Massager with Heat is designed to provide targeted massage and heat therapy to the neck area. This device features multiple massage modes and intensity levels, along with a constant temperature heat function, to help alleviate muscle tension and promote relaxation. Its cordless and portable design allows for convenient use in various settings.

## 2. PRODUCT FEATURES

---

- **5 Massage Modes:** Offers a variety of massage techniques including Scraping, Acupuncture, The massage, Pulse, Cupping, and Taps.
- **15 Intensity Levels:** Adjustable intensity to suit individual preferences.
- **Constant Temperature Heat:** Provides a soothing 107°F (41.6°C) heat fomentation for enhanced muscle relaxation.
- **Cordless & Portable Design:** Lightweight and U-shaped for comfortable wear and use on the go.
- **360° Rotary Fitting Arm:** Ensures optimal contact and fit around the neck.
- **304 Stainless Steel Contact Sheets:** Durable and skin-friendly electrodes.
- **Remote Control:** Easy adjustment of modes, intensity, and heat settings.


# 5 MODES

## Simulate 6 Traditional Massages

### Remote Control



Figure 2.1: Overview of the realink Neck Massager highlighting its 5 massage modes: Scraping, Acupuncture, The massage, Pulse, Cupping, and Taps. The device also features remote control.

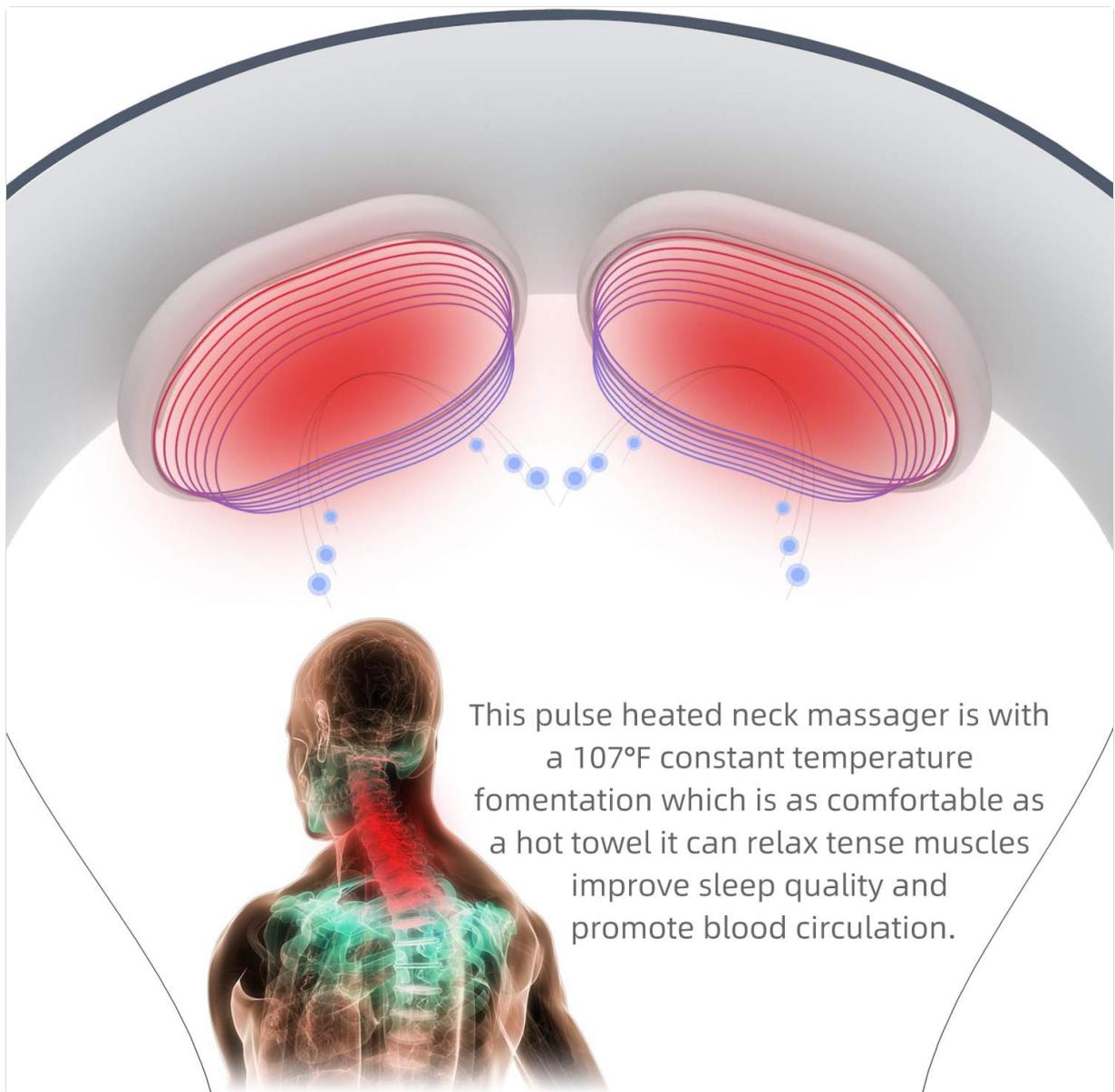

This pulse heated neck massager is with a 107°F constant temperature fomentation which is as comfortable as a hot towel it can relax tense muscles improve sleep quality and promote blood circulation.

Figure 2.2: Diagram illustrating the pulse heated neck massager's 107°F constant temperature fomentation, designed to relax tense muscles, improve sleep quality, and promote blood circulation.



Figure 2.3: Detailed view of the realink Neck Massager components, including Neck Protection, 360° Rotary Fitting Arm, Strengthened Elastic Arm, Thermostat Module, 304 Stainless Steel Contact Sheet, Micro-USB Charging Line, Remote Controller, and Fuselage Switch.



Figure 2.4: A person wearing the realink Neck Massager while working at a computer, demonstrating its wire-free and portable design for pulse massage.



Figure 2.5: A person wearing the realink Neck Massager while exercising outdoors, showcasing its breathable and fashionable design.

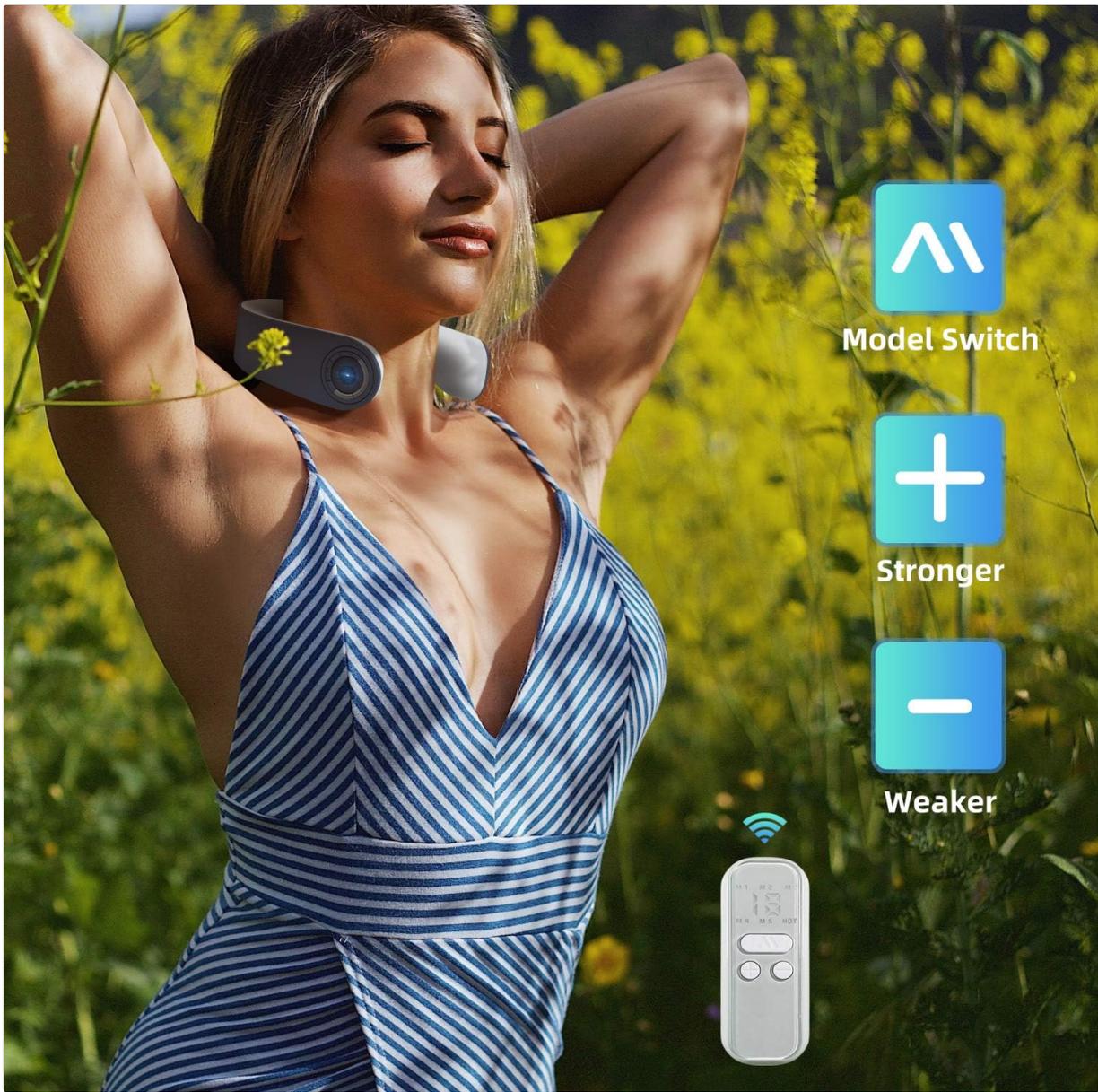


Figure 2.6: A person relaxing with the realink Neck Massager, with an overlay showing the remote control buttons for Model Switch, Stronger intensity (+), and Weaker intensity (-).

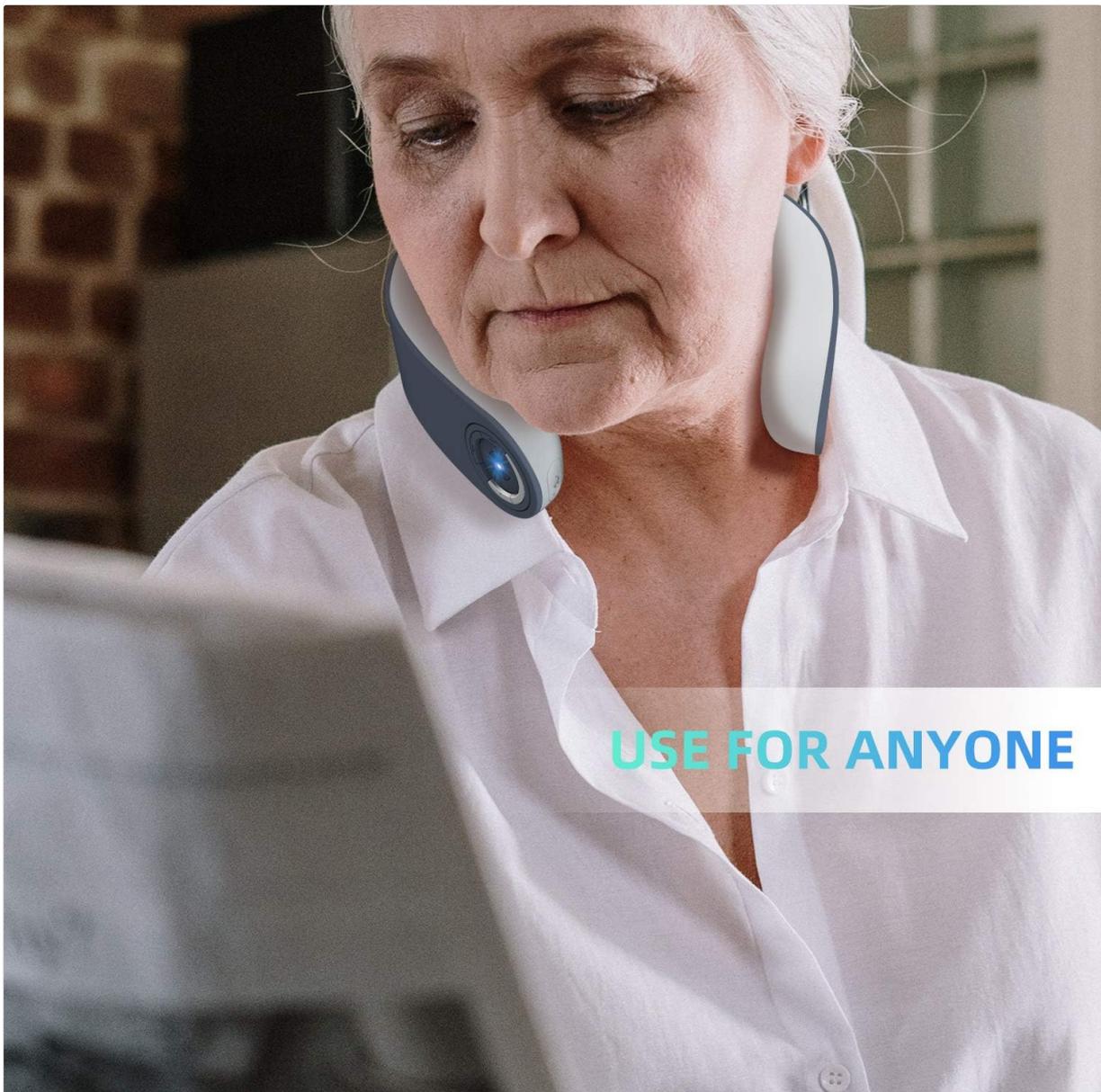


Figure 2.7: An elderly woman wearing the realink Neck Massager, illustrating its suitability for a wide range of users.

### 3. SETUP

---

1. **Charge the Device:** Before first use, fully charge the massager using the provided Micro-USB charging cable. The charging indicator will show the charging status.
2. **Prepare Your Neck:** For optimal contact and sensation, it is recommended to wipe your neck with a damp cloth or apply a small amount of moisturizer to the skin where the electrodes will make contact. This helps improve conductivity.

### 4. OPERATING INSTRUCTIONS

---

1. **Wear the Massager:** Place the massager around your neck, ensuring the 304 stainless steel contact sheets are firmly against your skin. Adjust the 360° rotary fitting arm for a comfortable and secure fit.
2. **Power On:** Press and hold the Fuselage Switch button on the device to turn it on. An indicator light will illuminate.
3. **Use the Remote Control:** The remote control allows you to adjust settings conveniently.
  - **Mode Selection (M):** Press the 'M' button to cycle through the 5 available massage modes.
  - **Intensity Adjustment (+/-):** Use the '+' and '-' buttons to increase or decrease the massage intensity

(15 levels). Start with a low intensity and gradually increase to a comfortable level.

- **Heat Function (HOT):** Press the 'HOT' button to activate or deactivate the heat function. The heat will reach a constant temperature of approximately 107°F (41.6°C).

4. **Session Duration:** The device typically has an automatic shut-off feature after a set period (e.g., 15 minutes) to prevent overuse. Refer to the device's specific timing if available.

5. **Power Off:** Press and hold the Fuselage Switch button on the device to turn it off when finished.

## 5. MAINTENANCE

- **Cleaning:** After each use, gently wipe the contact sheets and the device body with a soft, slightly damp cloth. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the massager in a cool, dry place, away from direct sunlight and extreme temperatures. Keep it in its original packaging or a protective case to prevent damage.
- **Battery Care:** To prolong battery life, avoid fully discharging the device frequently. Charge it regularly, even if not in use for extended periods.

## 6. TROUBLESHOOTING

- **Device Not Turning On:** Ensure the device is fully charged. Press and hold the power button firmly. If the issue persists, contact customer support.
- **No Sensation or Weak Sensation:**
  - Ensure the contact sheets are making firm and even contact with your skin.
  - Increase the intensity level using the remote control.
  - Wipe your neck with a damp cloth or apply a small amount of moisturizer to improve conductivity.
- **Uncomfortable Sensation (e.g., 'Electro Shock' feeling):**
  - Reduce the intensity level immediately.
  - Ensure the device is properly positioned and not pressing too hard on any specific point.
  - Ensure your skin is adequately moisturized.
  - If discomfort persists, discontinue use and consult a healthcare professional.
- **Remote Control Not Working:** Check the remote control's battery. Ensure there are no obstructions between the remote and the massager.

## 7. SPECIFICATIONS

Feature	Detail
Product Dimensions	5.59 x 5.91 x 1.97 inches
Weight	1 Pound
Material	Metal, Silicone
Use For	Neck
Special Feature	Heat
Product Benefits	Blood Flow Control
Manufacturer	realink

ASIN	B08B39K9LS
------	------------

## 8. SAFETY INFORMATION

---

Please read all safety warnings and instructions before using this product to avoid injury or damage.

- Consult your doctor before using this device if you are pregnant, have a pacemaker, suffer from heart disease, epilepsy, or any other medical condition.
- Do not use on broken, irritated, or sensitive skin.
- Keep out of reach of children.
- Do not use while sleeping or in situations where you might fall asleep.
- Avoid prolonged use on a single area to prevent skin irritation.
- Do not use if the device is damaged or malfunctioning.
- This device is not intended for medical use. It is intended only to provide a soothing massage.

## 9. WARRANTY AND SUPPORT

---

For warranty information, technical support, or customer service inquiries, please refer to the product packaging or contact the manufacturer directly through their official website or customer service channels. Please retain your proof of purchase for any warranty claims.