

MOBI 70113

MOBI Non-Contact Forehead Digital Thermometer User Manual

Model: 70113

1. INTRODUCTION

The MOBI Non-Contact Forehead Digital Thermometer is designed for safe, hygienic, and accurate temperature measurement. It utilizes infrared technology to provide readings without physical contact, making it suitable for all ages. This device can measure both body temperature and object/room temperature, offering versatility for various household needs.



Image: The MOBI Non-Contact Forehead Digital Thermometer alongside its retail packaging, highlighting its non-contact feature and digital display.

2. PRODUCT COMPONENTS AND FEATURES

Familiarize yourself with the main components and features of your thermometer:

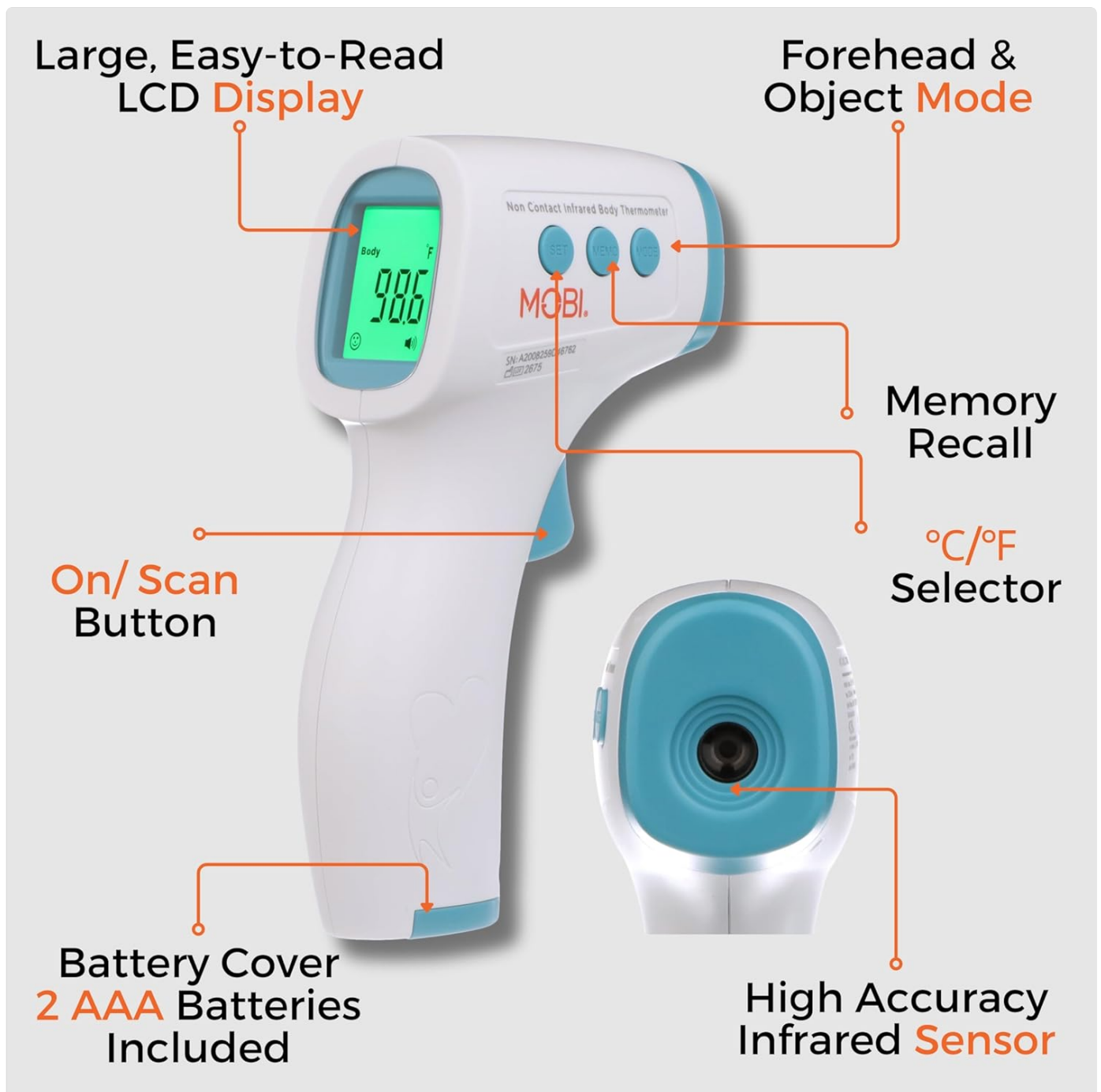


Image: A detailed diagram illustrating the various parts of the MOBI thermometer, including the LCD display, On/Scan button, Memory Recall, °C/°F Selector, Battery Cover, and High Accuracy Infrared Sensor.

- **Large, Easy-to-Read LCD Display:** Shows temperature readings clearly.
- **On/Scan Button:** Activates the device and initiates temperature measurement.
- **SET Button:** Used for accessing settings.
- **MEMO Button:** Accesses memory recall for previous readings.
- **MODE Button:** Switches between Body and Object temperature modes.
- **°C/°F Selector:** Allows switching between Celsius and Fahrenheit units.
- **High Accuracy Infrared Sensor:** Measures temperature from a distance.
- **Battery Cover:** Protects the battery compartment.

3. SETUP

3.1. Installing Batteries

Your thermometer requires 2 AAA batteries, which are included. To install:

1. Locate the battery cover at the bottom of the handle.
2. Slide the battery cover open.
3. Insert the 2 AAA batteries, ensuring correct polarity (+/-).
4. Close the battery cover securely.

3.2. Selecting Temperature Unit (°C/°F)

To switch between Fahrenheit (°F) and Celsius (°C):

- With the thermometer powered on, press the dedicated °C/°F Selector button (often labeled with the unit symbols or integrated with the MODE button).
- The display will change to reflect the selected unit.

4. OPERATING INSTRUCTIONS

4.1. Taking Body Temperature

Ensure the thermometer is in 'Body' mode (indicated on the display).

1. Point the thermometer sensor towards the center of the forehead, maintaining a distance of 1 to 2 inches (2.5 to 5 cm).
2. Press the On/Scan button.
3. The temperature reading will appear on the LCD display within one second.
4. The display backlight color will indicate the temperature status: Green for normal, Orange for elevated, and Red for high fever.

With Simple and Advanced Features

Infrared Sensor



Memory Recall of last 50 Reading



Fever Indicators
Green/Orange/Red



Image: A person demonstrating the use of the MOBI Non-Contact Forehead Digital Thermometer to measure another person's forehead temperature, illustrating the non-contact method.



Image: The MOBI thermometer's display showing different color backlights (green, orange, red) corresponding to normal, high, and fever temperature ranges, providing a visual fever indicator.

4.2. Taking Object/Room Temperature

Ensure the thermometer is in 'Object' mode (indicated on the display).

1. Point the thermometer sensor towards the object or area you wish to measure.
2. Press the On/Scan button.
3. The temperature reading will appear on the LCD display.

Check the temperature of either the body or any object from a safe distance in just 1 second

Food



Room



Object Reading



Image: The MOBI thermometer being used to measure the temperature of food, a room, and a baby bottle, demonstrating its versatility for object and room temperature readings.

4.3. Memory Recall

The thermometer stores the last 50 temperature readings. To access them:

- With the thermometer powered on, press the MEMO button.
- Each press will cycle through the stored readings, from most recent to oldest.

5. MAINTENANCE

5.1. Cleaning the Thermometer

To ensure accurate readings and hygiene, clean the thermometer regularly:

- Wipe the thermometer body with a soft, damp cloth.
- For the sensor lens, use a cotton swab lightly moistened with alcohol. Allow to dry completely before use.

- Do not immerse the thermometer in water or any other liquid.
- Do not use abrasive cleaning agents or thinners.

5.2. Storage

Store the thermometer in a cool, dry place, away from direct sunlight, extreme temperatures, and humidity. Keep it out of reach of children.

6. TROUBLESHOOTING

| Problem | Possible Cause | Solution |
|--------------------------------------|---|--|
| No display/Thermometer won't turn on | Batteries are dead or incorrectly inserted. | Replace batteries with new AAA batteries, ensuring correct polarity. |
| Inaccurate readings | Sensor is dirty; incorrect measurement distance; environmental factors (drafts, direct sunlight). | Clean the sensor lens; ensure proper distance (1-2 inches); measure in a stable environment. |
| Display shows 'Lo' or 'Hi' | Temperature is outside the measurable range for the selected mode. | Ensure you are in the correct mode (Body/Object) for the temperature being measured. |
| Memory recall not working | No readings stored or memory is full (older readings are overwritten). | Take new readings; the device stores the last 50 readings. |

7. SPECIFICATIONS







- **Model Number:** 70113
- **Product Dimensions:** 2 x 0.01 x 0.01 inches
- **Weight:** 5.29 ounces
- **Batteries:** 2 AAA batteries (included)
- **Display Type:** LCD
- **Special Features:** Fast Reading, Non-Contact, Accuracy, Large LCD, Auto Off, °F/°C, Room/Food Mode
- **Color:** Blue, White
- **Outer Material:** Plastic


8. WARRANTY AND SUPPORT

For warranty information or technical support, please contact MOBI customer service. Refer to the product packaging or the official MOBI website for the most current contact details.

MOBI Technologies, Inc.
Website: [MOBI Official Store on Amazon](#)

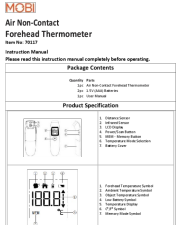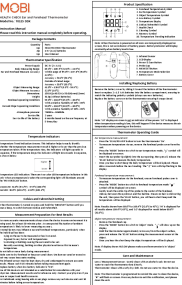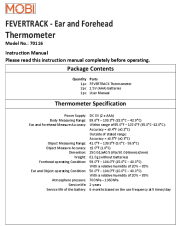Related Documents - 70113

| | |
|---|--|
|  | <p>MOBI Air Non-Contact Forehead Thermometer (Model 70117) - Instruction Manual</p> <p>Comprehensive instruction manual for the MOBI Air Non-Contact Forehead Thermometer (Item No: 70117), covering setup, operation, troubleshooting, specifications, and safety guidelines.</p> |
|  | <p>MOBI HEALTH CHECK Ear and Forehead Thermometer Model 70125 DDX Instruction Manual</p> <p>Comprehensive instruction manual for the MOBI HEALTH CHECK Ear and Forehead Thermometer, Model 70125 DDX. Learn how to use, maintain, and troubleshoot your digital thermometer for accurate temperature readings.</p> |
|  | <p>MOBI FEVERTRACK Thermometer (Model 70116) - User Manual & Specifications</p> <p>Comprehensive instruction manual and specifications for the MOBI FEVERTRACK Ear and Forehead Thermometer (Model 70116). Learn how to use, maintain, and troubleshoot your digital thermometer.</p> |
|  | <p>MOBI Smart Bluetooth Ear + Forehead DualScan Thermometer User Manual</p> <p>User manual for the MOBI Smart Bluetooth Ear + Forehead DualScan Thermometer (Model 700026). Provides instructions on usage, features, safety precautions, and troubleshooting for accurate temperature readings.</p> |
|  | <p>MOBI Digital Otoscope User Guide and Setup Instructions</p> <p>Comprehensive guide to setting up and using the MOBI Digital Otoscope in EZ Mode and AP Mode, including troubleshooting tips and FCC compliance information.</p> |
|  | <p>Turning Technologies Hardware Compatibility Charts</p> <p>Comprehensive compatibility charts for Turning Technologies' range of interactive learning hardware, including receivers, clickers, and whiteboards, detailing compatibility with various software and systems like Device Manager, ExamView, and Mobi.</p> |

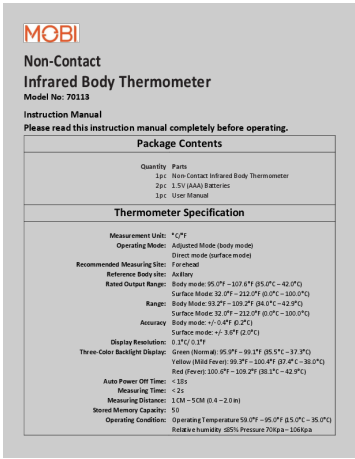


[pdf] User Manual Specifications

Kiumars Haziq MOBI Digital Non Contact ThermometerItem 70113 W Instruction Manuals Technologies

Inc Instructionmanual 09182020 v 1600470173 cdn shopify s files 1 1068 4858 |||

Non-Contact Infrared Body Thermometer Model No:**70113** Instruction Manual Please read this instruction manual completely before operating. Package Contents Quantity 1pc 2pc 1pc Parts Non-Contact Infrared Body Thermometer 1.5V AAA Batteries User Manual Thermometer Specification Measurement Unit:...

lang:en score:41 filesize: 302.57 K page_count: 8 document date: 2020-09-18

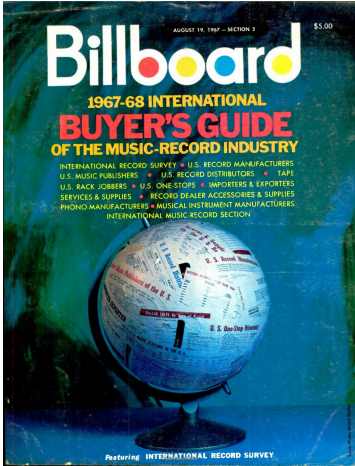

[pdf] Guide Buyer's Guide Accessories

BUYER S GUIDE Billboard 1967 08 19 II International Buyers Guide worldradiohistory Archive All Music 60s |||


BUYER S GUIDE 1967-68 INTERNATIONAL OF THE MUSIC-RECORD INDUSTRY INTERNATIONAL RECORD SURVEY U.S. R ... : Hartco, Ebb-Tide II, Goiden Prod. NEW ORLEANS Area Code 504 All South, 1924 Lafayette St. Zip: **70113**. Tel: 522-1157. Pres. Mgr.: Henry J. Hildebrand Jr.; Promo.: Stan Chaisson. Labels: A M, Aft... lang:en score:19 filesize: 36.8 M page_count: 296 document date: 2024-06-25



Billboard 1982-83 International Buyer's Guide: Music Industry Directory

The Billboard 1982-83 International Buyer's Guide is a definitive industry resource, featuring comprehensive listings of record companies, music publishers, wholesalers, services, and detailed statistics from over 70 countries. Essential for music business professionals navigating the early 1980s market.

lang:en score:16 filesize: 70.05 M page_count: 582 document date: 2019-10-07

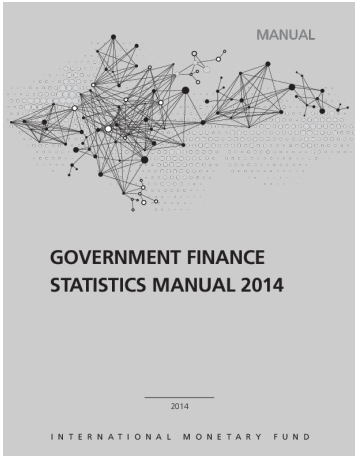

[pdf] User Manual Catalog

Government Finance Statistics Manual 2014 IMF

MANUAL GOVERNMENT FINANCE STATISTICS MANUAL 2014 2014

INTERNATIONAL MONETARY FUND 2014 Internatio ... ges and interest payments on government loans 70170 ; supervision of the banking industry 70411 . **70113** External affairs CS Administration of external affairs and services Operation of the minis...

lang:en score:12 filesize: 4.54 M page_count: 470 document date: 2017-08-17

