

Oclean OC-X PRO-BLUE

Oclean X Pro Sonic Electric Toothbrush Instruction Manual

Model: OC-X PRO-BLUE

INTRODUCTION

The Oclean X Pro Electric Toothbrush is an advanced sonic toothbrush designed to provide a superior cleaning experience. Featuring a powerful motor, smart touchscreen, and customizable brushing plans, it aims to enhance your daily oral hygiene routine. This manual provides essential information for setting up, operating, maintaining, and troubleshooting your Oclean X Pro toothbrush.

Key Features:

- **Sonic Cleaning:** Utilizes sonic waves for effective teeth and gum cleaning.
- **Long Battery Life:** Up to 30-day battery life with a 2-hour fast charge.
- **Advanced Features:** 3 modes and 32 intensities for personalized oral care.
- **Smart Timer:** Built-in timer ensures proper brushing duration.
- **Durable Construction:** Sturdy plastic body with a DuPont brush head.

SETUP

1. Unboxing and Initial Inspection

Carefully remove all components from the packaging. Ensure the toothbrush handle, brush head, charging base, and USB cable are present and undamaged.

SMART SCREEN

see your brushing results
and blind areas instantly



Image: The Oclean X Pro toothbrush handle and a separate brush head, ready for assembly.

2. Attaching the Brush Head

Align the brush head with the metal shaft on the toothbrush handle. Gently push the brush head down until it clicks securely into place. Do not twist the brush head.

3. Initial Charging

Connect the USB charging cable to the charging base and plug it into a compatible USB power adapter (not included). Place the toothbrush handle onto the charging base. The smart screen will display the charging status. A full charge takes approximately 2 hours and provides up to 30 days of use.

LONG BATTERY LIFE

30-day battery life
2-hour fast charge



Image: The Oclean X Pro toothbrush standing upright on its charging base, connected via a white USB cable.

OPERATING INSTRUCTIONS

1. Powering On/Off and Mode Selection

Press the power button to turn the toothbrush on. The smart screen will display the current mode and intensity. Swipe up or down on the touchscreen to cycle through the available brushing modes (Clean, Sensitive, Massage, White) and adjust the intensity (32 levels). Press the power button again to confirm your selection or to turn the toothbrush off.



Image: Close-up of the Oclean X Pro's smart screen showing brushing results and blind areas.

SMOOTH CLEAN FEELING

patented maglev motor,
up to

84,000

movements/min

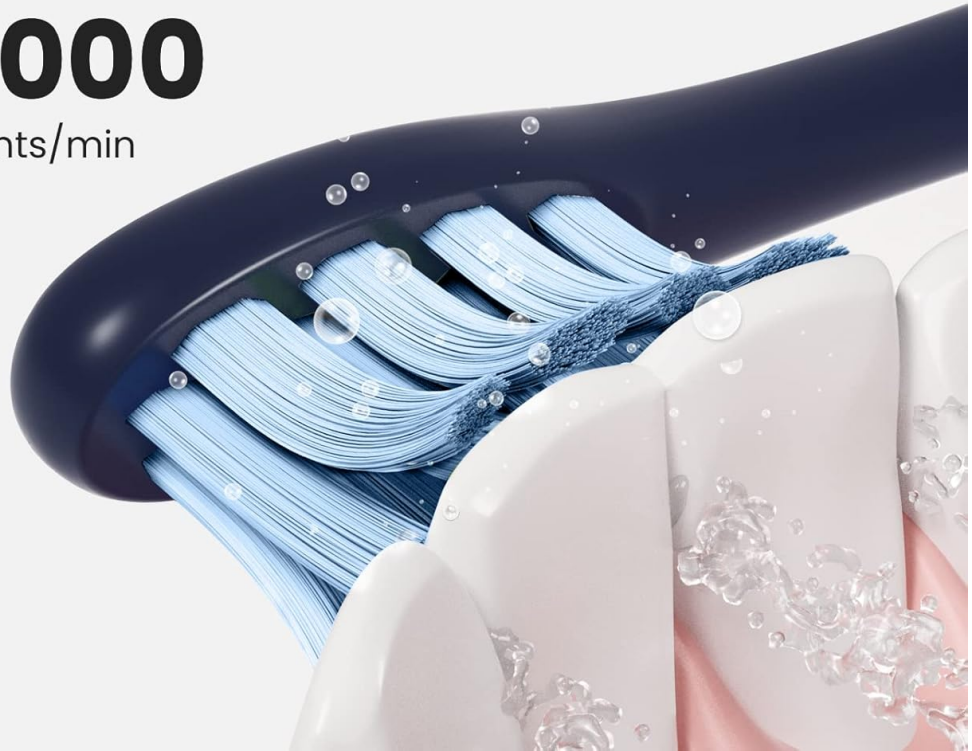


Image: The Oclean X Pro's smart screen showing options for brushing modes (Clean, White, Sensitive) and intensity levels.

2. Brushing Technique

Apply a small amount of toothpaste to the brush head. Place the bristles against your teeth at a slight angle towards the gum line. Move the brush slowly across each tooth, allowing the sonic vibrations to clean. The smart timer will guide you through different sections of your mouth, ensuring thorough cleaning.

3. Using the Oclean App

Download the Oclean app from your smartphone's app store. Connect your toothbrush via Bluetooth. The app allows for personalized brushing plans, real-time feedback on brushing coverage and pressure, and detailed reports on your oral hygiene habits. This helps identify missed spots and improve your brushing technique.

Your browser does not support the video tag.

Video: An official Oclean video introducing the Oclean X Pro, highlighting its features and app connectivity.

MAINTENANCE

1. Cleaning the Toothbrush

After each use, rinse the brush head thoroughly under running water. Wipe the handle with a damp cloth. Do not

immerse the entire handle in water for extended periods, although it is IPX7 waterproof for splashes and brief immersion.

2. Replacing the Brush Head

Dentists recommend replacing your brush head every 3 months or sooner if the bristles appear frayed. To replace, simply pull the old brush head straight off the handle and push a new Oclean DuPont brush head onto the shaft until it clicks.



Image: A close-up view of the Oclean X Pro's 'W' shape DuPont Tynex bristles, designed for gum and enamel protection.

TROUBLESHOOTING

- **Toothbrush not turning on:** Ensure the toothbrush is fully charged. If the issue persists, contact customer support.
- **App connectivity issues:** Ensure Bluetooth is enabled on your smartphone and the toothbrush is within range. Restart both the app and the toothbrush if necessary.
- **Inadequate cleaning results:** Check your brushing technique using the app's feedback. Adjust brushing mode and intensity as needed. Ensure the brush head is not worn out.
- **Battery draining quickly:** Higher intensity levels will shorten battery life. Ensure the toothbrush is not left on the charging base for extended periods when not actively charging.

SPECIFICATIONS

Feature	Detail
Product Dimensions	9.57 x 0.94 x 0.94 inches
Item Weight	3.49 ounces (99 Grams)
Model Number	OC-X PRO-BLUE
Power Source	Battery Powered (1 Lithium Ion battery required)
Battery Life	Up to 30 days
Charge Time	2 hours
Operation Mode	Semi-Automatic
Number of Settings	3 modes, 32 intensities
Bristle Material	DuPont Diamond + DuPont Classic
Waterproof Rating	IPX7

WARRANTY AND SUPPORT

Oclean provides a **2-year hassle-free replace/return warranty** for the Oclean X Pro Electric Toothbrush. For any product inquiries, technical assistance, or warranty claims, please contact Oclean customer support. They offer **7*24h customer support** to assist you.

Oclean is committed to improving oral health through technology and innovation, ensuring a complete and efficient oral care routine for its users.