

Sharper Image 843479126969

Sharper Image Acupressure Shiatsu Foot Massager Instruction Manual - Model 843479126969

INTRODUCTION

Thank you for choosing the Sharper Image Acupressure Shiatsu Foot Massager. This manual provides essential information for the safe and effective operation of your device. Please read it thoroughly before use and retain it for future reference.



Image: The Sharper Image Acupressure Shiatsu Foot Massager, showcasing its sleek design and foot inserts.

Safety Information

For dry use only. Use with caution near power sources. Not suitable for users with certain medical conditions. Always consult a physician before using any new massage device, especially if you have pre-existing medical conditions, are pregnant, or have a pacemaker. Do not use the device if it is damaged or malfunctioning. Keep out of reach of children.

Package Contents


Your Sharper Image Acupressure Shiatsu Foot Massager package includes:

- 1 x Acupressure Foot Massager Unit
- 1 x DC Adapter
- User Manual (this document)

WHAT YOU GET



Acupressure
Foot Massager

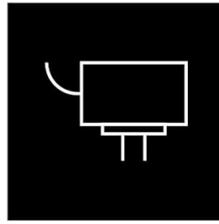

DC Adapter



WHAT YOU NEED



Feet



Image: Visual representation of the foot massager unit and its accompanying DC adapter, indicating what is included in the box.

PRODUCT FEATURES

The Sharper Image Acupressure Shiatsu Foot Massager is designed to provide a comprehensive foot massage experience with multiple features:

- **Acupressure Benefits:** Engages multiple acupressure points using reflexology micro nodes to rejuvenate and soothe tired feet.
- **Air Compression Comfort:** Immerses your feet in wraparound comfort, effectively soothing sore muscles and enhancing blood flow to improve circulation.
- **Soothing Heat Therapy:** Alleviates foot pain with soothing warmth, complementing acupressure and compression for complete muscle relaxation.
- **Deep Kneading Shiatsu:** Targets relaxation precisely with rolling shiatsu nodes across the arch of the foot.
- **Variable Strength:** Allows selection of gentle or more intense massage for customized relief.

SHARPER IMAGE®

ACUPRESSURE FOOT

MULTIPOINT SHIATSU MASSAGER



PINPOINT MASSAGE



AIR COMPRESSION



SHIATSU



HEATED



Image: A woman seated, using the foot massager, with icons illustrating the key features: pinpoint massage, air compression, shiatsu, and heated therapy.



WHERE LEADING-
EDGE TECHNOLOGY
MEETS EASE OF USE
WITH PRODUCTS
THAT ARE
ACCESSIBLE TO ALL.

Image: A diagram illustrating the various massage elements within the device, including air compression, heated zones, acupressure micro nodes, and shiatsu rolling nodes.

SETUP

1. **Unpack:** Carefully remove the foot massager and all accessories from the packaging.
2. **Placement:** Place the massager on a flat, stable surface in a comfortable seating area.
3. **Power Connection:** Connect the DC adapter to the massager's power input port, then plug the adapter into a standard electrical outlet.
4. **Hygiene:** The internal foot sleeves are removable for cleaning. Ensure they are properly inserted before use.



Image: The foot massager unit with its power cord plugged into an outlet, ready for use.

OPERATING INSTRUCTIONS

The foot massager can be operated using the integrated control panel or the included remote control.

Control Panel Overview



Image: A detailed view of the control panel, showing buttons for power, time, kneading intensity, air intensity, and heat.

- **Power Button:** Press to turn the device ON/OFF.
- **Time Button:** Cycle through 15, 20, or 30-minute massage sessions.
- **Knead Intensity Button:** Adjust the intensity of the kneading massage (Low, Medium, High).
- **Air Intensity Button:** Adjust the intensity of the air compression massage (Low, Medium, High).
- **Heat Button:** Turn the heat function ON/OFF.

Getting Started

1. Sit comfortably and place your feet into the foot inserts.
2. Press the **Power Button** on the control panel or remote to turn on the massager. The device will start with default settings.
3. Use the respective buttons to adjust **Knead Intensity**, **Air Intensity**, **Heat**, and **Time** according to your preference.
4. To end the session early, press the **Power Button** again. The device will automatically shut off after the selected time duration.

MAINTENANCE

- **Cleaning:** The internal foot sleeves are removable and washable. Unzip the covers, remove them, and hand wash with mild soap and water. Allow to air dry completely before reattaching.
- **Exterior:** Wipe the exterior of the massager with a soft, damp cloth. Do not use abrasive cleaners or immerse the unit in water.
- **Storage:** Store the massager in a cool, dry place away from direct sunlight when not in use.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Power cord not properly connected.	Ensure the power cord is securely plugged into the massager and a working electrical outlet.
Massage feels too weak/strong.	Intensity settings are not adjusted.	Adjust the Knead Intensity and/or Air Intensity settings using the control panel or remote.
Heat function is not working.	Heat function is not activated.	Press the Heat Button to turn on the heat function. Allow a few minutes for the heat to become noticeable.

SPECIFICATIONS

- **Product Dimensions:** 7 x 14 x 15 inches
- **Item Weight:** 1 Pound
- **Item Model Number:** 843479126969
- **Power Source:** Corded Electric
- **Material:** Plastic, Rubber
- **Manufacturer:** MerchSource
- **Date First Available:** July 15, 2020




Image: The foot massager with its dimensions (12.5 inches by 14 inches) indicated.

WARRANTY AND SUPPORT

For warranty information or customer support, please refer to the product packaging or contact Sharper Image customer service directly. Keep your purchase receipt for any warranty claims.

Related Documents - 843479126969

	<p>Sharper Image Shiatsu Dome Foot Massager Owner's Manual</p> <p>Official owner's manual for the Sharper Image Shiatsu Dome Foot Massager (Item No. 205360), detailing features, functions, operation, safety warnings, and warranty information for this Shiatsu, air compression, and heat foot massager.</p>
---	--



[Sharper Image Shiatsu Full Body Massager User Manual and Safety Guide](#)

Comprehensive guide for the Sharper Image Shiatsu Full Body Massager, including setup, features, care instructions, and important safety warnings for optimal and safe use.



[Sharper Image Shiatsu Foot + Calf Massager User Manual and Safety Guide](#)

Official user manual and safety instructions for the Sharper Image Shiatsu Foot + Calf Massager, detailing features, operation, care, and FCC compliance.




[Shiatsu Kneading and Rolling Foot Massager User Guide | Sharper Image 206997](#)

User guide for the Sharper Image Shiatsu Kneading and Rolling Foot Massager (Item No. 206997). Learn about features, basic operation, control pad functions, safety warnings, and warranty information.



[Sharper Image RealTouch Massage User Guide](#)

Instructions and safety warnings for the Sharper Image RealTouch Wireless Neck and Back Shiatsu Massager with Heat.

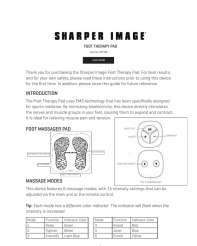

[Sharper Image Foot Therapy Pad User Guide - Model 207382](#)

Comprehensive user guide for the Sharper Image Foot Therapy Pad (Model 207382). Learn about its EMS technology, massage modes, operating instructions, practice modes, safety precautions, battery assembly, charging, maintenance, specifications, and warranty information.