

Globus G5944

Globus Mat 100 Memory Foam Magnetotherapy Mat User Manual

Model: G5944

1. INTRODUCTION

This manual provides essential instructions for the proper use, setup, and maintenance of your Globus Mat 100 Memory Foam Magnetotherapy Mat. The Mat 100 is designed as an accessory for pulsed electromagnetic field (PEMF) treatments, compatible with most Globus Magnum line magnetotherapy devices (excluding Magnum L). It is engineered to provide comfortable support during treatments, particularly for large body surfaces such as the spine.

2. PRODUCT OVERVIEW AND FEATURES

The Globus Mat 100 is a specialized memory foam mat equipped with four solenoids for magnetotherapy applications. Its design focuses on user comfort and effective treatment delivery.

- **PEMF Accessory:** Specifically designed for use with pulsed electromagnetic field (PEMF) therapy devices.
- **Targeted Treatment:** Ideal for treating large body surfaces, such as the spine, due to its generous dimensions.
- **Comfort Design:** Features a triple layer of padding, including memory foam, to comfortably accommodate the treated area, especially beneficial for inflamed or sore regions.
- **Compatibility:** Compatible with all Globus Magnum line magnetotherapy devices, with the exception of the Magnum L model.
- **Integrated Solenoids:** Contains 4 solenoids strategically placed for effective magnetic field distribution.
- **Washable Covers:** Equipped with a washable hypoallergenic outer cover and a waterproof, antibacterial inner lining for hygiene.



Figure 1: The Globus Mat 100 Memory Foam Magnetotherapy Mat. This image shows the overall appearance of the mat, highlighting its grey, textured outer cover and rectangular shape.



Figure 2: A cutaway view of the Globus Mat 100, revealing its internal structure. This image illustrates the multiple layers of padding, including memory foam, and indicates the placement of the four solenoids within the mat.

3. SETUP INSTRUCTIONS

1. **Unpack the Mat:** Carefully remove the Mat 100 from its packaging.
2. **Inspect for Damage:** Before first use, inspect the mat and its connectors for any signs of damage. Do not use if damaged.
3. **Position the Mat:** Place the Mat 100 on a flat, stable surface where you intend to conduct your magnetotherapy session. Ensure adequate space around the mat.
4. **Connect to Device:** Connect the Mat 100 to a compatible Globus Magnum series magnetotherapy device (e.g., Magnum 2500, Magnum XL Pro). Refer to your specific magnetotherapy device's manual for connection port details. The mat features connectors for its 4 solenoids.
5. **Ensure Compatibility:** Verify that your Globus Magnum device is compatible with accessories like the Mat 100. Note that the Mat 100 is **not** compatible with the Magnum L model.

4. OPERATING INSTRUCTIONS

The Mat 100 functions as an applicator for your Globus magnetotherapy device. Its operation is controlled

entirely by the connected device.

1. **Prepare for Treatment:** Once the Mat 100 is connected to your magnetotherapy device, select the desired treatment program on the device.
2. **Position Yourself:** Lie down or position the body part requiring treatment comfortably on the Mat 100. The mat's memory foam padding is designed to conform to your body, providing support for areas such as the spine.
3. **Start Treatment:** Initiate the treatment program on your magnetotherapy device according to its specific instructions. The solenoids within the Mat 100 will then generate the pulsed electromagnetic fields.
4. **During Treatment:** Remain still and relaxed during the treatment session. The mat's design is intended to maximize comfort and therapeutic contact.
5. **End Treatment:** Once the program concludes, turn off your magnetotherapy device. You may then carefully get up from the mat.

Important: Always refer to the user manual of your specific Globus Magnum magnetotherapy device for detailed program selection, intensity settings, and safety guidelines. The Mat 100 is an accessory and does not operate independently.

5. MAINTENANCE AND CARE

Proper maintenance ensures the longevity and hygiene of your Globus Mat 100.

- **Cleaning the Outer Cover:** The outer hypoallergenic fabric cover is washable. Remove it carefully and follow the washing instructions on the label, typically a gentle cycle with mild detergent. Allow it to air dry completely before reattaching.
- **Cleaning the Inner Lining:** The inner lining is made of washable, waterproof, and antibacterial fabric. This can be wiped clean with a damp cloth and mild disinfectant. Ensure it is dry before reassembling the mat.
- **Storage:** Store the mat in a clean, dry place away from direct sunlight and extreme temperatures. Avoid placing heavy objects on the mat to prevent deformation of the memory foam.
- **Inspection:** Periodically inspect the mat for any signs of wear, tear, or damage to the fabric or connectors.

6. TROUBLESHOOTING

The Globus Mat 100 is a passive accessory. Most operational issues will stem from the connected magnetotherapy device.

- **No Magnetic Field Sensation:**
 - Ensure the Mat 100 is securely connected to the magnetotherapy device.
 - Verify that the magnetotherapy device is powered on and the treatment program is active.
 - Check the device's error messages or indicators.
 - Confirm that your device is compatible with the Mat 100 (not Magnum L).
- **Physical Discomfort:**
 - Adjust your position on the mat.
 - Ensure the mat is placed on a flat, even surface.
 - If discomfort persists, discontinue use and consult a healthcare professional.

For issues related to the magnetotherapy device itself, refer to its dedicated user manual or contact Globus customer support.

7. TECHNICAL SPECIFICATIONS

Brand	Globus
Model Number	G5944
Dimensions (L x W x H)	70 x 44 x 6.5 cm
Item Weight	700 g
Material	Memory Foam
Number of Solenoids	4
Outer Cover	Washable hypoallergenic fabric with non-slip underside
Inner Lining	Washable waterproof and antibacterial fabric
Color	Grey

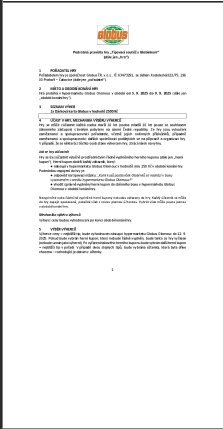
8. WARRANTY AND SUPPORT

Specific warranty information for the Globus Mat 100 is not provided in this manual. For details regarding warranty coverage, terms, and conditions, please refer to the documentation included with your purchase or visit the official Globus website.

For technical support, service, or any inquiries regarding your Globus Mat 100 or compatible magnetotherapy device, please contact Globus customer service through their official channels.

Manufacturer: Globus

Related Documents - G5944

	<p>Globus Olomouc 'Tipovací soutěž s Globínkem' Game Rules</p> <p>Detailed rules for the Globus Olomouc 'Tipovací soutěž s Globínkem' (Guessing Contest with Globínek) promotional game. Learn how to participate, prize details, winner selection, and personal data handling.</p>
---	---



Pravidla soutěže "Hledá se Globínek"
2024-2025

- 1. ÚČASTNÍCI**
Soutěž je určena pro všechny děti, které v době soutěže (od 1. 10. 2024 do 31. 12. 2024) mají věk od 6 do 14 let. Soutěž je určena pro děti, které v době soutěže (od 1. 10. 2024 do 31. 12. 2024) mají věk od 6 do 14 let.
- 2. CÍL SOUTĚŽE**
Cílem soutěže je seznámit děti s životem zvířat a naučit je, jak se chovat k přírodě a životním prostředí. Soutěž je určena pro děti, které v době soutěže (od 1. 10. 2024 do 31. 12. 2024) mají věk od 6 do 14 let.
- 3. ÚČASTNÍCI**
Soutěž je určena pro všechny děti, které v době soutěže (od 1. 10. 2024 do 31. 12. 2024) mají věk od 6 do 14 let. Soutěž je určena pro děti, které v době soutěže (od 1. 10. 2024 do 31. 12. 2024) mají věk od 6 do 14 let.
- 4. ÚČASTNÍCI**
Soutěž je určena pro všechny děti, které v době soutěže (od 1. 10. 2024 do 31. 12. 2024) mají věk od 6 do 14 let. Soutěž je určena pro děti, které v době soutěže (od 1. 10. 2024 do 31. 12. 2024) mají věk od 6 do 14 let.
- 5. ÚČASTNÍCI**
Soutěž je určena pro všechny děti, které v době soutěže (od 1. 10. 2024 do 31. 12. 2024) mají věk od 6 do 14 let. Soutěž je určena pro děti, které v době soutěže (od 1. 10. 2024 do 31. 12. 2024) mají věk od 6 do 14 let.

[Globus 'Hledá se Globínek' Game Rules - Plzeň-Chotíkov](#)

Official rules for the 'Hledá se Globínek' game organized by Globus at their Plzeň-Chotíkov hypermarket. Learn how to participate, win prizes, and understand personal data handling.