

## Medisana 40462

# Medisana PSD Mechanical Bathroom Scale (Model 40462) - User Manual

Your guide to setting up, operating, and maintaining your Medisana PSD Mechanical Bathroom Scale.

## 1. INTRODUCTION

Thank you for choosing the Medisana PSD Mechanical Bathroom Scale. This analogue scale combines a classic retro design with reliable mechanical operation, providing a straightforward way to monitor your body weight. This manual provides essential information for the correct setup, safe operation, and proper maintenance of your scale.



**Image Description:** A top-down view of the Medisana PSD Mechanical Bathroom Scale. The scale features a black weighing surface and a prominent circular analogue dial with weight markings and the Medisana logo.

## 2. IMPORTANT SAFETY INFORMATION

---

- Always place the scale on a firm, flat, and hard surface. Uneven or soft surfaces (like carpets) can affect accuracy.
- Step onto the scale barefoot for accurate readings.
- Do not use the scale on slippery surfaces, such as wet floors.
- The maximum weight capacity for this scale is 150 kg. Do not exceed this limit.
- Keep the scale clean and dry. Avoid exposing it to extreme temperatures or direct sunlight.
- This device is intended for personal use only and not for commercial or medical purposes.

## 3. PACKAGE CONTENTS

---

Upon unpacking, please ensure all components are present and undamaged:

- Medisana PSD Mechanical Bathroom Scale (Model 40462)
- User Manual

## 4. PRODUCT FEATURES


The Medisana PSD Mechanical Bathroom Scale offers the following key features:

- **Retro Design:** A classic, timeless aesthetic reminiscent of 1960s scales.
- **Mechanical Operation:** Battery-free and simple to use.
- **Large Analogue Display:** Clear weight readings in 500-gram steps.
- **Durable Construction:** Features a reinforced glass cover and a high-quality stainless steel weighing surface (Second Edition).
- **Non-Slip Surface:** Ensures stability during use.

Non-slip surface



Solid metal casing

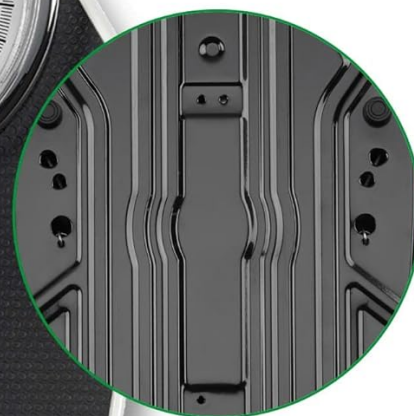

Trendy retro design



**Image Description:** An image highlighting the key features of the Medisana PSD scale. It shows a close-up of the non-slip textured surface, an internal view of the solid metal casing, and the overall retro design of the scale.

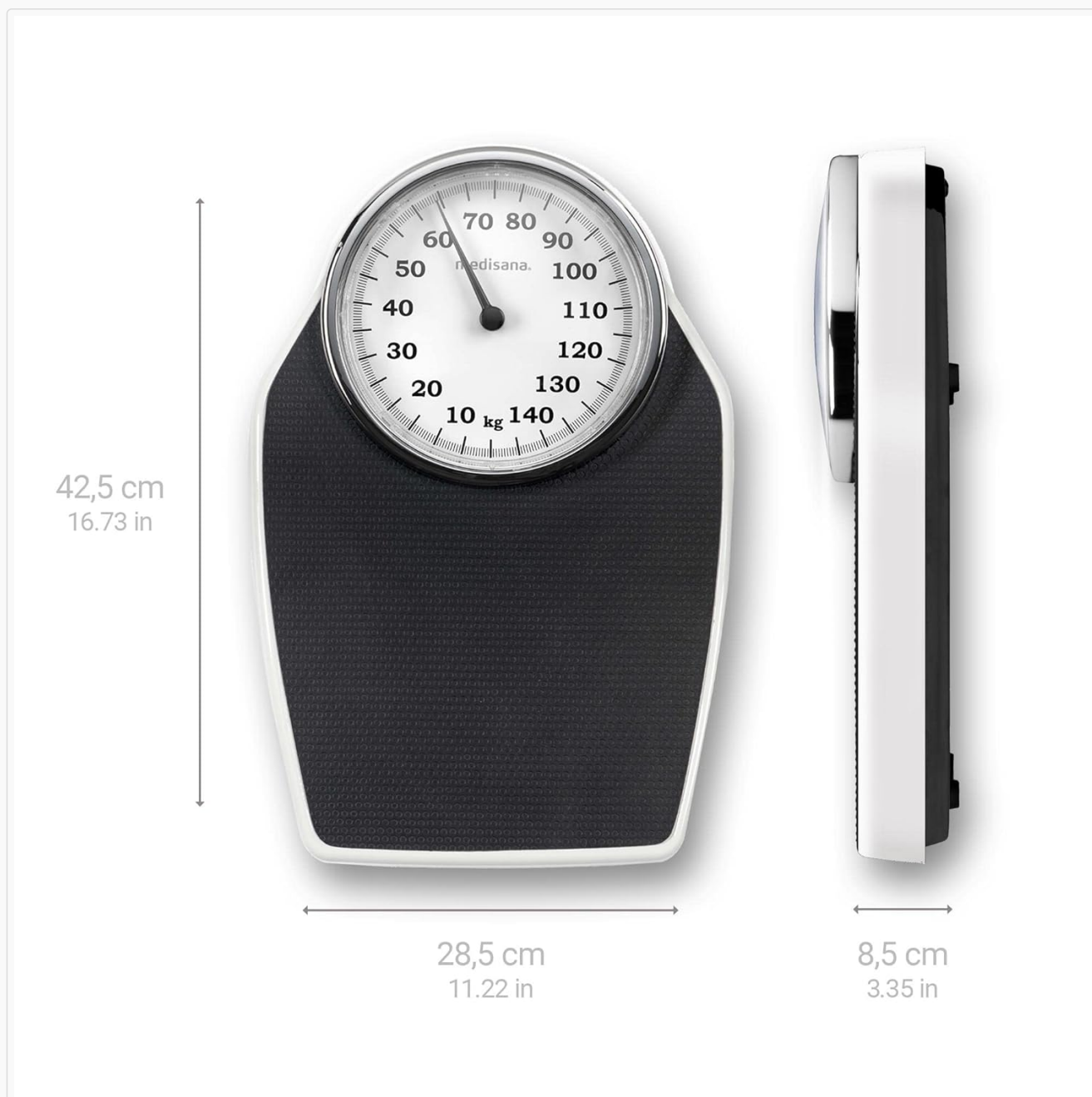


**Image Description:** A detailed close-up of the Medisana PSD scale's analogue dial. The dial clearly shows weight markings from 0 to 140 kg, with the Medisana logo visible. Text overlays indicate "Large analogue full view scale" and "Battery-free, mechanical operation."

## 5. SETUP AND INITIAL CALIBRATION

1. **Unpack the Scale:** Carefully remove the scale from its packaging.
2. **Placement:** Place the scale on a hard, flat, and stable surface. Avoid carpets, rugs, or uneven flooring, as these can affect accuracy.
3. **Initial Calibration:** If the scale has been transported or moved, it requires initial calibration.
  - Gently tap your foot on the surface of the scale.
  - Observe the pointer on the dial. If it does not point to 0.0 kg, locate the adjustment wheel (usually at the base or side of the scale).
  - Turn the adjustment wheel left or right until the pointer aligns precisely with the 0.0 kg mark.

- Repeat this process if necessary to ensure accurate zeroing.



**Image Description:** An image displaying the dimensions of the Medisana PSD scale. The scale is shown from a top-down view and a side profile, with measurements indicating a height of 42.5 cm (16.73 in), a width of 28.5 cm (11.22 in), and a depth of 8.5 cm (3.35 in).

## 6. OPERATING THE SCALE

1. Ensure the scale is placed on a hard, flat surface and the pointer is set to 0.0 kg.
2. Step onto the scale carefully and stand still. For best accuracy, step on barefoot.
3. Wait for the pointer to stabilize.
4. Read your weight from the analogue dial. The scale displays weight in 500-gram increments.
5. Step off the scale. The pointer will return to zero.



**Image Description:** A person's bare feet are shown standing on the Medisana PSD Mechanical Bathroom Scale, which is placed on a tiled floor. This illustrates the correct posture for weighing.

## 7. CARE AND MAINTENANCE

---

- **Cleaning:** Wipe the scale surface with a damp cloth and a mild detergent if necessary. Do not use abrasive cleaners or immerse the scale in water.
- **Storage:** Store the scale in a dry place at room temperature.
- **Avoid Impact:** Protect the scale from strong impacts or dropping, as this can damage the internal mechanism.

## 8. TROUBLESHOOTING

---

### **Inaccurate Readings:**

- Ensure the scale is on a hard, flat surface. Soft surfaces like carpets can cause inaccuracies.
- Verify that the pointer is set to 0.0 kg before each use. If not, recalibrate using the adjustment wheel as described in Section 5.

- Initial results may vary if the scale has just been moved or is new. Allow it to settle and perform the calibration steps.

**Pointer Not Returning to Zero:**

- This indicates the scale needs recalibration. Adjust the zero point using the adjustment wheel until the pointer rests exactly on 0.0 kg.

9. TECHNICAL SPECIFICATIONS

Attribute	Detail
Model Number	40462
Brand	Medisana
Type	Mechanical Bathroom Scale
Display Type	Analogue
Weight Limit	150 Kilograms
Read-out Accuracy	+/- 0.1 kg (100 grams)
Product Dimensions (L x W x H)	42.5 x 28.5 x 8.5 cm (approx. 16.73 x 11.22 x 3.35 inches)
Item Weight	3.1 Grams <i>(Note: This value appears to be a data entry error. Actual product weight for a mechanical scale is typically several kilograms.)</i>
Material	Stainless Steel, Glass (reinforced cover)
Special Features	Bodyweight measurement, Large Display, Battery-free

10. WARRANTY AND SUPPORT

For warranty information and customer support, please refer to the warranty card included with your product or visit the official Medisana website. Keep your purchase receipt as proof of purchase.  
Manufacturer: medisana



© 2023 Medisana. All rights reserved. This manual is for informational purposes only.  
For further assistance, please contact Medisana customer service.

