

Beurer 575.01

Beurer IPL Velvet Skin Pro User Manual

Model: 575.01

1. PRODUCT OVERVIEW

The Beurer IPL Velvet Skin Pro is an innovative device designed for long-lasting hair removal using advanced Intense Pulsed Light (IPL) technology. It provides a convenient and effective solution for achieving smooth skin from the comfort of your home. This epilator is suitable for use on various body areas, including legs, armpits, bikini line, chest, stomach, back, and face (for women, below the cheekbones).



Figure 1: Beurer IPL Velvet Skin Pro device and its packaging.

The device is cordless, powered by a rechargeable battery, offering flexibility during use. It features a large light surface for efficient treatment and includes a precision attachment for smaller, sensitive areas. Integrated safety features, such as a skin type and contact sensor, ensure safe operation.

2. KEY FEATURES

- **Long-lasting Hair Removal:** Utilizes advanced IPL technology for permanent hair reduction by deactivating hair roots.
- **Cordless Operation:** Rechargeable battery allows for comfortable, wire-free use.
- **Fast Treatment:** Approximately 21 minutes for a full body treatment.
- **Large Light Surface:** 4.5 cm² light surface for particularly fast application.
- **Adjustable Intensity:** 6 intensity levels to suit different skin types and treatment areas.
- **Precision Attachment:** Included for precise treatment of the face and bikini line.
- **Safety Sensors:** 2-in-1 skin type and contact sensor for safe application.
- **High Pulse Capacity:** Up to 600,000 light pulses for extended device lifespan.
- **Auto-Flash Mode:** For automatic triggering of light pulses, enhancing ease of use.

- **Visible Results:** Up to 50% hair reduction after just 3-4 treatments.
- **Integrated UV Filter:** Protects the skin during treatment.
- **App Integration:** Free 'beurer MyIPL' app guides you through the application and creates a personalized treatment program.

- **Technologie IPL** moderne pour une peau lisse et soyeuse
- **Des résultats visibles** réduction de la pilosité jusqu'à 50 % après 3 à 4 utilisations
- **Sans fil**
- **Avec application** en option
- **Détection automatique** des types de peau et de cheveux adaptés



Figure 2: The Beurer IPL Velvet Skin Pro, showcasing its modern IPL technology, visible results, cordless design, optional app integration, and automatic skin type detection.

3. UNDERSTANDING IPL TECHNOLOGY

Intense Pulsed Light (IPL) technology has been clinically proven safe and effective for long-term hair reduction for over 15 years. The IPL Velvet Skin Pro uses optical energy to deactivate hair growth. The melanin pigment in the hair absorbs the light from the flash emitted by the IPL device. This generates heat, which deactivates the hair follicle (the root of the hair). As a result, the hair falls out and does not regrow.

Comment fonctionne la méthode IPL ?



Figure 3: IPL Mechanism - Hair melanin absorbs light, generated heat deactivates the hair follicle, and hair falls out without regrowing.

IPL is most effective on certain hair and skin color combinations. The device is designed to work safely and effectively within these parameters.

Couleurs de peau

Couleur de poils	Lightest skin	Light skin	Medium skin	Dark skin	Very dark skin	Darkest skin
Dark hair	X	X	X	X	X	X
Light hair	✓	✓	✓	✓	✓	X

Figure 4: Skin tone and hair color compatibility for IPL treatment. IPL is generally effective on light to medium skin tones with dark hair, and not recommended for very dark skin or very light hair (white, grey, red).

4. SETUP

- Charge the Device:** Before using the device for the first time, ensure it is fully charged. Connect the provided power adapter to the device and a power outlet.
- Skin and Contact Sensor:** The device is equipped with a 2-in-1 skin type and skin contact sensor. This sensor ensures safe operation by only allowing flashes when the device is in full contact with

suitable skin.

3. **Power On and Select Intensity:** Press the power button to turn on the device. Select the appropriate intensity level based on your skin type and comfort. Refer to the skin compatibility chart for guidance.
4. **Attach Precision Head (if needed):** For smaller or sensitive areas like the face or bikini line, attach the precision head to the device.

5. OPERATING INSTRUCTIONS

For optimal and lasting results, regular use of the Beurer IPL Velvet Skin Pro is necessary. Follow these steps for effective treatment:

1. **Prepare the Skin:** Ensure the treatment area is clean, dry, and free of hair (shave prior to treatment). Do not wax or pluck hair before IPL treatment, as the hair follicle needs to be present for the light to be absorbed.
2. **Select Intensity Level:** Based on your skin tone and hair color, choose the appropriate intensity level. Start with a lower setting and gradually increase if comfortable.
3. **Position the Device:** Place the light surface of the device firmly against the skin. The skin contact sensor will ensure proper placement.
4. **Initiate Flash:** Press the flash button. For larger areas, use the 'Auto-Flash' mode by holding down the flash button and gliding the device over the skin.
5. **Treat All Areas:** Move the device systematically across the desired treatment area, ensuring all parts are covered without overlapping too much.
6. **Treatment Duration:** A full body treatment takes approximately 21 minutes.
7. **Post-Treatment Care:** After treatment, you may apply a soothing lotion if desired. Avoid direct sun exposure on treated areas immediately after.

Utilisation sur tout le corps

avec embout de précision pour les zones sensibles



Figure 5: The device can be used on various body parts, including armpits, with the precision attachment for sensitive areas.

Épilation longue durée


pour une peau douce et lisse



Figure 6: Experience long-lasting smooth skin with regular use of the IPL device.

Using the 'beurer MyIPL' App

The free 'beurer MyIPL' app provides a personalized treatment plan and guides you step-by-step throughout your hair removal journey. It offers reminders and helps you track your progress.

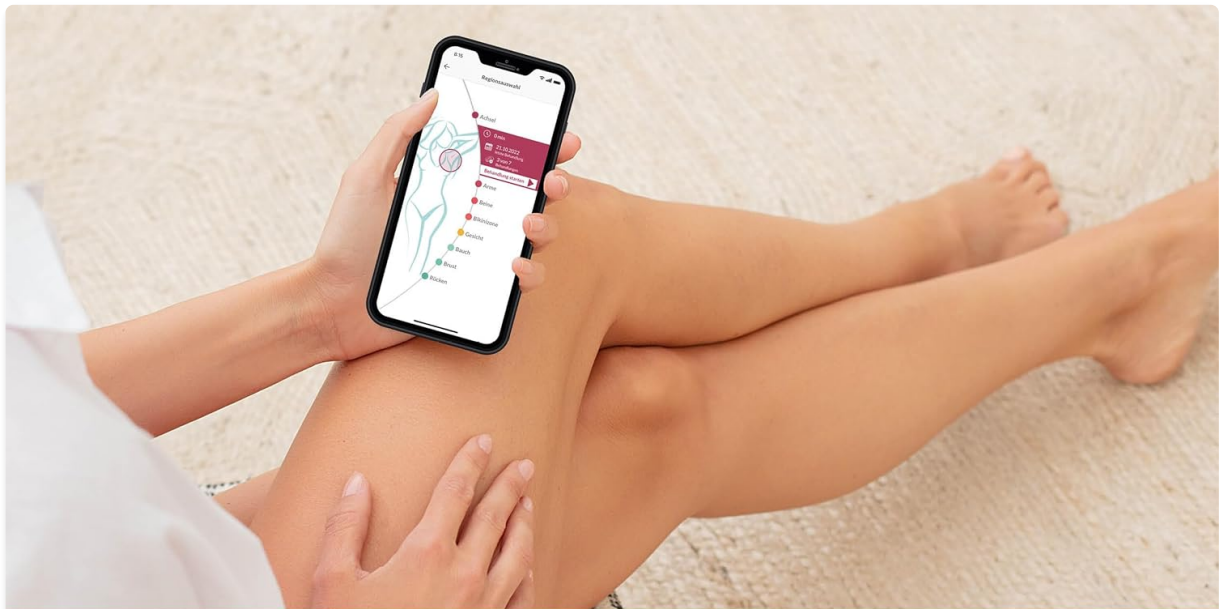

Utilisation optionnelle via l'application

Avec l'application « beurer MyIPL », vous recevez une recommandation adaptée à votre type de peau, un plan d'utilisation personnel avec des rappels pratiques et vous êtes guidé pas à pas tout au long du traitement.

Figure 7: The 'beurer MyIPL' app provides optional guidance and personalized treatment plans.

6. MAINTENANCE

To ensure the longevity and optimal performance of your Beurer IPL Velvet Skin Pro, follow these maintenance guidelines:

- **Cleaning:** After each use, gently wipe the light window and the device's surface with a soft, dry cloth. Do not use abrasive cleaners or solvents.
- **Storage:** Store the device in its original box or a protective case in a cool, dry place, away from direct sunlight and moisture.
- **Battery Care:** To maintain battery health, avoid completely discharging the battery frequently. Charge the device regularly, even if not in use for extended periods.

7. TROUBLESHOOTING

If you encounter issues with your Beurer IPL Velvet Skin Pro, consider the following common troubleshooting steps:

- **Device Not Turning On:** Ensure the device is fully charged. If it's charged, try pressing and holding

the power button for a few seconds.

- **No Flash Emitted:** Check if the device is in full contact with your skin. The skin contact sensor prevents flashes if contact is not sufficient. Also, ensure your skin tone is compatible with IPL treatment (refer to the skin compatibility chart).
- **Device Overheating:** If the device feels excessively warm, turn it off and allow it to cool down for a few minutes before resuming use. Ensure ventilation openings are not blocked.
- **App Connectivity Issues:** Ensure Bluetooth is enabled on your smartphone and the device is within range. Restart both the app and the device if connection problems persist.

If the problem persists after attempting these steps, please refer to the warranty and support section for further assistance.

8. SPECIFICATIONS

Feature	Specification
Model Number	575.01
Manufacturer	Beurer
Product Dimensions (L x W x H)	26.7 x 22.5 x 11.3 cm
Weight	680 grams
Power Source	Battery
Battery Type	2 AA (included)
Voltage	3.6 Volts (DC)
Charging Time	2 hours
Operating Time	21 Minutes (full body treatment)
Light Pulses	Up to 600,000
Light Surface Area	4.5 cm ²
Intensity Levels	6
Included Components	Precision attachment (for face and bikini)
Recommended Use	Body, face, and bikini hair removal

9. WARRANTY AND SUPPORT

The Beurer IPL Velvet Skin Pro comes with a **3-year warranty**, ensuring peace of mind regarding its quality and performance.

Utilisation sur tout le corps

avec embout de précision pour les zones sensibles



Figure 8: The Beurer IPL Velvet Skin Pro is backed by a 3-year warranty.

For detailed warranty terms, technical support, or service inquiries, please refer to the official Beurer website or contact their customer service department. Keep your proof of purchase for warranty claims.