

YA-MAN HRF40S

YA-MAN RF Facial Device Bright Lift White HRF40S User Manual

Model: HRF40S

1. PRODUCT OVERVIEW

The YA-MAN RF Facial Device Bright Lift White HRF40S is designed to enhance skin clarity and address common skin concerns. This device integrates Radio Frequency (RF), Electrical Muscle Stimulation (EMS), and LED light therapy to provide comprehensive skincare.

It operates through a three-step process: **Polish**, **Replenish**, and **Seal**, targeting various aspects of skin health.



Image: YA-MAN RF Facial Device Bright Lift White HRF40S. This image shows the device from the front, highlighting its sleek

Key Features:

- **RF (Radio Frequency) Technology:** Warms the skin to support effective skincare.
- **Ion Cleansing & Introduction:** Utilizes ion power to remove impurities and enhance the penetration of skincare ingredients.
- **EMS (Electrical Muscle Stimulation):** Stimulates facial muscles for a toning effect.
- **Blue LED Light:** Supports skin tightening and overall skin health.
- **Melano Barrier Pulse:** A new waveform in MOIST mode, designed to work with whitening cosmetics.
- **Cool-down Function:** Helps to soothe and tighten the skin after treatment.

2. SAFETY PRECAUTIONS

Please read all safety instructions carefully before using the device to prevent injury or damage.

Do Not Use If You Are:

- Using medical electrical equipment (e.g., implantable pacemakers, life support systems like artificial heart-lung machines, or wearable electrocardiographs).
- Suffering from heart disease, disability, or have a suspected condition.
- Unable to express your intentions or sensations (e.g., temperature or stimulation).
- Diagnosed with autoimmune diseases.

Consult a Doctor Before Use If You Are:

- Suffering from acute diseases or infectious diseases.
- Diagnosed with malignant tumors or febrile diseases.
- Experiencing high blood pressure or blood disorders.
- Having a fever.
- Experiencing sensory impairment due to severe peripheral circulatory dysfunction (e.g., diabetes).
- Diagnosed with neurological disorders.
- Suffering from skin diseases or atopic dermatitis.
- Currently undergoing medical treatment or taking medication.
- Prone to allergies, contact dermatitis, or have pre-existing skin conditions, or concerns about skin reactions during use.

Do Not Use On The Following Areas:

- Fractured areas.
- Areas where metal has been implanted in the body due to surgical procedures (e.g., plastic surgery, orthopedic surgery, cosmetic surgery).
- Areas with internal bleeding, swelling, or scabs after cosmetic medical procedures.

Note: For areas that have undergone cosmetic medical procedures (e.g., hyaluronic acid injection, silicone injection, Botox, thread lift, double eyelid surgery, HIFU), refrain from using the device for 2 weeks after the procedure. Consult your doctor, and if there are no issues, you may then use the device.

If you feel unwell or notice any skin abnormalities during use, immediately stop using the device. If the condition does not improve, consult a doctor.

3. SETUP

Unpacking and Initial Charge:

1. Carefully remove the device and all accessories from the packaging.
2. Place the device on the charging stand.
3. Connect the USB charging cable to the charging stand and then to the AC adapter.
4. Plug the AC adapter into a power outlet.
5. Charge the device for approximately 3 hours before first use. The indicator lights will show charging status.



Image: The YA-MAN RF Facial Device Bright Lift White HRF40S on its charging stand, ready for charging.

Attaching the Cotton Stopper (for CLEAN mode):

For the CLEAN mode, it is recommended to use a cotton pad. Attach the cotton stopper to secure a cotton pad onto the treatment head.

4. OPERATING INSTRUCTIONS

The device features three main modes: CLEAN (Polish), MOIST (Replenish), and COOL (Seal).

General Usage Steps:

1. Ensure your face is clean and free of makeup.
2. Apply appropriate skincare products as instructed for each mode.
3. Select the desired mode and intensity level.
4. Gently glide the device over your skin as directed.
5. After use, clean the device as per maintenance instructions.



Mode Details:

Step 1: CLEAN (Polish)

This mode uses ion power to remove impurities and old keratin from pores after washing your face, enhancing skin clarity and preparing it for subsequent skincare.

- Attach a cotton pad with the cotton stopper.
- Apply a suitable toner or cleansing lotion to the cotton pad.
- Select CLEAN mode and desired intensity.
- Gently glide the device over your face.




Image: Illustration of Ion Cleansing and Ion Introduction. The left panel shows positive ions attracting dirt, while the right panel shows negative ions aiding ingredient penetration.

Step 2: MOIST (Replenish)

This mode utilizes RF and the new Melano Barrier Pulse to help skincare ingredients penetrate the skin, promoting a brighter complexion. It is recommended for use with whitening cosmetics containing ingredients like Vitamin C derivatives or tranexamic acid.

- Apply your preferred serum or moisturizing cream to your face.
- Select MOIST mode and desired intensity.
- Gently glide the device over your face, allowing the RF warmth and pulse to work.

RF(ラジオ波)

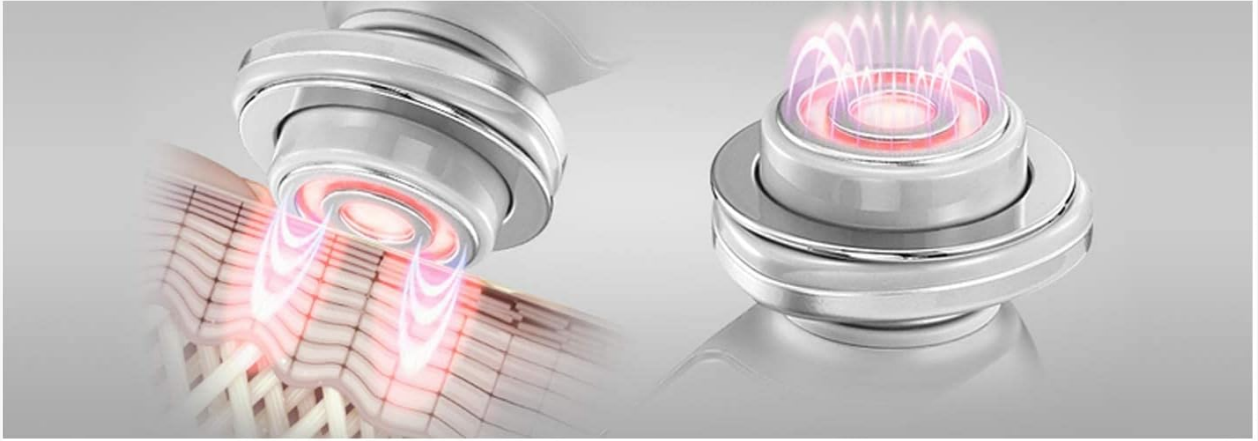


※角質層まで

高周波とも呼ばれる非常に細かい振動を肌に与えることで
角質層から温めるエステサロンや美容クリニック等でも人気の高いトリートメントです。
肌を温めることにより、効率のよいスキンケアをサポートします。

Image: Diagram illustrating RF (Radio Frequency) technology. It shows RF waves penetrating the skin layers, generating warmth.

EMS

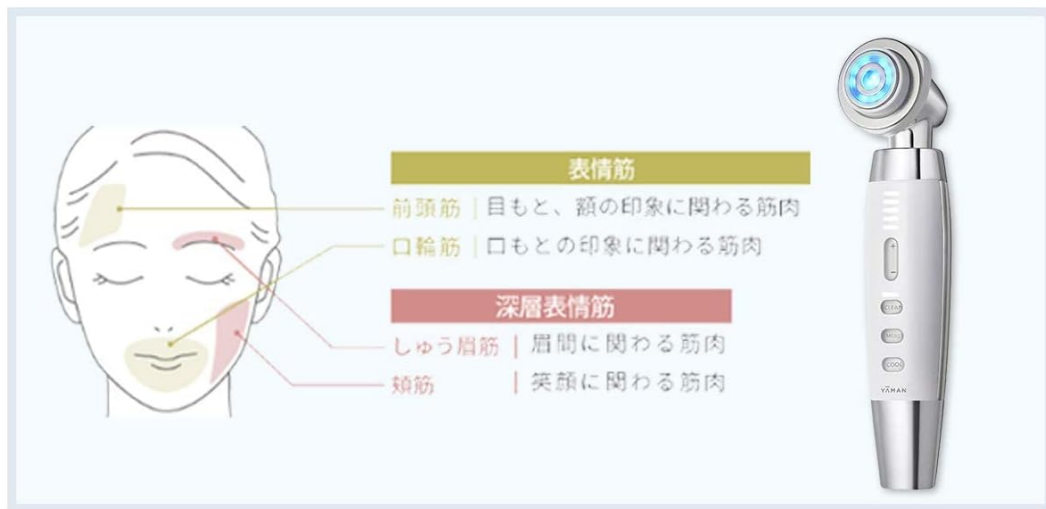

電気信号を表情筋に伝えることで、筋肉を刺激してトレーニングする機能を持つ美顔器です。
電極が接している面の筋肉を集中的に刺激するため、
ピンポイントなケアが可能です。



Image: Explanation of EMS (Electrical Muscle Stimulation). A facial diagram indicates various muscle groups, demonstrating how the device targets them.

Step 3: COOL (Seal)


The final step uses a cool-down function and blue LED light to help tighten and soothe the skin.

- After completing the MOIST mode, select COOL mode.
- Gently glide the cool head over your face.
- When using around the eyes, close your eyes and keep the electrode head at least 5cm away from the skin.

LED



まぶたや目の周りに使用する場合は、目を閉じ、
電極ヘッドをお肌から5cm以上離し
LED光を照射してください。

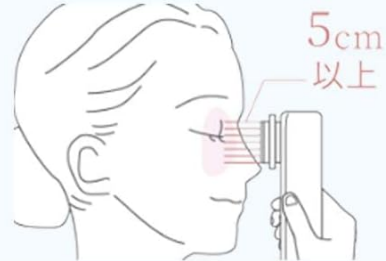

エステサロンでも話題の光を当てることでエイジングケア*1する技術LED（フォト）を搭載。

ブライトリフトはBLUE LEDを搭載しており、

まぶた*2やデコルテのケアも可能です。

※エイジングケア...年齢層にあったケア



Image: Blue LED light function and usage around the eyes. The image shows the device emitting blue light and a diagram indicating a minimum 5cm distance from the eye area during use.

5. MAINTENANCE

Cleaning the Device:

1. After each use, ensure the device is turned off and disconnected from the power source.
2. Wipe the electrode head and body of the device with a soft, damp cloth.
3. Do not use abrasive cleaners, alcohol, or harsh chemicals, as these may damage the device.
4. Ensure the device is completely dry before storing or recharging.

Storage:

Store the device in a cool, dry place, away from direct sunlight and extreme temperatures. Keep it out of reach of children.

6. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or not charged.	Charge the device for approximately 3 hours.

Problem	Possible Cause	Solution
No sensation during use.	Intensity level is too low or poor contact with skin.	Increase the intensity level. Ensure the electrode head is in full contact with your skin. Apply sufficient skincare product.
Skin irritation occurs.	Device used incorrectly or skin sensitivity.	Stop use immediately. Consult a doctor if irritation persists. Ensure proper application of skincare products and follow usage guidelines.

7. SPECIFICATIONS


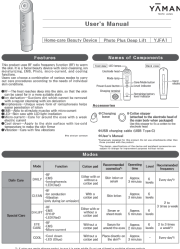




Feature	Detail
Model Number	HRF40S
Power Supply	DC5V 1.0A
Power Consumption	Approx. 4.5W
Battery	Lithium-ion rechargeable battery
Dimensions (with cotton stopper)	Approx. W45 × D66 × H211 (mm)
Weight	Approx. 230g
Charging Time	Approx. 3 hours
Operating Time	Approx. 40 minutes
Accessories	Charging stand, USB charging cable, cotton stopper, AC adapter
Color	Silver
Materials (Body)	ABS resin, polycarbonate resin
Materials (Electrode Head)	Polycarbonate resin / Stainless steel
Materials (Hand Electrode)	ABS resin (chrome plating)
Country of Manufacture	Japan

8. WARRANTY AND SUPPORT

For warranty information and customer support, please refer to the warranty card included with your product or visit the official YA-MAN website. Keep your proof of purchase for warranty claims.

Manufacturer Reference: 002-0011_22220

Related Documents - HRF40S

	<p>YA-MAN Photo PLUS Deep Lift User Manual</p> <p>Comprehensive user manual for the YA-MAN Photo PLUS Deep Lift home beauty device, detailing features, usage, safety precautions, and maintenance to ensure optimal and safe operation.</p>
	<p>YA-MAN Photo Plus Deep Lift YJFA1 User Manual - Home Care Beauty Device</p> <p>Comprehensive user manual for the YA-MAN Photo Plus Deep Lift YJFA1, a home-care beauty device featuring RF, EMS, Iontophoresis, LED, and Micro-current technologies for facial skin treatment.</p>
	<p>YA-MAN Liftlogy SP User Manual</p> <p>User manual for the YA-MAN Liftlogy SP home-care beauty device, providing instructions on usage, safety precautions, maintenance, and troubleshooting.</p>
	<p>YA-MAN YJMF4 Blue Green Mask Lift User Manual</p> <p>Comprehensive user manual for the YA-MAN YJMF4 Blue Green Mask Lift home-care beauty device, covering setup, usage instructions, safety precautions, troubleshooting, and product specifications.</p>
	<p>YA-MAN YJMF4 Blue Green Mask Lift User Manual</p> <p>Comprehensive user manual for the YA-MAN YJMF4 Blue Green Mask Lift, providing detailed instructions on usage, safety precautions, maintenance, and troubleshooting for this home-care beauty device.</p>
	<p>YA-MAN Photo Plus Deep Lift YJFA1 User Manual</p> <p>User manual for the YA-MAN Photo Plus Deep Lift YJFA1, a home beauty device featuring RF, Ion cleansing, Ion moisturizing, EMS, Photo (LED), microcurrent, and cool modes for facial care.</p>