

Microlife NC200

Microlife NC200 Non-Contact Forehead Thermometer User Manual

Model: NC200

1. INTRODUCTION

This manual provides detailed instructions for the safe and effective use of your Microlife NC200 Non-Contact Forehead Thermometer. Please read this manual thoroughly before first use and keep it for future reference. This device is designed for measuring human body temperature from the forehead and can also measure the surface temperature of objects and ambient room temperature.



Image 1.1: The Microlife NC200 Non-Contact Forehead Thermometer, showing its design and examples of its backlit display indicating normal (green) and high (red) temperatures.

2. SAFETY INFORMATION

- Do not use the device for purposes other than those described in this manual.
- Keep the device away from direct sunlight, high temperatures, humidity, and dust.
- Do not drop or subject the thermometer to strong shocks.
- Do not immerse the thermometer in water or other liquids.
- Consult a physician if you suspect a fever, even if the thermometer indicates a normal temperature.
- Keep out of reach of children. Small parts may be swallowed.
- The thermometer is not waterproof.

3. PRODUCT FEATURES

The Microlife NC200 thermometer offers several features designed for accurate and convenient temperature measurement:

- **Advanced Smart Accuracy:** Ensures measurements are taken within the optimal range (approximately 2 inches from the forehead/object) for accuracy.
- **Fast and Easy Non-Contact:** Provides quick readings for all ages without physical contact.
- **Self-Measurement Ready:** A visual green light on the handle indicates correct distance for self-measurement.
- **FeverGlow Technology:** The extra-large LCD screen illuminates green for normal temperatures and red for high temperatures (fever).
- **Sound On/Off:** Allows for silent operation, useful for not disturbing sleeping individuals.
- **Versatile Modes (3-in-1):** Measures body temperature, surface temperature of objects (e.g., baby's milk, bathwater), and ambient room temperature.
- **Memory Function:** Stores up to 30 readings.
- **Unit Selection:** Easily switch between Fahrenheit (°F) and Celsius (°C).

microlife

Non-Contact
Forehead
thermometer

with distance
control



3-second reading



Stores 30 readings in memory



Extra-large illuminated display
with date and time



High temperature (fever) alarm



Auto measurement
with distance control



Suitable for the entire family



FeverGlow™
Technology

NC200 Features

Copyright © 2020 Microlife corporation

Image 3.1: Overview of the Microlife NC200's key features, highlighting its speed, memory, display, and versatile use.

4. SETUP

4.1. Unpacking and Battery Installation

Upon opening the package, ensure all components are present:

- Microlife NC200 Infrared Forehead Thermometer
- Sensor Protection Cover
- 2 x AAA Batteries
- English / Spanish Instruction Manual (this document)

To install the batteries:

1. Open the battery compartment cover located on the back of the device.
2. Insert the two AAA batteries, ensuring correct polarity (+ and -).
3. Close the battery compartment cover securely.

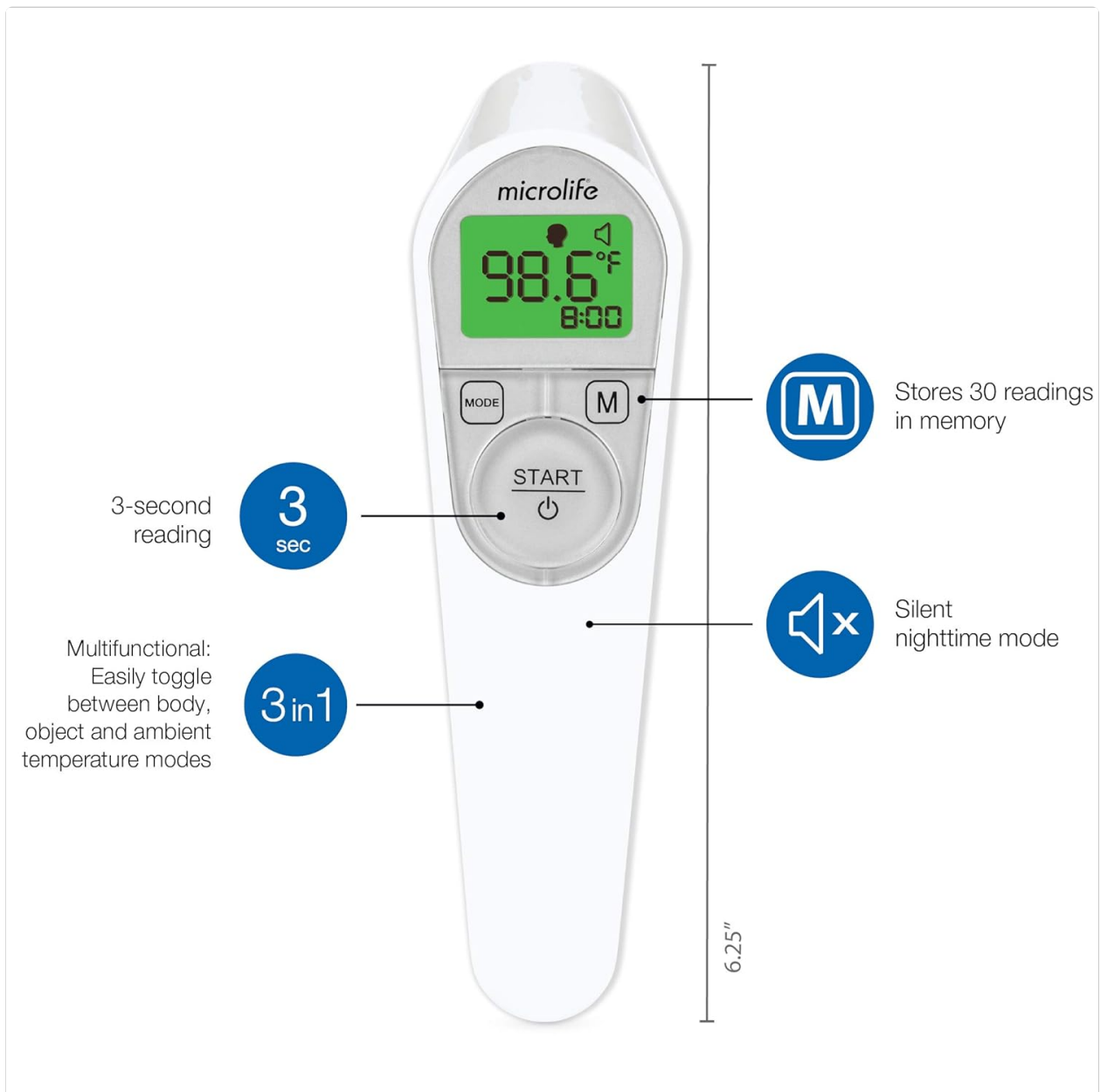


Image 4.1: The Microlife NC200 thermometer, its packaging, user manual, and two AAA batteries, as typically included in the box.

4.2. Initial Settings

After battery installation, the device may prompt for initial settings such as date, time, and temperature unit

(°F/°C). Follow the on-screen instructions to configure these settings.

5. OPERATING INSTRUCTIONS

5.1. Taking a Forehead Temperature Measurement

1. Ensure the thermometer is in body temperature mode (indicated by a small person icon on the display). If not, press the "MODE" button to switch.
2. Press the "START" button to activate the thermometer.
3. Aim the blue tracking light at the center of the forehead, approximately 2 inches away.
4. Move the unit closer to the forehead until the blue light stops blinking and a green light on the handle glows, indicating the correct distance.
5. A long beep will sound when the measurement is complete, and the temperature will be displayed on the screen. The screen will illuminate green for normal temperature or red for fever.

Auto Measurement = Accuracy



**automatically takes temperature
when held 2 inches away**

Image 5.1: Demonstrates the non-contact measurement process, highlighting the automatic distance control for accurate readings.

5.2. Measuring Object and Ambient Temperatures

The NC200 can also measure surface temperatures of objects and the ambient room temperature:

1. Press the "MODE" button to cycle through the available modes until the desired mode (object or ambient) is selected.
2. For object temperature, aim the thermometer at the surface of the object (e.g., baby bottle, bathwater) from a short distance.
3. For ambient temperature, simply place the thermometer in the room and allow it to stabilize.
4. Press the "START" button to take the measurement.



Body temperature



Object temperature



Ambient temperature

3 in 1 Accurate measurement

Copyright © 2020 Microlife corporation

Image 5.2: Visual representation of the thermometer's 3-in-1 functionality for body, object, and ambient temperature measurements.

5.3. Memory Recall

To view previous measurements, press the "M" (Memory) button. The thermometer stores up to 30 readings with date and time stamps.

5.4. Silent Mode

To enable or disable the sound, refer to the specific instructions in the full user manual. This feature is useful for taking measurements without disturbing sleeping individuals.

microlife®

Non-Contact Thermometer

Forehead thermometer

with distance control



Clinically Tested

Copyright © 2020 Microlife corporation

Image 5.3: Detailed view of the thermometer's interface, pointing out the memory button, mode button, and silent operation feature.

6. MAINTENANCE

6.1. Cleaning the Thermometer

To ensure accurate readings and hygiene, clean the thermometer after each use:

- Wipe the probe and the body of the thermometer with a soft cloth moistened with 70% isopropyl alcohol or an alcohol swab.
- Do not use abrasive cleaning agents, thinners, or benzene for cleaning.
- Ensure no liquid enters the interior of the device.
- Allow the thermometer to air dry completely before storing.

6.2. Storage

Store the thermometer in a clean, dry place, away from direct sunlight, extreme temperatures, and humidity. Keep the sensor protection cover on when not in use.

6.3. Battery Replacement

When the low battery indicator appears on the display, replace both AAA batteries with new ones. Follow the battery installation steps in Section 4.1.

7. TROUBLESHOOTING

Problem	Possible Cause	Solution
No display after pressing START.	Batteries are dead or incorrectly inserted.	Check battery polarity or replace batteries.
Inaccurate readings.	Incorrect measurement distance; dirty sensor; environmental factors (drafts, direct sunlight).	Ensure correct distance (approx. 2 inches); clean the sensor; measure in a stable environment. Allow the thermometer to acclimatize to room temperature for 30 minutes if moved from a different environment.
"Lo" or "Hi" displayed.	Temperature is outside the measurable range.	Ensure the thermometer is used within its specified operating temperature range.
Display shows an error message (e.g., "ERR").	Device malfunction or internal error.	Remove batteries, wait 1 minute, then reinsert. If the error persists, contact customer support.

8. SPECIFICATIONS

- **Model:** NC200
- **Measurement Range:**
 - Body Mode: 34.0 - 43.0 °C (93.2 - 109.4 °F)
 - Object Mode: 0 - 100.0 °C (32 - 212.0 °F)
- **Accuracy:**
 - Body Mode: ±0.2 °C (±0.4 °F) between 35.0 - 42.0 °C (95.0 - 107.6 °F)
 - Object Mode: ±1.0 °C (±2.0 °F)
- **Display:** LCD with backlight (FeverGlow Technology)
- **Memory:** 30 readings with date and time
- **Power Source:** 2 x 1.5V AAA batteries
- **Battery Life:** Approximately 2000 measurements
- **Dimensions:** 1.83 x 1.69 x 6.17 inches (46.5 x 42.9 x 156.7 mm)
- **Weight:** 4 ounces (113.4 g) (with batteries)
- **Operating Conditions:** 10 - 40 °C (50 - 104 °F), 15 - 95% relative humidity (non-condensing)
- **Storage Conditions:** -20 - +50 °C (-4 - +122 °F), 15 - 95% relative humidity (non-condensing)
- **Connectivity Technology:** Infrared

- **Special Feature:** High Accuracy, Non-Contact, Fast Read, Auto Measurement with Distance Control
- **Material:** Plastic

9. WARRANTY AND SUPPORT

9.1. Limited Warranty

Microlife provides a limited warranty for this device. For detailed warranty terms and conditions, please refer to the warranty card included with your product or visit the official Microlife website. The warranty typically covers manufacturing defects under normal use.

9.2. Customer Support


If you have any questions, encounter issues with your device, or require technical assistance, please contact Microlife customer support. Contact information can usually be found on the product packaging, the official Microlife website, or the included warranty documentation.

For additional resources, you may refer to the official Microlife website: www.microlifeusa.com

© 2023 Microlife. All rights reserved.

Related Documents - NC200

	Microlife NC200 Non-Contact Infrared Thermometer User Manual Comprehensive user manual for the Microlife NC200 Non-Contact Infrared Thermometer, detailing its features, operation, safety guidelines, technical specifications, and troubleshooting. Learn how to accurately measure body and surface temperatures.
	Microlife NC 400 Non-Contact Thermometer User Manual User manual for the Microlife NC 400 non-contact infrared thermometer, providing instructions on usage, features, safety, and technical specifications.
	Microlife NC150 BT Bluetooth Infrared Thermometer User Manual Comprehensive user manual for the Microlife NC150 BT Bluetooth infrared thermometer, detailing features, usage instructions, safety precautions, technical specifications, and troubleshooting.
	Microlife NC 150 BT Non-Contact Infrared Thermometer - User Manual & Guide Detailed user manual for the Microlife NC 150 BT Non-Contact Infrared Thermometer. This guide covers features like Bluetooth connectivity with the Connected Health+ app, operation, safety instructions, technical specifications, and troubleshooting for accurate and hygienic temperature measurements.
	Microlife NC150 BT Non Contact Thermometer - Accurate and Hygienic Temperature Measurement Discover the Microlife NC150 BT, a state-of-the-art non-contact thermometer offering fast, accurate, and hygienic temperature readings. Features Bluetooth connectivity for health tracking via the Microlife Connected Health+ app.

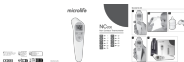	Microlife NC 200 Non-Contact Thermometer User Manual and Instructions User manual for the Microlife NC 200 non-contact infrared thermometer. This guide provides detailed instructions on how to use the device, safety precautions, troubleshooting, and technical specifications for accurate and hygienic temperature measurement.
---	---

Documents - Microlife – NC200



[\[pdf\]](#) User Manual Instructions

Layout 1 Jennifer McRae Forehead thermometer Microlife instruction manual Measurement in a matter of seconds The innovative infrared technology can take measurement without touching the object This guarantees safe and sanitary measurements within 3 Auto with distance control automatically when Be aware that basic physiological effect called vasoconstriction occur early stages fever resulting cool skin measured temperatures using this during NC200 Instruction Manual 041720 W microlifeusa 2020 11 User Manual Forehead thermometer with distance control **NC200** Non-Contact Thermometer Important Product and Safety Information Follow Instructions for Use. This document provides important product information regarding this thermometer. Please read this operation and safety information thoroughl...

lang:en score:44 filesize: 372.88 K page_count: 53 document date: 2020-11-10

microlife

Forehead
thermometer
with distance
control

User Manual



NC200
Non-Contact
Thermometer

[\[pdf\]](#) User Manual Instructions

var www.html.sideloads.microlife.nc200.manual.original.internal. |||

User Manual Forehead thermometer with distance control **NC200** Non-Contact Thermometer Important Product and Safety Information Follow Instructions for Use. This document provides important product information regarding this thermometer. Please read this operation and safety information thoroughl...

lang:en score:42 filesize: 372.88 K page_count: document date:

microlife

Forehead
thermometer
with distance
control

User Manual



NC₂₀₀
Non-Contact
Thermometer

[\[pdf\]](#) User Manual Instructions

var www.html.sideloads.microlife.nc200.manual.optimized.internal.iii

User Manual Forehead thermometer with distance control **NC200** Non-Contact Thermometer Important Product and Safety Information Follow Instructions for Use. This document provides important product information regarding this thermometer. Please read this operation and safety information thoroughl...

lang:en score:42 filesize: 378.15 K page_count: document date:



March 1, 2022

Microlife Intellectual Property GmbH
% Vaibhav Rajal
Official Correspondent for Microlife Intellectual Property GmbH
reel Consultants, Inc.
25 Northern Blvd., Suite 200
Great Neck, New York 11021

Re: K211776
Trade/Device Name: Microlife Non-Contact Infrared Forehead Thermometer, Model FR1MF1-B
(NC150 BT)
Regulation Number: 21 CFR 880.2910
Regulation Name: Class II Electronic Thermometer
Regulatory Class: Class II
Product Code: FLL
Date: January 27, 2022
Received: January 26, 2022

Dear Vaibhav Rajal:

We have reviewed your Section 510(k) (premarket notification of intent to market the device referenced above) and have determined the device is substantially equivalent (for the indications for use stated in the notification) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976; the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. Although this letter refers to your product as a device, please be aware that some cleared products may instead be combination products. The 510(k) Premarket Notification Database located at <https://www.accessdata.fda.gov/cdrh/devices/fda/510knew.htm> identifies combination product submissions. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration. Please note: CDRL does not evaluate information related to contract liability warranties. We cannot, you, however, that device labeling must be truthful and not misleading.

If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 900. In addition, FDA may publish further announcements concerning your device in the [Federal Register](#).

Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal

U.S. Food & Drug Administration
10801 New Hampshire Avenue
Silver Spring, MD 20910
FDA@FDA.GOV
2022-03-01

[\[pdf\]](#) User Manual

March 1 2022 Microlife Intellectual Property GmbH Vaibhav mar — Trade Device Name Non Contact Infrared Forehead Thermometer As states in Draft Instruction Manual for FR1MF1 B NC150 BT Page 12 K211776 accessdata.fda.gov/cdrh/docs/21.iii

March 1, 2022 Microlife Intellectual Property GmbH Vaibhav Rajal Official Correspondent for Microlife ... R 880.2910 FLL Class II K191829, Microlife Non-Contact Infrared Forehead Thermometer, Model FR1DG1 **NC200** Device Description and Mode of Operation The Microlife Non-Contact Infrared Forehead Thermom... lang:en score:42 filesize: 400.07 K page_count: 11 document date: 0000-00-00

[Microlife NC 200 Non-Contact Thermometer: Accurate Temperature Measurement](#)

Discover the Microlife NC 200 Non-Contact Thermometer for fast, hygienic, and accurate temperature readings. This guide covers its key features, usage, and benefits for home health monitoring.

lang:en score:40 filesize: 1.02 M page_count: 64 document date: 2020-11-19





[\[pdf\]](#) User Manual Instructions

Instruction Manual Jennifer McRae Microlife USA Forehead thermometer Axillar 94 5 99 1 °F 34 7 37 3 °C por adultos niños y bebés Comuníquese con su médico en caso de fiebre prolongada Follow Instructions for Use This document provides important product information regard ing this Please read operation and safety thoroughly before using the device keep future reference digital infrared noncontact forehead is a clinical an sensor to NC200 MLU IB FINAL 041720 v2 W microlifeusa 2021 03 User Manual Forehead thermometer with distance control **NC200** Non-Contact Thermometer Important Product and Safety Information Follow Instructions for Use. This document provides important product information regarding this thermometer. Please read this operation and safety information thoroughly...

lang:en **score:39** filesize: 372.56 K page_count: 52 document date: 2021-03-16



[\[pdf\]](#) User Manual Instructions Recall

Untitled Instruction manual NC 200 0 9 MB Microlife AG ib nc n v10 0222 microlife media 2403 IB N V10 v 11 |||

Microlife AG Espenstrasse 139 9443 Widnau / Switzerland www.microlife.com Microlife UAB P. Luksio g. 32, 08222 Vilnius / Lithuania IB NC 200 N-V10 0222 Revision Date: 2021-12-15 Microlife NC 200 1 2a 2b **NC200** Non Contact Thermometer EN 2 LV 32 SV 8 LT 38 3 4 FI 14 EE 44 DA 20 RU ...

lang:en **score:37** filesize: 899.27 K page_count: 63 document date: 2022-01-10



[\[pdf\]](#) User Manual Instructions Recall

IB NC 200 N V10 2522 joerg hezel Bedienungsanleitung 2 MB Instruction manual Microlife AG ib nc n v10 microlife media 2403 v 12 |||

Microlife AG Espenstrasse 139 9443 Widnau Switzerland www.microlife.com Microlife UAB P. Luksio g. 32 08222 Vilnius Lithuania IB NC 200 N-V10 2522 Revision Date: 2022-05-31 Microlife NC 200 1 2a 2b **NC200** Non Contact Thermometer EN 2 LV 32 SV 8 LT 38 3 4 FI 14 ET 44 DA 20 RU 50 ...

lang:en **score:37** filesize: 2.12 M page_count: 62 document date: 2022-07-21



[Microlife NC 200 Non-Contact Infrared Thermometer User Manual and Specifications](#)

Comprehensive user guide for the Microlife NC 200 non-contact infrared thermometer, covering features, safety instructions, operation, error messages, and technical specifications. Learn how to accurately measure body and object temperatures.


lang:en **score:36** filesize: 723.66 K page_count: 16 document date: 2023-02-23



[\[pdf\]](#) Borchure

Microlife NC200 Non Contact Thermometer ENG NC 200 brochure medidyne dk Brochure thermometer berøringsfri Medidyne |||

3sec Automatic measurement **NC200** Non Contact Thermometer Key features
 Clinically tested Several studies conducted in Austria prove the high accuracy and reliability of Microlife s infrared thermometers. Auto measurement with distance control CN Pat. ZL 201410095361.7; Others are patent pending T...
 lang:en score:30 filesize: 723.73 K page_count: 1 document date: 2018-08-30

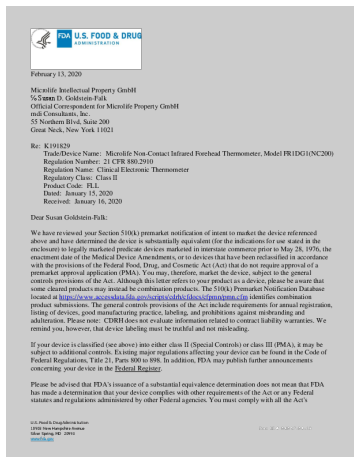

[\[pdf\]](#)

February 13 2020 accesdata fda gov Thermometer Model FR1DZ1 Device Description and Mode of Operation The Microlife Non Contact Infrared Forehead FR1DG1 NC200 is an electronic thermometer using infrared sensor to measure energy radiated from the forehead This collected through lens converted a temperature value K191829 cdrh docs 19 |||
 February 13, 2020 Microlife Intellectual Property GmbH Susan D. Goldstein-Falk
 Official Corresponde ... 21 Re: K191829 Trade/Device Name: Microlife Non-Contact Infrared Forehead Thermometer, Model FR1DG1 **NC200** Regulation Number: 21 CFR 880.2910 Regulation Name: Clinical Electronic Thermometer Regulator...
 lang:en score:29 filesize: 451.19 K page_count: 11 document date: 0000-00-00

[\[pdf\]](#) User Manual Instructions Recall

OBJ BUCH 448 001 book anita steck manual microlife ib nc 200 en pl 3720 media 4563 |||
 Microlife AG Espenstrasse 139 9443 Widnau / Switzerland www.microlife.com CHDE
 Polska S.A. Biesiadna 7, 35-304 Rzeszów Tel. 0-17 229-37-89 Fax. 0-17 230-21-14 E-mail: biuro.rzeszow chde.pl BDO 000046780 IB NC 200 EN-PL 3720 Revision Date: 2020-09-15 Microlife NC 200 1 2a 2b **NC200** Termometr...
 lang:en score:27 filesize: 1.51 M page_count: 13 document date: 2020-09-15

[\[pdf\]](#) User Manual Instructions Recall

OBJ BUCH 448 001 book joerg hezel manual microlife ib nc 200 en pl 2020 media 2911 |||
 Microlife AG Espenstrasse 139 9443 Widnau / Switzerland www.microlife.com CHDE
 Polska S.A. Biesiadna 7, 35-304 Rzeszów Tel. 0-17 229-37-89 Fax. 0-17 230-21-14 E-mail: biuro.rzeszow chde.pl BDO 000046780 IB NC 200 EN-PL 2020 Revision Date: 2020-05-12 Microlife NC 200 1 2a 2b **NC200** Termometr...
 lang:en score:27 filesize: 1.26 M page_count: 13 document date: 2020-05-12



[\[pdf\]](#) User Manual Instructions Recall

OBJ BUCH 525 001 book anita steck manual microlife ib nc 200 sk v3 3320 media 4195 |||

Microlife AG Espenstrasse 139 9443 Widnau / Switzerland www.microlife.com IB NC 200 SK-V3 3320 Revision Date: 2020-08-28 Microlife NC 200 1 2a 2b **NC200** Non Contact Thermometer EN 2 SK 8 CZ 14 3 4 Microlife NC 200 2 1 3 4 6 8 9 AQ AT AR AK AS AL BT AM BK AN BL AO BM AP...

lang:en **score:26** filesize: 1.97 M page_count: 19 document date: 2020-09-02





[\[pdf\]](#)

3sec MSA Microlife AG Espenstrasse 139 9443 Widnau Switzerland microlife Auto measurement with distance control This thermometer takes automatically a GB NC 200 N V10 2017 msa sm ee ctrl en Fail laadi alla 64509 |||

3sec Automatic measurement **NC200** Non Contact Thermometer EN Non Contact Thermometer Touch-free measurement Gentle and hygienic Auto measurement with distance control SV Berringsfri termometer Mter temperaturer utan kontakt Anvndarvnlig och hygienisk Automatisk mtning med avstnd kontroll F...

lang:en **score:25** filesize: 1.05 M page_count: 1 document date: 2019-04-10

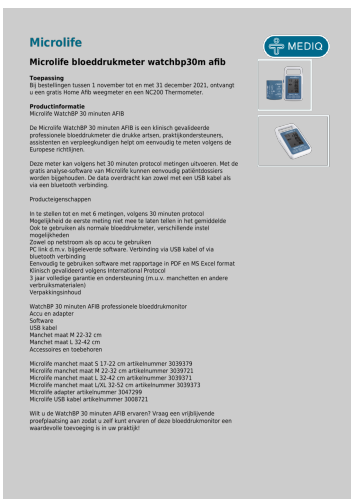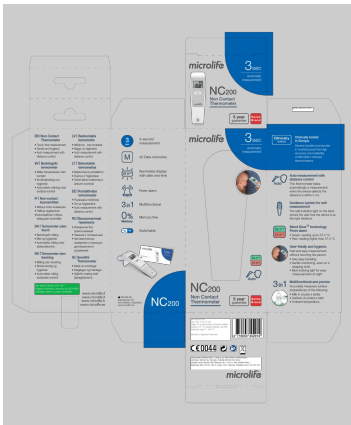

[\[pdf\]](#) Catalog

Download pagina Microlife bloeddrukmeter watchbp30m afib Catalogus Pagina medbis nl catalog home product 6851 html |||

Microlife Microlife bloeddrukmeter watchbp30m afib Toepassing Bij bestellingen tussen 1 november tot en met 31 december 2021, ontvangt u een gratis Home Afib weegmeter en een **NC200** Thermometer. Productinformatie Microlife WatchBP 30 minuten AFIB De Microlife WatchBP 30 minuten AFIB is een klinisch g...

lang:nl **score:23** filesize: 1.22 M page_count: 2 document date: 2023-06-03

[\[pdf\]](#) Instructions Recall

OBJ BUCH 448 001 book DSO instrukcja obsługi microlife nc 200 termometr bezdotykowy Read the instructions carefully before using this device Type BF applied part This Microlife thermometer is a high quality product incorporating metler pl img cms Termometr bezkontaktowy Metler AG Espenstrasse 139 9443 Widnau Switzerland NC 1 2a 2b 3 4 EN 2 PL 8 Guarantee Card 9 AT AK AL AM AN AO AP AQ AR AS BT BK BL BO BP BR BS BQ BM BN 6 5 7 Name of Purchaser Imię i nazwisko nabywcy Serial Number Numer seryjny |||

Microlife AG Espenstrasse 139 9443 Widnau / Switzerland www.microlife.com IB NC 200 EN-PL 4917 Microlife NC 200 1 2a 2b **NC200** Termometr bezkontaktowy EN 2 PL 8 3 4 Microlife NC 200 2 1 4 6 3 8 9 AQ AT AR AK AS AL BT AM BK AN BL AO BM AP BN 5 7 BO BP BQ BR BS Microlife ...

lang:en **score:23** filesize: 787.86 K page_count: 13 document date: 2017-12-04



[\[pdf\]](#) Instructions Recall

OBJ BUCH 620 001 book anita steck Термометр медичний безконтактний Microlife ib nc 200 ua v3

4120 microlife media 4274 IB NC UA V3 v 1 |||

: TOB , . a, . 17, . , 49044, , .: 0562 36-77-58, office vega-healthlife.com UA.TR.121 :

TOB ... om, www.microlife.ua IB NC 200 UA-V3 4120 Revision Date: 2020-10-27

Microlife NC 200 1 2a 2b **NC200** Microlife EN 2 RU 10 UA 18 3 4 Microlife NC 200 2 1 3 4 6 8 9 AQ AT AR AK ...

lang:en **score:21** filesize: 2.9 M page_count: 25 document date: 2020-10-27



[\[pdf\]](#) Instructions

IB NC 200 e SR V2 MDR 5124 PREUZMITE UPUSTVOBeskontaktni toplomer R B Medical

Companybeskontaktni nc200 0425rbmedical rs 2025 beskontaktni 0425 |||

Microlife Corporation 9F, No. 431, RuiGuang Road, NeiHu Taipei 114, Taiwan, China

www.microlife.com Microlife UAB P. Luksio g. 32 08222 Vilnius Lithuania IB NC 200-e

SR-V3 5124 Revision Date: 2024-11-19 Microlife NC 200 1 2a 2b **NC200** Non

Contact Thermometer SR 3 4 Microlife NC 200 2 1 3 ...

lang:af **score:20** filesize: 668.71 K page_count: 10 document date: 2025-04-11



[\[pdf\]](#) Instructions Recall

OBJ BUCH 445 001 book joerg hezel NC200 Non Contact Thermometer 18 nov 2020 — IB NC 200 e E

V11 5020 ib nc v11 microlife media 2407 v 5

Microlife AG Espenstrasse 139 9443 Widnau / Switzerland www.microlife.com IB NC

200-e E-V11 5020 Revision Date: 2020-11-18 Microlife NC 200 1 2a 2b **NC200** Non

Contact Thermometer EN 2 SL 42 RU 8 SR 48 3 4 BG 16 HU 54 RO 24 HR 60 CZ 30

PL 66 SK 36 Microlife NC 200 2 1 3 4 6...

lang:en **score:20** filesize: 6.16 M page_count: 71 document date: 2020-11-18



[\[pdf\]](#) Instructions Recall

OBJ BUCH 456 001 book anita steck Non Contact Thermometer c86184308eec6171f18a80512cd2e89e

microlife shop by upload iblock c86 IB NC 200 BY V3 3618 Microlife AG Espenstrasse 139 9443 Widnau

Switzerland 1 2a 2b 3 4 |||

Microlife AG Espenstrasse 139 9443 Widnau / Switzerland www.microlife.com IB NC


200 BY-V3 3618 Microlife NC 200 1 2a 2b **NC200** Non Contact Thermometer EN 2

RU 8 BG 14 3 4 www.microlife.by www.microlife.bg www.microlife.com Microlife NC

200 2 1 3 4 6 8 9 AQ AT AR AK AS AL BT ...

lang:en **score:20** filesize: 2 M page_count: 19 document date: 2019-02-07



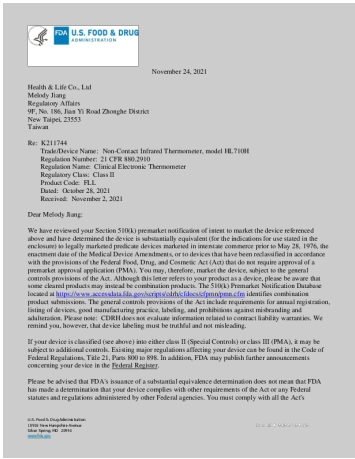

[pdf]

K211744 accessdata fda gov cdrh docs 21 |||

November 24, 2021 Health Life Co., Ltd Melody Jiang Regulatory Affairs 9F, No. 186, Jian Yi Road Z ... eared Device Predicate Device : Microlife Non-Contact Infrared Forehead Thermometer, Model: FR1DGI **NC200** KI 91829 1 4. Device Description: The Non-Contact Infrared Forehead Thermometer, Model HL...


lang:en score:19 filesize: 914.3 K page_count: 14 document date: 0000-00-00



[Microlife NC200 Non-Contact Infrared Thermometer User Manual](#)

Comprehensive user manual for the Microlife NC200 Non-Contact Infrared Thermometer, detailing its features, operation, safety guidelines, technical specifications, and troubleshooting. Learn how to accurately measure body and surface temperatures.


lang:en score:19 filesize: 750.4 K page_count: 16 document date: 2023-03-01

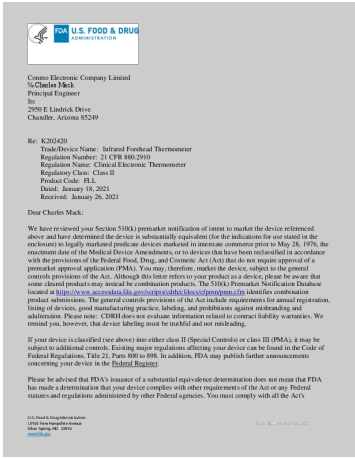

[pdf]

Conmo Electronic Company Limited % Charles Mack Principal K202420 accessdata fda gov cdrh docs 20 Conmo Electronic Company Limited Charles Mack Principal Engineer Irc 2950 E Lindrick Drive Chandler ... Property GmbH Predicate Device Microlife Non-Contact Infrared Forehead Thermometer, Model FR1DGI **NC200** 510 k Number K191829 4. Device Description: The infrared forehead thermometer, Model DPT-I...

lang:en score:19 filesize: 463.3 K page_count: 12 document date: 0000-00-00

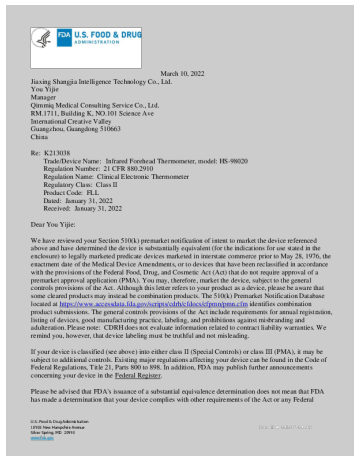

[pdf]

K213038 March 2022 510 K Clearances FDA accessdata fda gov cdrh docs 21 |||


March 10, 2022 Jiaxing Shangjia Intelligence Technology Co., Ltd. You Yijie Manager
Qimmiq Medical C ... 91829 Trade Name: Microlife Non-Contact Infrared Forehead
Thermometer Page 1/7 Model: FR1DG1 **NC200** Classification name: Clinical
electronic thermometer Review panel: General Hospital Produc...
lang:en **score:18** filesize: 521.44 K page_count: 11 document date: 0000-00-00

[pdf] Instructions

OBJ BUCH 453 001 book Joerg hezel NC 200 Microlife AG ib nc s v11 0520 microlife media 2416 IB S
V11 v 4 |||

Microlife AG Espenstrasse 139 9443 Widnau / Switzerland www.microlife.com IB NC
200 S-V11 0520 Microlife NC 200 1 2a 2b **NC200** Non Contact Thermometer EN 2 NL
38 ES 8 GR 44 3 4 FR 14 TR 50 IT 20 AR 56 DE 26 FA 62 PT 32 Microlife NC 200 2 1
3 4 6 8 9 AQ AT AR AK AS A...

lang:en **score:17** filesize: 6.4 M page_count: 70 document date: 2020-01-24

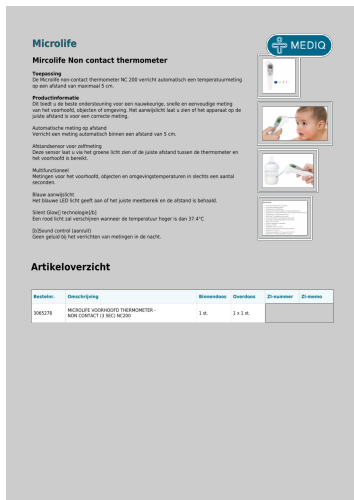

[pdf] Catalog

Download pagina Mircolife Non contact thermometer Catalogus Pagina m medbis nl catalog home
product 7033 html |||

Microlife Mircolife Non contact thermometer Toepassing De Microlife non-contact
thermometer NC 200 v ... tikeloverzicht Bestelnr. 3065278 Omschrijving MICROLIFE
VOORHOOFD THERMOMETER NON CONTACT 3 SEC **NC200** Binnendoos
Overdoos ZI-nummer ZI-memo 1 st. 1 x 1 st. ...

lang:nl **score:17** filesize: 2 M page_count: 1 document date: 2023-06-06

[pdf] Instructions Recall

OBJ BUCH 453 001 book DSO IB NC 200 S V11 3217 Microlife AG Espenstrasse 139 9443 Widnau
Switzerland microlife Non Contact Thermometer 1 2a 2b 3 4 Handleiding bloeddrukmeterswebshop nl
amfilerating file id 89 |||

Microlife NC 200 1 2a 2b 3 4 Microlife AG Espenstrasse 139 9443 Widnau /
Switzerland www.microlife.com IB NC 200 S-V11 3217 **NC200** Non Contact
Thermometer EN 2 ES 8 FR 14 IT 20 DE 26 PT 32 NL 38 GR 44 TR 50 AR 56 FA 62
Microlife NC 200 2 1 3 4 6 8 9 AQ AT AR AK AS AL BT ...

lang:en **score:17** filesize: 1.45 M page_count: 67 document date: 2017-09-04



[\[pdf\]](#) Instructions

OBJ BUCH 453 001 book DSO Karabinis Medical IB NC 200 S V11 5218 Microlife AG Espenstrasse 139 9443 Widnau Switzerland microlife Non Contact Thermometer 1 2a 2b 3 4 karabinismedical gr gggg 2019 01 |||



Microlife AG Espenstrasse 139 9443 Widnau / Switzerland www.microlife.com IB NC 200 S-V11 5218 Microlife NC 200 1 2a 2b **NC200** Non Contact Thermometer EN 2 ES 8 FR 14 IT 20 DE 26 PT 32 NL 38 GR 44 TR 50 AR 56 FA 62 3 4 Microlife NC 200 2 1 4 6 3 8 9 AQ AT AR AK AS AL BT ...

lang:en **score:17** filesize: 1.53 M page_count: 68 document date: 2019-01-29

[\[pdf\]](#) Instructions

OBJ BUCH 449 001 book DSO NC 200 Microlife AG ib nc n v10 2419 microlife media 2403 IB N V10 v 5 |||



Microlife AG Espenstrasse 139 9443 Widnau / Switzerland www.microlife.com IB NC 200 N-V10 2419 Microlife NC 200 1 2a 2b **NC200** Non Contact Thermometer EN 1 LV 32 SV 8 LT 38 3 4 FI 14 EE 44 DA 20 RU 50 NO 26 IS 56 Microlife NC 200 2 1 3 4 6 8 9 AQ AT AR AK AS AL BT A...

lang:en **score:17** filesize: 1.45 M page_count: 62 document date: 2019-06-21

[\[pdf\]](#) Instructions Recall

OBJ BUCH 449 001 book joerg hezel NC 200 Microlife AG ib nc n v10 2020 microlife media 2403 IB N V10 v 7 |||



Microlife AG Espenstrasse 139 9443 Widnau / Switzerland www.microlife.com IB NC 200 N-V10 2020 Revision Date: 2020-05-07 Microlife NC 200 1 2a 2b **NC200** Non Contact Thermometer EN 2 LV 32 SV 8 LT 38 3 4 FI 14 EE 44 DA 20 RU 50 NO 26 IS 58 Microlife NC 200 2 1 3 4 6 8 9 AQ ...


lang:en **score:16** filesize: 5.43 M page_count: 63 document date: 2020-05-07

[\[pdf\]](#)

K211752 January 2022 510 K Clearances FDA accessdata fda gov cdrh docs 21 |||

January 21, 2022 Beijing Choice Electronic Technology Co., Ltd. Haiying Zhao
Official Correspondent ... r Device Trade/Proprietary Name: Microlife Non-Contact
Infrared Forehead Thermometer Model: FR1DG1 **NC200** Product Code: FLL
Regulation Number: 21 CFR 880.2910 Device Class: Panel: General Hospital M...

lang:en **score:16** filesize: 499.35 K page_count: 18 document date: 0000-00-00

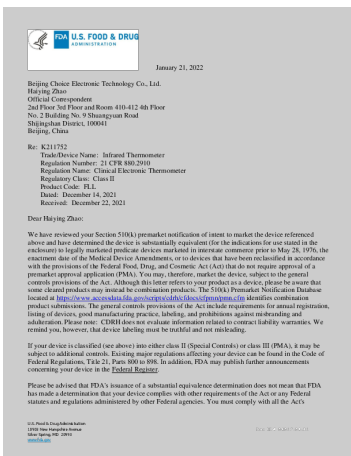


[\[pdf\]](#) Instructions Recall

Untitled Инструкция NC 200 1 5 MB инструкция Microlife Микролайф ib nc e s v11 1223 filename* utf 8 IB S V11 microlife media 2416 v 10 |||



Microlife AG Espenstrasse 139 9443 Widnau Switzerland www.microlife.com Microlife UAB P. Luksio g. 32 08222 Vilnius Lithuania IB NC 200-e S-V11 1223 Revision Date: 2023-02-24 Microlife NC 200 1 2a 2b **NC200** Non Contact Thermometer EN 2 NL 38 ES 8 GR 44 3 4 FR 14 TR 50 IT 20 AR 56...

lang:en **score:16** filesize: 1.51 M page_count: 67 document date: 2023-02-24

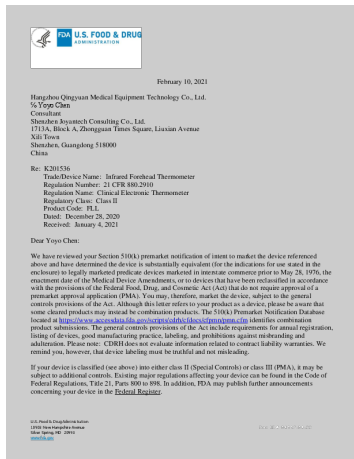

[pdf]

K201536 February 2021 510 K Clearances FDA accessdata fda gov cdrh docs 20 |||

February 10, 2021 Hangzhou Qingyuan Medical Equipment Technology Co., Ltd.
Yoyo Chen Consultant She ... lectual Property GmbH Device name Microlife Non-Contact Infrared Forehead Thermometer Model FR1DG1 **NC200** 510 K Number: K191829 Product Code FLL 4. Device Description The infrared forehead thermomet...
lang:en **score:16** filesize: 520.73 K page_count: 10 document date: 0000-00-00

[pdf] Instructions

OBJ BUCH 453 001 book DSO ib nc 200 s v11 2619 microlife media 2416 IB NC S V11 v 3

Microlife AG Espenstrasse 139 9443 Widnau / Switzerland www.microlife.com IB NC 200 S-V11 2619 Microlife NC 200 1 2a 2b **NC200** Non Contact Thermometer EN 2 NL 38 ES 8 GR 44 3 4 FR 14 TR 50 IT 20 AR 56 DE 26 FA 62 PT 32 Microlife NC 200 2 1 3 4 6 8 9 AQ AT AR AK AS A...
lang:en **score:16** filesize: 1.17 M page_count: 68 document date: 2019-06-25



[pdf]

Confirmation Letter Reference Mar 27 2024 · This letter confirms that SGS Belgium NV a Notified Body NB designated against Regulation EU 2017 745 MDR and identified by the c32868 gimaitaly DocumentiGIMA CE |||
Microlife Corporation 9F, No. 431, RuiGuang Road, NeiHu, Taipei, 114, Taiwan, R.O.C. 27Th Marc ... 0404 882 750 Belfius 550-3560000-93 ONBO Electronic Shenzhen Co., Ltd.: 697038816IR3D FR1DG1 **NC200** FR1DS1 NC400 FR1MA1 NC100 FR1MF1 NC150 FR1MF1-B NC150 BT FR1MY1 NC300 IR1DF1-1 IR...
lang:en **score:16** filesize: 2.87 M page_count: 12 document date: 2024-05-16

[pdf] Instructions Recall

Untitled ib nc 200 e v11 0222 microlife ru media 2407 IB NC E V11 |||


Microlife AG Espenstrasse 139 9443 Widnau / Switzerland www.microlife.com
Microlife UAB P. Luksio g. 32, 08222 Vilnius / Lithuania IB NC 200-e E-V11 0222
Revision Date: 2022-01-18 Microlife NC 200 1 2a 2b **NC200** Non Contact Thermometer EN 2 SL 42 RU 8 SR 48 3 4 BG 16 HU 54 RO 24 H...
lang:en **score:16** filesize: 1.05 M page_count: 71 document date: 2022-01-18

[pdf] Instructions Recall

OBJ BUCH-453-001 book joerg hezel ib-nc-200-e-s-v11-2820 https microlife media 2416 download IB NC 200-e S-V11 2820 v 6

Microlife AG Espenstrasse 139 9443 Widnau / Switzerland www.microlife.com IB NC 200-e S-V11 2820 Revision Date: 2020-07-28 Microlife NC 200 1 2a 2b **NC200** Non Contact Thermometer EN 2 NL 40 ES 8 GR 46 3 4 FR 14 TR 52 IT 20 AR 58 DE 26 FA 64 PT 32 Microlife NC 200 2 1 3 4 6...
lang:en **score:16** filesize: 5.92 M page_count: 69 document date: 2020-07-28




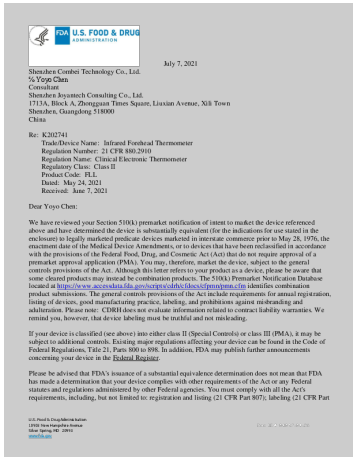

[pdf]

K202741 July 2021 510 K Clearances FDA accessdata fda gov cdrh docs 20 |||

July 7, 2021 Shenzhen Combei Technology Co., Ltd. Yoyo Chen Consultant
Shenzhen Joyantech Consultin ... e Microlife Intellectual Property GmbH Microlife
Non-Contact Infrared Forehead Thermometer FR1DG1 **NC200** K191829 880.2910
Class II FLL 4. Device Description The Infrared Forehead Thermometer, mode...
lang:en **score:15** filesize: 447.28 K page_count: 9 document date: 0000-00-00

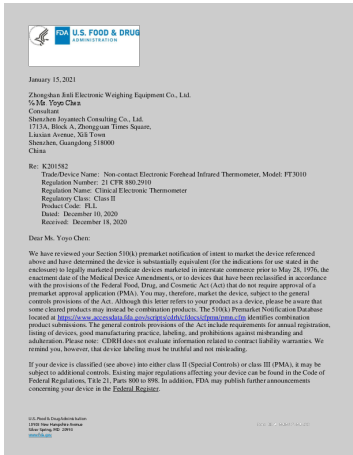

[pdf]

K201582 January 2021 510 K Clearances FDA accessdata fda gov cdrh docs 20 |||

January 15, 2021 Zhongshan Jinli Electronic Weighing Equipment Co., Ltd. Ms. Yoyo
Chen Consultant S ... lectual Property GmbH Device name Microlife Non-Contact
Infrared Forehead Thermometer Model FR1DG1 **NC200** 510 K Number: K191829
Product Code FLL 4. Device Description Non-contact Electronic Forehead...
lang:en **score:15** filesize: 800.58 K page_count: 11 document date: 0000-00-00

[pdf] Instructions Recall

sabine sagmeister ib nc 200 n v10 2523 filename* utf 8 IB NC N V10 microlife lv media 5277 v 2 |||

Microlife Corporation 9F, 431, RuiGuang Road, NeiHu Taipei 11492, Taiwan, China
www.microlife.com Microlife UAB P. Luksio g. 32 08222 Vilnius Lithuania IB NC 200
N-V10 2524 Revision Date: 2023-05-25 Microlife NC 200 1 2a 2b **NC200** Non Contact
Thermometer EN 2 LV 32 SV 8 LT 38 3 4 FI ...

lang:en **score:15** filesize: 1.78 M page_count: 63 document date: 2023-05-26

[Microlife NC 200 Non-Contact Thermometer User Manual](#)

User manual for the Microlife NC 200 Non-Contact Infrared Thermometer. Learn how
to accurately measure temperature, understand safety guidelines, and utilize features
like the fever alarm and memory recall. Visit Microlife for support.

lang:en **score:15** filesize: 2.04 M page_count: 63 document date: 2024-08-22

[pdf] Instructions Recall

IB NC 200 e S V11 3923 sabine sagmeister ib nc s v11 filename* utf 8 microlife media 2416 1223 |||

Microlife Corporation 9F, 431, RuiGuang Road, NeiHu Taipei 11492, Taiwan, China
www.microlife.com Microlife UAB P. Luksio g. 32 08222 Vilnius Lithuania IB NC 200-e
S-V11 3923 Revision Date: 2023-08-31 Microlife NC 200 1 2a 2b **NC200** Non Contact
Thermometer EN 1 PT 42 ES 7 NL 49 3 4 FR...

lang:en **score:15** filesize: 2.39 M page_count: 75 document date: 2023-08-31



[\[pdf\]](#) Instructions

IB NC 200 e E V12 4424 kurt steinebrunner Non Contact Thermometer ib nc v12 filename* utf 8 microlife
media 2407 v |||



Microlife UAB P. Luksio g. 32, 08222 Vilnius Lithuania Microlife Corporation 9F, 431, RuiGuang Road, NeiHu Taipei 11492, Taiwan, China www.microlife.com IB NC 200-e E-V12 4424 Revision Date: 2024-10-17 Microlife NC 200 1 2a 2b **NC200** Non Contact Thermometer EN 2 SL 40 RU 8 SR 46 3 4 BG...

lang:en score:15 filesize: 1.56 M page_count: 77 document date: 2024-10-17

[\[pdf\]](#) Instructions Recall

IB NC 200 e S V11 3824 sabine sagmeister Termómetro sin contacto con auto medición Microlife
AGmicrolife es media 2416 v 12 |||



Microlife Corporation 9F, 431, RuiGuang Road, NeiHu Taipei 11492, Taiwan, China
www.microlife.com Microlife UAB P. Luksio g. 32 08222 Vilnius Lithuania Microlife AG
Eспенstrasse 139 9443 Widnau / Switzerland IB NC 200-e S-V11 3824 Revision
Date: 2024-08-21 Microlife NC 200 1 2a 2b **NC200** Non ...


lang:en score:15 filesize: 2.04 M page_count: 63 document date: 2024-08-22

Microlife NC 200 Non-Contact Thermometer User Manual and Instructions

User manual for the Microlife NC 200 non-contact infrared thermometer. This guide provides detailed instructions on how to use the device, safety precautions, troubleshooting, and technical specifications for accurate and hygienic temperature measurement.



lang:en score:15 filesize: 2.68 M page_count: 75 document date: 2025-09-10

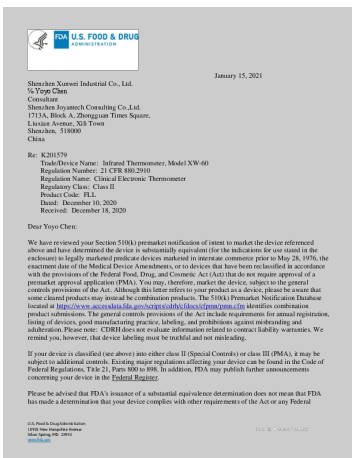

[\[pdf\]](#)

K201579 January 2021 510 K Clearances FDA accesdata fda gov cdrh docs 20 |||

Shenzhen Xunwei Industrial Co., Ltd. Yoyo Chen Consultant Shenzhen Joyantech Consulting Co.,Ltd. 17 ... Intellectual Property GmbH Device name Microlife Non-Contact Infrared Forehead Thermometer Model FR1DG1 **NC200** 510 K Number: K191829 Product Code FLL 4. Device Description The device is intended for the...

lang:en score:15 filesize: 989.42 K page_count: 11 document date: 0000-00-00

[\[pdf\]](#) Instructions Recall

Untitled Non Contact Thermometer ib nc 200 e s v11 0222 microlife pt media 2416 IB NC S V11 v 9 |||



Microlife AG Espenstrasse 139 9443 Widnau / Switzerland www.microlife.com EC
 REP Microlife UAB P. Luksio g. 32 08222 Vilnius / Lithuania IB NC 200-e S-V11 0222
 Revision Date: 2021-12-15 Microlife NC 200 1 2a 2b **NC200** Non Contact
 Thermometer EN 2 NL 38 ES 8 GR 44 3 4 FR 14 TR 50 IT...


lang:en score:15 filesize: 1022.35 K page count: 67 document date: 2022-01-11



[Microlife NC 200 Non-Contact Thermometer User Manual](#)

User manual for the Microlife NC 200 Non-Contact Thermometer, providing instructions, safety information, and technical specifications for accurate temperature measurement. Learn about its features, usage, and maintenance.

lang:en score:15 filesize: 5.46 M page_count: 63 document date: 2020-08-18

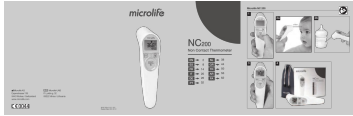

[Microlife NC 200 Non Contact Thermometer User Manual](#)

Comprehensive user manual for the Microlife NC 200 Non Contact Thermometer, detailing features, usage instructions, safety precautions, and technical specifications for accurate temperature measurement.

lang:en score:15 filesize: 1005.93 K page_count: 67 document date: 2021-02-12