

usmile P1

usmile P1 Sonic Electric Toothbrush User Manual

Model: P1

[Setup](#)

[Operation](#)

[Maintenance](#)

[Troubleshooting](#)

[Specifications](#)

[Warranty & Support](#)

1. SETUP

1.1 Unboxing and Components

Carefully unpack all items from the box. Ensure all components are present:

- Brush Handle (x1)
- Soft Brush Head (x1)
- Pro Brush Head (x1)
- User Manual (x1)
- Type-C Cable (x1)

EXPERIENCE AN AWARD-WINNING ORAL CARE

17+

International Awards

450+

Technical Patents

22 million+

Global Users



Japan Good
Design Award



Contemporary
Good Design



China Red Star
Design Award



Golden Pin
Design Award



iF Design Award



reddot award
Winner



CMF Design Award



CHEAA

Image: Package Contents. This image displays the usmile P1 Sonic Electric Toothbrush, including the brush handle, two brush heads (soft and pro), the user manual, and a USB Type-C charging cable, all neatly arranged in their packaging.

1.2 Initial Charging

Before first use, fully charge the toothbrush. Connect the USB Type-C cable to the charging port on the bottom of the toothbrush handle and plug the other end into a compatible USB power adapter (not included).

- A full charge takes approximately 4 hours.
- The battery indicator will show charging status and full charge.
- A single full charge can provide up to 6 months of use (based on 2 minutes, twice daily brushing).

1 Full Charge = 180 Days *



* at a using frequency of 2 times/day,
2 minutes/time

Image: Extended Battery Life. This image shows the usmile P1 electric toothbrush alongside a clock graphic, highlighting its impressive 180-day battery life on a single full charge, based on typical usage.

1.3 Attaching the Brush Head

Align the brush head with the metal shaft on the toothbrush handle and push it down firmly until it clicks into place. Ensure the brush head is securely attached before use.

2. OPERATION

2.1 Powering On/Off and Mode Selection

Press the power button to turn the toothbrush on. Press the button again to cycle through the available cleaning modes. Press and hold the power button to turn the toothbrush off.



Image: Cleaning Modes. This graphic illustrates the three distinct cleaning modes available on the usmile P1 toothbrush: 'Clean' for daily use, 'White' for teeth whitening, and 'Soft' suitable for beginners or sensitive gums.

2.2 Cleaning Modes

The usmile P1 offers 3 cleaning modes, which can be combined with 2 different brush heads for a total of 6 brushing experiences:

- **Clean:** For daily plaque removal and general cleaning.
- **White:** Designed for effective stain removal and teeth whitening.
- **Soft:** Gentle cleaning for sensitive teeth and gums.

Healthier Oral in 4 Weeks*



Image: Sonic Whitening Cleaning. This close-up image demonstrates the sonic cleaning action of the usmile P1 toothbrush, with water particles indicating high-frequency vibrations effectively cleaning a tooth surface.

2.3 Brushing Technique and Smart Timer

Apply a small amount of toothpaste to the wet brush head. Place the brush head against your teeth at a slight angle towards the gum line. Move the brush slowly across your teeth, allowing the sonic vibrations to do the work.

- The toothbrush features a **2-minute smart timer** to ensure you brush for the dentist-recommended duration.
- A **30-second zone change reminder** will gently pause the vibrations, prompting you to move to the next quadrant of your mouth.
- The toothbrush will automatically turn off after 4 minutes.

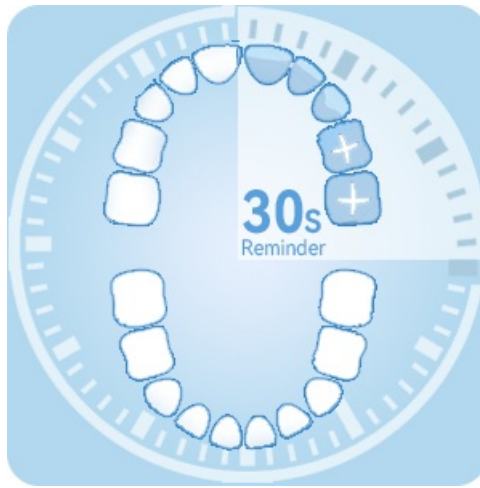


Image: Smart Timer Function. This diagram illustrates the 30-second zone change reminder feature, showing a mouth divided into quadrants to guide users in brushing each section evenly.

2.4 Official Product Videos

New Sonic Electric Toothbrush Extra Soft for adults

Video: This video demonstrates the usmile P1 Sonic Electric Toothbrush in action, highlighting its extra soft bristles suitable for adults. It showcases the brushing process and the gentle yet effective cleaning.

Extra Soft Sonic Toothbrush for Sensitive Gum and Oral Care

Video: A short demonstration focusing on the extra soft bristles of the usmile P1, emphasizing its suitability for sensitive gums and overall oral care. It shows the brush gently cleaning surfaces.

usmile Sonic Electric Toothbrush Battery Life Miracle

Video: This video highlights the exceptional battery life of the usmile P1 Sonic Electric Toothbrush, showcasing its long-lasting performance on a single charge.

Super Whitening of usmile electric toothbrush

Video: A brief video demonstrating the whitening capabilities of the usmile P1 electric toothbrush, showing its effectiveness in removing stains and enhancing tooth brightness.

3. MAINTENANCE

3.1 Cleaning the Toothbrush

Regular cleaning ensures optimal performance and hygiene.

- After each use, rinse the brush head thoroughly under running water.
- Remove the brush head and rinse the metal shaft area.

- Wipe the handle with a damp cloth. The usmile P1 is IPX7 waterproof, allowing for easy rinsing.

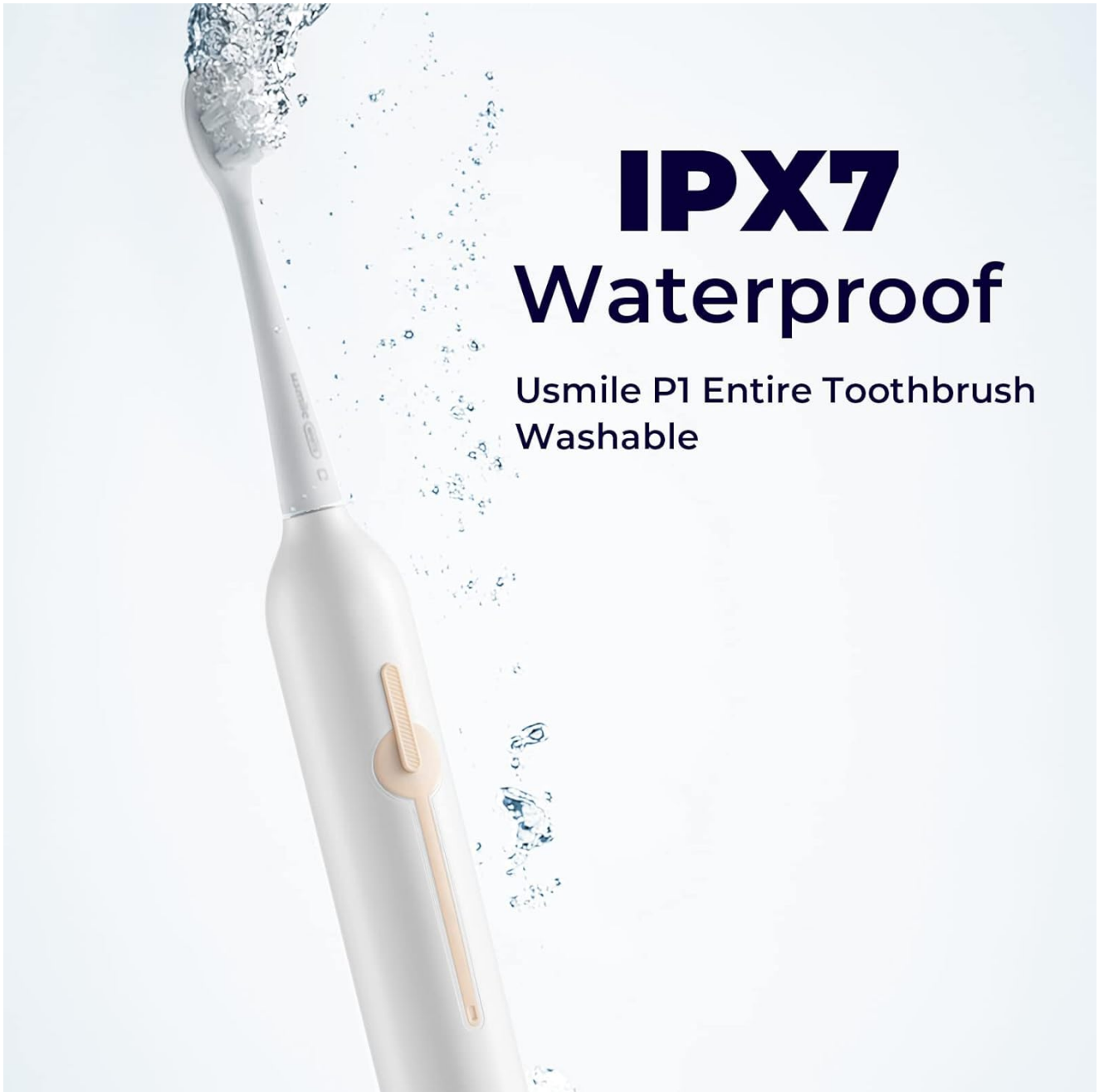


Image: IPX7 Waterproof Rating. This image shows the usmile P1 electric toothbrush being submerged in water, illustrating its IPX7 waterproof design which allows for easy and safe cleaning under running water.

3.2 Brush Head Replacement

Dentists recommend replacing your toothbrush head every 3 months, or sooner if the bristles become frayed. Replacement brush heads are available from usmile.

3.3 Storage

Store the toothbrush in a dry place when not in use. For travel, utilize the provided travel case to protect the brush head and handle.

TRAVEL READY

Travel Case



Image: Travel Ready Design. This image features the usmile P1 electric toothbrush and its travel case alongside travel accessories like a hat and sunglasses, emphasizing its portability and convenience for travel.

4. TROUBLESHOOTING

If you encounter any issues with your usmile P1 Sonic Electric Toothbrush, please refer to the following common solutions:

Problem	Possible Cause	Solution
Toothbrush does not turn on.	Low battery or no charge.	Charge the toothbrush for at least 4 hours. Ensure the charging cable is securely connected.
Reduced vibration intensity.	Worn-out brush head or low battery.	Replace the brush head. Recharge the toothbrush.

Problem	Possible Cause	Solution
Toothbrush stops during use.	Automatic 2-minute timer or 4-minute auto-off feature activated.	This is normal operation. Restart if needed for further cleaning.
Charging indicator not working.	Faulty cable or power adapter.	Try a different USB Type-C cable or power adapter. Ensure the charging port is clean.

5. SPECIFICATIONS

Key technical specifications for the usmile P1 Sonic Electric Toothbrush:

- **Model Name:** P1
- **Product Dimensions:** 9.92 x 1.14 x 1.18 inches
- **Weight:** 9.88 ounces
- **Power Source:** Battery Powered (1 Lithium Ion battery required)
- **Charging:** USB Type-C, 4-hour fast charge for up to 6 months use
- **Vibrations:** Up to 38,000 sonic vibrations per minute
- **Cleaning Modes:** 3 (Clean, White, Soft)
- **Timer:** 2-minute smart timer, 30-second zone change reminder, 4-minute auto-off
- **Waterproof Rating:** IPX7
- **Bristle Material:** Nylon (DuPont)
- **Handle Material:** Plastic
- **Color:** White

6. WARRANTY & SUPPORT

6.1 Warranty Information

usmile provides a **30-day money-back guarantee** and a **1-year warranty** for quality-related issues from the date of purchase. Please retain your proof of purchase for warranty claims.

6.2 Customer Support

For any questions, concerns, or warranty claims, please contact usmile customer service. You can find contact information on the official usmile website or through your retailer.

Please note: The warranty does not cover damage caused by misuse, accident, unauthorized modification, or normal wear and tear.

