

## Oral-B Pro 100

# Oral-B Pro 100 Charcoal Battery Powered Electric Toothbrush User Manual

BRAND: ORAL-B | MODEL: PRO 100

## 1. INTRODUCTION

The Oral-B Pro 100 Charcoal Battery Powered Electric Toothbrush is designed to provide a superior clean compared to a manual toothbrush. It features charcoal-infused bristles for effective stain removal and a rotating action for comprehensive plaque removal. This manual provides essential information for the proper setup, operation, and maintenance of your new toothbrush.

### Key Features:

- End-rounded, charcoal-infused bristles for effective cleaning.
- Cup-shaped bristles designed to remove surface stains from tea and coffee.
- Consistent battery performance with two AA batteries included.
- Compatible with various Oral-B brush heads including CrossAction, Sensitive Gum Care, FlossAction, and 3D White.
- Equipped with a 2-minute timer to encourage dentist-recommended brushing duration.





Figure 1: The Oral-B Pro 100 Charcoal Electric Toothbrush as packaged, highlighting its key features and battery-powered nature.

## 2. SETUP

### 2.1 Unboxing and Initial Inspection

Carefully remove all components from the packaging. Ensure that the toothbrush handle, one brush head, and two AA batteries are present and undamaged.

## 2.2 Installing Batteries

1. To open the battery compartment, firmly press the two small indentations on either side of the bottom cap and pull downwards.
2. Insert two AA batteries, ensuring correct polarity (+/-) as indicated inside the compartment.
3. Replace the cap by aligning it and pushing upwards until it clicks securely into place.



Figure 2: Illustration of how to open the battery compartment and insert the two AA batteries into the Oral-B Pro 100 handle.

## 2.3 Attaching the Brush Head

4. Align the brush head with the metal shaft on the toothbrush handle.
5. Push the brush head down firmly until it locks into place. Ensure the brush head is facing the same direction as the power button.





Figure 3: The Oral-B Pro 100 electric toothbrush with its charcoal-infused brush head securely attached to the handle.

### 3. OPERATING INSTRUCTIONS

---

#### 3.1 Brushing Your Teeth

5. Wet the brush head and apply a small amount of toothpaste.
6. Place the brush head in your mouth before turning it on to prevent splashing.
7. Press the power button to activate the toothbrush.
8. Gently guide the brush head from tooth to tooth, allowing the rotating bristles to clean each surface. Do not scrub aggressively.

# CONSTANT POWER

STRONGER BATTERY  
THROUGHOUT THE  
LIFE OF THE PRODUCT



Figure 4: Proper technique for using the Oral-B Pro 100 electric toothbrush, guiding the brush head gently across teeth surfaces.

## 3.2 Timer Function

The Oral-B Pro 100 includes a built-in 2-minute timer. The toothbrush will briefly pause its motor three times during your brushing session, indicating that you have brushed for the dentist-recommended duration. The toothbrush will not automatically turn off after 2 minutes, allowing you to continue brushing if desired.

## 4. MAINTENANCE

### 4.1 Cleaning Your Toothbrush

7. After each use, rinse the brush head thoroughly under running water with the handle still turned on.
8. Turn off the toothbrush and remove the brush head.
9. Rinse both the brush head and the handle separately under running water. Wipe the handle dry with a cloth.





**PRO 100**

**ROUND CLEANS  
BETTER**

Dentist-inspired round  
brush head delivers  
superior clean\*



**PROCORE  
POWER**

ProCore constant  
power delivers optimal  
performance every  
time

**2-MINUTE  
COACHING  
TIMER**

Helps keep your  
brushing on track

timer

**PREMIUM ORAL-B  
PRO 100 POWER**

Delivers 9,000  
revolutions per minute  
for superior plaque  
removal\*



\*vs. a regular manual  
toothbrush







Figure 5: Cleaning the Oral-B Pro 100 electric toothbrush by rinsing the brush head and handle after use.

## 4.2 Replacing the Brush Head

8. Dentists recommend replacing your brush head every 3 months, or sooner if the bristles become splayed or worn.
9. To remove the old brush head, simply pull it straight off the metal shaft.
10. Attach a new compatible Oral-B brush head as described in Section 2.3.



Figure 6: A selection of compatible Oral-B replacement brush heads, offering different cleaning benefits.



## 5. TROUBLESHOOTING

---

### 5.1 Reduced Brush Performance

If the brush head's rotation slows down significantly or the motor seems weaker, it indicates that the batteries are running low. Replace both AA batteries with new ones following the instructions in Section 2.2.

### 5.2 Difficulty Opening Battery Compartment

If you experience difficulty opening the battery compartment, ensure you are pressing both indentations simultaneously and firmly. Regular cleaning of the area around the battery cap and the rubber gasket can help prevent buildup that might hinder opening.

## 6. SPECIFICATIONS

---

Product Dimensions	2.17 x 1.5 x 9.45 inches
Item Weight	4.91 ounces
Item Model Number	10069055130510
Power Source	Battery Powered (2 AA batteries)
Special Feature	Charcoal Infused Bristle
Item Firmness Description	Extra Soft
Manufacturer	Procter & Gamble

## 7. SAFETY INFORMATION

---

For safety information regarding this product, please refer to the product packaging or the official Oral-B website. Always follow general electrical safety guidelines when handling battery-powered devices.

## 8. WARRANTY AND SUPPORT

---

For warranty details and customer support, please refer to the information provided on the product packaging or visit the official Oral-B website. Keep your purchase receipt as proof of purchase for any warranty claims.

## 9. LEGAL DISCLAIMER


---

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This product is intended for oral hygiene purposes only.



### [Oral-B Vitality Electric Toothbrush User Manual](#)

Comprehensive user manual for the Oral-B Vitality electric toothbrush, covering features, usage instructions, safety guidelines, brush head information, cleaning, and warranty. Includes model information for Charger Type 3757 and Handle Type 3710.

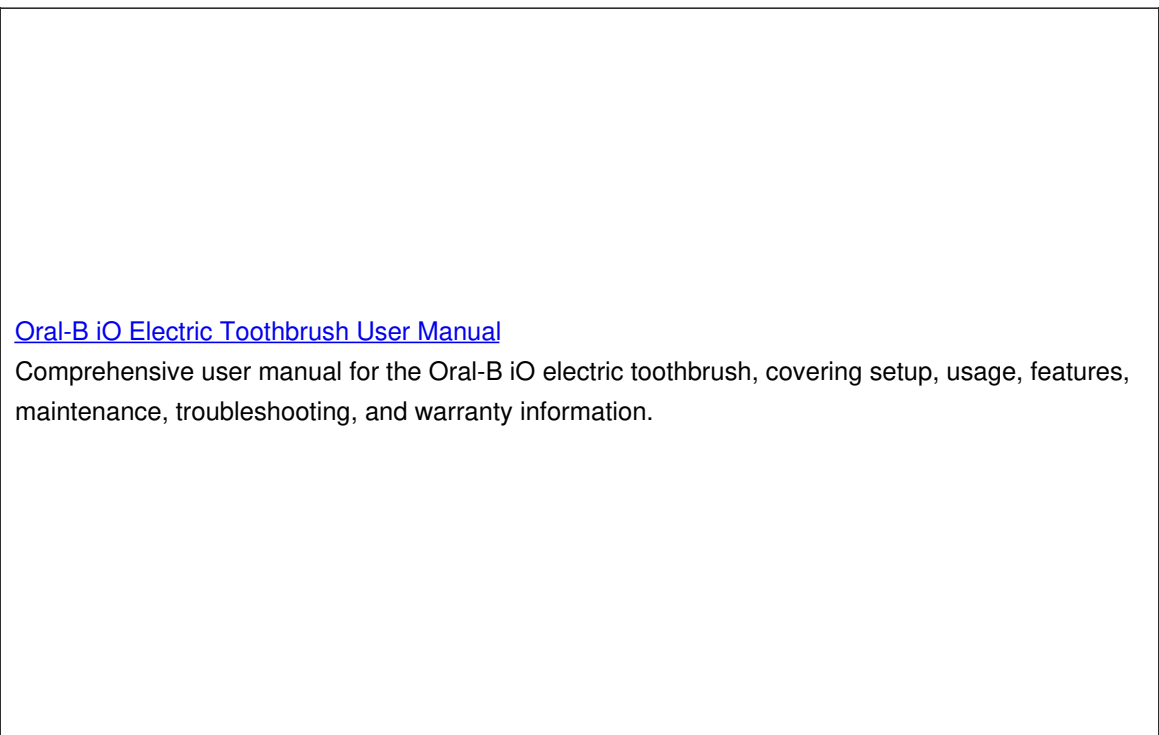

Comprehensive user manual for the Oral-B iO electric toothbrush, covering setup, usage, features, maintenance, troubleshooting, and warranty information.

Comprehensive user manual for the Oral-B iO electric toothbrush, covering setup, usage, features, maintenance, troubleshooting, and warranty information.



- Customer Service**
- DE 00 800 27 28 64 63
  - GB 08 44 - 88 40 10
  - FR 0 800 731 1792
  - ES 1 800 003 448
  - IT 0 800 844 803 (service at Apple points)
  - SE 0 800 14 592
  - RU 800 21172
  - BE 900 814 208
  - PT 800 20 00 33
  - FI 800 124 600
  - NL 0800 445 53 88
  - DK 70 15 00 13
  - NO 22 63 00 93
  - SE 020 - 21 33 21
  - FI 020 377 877

**Internet:**  
www.oralb.com  
www.braun.com  
www.service.braun.com

Braun GmbH  
Friedrichstraße 145  
54756 Krefeld, Germany

DELUXEFLUXX FÜR KINDER  
1070046-01-15

## [Oral-B Kids Electric Toothbrush - User Manual](#)

Comprehensive user manual for the Oral-B Kids electric toothbrush, powered by Braun. Includes safety instructions, charging, usage, modes, cleaning, and warranty information for children's oral hygiene.



© 2015 P&G



1-800-775-1752  
1-800-200-6466  
1-800-200-20-20

InterWeb  
www.oralb.com  
www.oralb.com  
www.oralb.com

Brush Doctor  
P.O. Box 1000  
20150 Northridge, CA 91329

Charger Type 3757  
Handle Type 3767

© 2015 P&G

### [Oral-B SMART/JUNIOR Electric Toothbrush User Manual and Guide](#)

Comprehensive user manual for the Oral-B SMART/JUNIOR electric toothbrush, covering setup, operation, features, safety, troubleshooting, and warranty information. Learn how to use your smart toothbrush effectively.



© 2015 P&G



1-800-775-1752  
1-800-200-6466  
1-800-200-20-20  
1-800-200-6466  
1-800-200-6466

InterWeb  
www.oralb.com  
www.oralb.com  
www.oralb.com

Brush Doctor  
P.O. Box 1000  
20150 Northridge, CA 91329

Charger Type 3757  
Handle Type 3767

© 2015 P&G

### [Oral-B SMART/Teen Electric Toothbrush User Manual and Guide](#)

Comprehensive user manual for the Oral-B SMART/Teen electric toothbrush (Charger Type 3757, Handle Type 3767). Learn about features, usage, charging, cleaning, troubleshooting, and warranty.



### [Oral-B GENIUS 8000-9000 Electric Toothbrush User Manual](#)

User manual for the Oral-B GENIUS 8000-9000 electric toothbrush, covering operation, features, app connectivity, maintenance, and troubleshooting.