

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [Oral-B](#) /
- › [Oral-B Smart 1500 Electric Power Rechargeable Toothbrush User Manual](#)

Oral-B Smart 1500

Oral-B Smart 1500 Electric Power Rechargeable Toothbrush User Manual

Model: Smart 1500 | Brand: Oral-B

INTRODUCTION

Welcome to your Oral-B Smart 1500 Electric Power Rechargeable Toothbrush. This manual provides essential information for the proper use, maintenance, and care of your new toothbrush. The Oral-B Smart 1500 is designed to remove significantly more plaque than a manual toothbrush, promoting healthier gums. Its advanced 3D cleaning action, combined with a visible pressure sensor and quadrant timer, ensures an effective and safe brushing experience.

IMPORTANT SAFETY INFORMATION

- Keep out of reach of children.
- This product contains small parts that may pose a choking hazard if swallowed.
- Always supervise children under 3 years of age while brushing.
- Replace toothbrush head every 3 months or sooner if bristles become frayed.
- Do not use for purposes other than oral hygiene.
- Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

PACKAGE CONTENTS

Your Oral-B Smart 1500 package typically includes:

- 1 Oral-B Smart 1500 Rechargeable Toothbrush Handle
- 1 Oral-B CrossAction Brush Head
- 1 Charging Station

SETUP

Initial Charging

Before first use, charge your toothbrush for at least 12 hours.

1. Plug the charging station into an electrical outlet.
2. Place the toothbrush handle onto the charging station. The battery indicator light will show charging status.



Image: Oral-B Smart 1500 toothbrush handle placed on its charging station for initial charging.

Attaching the Brush Head

1. Align the brush head with the shaft on the toothbrush handle.
2. Push the brush head down firmly until it clicks into place.

OPERATING INSTRUCTIONS

Brushing Technique

Dentists recommend brushing for 2 minutes, twice a day. The Oral-B Smart 1500 features a built-in timer to assist you.


1. Wet the brush head and apply a small amount of toothpaste.
2. Place the brush head against your teeth at a slight angle towards the gumline.
3. Press the power button to turn on the toothbrush.
4. Guide the brush head slowly from tooth to tooth, allowing it to clean each surface for a few seconds. Do not scrub.
5. The 360-degree visible pressure sensor will light up red if you are brushing too hard. Reduce pressure immediately to protect your gums.
6. The in-handle quadrant timer will pulse every 30 seconds to indicate when to move to a different quadrant of your mouth. A longer pulse indicates the completion of the 2-minute brushing cycle.
7. After brushing, press the power button to turn off the toothbrush.



Image: Close-up view of the brush head demonstrating 360-degree cupping technology surrounding a tooth.



Image: The Oral-B Smart 1500 toothbrush highlighting its 2-minute quadrant timer feature.

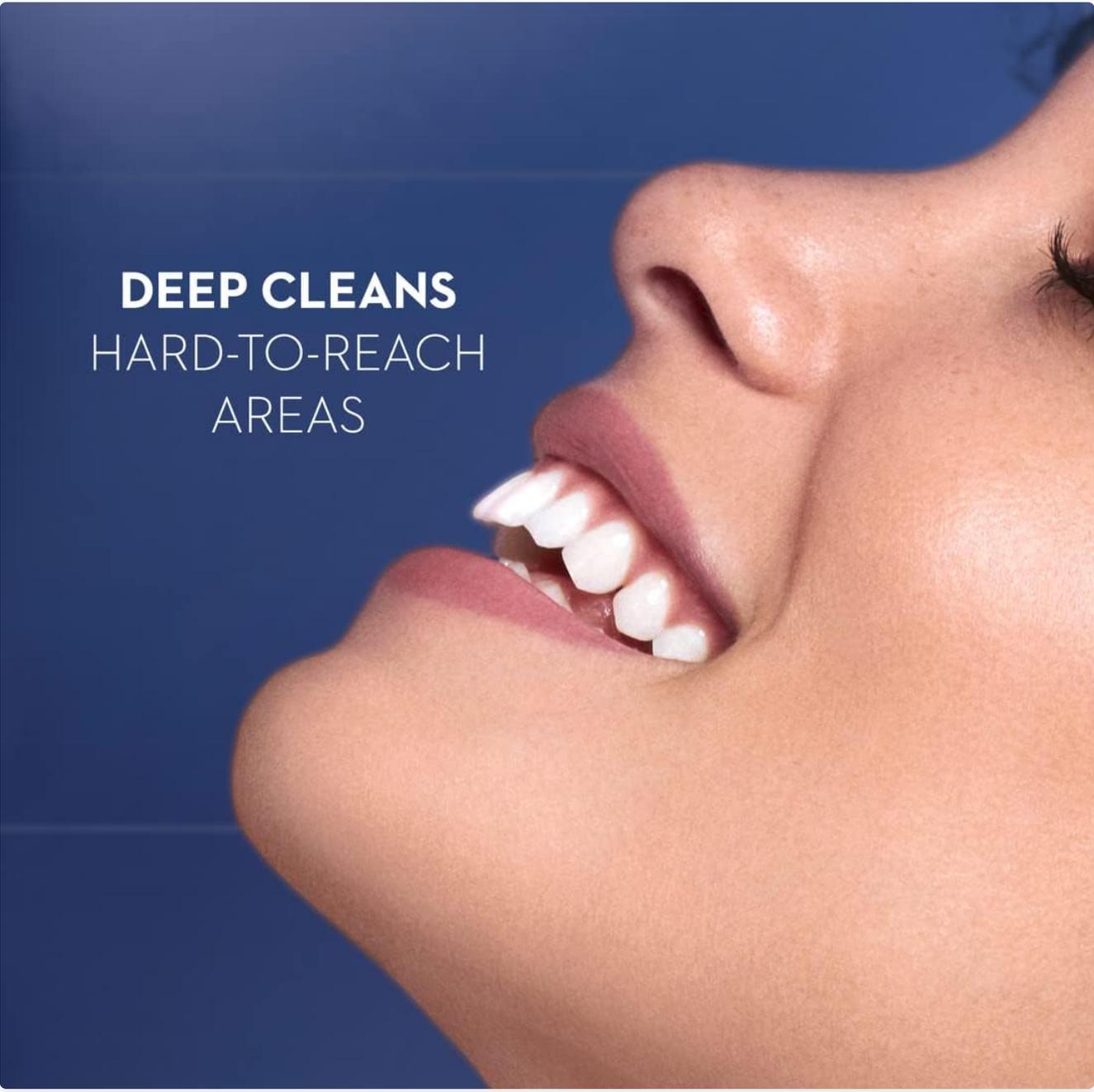

DEEP CLEANS HARD-TO-REACH AREAS

Image: A person demonstrating the use of the Oral-B toothbrush for deep cleaning, reaching hard-to-reach areas.

Cleaning Modes

The Oral-B Smart 1500 offers multiple cleaning modes. Press the power button repeatedly to cycle through the available modes:

- **Daily Clean:** Standard mode for daily plaque removal.
- **Sensitive:** Gentle speed for sensitive areas.
- **Whitening:** Alternating speeds to polish away surface stains.

MAINTENANCE

Cleaning the Toothbrush

1. After each use, rinse the brush head thoroughly under running water with the handle switched on.
2. Remove the brush head from the handle and clean both parts separately under running water.

3. Wipe the handle with a damp cloth. Do not immerse the handle in water.
4. Clean the charging unit regularly with a damp cloth.

Replacing the Brush Head

Dentists recommend replacing your brush head every 3 months, or sooner if the bristles are splayed or worn. Indicator bristles will fade halfway to signal replacement is needed.

1. Pull the old brush head straight off the handle.
2. Push a new Oral-B replacement brush head onto the handle until it clicks.



Image: A close-up view of a clean Oral-B brush head, showing its distinct bristle pattern.

TROUBLESHOOTING

Toothbrush Not Turning On

- Ensure the toothbrush is fully charged. Place it on the charger for at least 12 hours.

- Check if the charging station is properly plugged into a working outlet.

Toothbrush Not Charging

- Verify the charging station is connected to a live power source.
- Ensure the toothbrush handle is correctly seated on the charging station.
- If the issue persists, the charger or battery may need replacement.

Pressure Sensor Lights Up Frequently

- You may be brushing too hard. Apply less pressure during brushing. The toothbrush is designed to do the work for you.

Reduced Performance

- The brush head may be worn. Replace it with a new Oral-B brush head.
- The battery may be low. Recharge the toothbrush fully.

SPECIFICATIONS

Feature	Detail
Brand	Oral-B
Model Number	069055129029
Power Source	Battery Powered, Corded Electric
Special Features	Pressure Sensor, Rechargeable
Product Dimensions	7.01 x 3.94 x 9.96 inches
Item Weight	9.46 ounces
Age Range	Adult
Manufacturer	Procter & Gamble

WARRANTY AND SUPPORT

For warranty information and customer support, please refer to the documentation included with your product or visit the official Oral-B website. Keep your proof of purchase for warranty claims.



