

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [Hairmax](#) /

› [Hairmax LaserBand 82 ComfortFlex Hair Growth Laser Light Device Instruction Manual](#)

Hairmax 1030016

Hairmax LaserBand 82 ComfortFlex Instruction Manual

Hair Growth Laser Light Device

1. PRODUCT OVERVIEW

The Hairmax LaserBand 82 ComfortFlex is a home-use medical device designed to treat hair loss and promote hair growth. It utilizes therapeutic laser light energy to stimulate hair follicles, restore the natural hair growth cycle, and increase hair density and fullness. This device is clinically proven and FDA cleared for men and women.



Figure 1: Hairmax LaserBand 82 ComfortFlex in charging cradle.

Key Features:

- **Clinically Validated:** Technology supported by 7 clinical studies, 14 international medical device licenses, and 6 published articles.
- **Effective Hair Regrowth:** Proven to regrow hair with an average increase of 129 additional new hairs per square inch after six months.
- **Medical-Grade Lasers:** Features 82 medical-grade lasers (no LEDs) providing 1230mW of therapeutic laser light.
- **Patented Hair Parting Teeth:** Designed to part hair for optimal laser light delivery directly to hair follicles.
- **Fast Treatment:** Requires only 90 seconds, three times a week.
- **ComfortFlex Design:** Flexible band and soft-touch teeth ensure comfortable use.

2. SAFETY INFORMATION

Please read and understand all safety warnings before using the Hairmax LaserBand 82 ComfortFlex.

- **External Use Only:** Do not ingest or allow the product to come into contact with eyes, mouth, or other mucous membranes.
- **Allergic Reactions:** Discontinue use if you experience scalp irritation, redness, or rash. Consult a healthcare professional if symptoms persist.
- **Pregnancy and Breastfeeding:** Not recommended for use by pregnant or nursing women.
- **Children:** Keep out of reach of children. Not for use by individuals under 18 years old.
- **Scalp Conditions:** Do not use on damaged, sunburned, or irritated scalp. Consult a healthcare professional if you are using other medications, especially those applied to the scalp.
- **General Precaution:** Always follow the instructions for use carefully.

3. PRODUCT COMPONENTS

The Hairmax LaserBand 82 ComfortFlex package includes the following items:

- Hairmax LaserBand 82 ComfortFlex Device
- Charging Base/Cradle
- Power Supply & Cord
- Quick Start Guide
- User Manual

LASER DEVICES

LASERBAND 82 COMFORTFLEX

- Delivers nourishing laser light to stimulate hair growth & treat hair loss
- Increases density & fullness
- Reverses the thinning process
- Revitalizes damaged hair
- Drug-free, no adverse side effects
- Fastest treatment time in the market



**Hairmax laser devices are indicated to promote hair growth in males with Androgenetic alopecia who have Norwood-Hamilton Classifications of IIa - V, and in females with Androgenetic alopecia who have Ludwig (Savin) Classifications I - II or frontal patterns of hair loss and who both have Fitzpatrick Skin Types I - IV.*



Figure 2: Hairmax LaserBand 82 ComfortFlex and accessories.

4. SETUP

Before first use, ensure the device is fully charged.

1. Unpack the Hairmax LaserBand 82 ComfortFlex device and its charging base.
2. Connect the power supply cord to the charging base.
3. Place the LaserBand device into the charging base. The device will indicate charging status.
4. Allow the device to charge completely before its initial use.


LASERBAND 82 COMFORTFLEX

PATENTED HAIR PARTING TEETH

PARTS HAIR FOR OPTIMAL LASER LIGHT DELIVERY & FOLLICLE STIMULATION

FAST TREATMENT TIME

90 SECONDS, 3 TIMES A WEEK

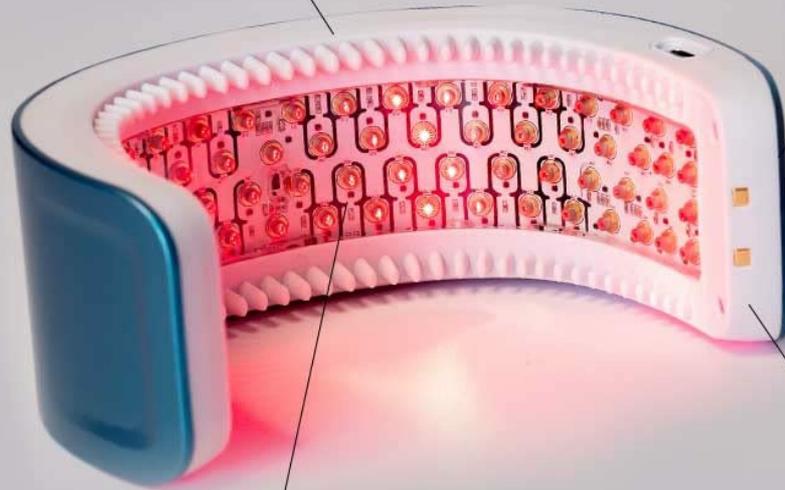

82 MEDICAL GRADE LASERS

STIMULATE HAIR FOLLICLES & EXTEND HAIR GROWTH CYCLES

FLEXIBLE BAND DESIGN

FOR COMFORTABLE TREATMENT & MOVEMENT AT HOME OR TRAVELS

Figure 3: Charging the Hairmax LaserBand 82 ComfortFlex.

5. OPERATING INSTRUCTIONS

Follow these steps for an effective hair growth treatment:

1. Ensure your hair is clean and dry.
2. To start your treatment, press the power button once.
3. Position the device on your scalp at the hairline (Position 1), ensuring the patented hair parting teeth make contact with your hair.
4. Leave it in place until it vibrates (approximately 30 seconds).
5. When you feel the vibration, slide your Hairmax backwards about 3 inches to Position 2 (the middle portion of your scalp).
6. Repeat the 30-second procedure for Position 3 (the back of your scalp).
7. The total treatment time is 90 seconds. The device will automatically turn off after completing the full treatment.
8. Use the device three times a week for best results.

Treatment instructions.



01

To start your treatment, press the power button once. Place your Hairmax® at your hairline and leave it in place until it vibrates. This is position 1.



02

When you feel the vibration, slide your Hairmax® backwards about 3 inches to position 2, the middle portion of your scalp.



03

Use light pressure so the teeth can separate your hair. Leave it in position 2 until it vibrates again. Repeat the 30 second procedure on position 3.

Figure 4: Hairmax LaserBand 82 treatment instructions for scalp coverage.

Official Product Videos:

Your browser does not support the video tag.

Video 1: Hair Regrowth by Hairmax Laser Band 82. This video demonstrates the proper technique for using the Hairmax LaserBand 82 for optimal hair regrowth.

Your browser does not support the video tag.

Video 2: HairMax LaserBand 82 ComfortFlex for Hair Growth. This video highlights the comfort and flexibility of the Hairmax LaserBand 82 during use.

6. UNDERSTANDING THE TECHNOLOGY

The Hairmax LaserBand 82 utilizes Low-Level Laser Therapy (LLLT) to stimulate hair follicles. Unlike LED devices, Hairmax uses medical-grade lasers that deliver a precise, collimated light directly to the hair follicle root, which is crucial for stimulating hair growth.



Figure 5: Hairmax Medical Grade Laser vs. LED Devices.

The patented hair parting teeth ensure that the laser light reaches the scalp effectively by moving hair aside, maximizing the energy delivered to the follicles.

Exclusive hair parting teeth *maximize effectiveness.*

On Hairmax® devices, the teeth part hair for optimal laser light delivery and follicle stimulation.

On devices without teeth, hair blocks laser light from reaching your scalp.

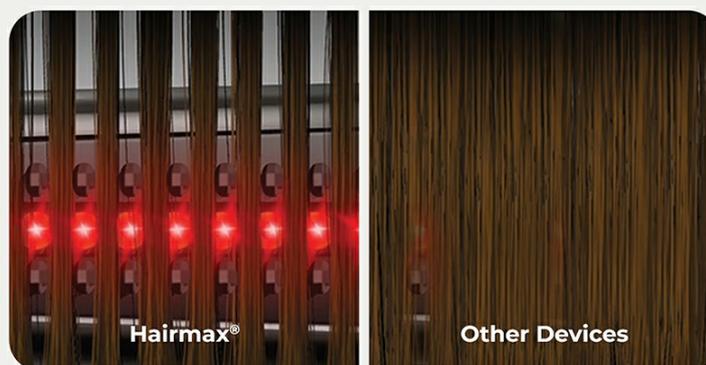


Figure 6: Exclusive hair parting teeth maximize effectiveness.

7. MAINTENANCE

To ensure the longevity and optimal performance of your Hairmax LaserBand 82 ComfortFlex, follow these maintenance guidelines:

- **Cleaning:** Gently wipe the device with a soft, damp cloth after each use. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the device in its charging cradle or a safe, dry place away from direct sunlight and extreme temperatures.
- **Battery Care:** Fully charge the device regularly, even if not in frequent use, to maintain battery health.
- **Avoid Damage:** Handle the device with care. Avoid dropping it or exposing it to harsh impacts.

8. TROUBLESHOOTING

If you encounter any issues with your Hairmax LaserBand 82 ComfortFlex, please refer to the following common troubleshooting tips:

- **Device Not Turning On:** Ensure the device is fully charged. Press and hold the power button for a few seconds.
- **Lasers Not Emitting:** The lasers are designed to activate only when the device is in proper contact with the scalp. Ensure the device is correctly positioned.
- **Charging Issues:** Verify that the power cord is securely connected to both the charging base and a working power outlet. Ensure the device is correctly seated in the charging cradle.
- **Irregular Vibrations:** If the device vibrates inconsistently, ensure it is held firmly against the scalp and moved smoothly between positions.
- **No Results After Use:** Hair growth is a gradual process. Consistent use over 3-6 months is typically required to observe visible results. Ensure you are following the recommended treatment schedule.

If the problem persists, please contact Hairmax Customer Support for further assistance.

9. SPECIFICATIONS

Feature	Detail
Product Dimensions	11.5 x 9.5 x 4.75 inches

Item Model Number	1030016
Batteries	2 Lithium Polymer batteries required (included)
Brand	Hairmax
Item Weight	3.2 Pounds
Item Form	Headband
Age Range (Description)	Adult

10. WARRANTY & SUPPORT

For warranty information, product registration, or technical support, please visit the official Hairmax website or contact their customer service department. Keep your proof of purchase for warranty claims.

Hairmax Official Website: [Visit the Hairmax Store on Amazon](#)