

SKG K6E

SKG K6E Neck Massager User Manual

1. INTRODUCTION

This manual provides detailed instructions for the safe and effective use of your SKG K6E Neck Massager. Please read this manual thoroughly before operating the device and retain it for future reference. The SKG K6E is designed to provide targeted neck massage with integrated heating, aiming to relieve muscle tension and promote relaxation.

2. PRODUCT COMPONENTS AND CONTROLS

Familiarize yourself with the various parts and controls of your SKG K6E Neck Massager as illustrated below.

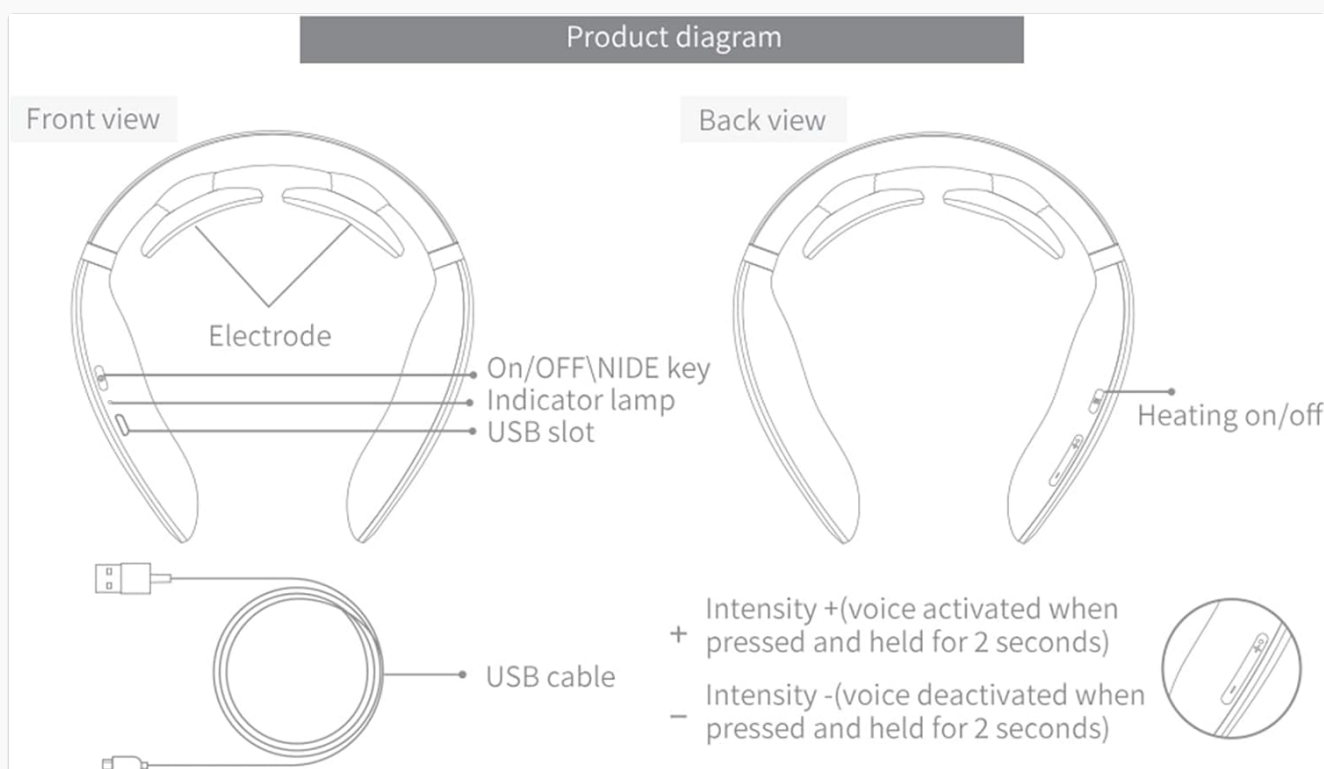


Image: Front and back view of the SKG K6E Neck Massager, detailing the electrodes, On/Off/Mode key, indicator lamp, USB slot, heating on/off button, and intensity adjustment buttons.

- **Electrodes:** Gold-plated conducting strips for massage impulses.
 - **On/Off/Mode Key:** Powers the device on/off and cycles through massage modes.
 - **Indicator Lamp:** Displays device status (e.g., charging, operating mode).
 - **USB Slot:** For charging the device using the provided USB cable.
 - **Heating On/Off Button:** Activates and adjusts the heating function.
 - **Intensity + Button:** Increases massage intensity. Press and hold for 2 seconds to activate/deactivate voice prompts.
 - **Intensity - Button:** Decreases massage intensity. Press and hold for 2 seconds to activate/deactivate voice prompts.
 - **USB Cable:** For charging the massager.
-

3. SETUP

3.1 Charging the Device

Before first use, fully charge the device. Connect the provided USB cable to the USB slot on the massager and to a compatible USB power adapter (not included). The indicator lamp will show charging status. A full charge provides up to 14 15-minute massage sessions.

3.2 Device Placement

Position the massager around your neck with the electrodes making direct contact with your skin. Ensure a comfortable and secure fit. The ergonomic U-shaped design is suitable for various neck sizes. For optimal conductivity, you may moisten your neck slightly with water or a damp cloth before use.

4. OPERATING INSTRUCTIONS

4.1 Power On/Off

Press and hold the **On/Off/Mode key** for approximately 2 seconds to power on or power off the device. A voice prompt will confirm the action.

4.2 Select Massage Mode

With the device powered on, short press the **On/Off/Mode key** to cycle through the 4 available massage modes: Daily, Leisure, Exercise, and Therapy. The voice prompt will announce the selected mode.

4.3 Adjust Intensity

Press the **'+' button** to increase massage intensity (15 levels available). Press the **'-' button** to decrease intensity. The voice prompt will confirm the intensity level. To activate or deactivate voice prompts, press and hold either the **'+' or '-' button** for 2 seconds.

4.4 Activate Heating Function

Press the **'Heating on/off' button** located on the back of the device to activate the heating function. Short press to cycle through the 3 heating levels: Low (96.8°F-102.2°F), Middle (102.2°F-105.8°F), and High (105.8°F-111.2°F). Press again to turn off heating.

3 LEVEL HEATING COMPRESS

×

Further reduce cervical pressure



Image: Illustration of the SKG K6E Neck Massager's 3-level heating compress, showing heat distribution from the electrodes.

4.5 Automatic Shut-off

For safety and optimal use, the device is programmed to automatically shut off after 15 minutes of continuous operation.

5. MAINTENANCE

5.1 Cleaning

Wipe the electrodes and the silicone surface with a soft, damp cloth after each use. Do not use abrasive cleaners or immerse the device in water.

5.2 Storage

Store the massager in a cool, dry place, away from direct sunlight and extreme temperatures. Keep out of reach of children.

6. TROUBLESHOOTING

- **Device does not turn on:** Ensure the device is fully charged. Press and hold the On/Off/Mode key for at least 2 seconds.
- **No massage sensation:** Ensure the electrodes are in direct contact with your skin. You may moisten your neck slightly with water or a damp cloth to improve conductivity. Increase the intensity level using the '+' button.
- **Heating function not working:** Ensure the heating function is activated by pressing the 'Heating on/off' button. Check if the device is adequately charged.
- **Voice prompt issues:** If voice prompts are not desired, press and hold the '+' or '-' button for 2 seconds to deactivate them. Repeat to reactivate.

7. SPECIFICATIONS

Feature	Specification
Model	SKG K6E
Dimensions	14.2 x 13.9 x 3.7 cm
Weight	140 g (0.14 Kilograms)
Material	Plastic, Silicone
Power Source	Battery Powered
Battery Capacity	1500 mAh
Usage Time	Up to 14 sessions (15 minutes each) per charge
Massage Modes	4 (Daily, Leisure, Exercise, Therapy)
Intensity Levels	15
Heating Levels	3 (Low: 96.8°F-102.2°F, Mid: 102.2°F-105.8°F, High: 105.8°F-111.2°F)
Automatic Timer	15 minutes

8. SAFETY INFORMATION

Please observe the following safety precautions to ensure safe operation and prevent injury or damage to the device:

- Consult a physician before using this device if you have a pacemaker, are pregnant, have skin conditions, or any other medical conditions.

- Do not use while driving or operating heavy machinery.
- Keep out of reach of children.
- Do not disassemble or modify the device.
- Avoid using the device in high humidity environments or near water.
- Discontinue use immediately if you experience any discomfort or adverse reactions.




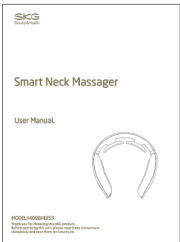

9. WARRANTY AND SUPPORT


For warranty information and customer support, please refer to the official SKG website or contact your point of purchase. You can also visit the SKG store online for additional resources and product information:

[SKG Official Store](#)

© 2023 SKG. All rights reserved.

Related Documents - K6E

	<p>SKG K6E Neck Massager User Manual: Operation, Safety, and Maintenance</p> <p>Comprehensive user manual for the SKG K6E Neck Massager. Learn about its features, how to wear, charge, operate, maintain, and important safety precautions for effective and safe use.</p>
	<p>SKG K6E Neck Massager User Manual</p> <p>Comprehensive user manual for the SKG K6E Neck Massager, detailing features, operation, safety precautions, troubleshooting, and specifications.</p>
	<p>SKG K6-EN Neck Massager User Manual: Operation Guide and Safety Precautions</p> <p>This user manual for the SKG K6-EN Neck Massager provides detailed instructions on operation, maintenance, safety guidelines, and troubleshooting. Learn how to use your SKG neck massager effectively and safely.</p>
	<p>SKG Smart Neck Massager User Manual (Models 4098/4353)</p> <p>User manual for the SKG Smart Neck Massager (Models 4098 and 4353), providing instructions on operation, precautions, maintenance, and troubleshooting. Features include heating and adjustable intensity for neck pain relief.</p>
	<p>SKG T5-1 Shoulder and Neck Massager User's Guide</p> <p>User's guide for the SKG T5-1 Shoulder and Neck Massager, detailing usage instructions, technical parameters, troubleshooting, and safety precautions for effective pain relief.</p>

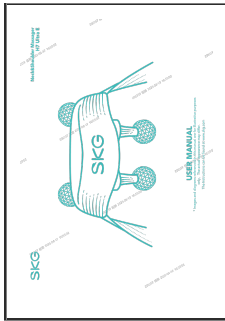

[SKG H7 Ultra E Neck & Shoulder Massager User Manual](#)

Comprehensive user manual for the SKG H7 Ultra E Neck & Shoulder Massager, covering features, operation, safety precautions, troubleshooting, and technical specifications.