

## DYRABREST WL-JMLB1218

# DYRABREST Pressotherapy Air Pressure Machine Slimming Suit User Manual

MODEL: WL-JMLB1218

## 1. Product Overview

The DYRABREST Pressotherapy Air Pressure Machine is designed to provide a body slimming and beauty massage experience through air pressure. It features a super silent compressor for quiet operation and adjustable settings for personalized use. This device aims to alleviate various conditions related to blood circulation and nerve function.

- **Super Silent Compressor:** Engineered for quiet operation, minimizing noise during use.
- **Uniform Power Release:** Delivers consistent air pressure for a relaxed and comfortable experience.
- **Time Setting Function:** Allows users to set massage duration from 0 to 30 minutes.
- **Adjustable Pressure:** Pressure levels can be customized to individual user needs.
- **Therapeutic Benefits:** May help alleviate symptoms of chronic venous insufficiency, blood circulation disorders, neuralgia, numbness, cold sensations, insomnia, and post-stroke sequelae.



Figure 1: DYRABREST Pressotherapy Air Pressure Machine with full body suit components. The image shows the main control unit and the purple air pressure suit, including arm, waist, and leg components, laid out on a white surface.

# Air Pressure Pressotherapy Suit

Prevent Phlebothrombosis



Figure 2: User wearing the DYRABREST Pressotherapy suit for full body treatment. The image displays a person wearing the purple air pressure suit, covering arms, torso, and legs, connected to the control unit.


## 2. Safety Information

Please read and understand all safety precautions before using the device. Failure to do so may result in injury or device malfunction.

- Do not use if you have severe limb infections that are not effectively controlled.
- Avoid use if you have recent deep venous thrombosis of the lower extremities.
- Do not use on areas with large ulcerative rashes.
- Exercise caution if you have a tendency to bleed.
- Not recommended for pregnant women.
- Avoid use if you have skin diseases such as scalds or reddening.
- Women in the menstrual period should avoid using the device.
- Consult a physician before use if you have any pre-existing medical conditions.

### 3. Setup Instructions

1. **Unpack the Massager:** Carefully remove all components from the packaging. Ensure all parts are present and undamaged.
2. **Connect to Power:** Plug the massager's power adapter into a suitable electrical outlet.
3. **Connect Air Hoses:** Identify the air hoses for the arm, waist, and leg components. Connect each hose securely to the corresponding port on the main control unit. Ensure all connections are firm to prevent air leakage.
4. **Wear the Suit:** Put on the arm, waist, and leg components of the pressotherapy suit. Adjust the straps to ensure a snug but comfortable fit. The suit should not be excessively tight.

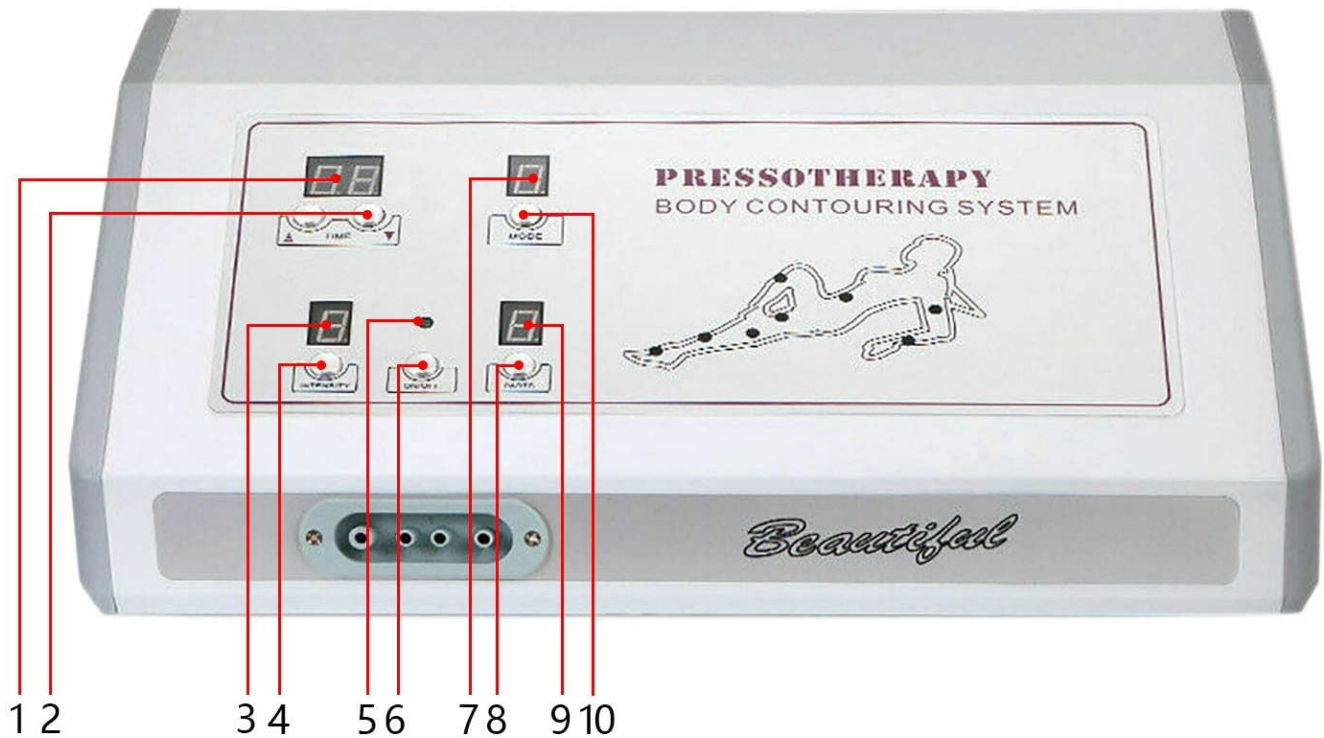

- |                                 |                              |
|---------------------------------|------------------------------|
| 1:Timer                         | 6: ON/OFF Button             |
| 2:Time UP/DOWN                  | 7: MODE of Airpressure       |
| 3: The INTENSITY of Airpressure | 8: Selection Button of PARTS |
| 4:SelectionButtonofIntensity    | 9: Massage PARTS             |
| 5: ON/OFF Indicator Light       | 10: Selection Button of MODE |

Figure 3: Detailed view of the DYRABREST Pressotherapy control panel with numbered components. This image highlights the timer, mode, intensity, ON/OFF, and parts selection buttons, along with the air pressure connection ports.

Video 1: This video demonstrates the connection of the air hoses to the DYRABREST Pressotherapy machine and how to wear the suit components, ensuring proper setup for operation.

### 4. Operating Instructions

1. **Power On:** Press the 'ON/OFF' button on the control unit to turn on the machine. The indicator light will

illuminate.

2. **Set Massage Time:** Use the 'TIME UP/DOWN' buttons to select the desired massage duration (0-30 minutes).
3. **Select Mode:** Press the 'MODE' button to cycle through available air pressure massage modes. Refer to the control panel diagram (Figure 3) for mode indicators.
4. **Adjust Pressure Intensity:** Use the 'INTENSITY' buttons to increase or decrease the air pressure level according to your comfort.
5. **Select Massage Parts:** If applicable, use the 'PARTS' selection button to target specific areas (arms, waist, legs, or feet) for the air pressure massage.
6. **Begin Massage:** The machine will start the selected massage program.
7. **End Session:** The machine will automatically shut off after the set time. Alternatively, press the 'ON/OFF' button to stop the massage manually.
8. **Disconnect:** After use, disconnect the massager from the power source and carefully remove the suit components.

Video 2: This video illustrates the general usage and massage patterns of the DYRABREST Pressotherapy machine, highlighting its powerful circulation and hand massage type modes for lower body treatment.

## 5. Maintenance

- **Cleaning:** Ensure the device is unplugged before cleaning. Wipe the control unit and suit components with a soft, damp cloth. Do not use abrasive cleaners or solvents.
- **Storage:** Store the device and suit in a cool, dry place away from direct sunlight and extreme temperatures. Keep out of reach of children.
- **Inspection:** Regularly check the air hoses and connections for any signs of wear or damage. Replace any damaged parts immediately.

## 6. Troubleshooting

If you encounter any issues with your DYRABREST Pressotherapy machine, please refer to the following common problems and solutions:

- **No Power:** Ensure the power adapter is securely connected to both the wall outlet and the control unit. Check if the power outlet is functional.
- **No Air Pressure:** Verify that all air hoses are correctly and firmly inserted into their respective ports on the control unit and suit components. Ensure there are no kinks or blockages in the hoses.
- **Uneven Pressure:** Check that the suit components are worn correctly and adjusted snugly. Ensure no air leaks are present from the connections or the suit itself.
- **Machine Overheating:** If the control unit feels excessively warm, turn off the device and allow it to cool down for at least 10 minutes before resuming use. Prolonged use at maximum intensity may cause temporary heating.

## 7. Specifications

Feature	Detail
Product Dimensions	15.75 x 10.63 x 5.12 inches

Item Model Number	WL-JMLB1218
Item Weight	26.5 Pounds
Brand	DYRABREST
Power Source	Corded Electric
Material	Acrylonitrile Butadiene Styrene, Metal
Manufacturer	Shibel store



Figure 4: Dimensions of the DYRABREST Pressotherapy control unit and suit components. This image provides measurements for the control unit and various parts of the purple air pressure suit.

## 8. Warranty and Support

---

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. For any product-related inquiries, technical support, or warranty claims, please contact the manufacturer, Shibel store, directly. Refer to your purchase documentation for specific warranty terms and contact information.