

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

› [COMFIER](#) /

› [COMFIER Hand Massager with Heat, Cordless Hand Massage Tool with Compression & Heating, 3 Modes & 3 Levels Pressure Point Hand Massager for Arthritis and Carpal Tunnel from Finger to Wrist](#)

COMFIER 4805

COMFIER Hand Massager with Heat - Instruction Manual

Model: 4805 | Brand: COMFIER

INTRODUCTION

The COMFIER Hand Massager with Heat is designed to provide ultimate relaxation and relief for your hands, from fingers to wrist. Utilizing air compression massage and optional heat therapy, it targets pressure points to soothe discomfort associated with arthritis, carpal tunnel, finger numbness, and general hand fatigue. This manual provides essential information for safe and effective operation of your device.



The COMFIER Hand Massager offers comprehensive palm and finger relaxation with air massage, adjustable intensity, multiple modes, and soothing heat.

SAFETY INFORMATION

Please read all safety warnings carefully before using this product.

- Avoid using the massager near water or in damp environments.

- Ensure the device is properly charged using the provided adapter and USB cable.
- Consult a healthcare professional before use if you have any medical conditions, especially arthritis, carpal tunnel, or circulatory issues.
- Do not use if the device is damaged or malfunctioning.
- Keep out of reach of children and supervise during use.
- Do not attempt to disassemble or repair the device yourself. Contact customer support for assistance.

PACKAGE CONTENTS

Verify that all items are included in your package:

- COMFIER Hand Massager Unit
- Handheld Finger Massager
- User Manual
- Decorative Stickers
- USB Cable
- Power Adapter



The package includes the hand massager unit, a handheld finger massager, user manual, decorative stickers, USB cable, and power adapter.

SETUP

Unboxing

Carefully remove all components from the packaging. Inspect the hand massager and accessories for any signs of damage. Retain the packaging for future storage or transport.

Charging the Device

The COMFIER Hand Massager is cordless and rechargeable. Before first use, fully charge the device. Connect the provided USB cable to the massager and plug the adapter into a standard wall outlet. The power indicator light will show charging status. The device cannot be operated while charging.



The massager can be charged using the provided USB cable and adapter, or via a laptop or power bank. A full charge provides 80-90 minutes of use.

OPERATING INSTRUCTIONS

Familiarize yourself with the control panel for optimal use.

HANDHELD FINGER MASSAGER

Increases circulation, reduces stiffness and soothes arthritis and joint pain in your hands



The control panel features intuitive buttons for Power, Timer, Mode, Heating, and Intensity.

1. **Power On/Off:** Press the **Power** button to turn the device on or off.
2. **Insert Hand:** Place your hand comfortably into the massager, ensuring your palm and fingers are positioned correctly for the air compression.
3. **Select Mode:** Press the **Mode** button to cycle through the three massage modes:
 - **Refresh:** Provides a general compression massage for overall hand revitalization.
 - **Health Care:** Focuses on pressure points and includes wrist massage for targeted relief.
 - **Relax:** Offers a gentle and soothing massage for deep relaxation.

COMPRESSION AND SOOTHING HEAT



Compression




3 Intensity



3 Modes



Heating

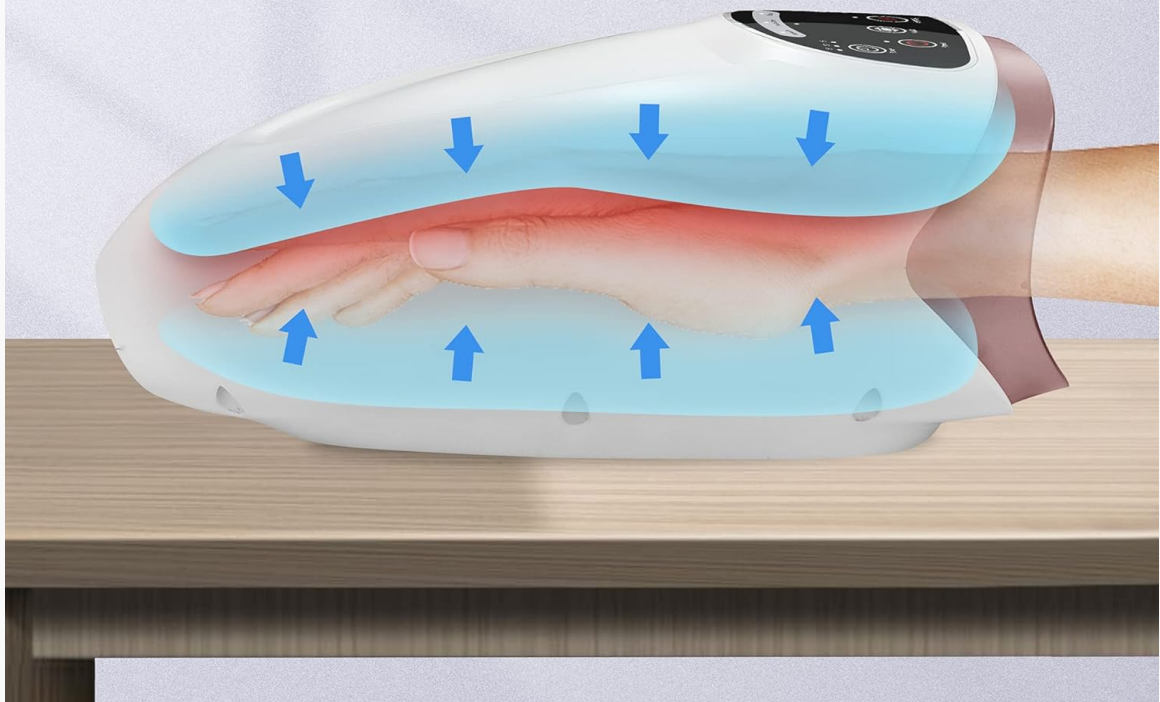

Experience comprehensive compression and soothing heat with adjustable intensity and multiple massage modes.

4. **Adjust Intensity:** Press the **Intensity** button to choose from three levels: Low (L), Medium (M), or High (H). Start with Low and adjust as desired.
5. **Activate Heating:** Press the **Heating** button to turn on the heat function. The heat provides gentle warmth to promote blood circulation and ease discomfort.
6. **Set Timer:** Press the **Timer** button to select a massage duration of 5, 10, or 15 minutes. The device will automatically shut off after the selected time.

Using the Handheld Finger Massager

The included handheld finger massager can be used independently to target specific fingers or areas of the hand.

- To massage fingers, place a finger between the rollers and gently move the massager up and down the finger.
- The rollerball on the opposite end can be used for targeted massage on palms or wrists.

WIRELESS AND RECHARGEABLE

It can work 80-90 minutes on one full charge



The handheld finger massager increases circulation, reduces stiffness, and soothes arthritis and joint pain.

MAINTENANCE

To ensure the longevity and proper functioning of your COMFIER Hand Massager, follow these maintenance guidelines:

- **Cleaning:** Wipe the exterior of the massager with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the massager in a cool, dry place away from direct sunlight and extreme temperatures.
- **Care:** Avoid placing heavy objects on the massager. Do not attempt to open the device casing.

TROUBLESHOOTING

If you encounter any issues with your hand massager, refer to the table below for common problems and solutions.

Trouble	Reason/Solution
---------	-----------------

Trouble	Reason/Solution
Cannot start the massager	Battery Low / Please charge it before using.
The massager stopped abnormally	Battery runs out of power / Charge the battery or Self-timed 15 mins is out, please restart the product.
Massager slight noise	It is normal working noise of the internal mechanism.

Troubleshooting	
Trouble	Reason/Solution
Can not start the massager	Battery Low/ Please charge it before using
The massager stoped abnormally	Battery runs out of power/charge the battery or Self-timed 15 mins is out, restart the product
Massager slight noise	It is normal working noise of the internal mechanism

Refer to this troubleshooting guide for quick solutions to common operational issues.

SPECIFICATIONS

Feature	Detail
Product Dimensions	7.5 x 5.5 x 11.8 inches
Item Model Number	4805
Item Weight	2.29 Pounds (1.04 Kilograms)
Brand	COMFIER
Power Source	Battery Powered
Special Feature	Rechargeable
Use for	Whole Body (specifically hands)















Key dimensions of the COMFIER Hand Massager.

WARRANTY & SUPPORT

For detailed warranty information and customer support, please refer to the official user manual included in your package or visit the COMFIER website. Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

© 2023 COMFIER. All rights reserved.

Related Documents - 4805

<div><div>CONFIER</div><div>Air Compression Leg Massager</div></div> <div><div>CE</div><div>User Manual</div><div>FE-7204B</div></div>	<p>COMFIER FE-7204B Air Compression Leg Massager User Manual</p> <p>Comprehensive user manual for the COMFIER FE-7204B Air Compression Leg Massager, detailing safety instructions, operating procedures, features, technical specifications, and troubleshooting for optimal use and care.</p>
<div><div>CONFIER</div><div>Hand Massager</div></div> <div><div>User Manual</div><div>CF-4805B</div></div>	<p>COMFIER CF-4805B Hand Massager User Manual - Features, Operation, and Warranty</p> <p>Comprehensive user manual for the COMFIER CF-4805B Hand Massager. Learn about its features, technical specifications, operating instructions, charging, maintenance, troubleshooting, and warranty information.</p>
<div><div>CONFIER</div><div>Cordless Neck Massager</div></div> <div><div>User Manual</div><div>CF-6026</div></div>	<p>Comfier Cordless Neck Massager CF-6026 User Manual</p> <p>Comprehensive user manual for the Comfier Cordless Neck Massager (Model CF-6026), detailing safety instructions, features, technical specifications, usage guidelines, troubleshooting tips, care and maintenance procedures, and warranty information.</p>
<div><div>1. Sicherheitshinweise</div><div>WICHTIGE HINWEISE ! UNBEDINGT AUFBEWAHREN ! <small>Lesen Sie die Gebrauchsanweisung, insbesondere die Sicherheitsanweisungen, sorgfältig durch, bevor Sie das Gerät benutzen und beachten Sie die Sicherheitsanweisungen, die für dieses Modell gelten. Wenn Sie die Geräte in Dritte weitergeben, geben Sie unbedingt diese Sicherheitsanweisungen weiter.</small></div><div>Zweckbestimmung <small>Dieses Massagegerät gehört zu diesem Gerät. Es enthält wichtige Informationen zur Identifizierung und Handhabung. Lesen Sie diese Sicherheitsanweisungen vollständig. Die Nichtbeachtung dieser Angaben kann zu schweren Verletzungen oder Schäden am Gerät führen.</small></div><div>WICHTIG! <small>Diese Massagegeräte müssen eingereguliert werden, um mögliche Verletzungen des Benutzers zu vermeiden.</small></div><div>ACHTUNG! <small>Diese Geräte müssen sorgfältig eingesetzt werden, um mögliche Beschädigungen am Gerät zu vermeiden.</small></div><div>WICHTIG! <small>Diese Hinweise geben Ihnen wichtige Informationen aufzufassen oder zum Lesen.</small></div><div>Schaltknoten 3</div><div>LOT <small>Lot Nummer</small></div><div>Hersteller</div><div>3</div></div>	<p>COMFIER FE-7204B Compression Leg Massager User Manual and Safety Guide</p> <p>Comprehensive user manual and safety guide for the COMFIER FE-7204B Compression Leg Massager, covering features, operation, maintenance, troubleshooting, and warranty information.</p>
<div><div>CONFIER</div><div>Shiatsu Neck and Back Massager with Heat</div></div> <div><div>User Manual</div><div>CF-2303-KR</div></div>	<p>COMFIER CF-2303-KR Shiatsu Neck and Back Massager with Heat User Manual</p> <p>Comprehensive user manual for the COMFIER CF-2303-KR Shiatsu Neck and Back Massager with Heat, covering safety instructions, features, technical data, setup, operation, controller functions, care, cleaning, troubleshooting, and warranty information.</p>

COMFIER
Hand Massager



User Manual
CF-4805S

[Comfier CF-4805S Hand Massager User Manual](#)

User manual for the Comfier CF-4805S Hand Massager, providing detailed instructions on safety, features, operation, charging, troubleshooting, maintenance, and warranty information.