

[Manuals.plus](#) /

› [Xllent](#) /

› Xllent Hand Massager with Compression & Heating (Model NX-S168) - Instruction Manual

Xllent NX-S168

Xllent Hand Massager with Compression & Heating

Model: NX-S168

Brand: Xllent

1. INTRODUCTION AND OVERVIEW

The Xllent Hand Massager is designed to provide comprehensive hand care through advanced acupressure therapy, compression, and gentle heat. It features four strategically placed airbags that compress and stimulate key pressure points on the wrist, palm, and fingers. This process aims to enhance blood circulation, relieve hand fatigue, and promote deep relaxation. It is particularly beneficial for individuals experiencing discomfort from arthritis, carpal tunnel, or general daily hand stress. The device is built with high-quality, practical components to support a healthier lifestyle.



The Xllent Hand Massager is a white, compact device designed for hand therapy, featuring compression and heating functions. It is shown alongside its retail packaging and an included finger roller.

2. SAFETY INFORMATION

Please read the user manual for detailed instructions and precautions before use. This product is intended for personal relaxation and pain relief. Consult a healthcare professional if you have pre-existing medical conditions, especially related to circulation, nerve damage, or skin sensitivities. Do not use on broken or irritated skin. Keep out of reach of children. The device is equipped with a 15-minute auto shut-off function to ensure safety and prevent overheating.

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

3. SETUP

Before first use, ensure the device is fully charged. The Xllent hand massager is built with a 2200mAh rechargeable battery and provides 60-90 minutes of use per charge via the included USB cable. To charge, connect the provided USB cable to the Type-C interface on the device and a suitable power source.



This image provides a detailed look at the massager's construction, emphasizing its durable ABS shell, comfortable breathable lining, stable anti-skid pads, and convenient Type-C fast charging port. Dimensions are also provided.

To use the hand massager, place your hand comfortably inside the device. Ensure your hand is positioned correctly to allow the airbags to compress the wrist, palm, and fingers effectively. Check the product size before purchasing to ensure a perfect fit.

4. OPERATING INSTRUCTIONS

The device features an intuitive control panel for easy operation:

1. **Power On/Off:** Long press the power button (gray button in the center) for 2-3 seconds to turn the device on. The default setting is usually the first massage mode with a moderate force. To shut down, tap the power button.
2. **Adjust Massage Force (Velocity Key):** If the massage strength is not sufficient, use the Velocity Key (indicated by arrows) to adjust through 5 gears. The force can be cycled from bottom to top.
3. **Activate Heating/Vibration (Function Key):** To turn on the heating and vibration functions, press the Function Key (indicated by wavy lines for heat and wavy arrows for vibration) to select your preferred mode. The device offers 5 heat levels and multiple vibration frequencies.
4. **Auto Shut-off:** The massager is equipped with a 15-minute auto shut-off function for safety.

How to Use

01

Long press for 2-3s
to turn on



After start up, the force is 1 gear, the massage mode is the first, and the function is none.

02



If the massage strength is not large enough, you can adjust the Velocity button to find the right gear.

03



If you need to turn on the heating and vibration function, please press the function key to select.

04



If you need to shut down, please tap the power button.

This visual guide shows the basic operation: 1. Long press to turn on. 2. Adjust velocity. 3. Press function key for heating/vibration. 4. Tap power button to shut down.

BUTTON OPERATION SEQUENCE



Vibration Key — 1 gear

the default breathing state is adopted on the system;
If you need to turn on the vibration function,
press the vibration keys to open



Velocity Key — 5 gears

from bottom to top, can be cycled.
Controls the suction strength in the cup body.



Heating Key — 1 gear

from the bottom up, can be recycled, control
the heat in the cup body.

- When the product is charged,
the vibration key indicates the lights
flickering, and the indicator light is always
on when full.



The control panel features a central power button surrounded by controls for vibration, massage intensity (velocity), and heat. Indicator lights show the selected settings.

5. FEATURES AND MODES

The Xllent Hand Massager offers a customizable experience with various features:

- **Acupressure & Compression:** Utilizes 4 airbags to compress and stimulate key pressure points, enhancing blood circulation and relieving hand fatigue.
- **Heat Therapy:** Features gentle heat compression to keep hands warm, prevent dryness, and ease muscle soreness, numbness, and cold fingers. Offers 5 adjustable heat levels.
- **Adjustable Pressure Modes:** Provides 5 adjustable pressure modes to customize the massage intensity.
- **Vibration:** Multiple vibration frequencies are available to further enhance relaxation.
- **Portable Design:** Lightweight and ergonomic, allowing for use at home, in the office, or on the go.

Air Compression & Warm Temperature



0-5 Gear
Circulatio




Vibration
Key



Heating
Key



15 Mins
Timer

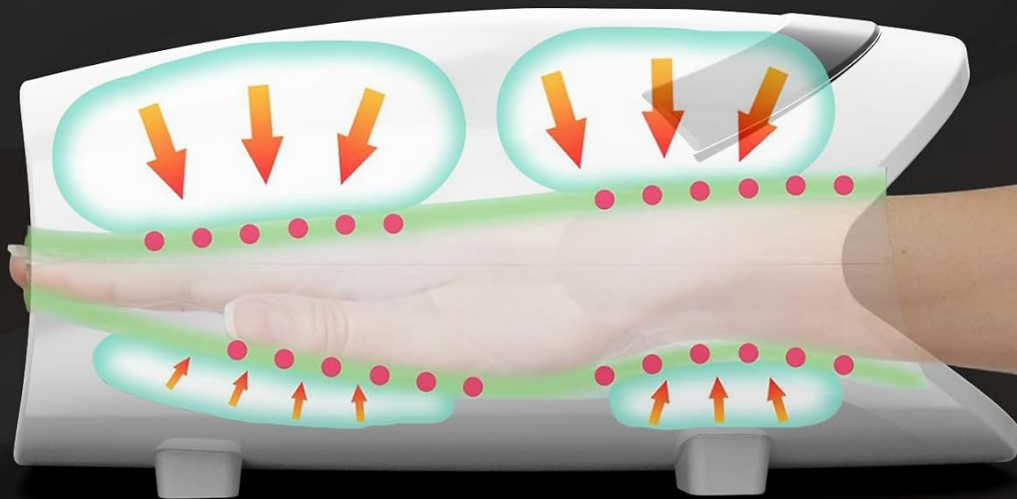

This diagram visually explains how the hand massager operates, with arrows indicating air compression and heat distribution around the hand for therapeutic effect.



This diagram explains the heat therapy feature, showing how to activate and adjust the temperature settings (35°C and 45°C) to promote blood circulation and relieve cold hands.

Helps a Lot in Relaxing Your Hands



This image highlights the diverse benefits of the hand massager for different hand conditions, suggesting relief for common issues like fatigue, stiffness, and discomfort.

6. MAINTENANCE

To ensure the longevity and optimal performance of your Xilent Hand Massager, follow these maintenance guidelines:

- **Cleaning:** Wipe the exterior of the device with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the massager in a cool, dry place away from direct sunlight and extreme temperatures.
- **Battery Care:** Charge the device regularly, even if not in frequent use, to maintain battery health. Avoid completely draining the battery before recharging.
- **Inspection:** Periodically check the device and charging cable for any signs of damage. If any damage is found,

discontinue use and contact customer support.

7. TROUBLESHOOTING

If you encounter any issues with your Xllent Hand Massager, please refer to the following common troubleshooting tips:

- **Device Not Turning On:** Ensure the device is fully charged. Long press the power button for 2-3 seconds.
- **Weak Massage Force:** Adjust the massage force using the Velocity Key. Ensure your hand is properly positioned inside the device.
- **No Heat/Vibration:** Press the Function Key to activate or adjust the heat and vibration settings.
- **Device Shuts Off Unexpectedly:** The device has a 15-minute auto shut-off feature. If it shuts off before 15 minutes, ensure the battery is charged.
- **Charging Issues:** Ensure the USB cable is securely connected to both the device and the power source. Try a different USB port or adapter if necessary.

If problems persist, please contact customer support for further assistance.

8. SPECIFICATIONS

Brand	Xllent
Model Number	NX-S168
Product Dimensions	6.69 x 5 x 10 inches
Item Weight	1.59 Pounds (720 Grams)
Power Source	Battery Powered (1 Lithium Ion battery required)
Material	Plastic
Auto Shut-off	15 minutes

9. WARRANTY AND SUPPORT

For any reason, if this heated hand massager does not meet your expectations, a free return is available within 30 days of purchase. Please refer to your purchase platform for specific return policies and procedures.

Legal Disclaimer: Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question