

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [EMP EZIMEDPIL](#) /
- › [EziMedPil Automatic Pill Dispenser User Manual](#)

EMP EZIMEDPIL B08RZ5PR98

EziMedPil Automatic Pill Dispenser User Manual

Model: B08RZ5PR98

PRODUCT OVERVIEW

The EziMedPil Automatic Pill Dispenser is designed to assist with medication adherence by providing timely audio and visual alerts. It supports programming for up to 6 alarms per day and includes 6 interchangeable dosage templates. The device features an easy-to-read LCD display and a secure locking mechanism with a key to prevent unauthorized access to medications.



Image: EziMedPil Automatic Pill Dispenser, showing both clear and frosted lid options.

This dispenser is ideal for individuals seeking to maintain independence in medication management, as well as for caregivers and home-care agencies. It helps ensure that the correct person takes the right dose of medication at the right time, improving overall medication adherence.

SETUP INSTRUCTIONS

1. Inserting Batteries

The dispenser requires 4 AA batteries (included). To install:

1. Slide the battery door latch on the bottom of the unit to the open position.
2. Remove the battery door.
3. Insert 4 fresh AA batteries, ensuring correct polarity as indicated inside the compartment.
4. Replace the battery door and close the latch securely.



Image: Bottom view of the dispenser with the battery compartment open, showing battery insertion.

2. Opening and Locking the Top Cover

To open the top cover:

1. If the top cover is locked, locate the cover lock on the bottom of the unit.
2. Insert the included key into the lock and turn clockwise to unlock.
3. Place the pill dispenser on a stable surface. Press the 'PRESS' button on the cover handle to release the cover latch.
4. Lift the cover handle to open the cover.



Image: The dispenser with the key inserted into the bottom lock for unlocking the cover.

3. Cleaning and Preparing the Removable Tray

Before first use, and periodically thereafter, clean the removable medication tray:

1. The removable tray lifts straight out of the pill dispenser.
2. Wash the removable tray with warm, soapy water.
3. Allow the tray to dry completely before re-inserting.

Important: DO NOT submerge the pill dispenser base in water. Do not clean any part of the pill dispenser in a dishwasher.



Image: The removable medication tray, ready for cleaning or filling.

4. Placing the Dosage Ring

The dispenser comes with 6 interchangeable dosage rings to match your medication schedule (1 to 6 doses per day). Each ring indicates the number of doses per day and the corresponding programmable time period (e.g., 1 dose for 28 days, 3 doses for 9 days).

1. Insert the removable tray back into the pill dispenser, aligning the arrows.
2. Select the dosage ring that aligns with your medication schedule.
3. Place the selected dosage ring in the removable tray with the tab resting in the slot, aligning with the arrows.
4. Fold the tab downwards to secure the ring.

POWERED BY 4 AA BATTERIES

WIRELESS OPERATION WITH
INCLUDED BATTERIES



FRESBEE DESIGN
EXPLOSION-PROOF
ANTI EXTRUSION

Image: A selection of the 6 included dosage rings, illustrating different dosing schedules.

5. Filling the Dispenser

Carefully place pills into the correct compartments according to your dosage ring and schedule. The dosage in the slot immediately to the right of the tab on the dosage ring will be your starting dose.

Important: Do not place pills in compartments corresponding to the red markings on the dosage ring, as these indicate non-dosing periods.

Once filled, close the cover securely. For extra security, you can lock the cover using the key.

AUTOMATICALLY ROTATES WHEN IT'S TIME FOR YOUR NEXT DOSE

ALARM AND LIGHTS
WILL SOUND FOR 30 MINUTES
OR UNTIL DISPENSER IS INVERTED



SETUP ALARM
TIME



CONTROL 3
LEVEL VOLUME



SELECT 2
ALARM TONES



Image: The dispenser with its tray being filled with various pills, showing the dosage ring in place.

Video: Setup and Operation Guide

Your browser does not support the video tag.


Video: A comprehensive guide demonstrating the setup process, including battery insertion, cover operation, tray preparation, dosage ring placement, and filling the dispenser. It also shows how to take medication and stop alarms.

OPERATING INSTRUCTIONS

Setting Time and Alarms

The EziMedPill dispenser features a large, easy-to-read LCD display and a simple 3-button interface (A, B, C) for programming. You can set the time, alarms, volume, and beep tones (12HR or 24HR format).

28 REMOVABLE SPACIOUS COMPARTMENTS

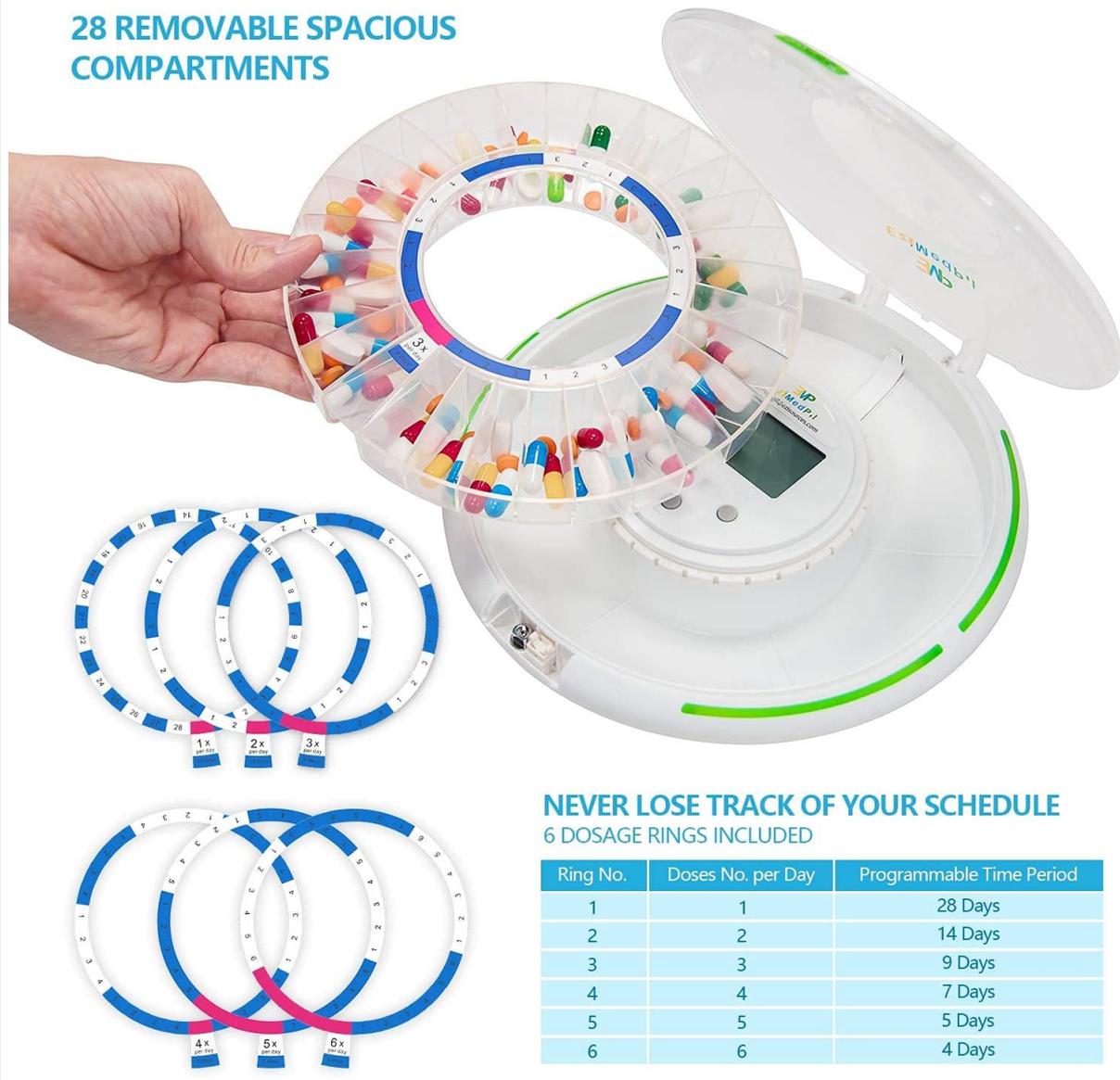

NEVER LOSE TRACK OF YOUR SCHEDULE 6 DOSAGE RINGS INCLUDED

Ring No.	Doses No. per Day	Programmable Time Period
1	1	28 Days
2	2	14 Days
3	3	9 Days
4	4	7 Days
5	5	5 Days
6	6	4 Days

Image: Close-up of the dispenser's LCD display, showing time, alarm settings, volume options, and the three control buttons (A, B, C).

Refer to the detailed instructions in the included user manual for specific steps on programming time and up to 6 daily alarms. Settings are saved even after changing batteries.

Taking Medication and Stopping Alarms

When it is time to take your medication, an alarm will sound, and the light on the dispenser will flash. The alarm will continue for 30 minutes or until deactivated.

1. Open the drug window on the cover by sliding it.
2. Turn the pill dispenser upside down to pour the medications into your hand or a suitable container.
3. This motion will deactivate the alarm and lights.

The pill dispenser is now ready to dispense the next dose on your schedule.



EZIMEDPIL AUTOMATIC PILL DISPENSER



Frosted Lid



Clear Lid

Available in Frosted or Clear Lids

Image: A hand receiving pills dispensed from the EziMedPill Automatic Pill Dispenser.

MAINTENANCE

Regular cleaning of the removable tray is recommended to ensure hygiene and proper functioning. As mentioned in the setup section, wash the removable tray with warm, soapy water and allow it to dry completely. Do not submerge the main unit in water or place any part of the dispenser in a dishwasher.

Periodically check the battery level on the LCD display. Replace batteries promptly when the low battery alert appears to ensure continuous operation.

TROUBLESHOOTING

- **Pills not dispensing or getting stuck:** This dispenser is designed for various pill sizes, but very small or very large pills may occasionally cause issues. Ensure pills are not overcrowded in compartments. If pills are consistently getting stuck, consider if the size is compatible with the dispenser's mechanism.
- **Pills accessible when locked:** The locking mechanism is designed to prevent unauthorized opening of the lid. However, if the unit is inverted, smaller pills might still be able to fall out through the dispensing slot. Always keep the dispenser upright to prevent accidental dispensing.
- **"Err" message on display:** An "Err" message typically indicates a malfunction. First, try replacing

the batteries with fresh ones. If the issue persists, attempt a manual reset (refer to the full instruction manual for reset procedures). If the error continues, contact customer support.

- **Alarm volume too low:** The dispenser offers 3 volume options. Ensure the volume is set to an appropriate level for the user. Refer to the programming instructions to adjust the alarm volume.

SPECIFICATIONS

Product Dimensions	8.86 x 8.86 x 2.68 inches; 1.9 Pounds
Batteries	4 AA batteries required (included)
Material	Acrylonitrile Butadiene Styrene (ABS)
Special Feature	Automatic Dispensing, Audio & Visual Alerts, Lockable
Included Components	User Manual, 2 Keys, Pill Tray, 4 AA Batteries, Dosage Rings

WARRANTY INFORMATION

Specific warranty details for the EziMedPil Automatic Pill Dispenser are not provided in the available product information. Please refer to the documentation included with your purchase or contact the manufacturer directly for warranty terms and conditions.

CUSTOMER SUPPORT

For further assistance, technical support, or inquiries regarding your EziMedPil Automatic Pill Dispenser, please contact the manufacturer:

- **Brand:** EMP EZIMEDPIL
- **Email:** ezimedpil@ezisources.com