



Manuals.plus /

› **Nicwell /**

› Nicwell Water Dental Flosser Teeth Pick - Cordless Portable 4 Modes Rechargeable Flossers Oral Irrigator IPX7 Waterproof Travel Gum Floss Cleaner Deep Cleaning Flossing Picks for Home Daily Shower

Nicwell

Nicwell Water Dental Flosser Instruction Manual

Model: F5025

PRODUCT OVERVIEW

The Nicwell Water Dental Flosser is a cordless, portable oral irrigator designed to provide a comprehensive and effective oral cleaning experience. It utilizes advanced water pulse technology to remove plaque, improve gum health, and ensure fresh breath. This manual provides essential information for the safe and efficient operation, maintenance, and troubleshooting of your device.



Figure 1: Nicwell Water Dental Flosser and included accessories.

150% More Effective



Figure 2: Illustration of the dual-thread water pulse technology for enhanced cleaning.

SETUP

1. Initial Charging

Before first use, fully charge the Nicwell Water Dental Flosser. A full charge typically takes approximately 4 hours and provides up to 30 days of battery life, depending on usage.

1. Connect the USB charging cable to the charging port on the device.
2. Plug the other end of the USB cable into a compatible power adapter (not included), laptop, car charger, or power bank.
3. The indicator light will show charging status. Ensure the waterproof plug is closed tightly after charging to maintain IPX7 rating.

IPX7 waterproof

suitable for use in wet environments
without the risk of damage



Figure 3: Charging options and battery life indication.

2. Attaching Nozzles

The flosser comes with multiple interchangeable jet tips for various oral care needs.

1. Select the desired nozzle (Classic, Tongue Scraper, or Orthodontic).
2. Insert the nozzle into the hole at the top of the flosser handle until it clicks securely into place.
3. To remove, press the nozzle release button located on the side of the handle and pull the nozzle out.

Lightweight Water Tank

Perfect water tank, Capacious to contain ample water yet compact to fit into bags



Figure 4: Overview of the five interchangeable jet tips and their functions.

3. Filling the Water Tank

The device features a 200ml detachable water tank for easy filling and cleaning.

1. Detach the water tank from the main unit by twisting or pulling it off (refer to product diagram for specific detachment method).
2. Fill the tank with lukewarm water. You may also add mouthwash as recommended by your dental professional.
3. Reattach the water tank securely to the main unit.


5 interchangeable jet tips included for your whole family




various nozzles & function



3 CLASSIC NOZZLES
Daily Cleaning

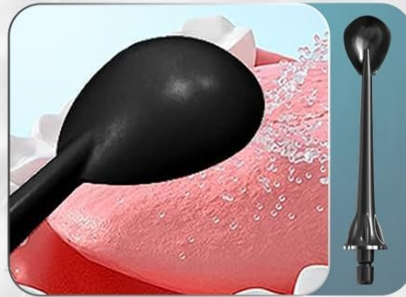

1 TONGUE SCRAPER
Removes bacteria build
up for fresher breath

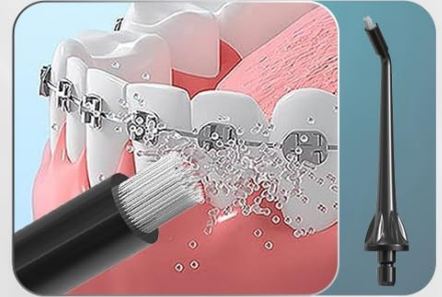

1 ORTHODONTIC NOZZLE
Remove food debris
in the brace

Figure 5: Detachable water tank for easy filling and cleaning.

OPERATING INSTRUCTIONS

1. Power On/Off and Mode Selection

The Nicwell Water Dental Flosser offers four distinct cleaning modes and 11 adjustable intensity levels to suit individual needs.

1. Press the 'ON/OFF' button to turn the device on.
2. Press the 'MODE' button to cycle through the available modes: Normal, Soft, Pulse, and Custom. The selected mode will be indicated by a light.
3. In Custom mode, you can manually adjust the pressure intensity to your preference. The device will remember your last used setting for convenience.
4. Press the 'ON/OFF' button again to turn the device off.



Figure 6: Control panel with mode selection options.

2. Proper Usage Technique

For optimal results and to prevent splashing, follow these steps:

1. Place the jet tip into your mouth before turning on the device.
2. Lean over a sink to allow water to drain.
3. Aim the tip at a 90-degree angle to your gumline.
4. Press the 'ON/OFF' button to start the water flow.
5. Guide the water stream along your gumline and between your teeth, pausing briefly between each tooth.
6. Continue until you have cleaned all areas of your mouth.

Video 1: A comprehensive overview of the Nicwell Water Dental Flosser, demonstrating its features and benefits.

Video 2: A brief demonstration highlighting the portability and lightweight design of the water flosser.

1. Cleaning the Unit and Nozzles

Regular cleaning ensures optimal performance and hygiene.

- After each use, empty any remaining water from the tank.
- Detach the water tank and rinse it thoroughly with clean water. The detachable design allows for easy access and cleaning.
- Rinse the nozzle under running water.
- Wipe the main unit with a soft, damp cloth. Do not immerse the main unit in water, although it is IPX7 waterproof for shower use, it should not be submerged for cleaning.
- Allow all parts to air dry before reassembling or storing.

The Best For The GO

lightweight and long life battery,
more suitable for travel.

30 Days Last for

200ml Water Tank

Filled with warm water or mouth wash is
beneficial for teeth

POWER ADAPTER

LAPTOP

POWER BANK

PORTABLE

LIGHTWEIGHT

Figure 7: The IPX7 waterproof rating allows for safe use in wet environments.

2. Storage

Store the flosser in a cool, dry place away from direct sunlight when not in use. Ensure the water tank is empty and dry to prevent mildew or bacterial growth.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on or charge.	Low battery; charging cable not properly connected; internal fault.	Ensure the device is fully charged. Check USB cable connection. Try a different charging source. If issues persist, contact customer support.
Weak water pressure.	Low battery; nozzle clogged; incorrect mode selected.	Recharge the device. Remove and clean the nozzle. Select a higher pressure mode.
Water leaks from the unit.	Water tank or nozzle not securely attached; damaged seal.	Ensure the water tank and nozzle are firmly attached. Check for any visible damage to seals. If leaking continues, contact customer support.
Excessive splashing during use.	Incorrect technique; water pressure too high.	Ensure the tip is placed in the mouth before turning on. Lean over the sink. Adjust to a lower pressure mode.

SPECIFICATIONS

- **Model Number:** F5025
- **Power Source:** Battery Powered (1 Lithium Ion battery required)
- **Water Tank Capacity:** 200ml
- **Water Pressure Pulses:** 1400-1800 times per minute
- **Waterproof Rating:** IPX7
- **Dimensions:** 8.03 x 4.02 x 2.72 inches
- **Weight:** 14.07 ounces
- **Manufacturer:** Nicwell
- **Country of Origin:** China

WARRANTY AND SUPPORT

Nicwell is committed to providing reliable products and excellent customer service.

- **Warranty:** The product comes with a 2-year warranty.
- **Customer Support:** 24-hour customer support is available.
- **Technical Support:** Lifetime technical support is provided.



Figure 8: Nicwell's commitment to customer service and support.

Important Safety Information

This product is not a toy. Keep out of reach of children.

Legal Disclaimer

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

