

Scienlodic c664352d-9198-43bb-9b41-d4681c9ae1dd

Scienlodic 24-Piece Electric Cupping Therapy Set User Manual

Brand: Scienlodic

Model: c664352d-9198-43bb-9b41-d4681c9ae1dd

INTRODUCTION

The Scienlodic 24-Piece Electric Cupping Therapy Set is designed for home use to provide cupping therapy, promote blood circulation, and aid in muscle relaxation and recovery. This portable and rechargeable device offers adjustable suction levels and includes various cup sizes for different body parts.

Key Features:

- **Rechargeable Electric Pump:** Provides consistent and adjustable suction without manual pumping.
- **5 Adjustable Suction Levels:** Customize intensity for various body areas and comfort levels.
- **24 Vacuum Suction Cups:** Includes multiple sizes and curved cups for versatile application on different body parts and joints.
- **Lightweight and Portable:** The pump weighs approximately 300g, making it easy to handle and transport.
- **Automatic Magnetic Needles:** Designed to self-adjust during suction to minimize skin discomfort.
- **Gua Sha Functionality:** Can be used for scraping therapy to promote blood flow.

SAFETY INFORMATION AND PRECAUTIONS

Before using the Scienlodic Electric Cupping Therapy Set, please read and understand the following safety guidelines:

- **Consult a Doctor:** Individuals with malignant tumors, acute diseases, severe hypertension, pacemakers, or electronic implants, as well as pregnant or menstruating women, those with varicose veins, or skin wounds, should consult a medical professional before use.
- **Skin Preparation:** Apply massage oils or lubricants to the skin before cupping to ensure smooth application and prevent skin damage.
- **Not Waterproof:** The electric pump and machine components are not waterproof. Do not immerse in water or use in wet environments. Clean the surface of the machine after each use with a dry or slightly damp cloth.

- **Proper Seal:** Ensure a good seal between the cupping cup and the machine/extension tube to achieve sufficient suction.
- **Adjust Suction:** Start with the lowest suction level and gradually increase as needed. Different body parts may require different suction intensities.
- **Avoid Over-Suction:** Do not leave cups on for excessive periods or apply excessively strong suction, as this may cause discomfort or skin marks.
- **Charging:** Ensure the device is fully charged before use for optimal performance.

PACKAGE CONTENTS

The Scienlodoc 24-Piece Electric Cupping Therapy Set includes the following items:

- 1x Electric Cupping Pump
- 24x Vacuum Suction Cups (various sizes, including curved cups)
- 1x Flexible Extension Tube
- 1x Charging Cable
- 1x Gua Sha Tool
- Magnetic Needles (quantity not specified, but mentioned in features)
- User Manual (this document)

Dimensions of 24 Cups



Image: Various sizes of the 24 vacuum suction cups included in the set, with their respective dimensions.



Image: Close-up of the durable, high-quality plastic cupping cups.



Image: The flexible extension tube for reaching difficult areas.



Image: The self-adjusting magnetic needles for enhanced therapy.



Image: The smooth Gua Sha tool for scraping therapy.

SETUP

1. **Charge the Device:** Before first use, fully charge the electric cupping pump using the provided charging cable. The indicator lights will show charging status.
2. **Prepare Skin:** Clean the area of skin where cupping will be applied. Apply a generous amount of massage oil or lubricant to the skin to ensure a good seal and allow for smooth movement if sliding cupping is desired.
3. **Attach Extension Tube (Optional):** For hard-to-reach areas like the back, attach one end of the flexible extension tube to the electric pump and the other end to the cupping cup. Ensure the arrowhead end of the extension tube is attached to the cup.
4. **Select Cup Size:** Choose an appropriate cup size for the body part you intend to treat. Smaller cups are suitable for joints and smaller muscle groups, while larger cups are for broader areas.



Image: Visual guide for setting up the electric cupping device, including attaching the extension tube and cups.

OPERATING INSTRUCTIONS

1. **Power On:** Long press the power button on the electric pump for 2 seconds to turn it on. The indicator lights will illuminate.
2. **Adjust Suction Level:** Click the power button repeatedly to cycle through the 5 available suction levels. Select a level appropriate for your comfort and the treatment area. Start with the lowest setting.
3. **Apply Cup:** Place the chosen cupping cup firmly onto the prepared skin, ensuring a tight seal.
4. **Start Suction:** Press the trigger button on the electric pump to begin the suction process. The pump will draw air out of the cup, creating a vacuum that pulls the skin and underlying tissue into the cup.
5. **Monitor Skin:** Observe the skin inside the cup. It should be pulled up, but not to the point of severe pain. Adjust suction if necessary.
6. **Release Cup:** To release the cup, twist the yellow air valve on top of the cup to break the vacuum. The cup will detach easily.
7. **Sliding Cupping (Gua Sha):** If performing sliding cupping or Gua Sha, ensure sufficient oil is applied to the skin. With the cup attached and suction applied, gently slide the cup across the skin in the desired direction. Alternatively, use the dedicated Gua Sha tool directly on oiled skin.



Efficient & Time-Saving


Strong Suction Electric Pump
for Easy Cupping

Image: A hand operating the electric pump to apply cupping cups to a person's back, demonstrating efficient use.

Ergonomic Curved Cups for Joint Relief



Image: Examples of curved cupping cups being applied to various joints like the elbow, knee, and shoulder for targeted relief.

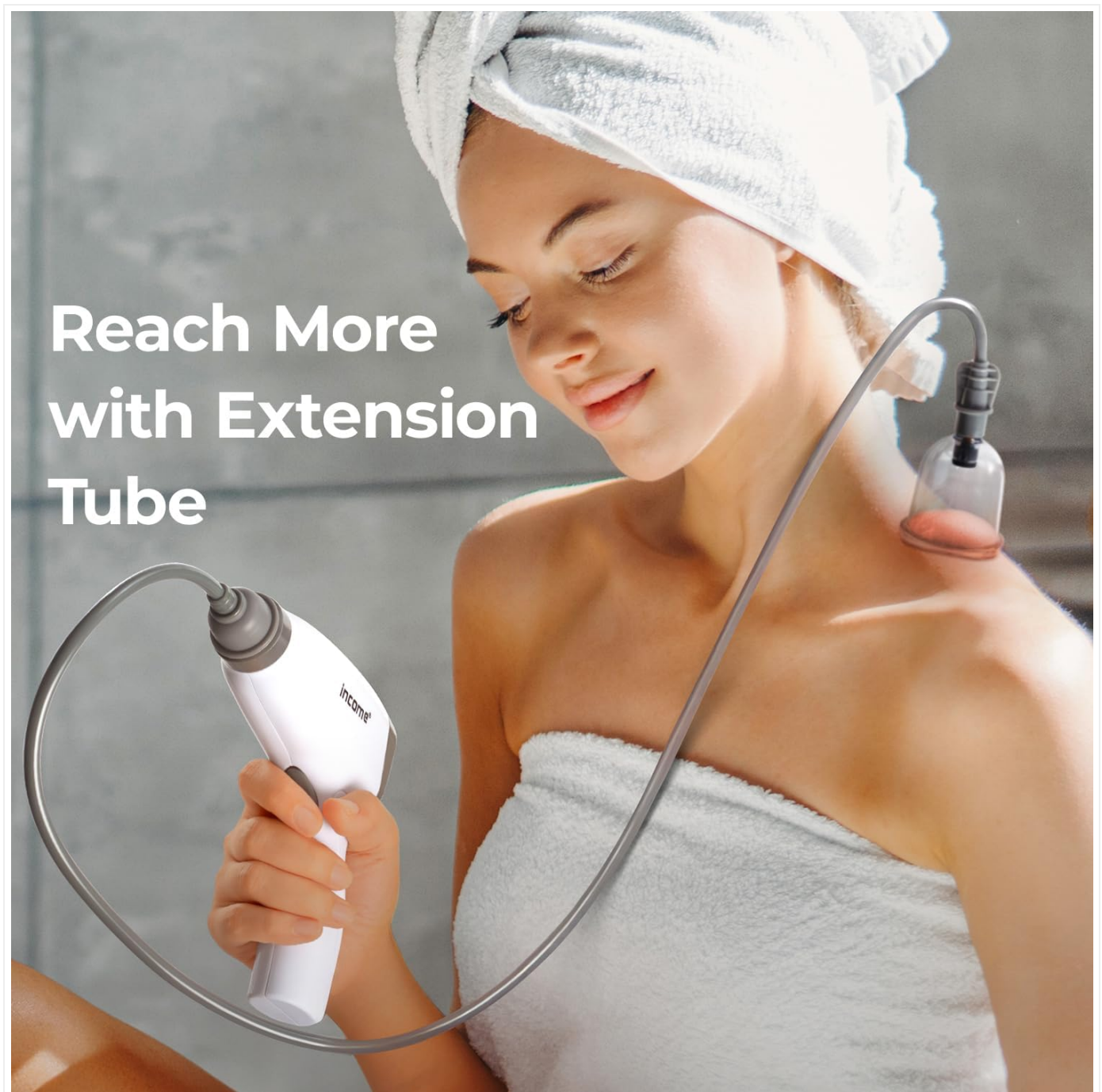

Reach More with Extension Tube

Image: A user demonstrating how to use the extension tube to reach and apply cupping to their own shoulder.

Your browser does not support the video tag.

Video: Demonstrates the use of the Portable Electric Cupping Set, including attaching cups, applying suction, and using the Gua Sha function.

MAINTENANCE

- **Cleaning Cups:** After each use, thoroughly clean the cupping cups with mild soap and water. Ensure all oil and skin residue are removed. Rinse well and allow to air dry completely before storage.
- **Cleaning Pump:** Wipe the surface of the electric pump with a soft, dry, or slightly damp cloth. Do not use abrasive cleaners or immerse the pump in water.
- **Storage:** Store the cupping set in a clean, dry place, away from direct sunlight and extreme temperatures. Keep all components organized in their original packaging or a suitable storage container to prevent damage.
- **Battery Care:** To prolong battery life, avoid fully discharging the device frequently. Charge it regularly, even if not in constant use.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Insufficient Suction	<ul style="list-style-type: none">◦ Poor seal between cup and skin.◦ Extension tube not properly attached.◦ Low battery.◦ Air valve on cup is open.	<ul style="list-style-type: none">◦ Apply more massage oil and ensure cup is flush against skin.◦ Check connections of the extension tube to both pump and cup.◦ Recharge the electric pump.◦ Ensure the air valve on the cup is twisted shut before applying suction.
Pump not turning on	<ul style="list-style-type: none">◦ Low or no battery charge.◦ Power button not pressed long enough.	<ul style="list-style-type: none">◦ Charge the device fully.◦ Press and hold the power button for 2 seconds until indicator lights appear.
Skin discomfort or pain	<ul style="list-style-type: none">◦ Suction level too high.◦ Cup left on for too long.◦ Skin not adequately lubricated.	<ul style="list-style-type: none">◦ Reduce the suction level.◦ Limit cupping sessions to recommended durations (e.g., 5-10 minutes per area).◦ Apply more massage oil to the skin.

SPECIFICATIONS

Feature	Detail
Product Dimensions	15.75 x 10.24 x 3.15 inches
Item Weight	2.29 Pounds
Model Number	c664352d-9198-43bb-9b41-d4681c9ae1dd
Power Source	Battery Powered (1 Lithium Polymer battery included)
Special Features	Rechargeable, Adjustable Suction (5 levels), Automatic Magnetic Needles
Massager Form	Massage Cup
Specific Uses	Electronic Cupping Device, Scraping (Gua Sha)
Product Benefits	Metabolism Improvement, Muscle Relaxation, Blood Circulation
Water Resistance Level	Not Water Resistant
Material Type	Silicone (for certain components, cups are plastic)
Color	White

WARRANTY AND SUPPORT

Scienlodic provides a free replacement warranty for any quality and logistics issues. For further assistance, troubleshooting, or warranty claims, please contact Scienlodic customer support through your purchase platform or the official Scienlodic website.

