

SINGCALL 6100

SINGCALL Wireless Caregiver Calling System Receiver User Manual

Brand: SINGCALL | Model: 6100

1. PRODUCT OVERVIEW

The SINGCALL Wireless Caregiver Calling System Mobile Watch Receiver (Model 6100) is designed to work with compatible call buttons/pagers to create a simple and efficient calling system. This receiver is ideal for various environments such as hospitals, care facilities, restaurants, and homes, enabling quick and discreet communication between individuals needing assistance and caregivers or service staff. Please note: This product is a wireless watch receiver only; call buttons/pagers are not included and must be purchased separately to form a complete system.

Key Features

- Mobile pager receiver with lanyard for convenient wearing around the neck.
- Integrated calendar, clock, and alarm clock functions.
- Rechargeable battery with extended standby time (2-3 days).
- Fast delete function to remove individual or all registered numbers with a single key.
- Encryption function to prevent accidental loss of information.
- Supports downloading registration information via software for easier setup.



Figure 1.1: SINGCALL Mobile Watch Receiver with Lanyard. This image shows the compact, circular receiver with its display and a black lanyard attached, ready for wear.



Figure 1.2: Front view of the SINGCALL Mobile Watch Receiver. The image highlights the clear digital display showing time and battery status, along with 'FUN' and 'ENT' buttons.



Figure 1.3: Close-up view of the watch receiver's display, showing the time (12:00) and date (08/08) clearly.



Figure 1.4: Various application areas for the SINGCALL system, including beadhouse, cadre's sanitarium, family settings, and hotels, illustrating its versatility.

2. SETUP AND PAIRING

The SINGCALL 6100 receiver requires pairing with compatible SINGCALL call buttons/pagers to function as a complete calling system. The receiver can store up to 500 call numbers. Registration information can be managed via dedicated software, simplifying the setup process.

System Integration

To set up your system, position the call buttons in areas where assistance may be needed. When a button is pressed, the receiver will display the corresponding number and provide an alert (vibration/sound).



Figure 2.1: Schematic diagram illustrating the usage of the SINGCALL restaurant calling system. It shows diners pressing pagers, which send signals to a screen receiver and a mobile receiver worn by service staff, enabling efficient service.

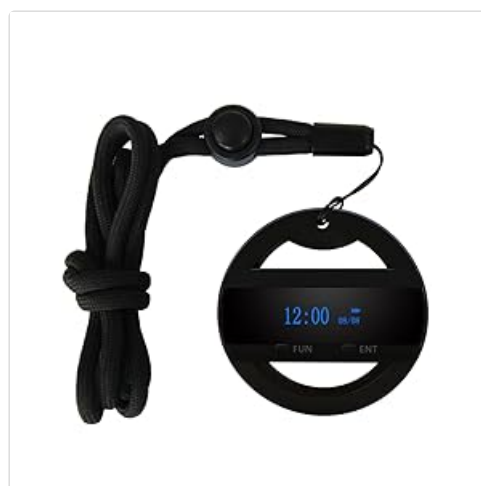


Figure 2.2: An overview of the SINGCALL watch receiver alongside several compatible call buttons, demonstrating the components of a complete calling system.

Video Demonstration: Receiver and Button Interaction

Video 2.1: This video demonstrates the interaction between the SINGCALL 6100 receiver and call buttons. It shows how pressing a call button triggers an alert on the receiver, displaying the corresponding call number.

3. OPERATING INSTRUCTIONS

Receiving Calls

When a paired call button is pressed, the receiver will alert the user through vibration and/or sound, and the corresponding call number will be displayed on the screen. This allows for immediate identification of the location or individual requiring assistance.

Managing Call Records

The receiver features a 'Fast Delete Function' which allows you to quickly clear call records. You can choose to delete a single registered number or all registered numbers with a simple operation, ensuring privacy and efficient management of alerts.

Time and Alarm Functions

The device includes a built-in calendar, clock, and alarm clock function, providing additional utility beyond its primary calling system role. Ensure the time is set correctly for optimal use of these features.



Figure 3.1: A waiter wearing the SINGCALL watch receiver while attending to a customer in a restaurant, demonstrating practical use in a service environment.

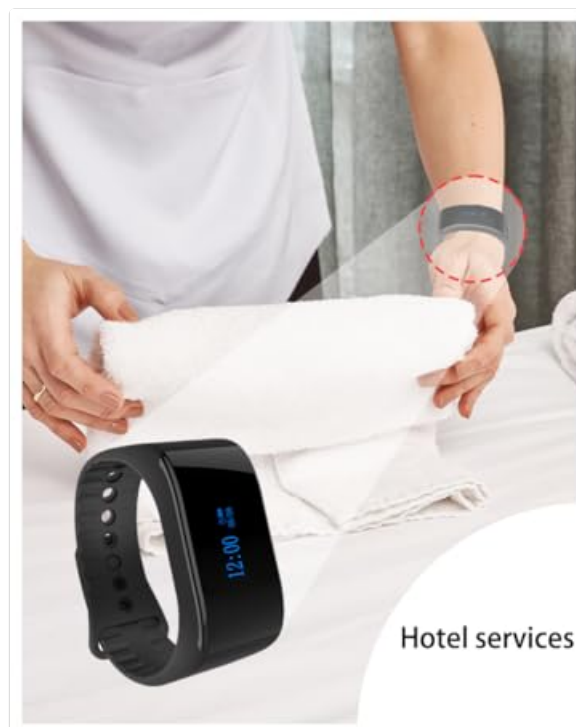


Figure 3.2: A caregiver or hotel staff member wearing the SINGCALL watch receiver, illustrating its use in hotel services for efficient communication.

4. MAINTENANCE

- **Cleaning:** Wipe the receiver with a soft, dry cloth. Do not use abrasive cleaners or solvents.
- **Charging:** Recharge the battery regularly to ensure optimal performance and standby time. Avoid

overcharging or completely draining the battery for extended periods.

- **Storage:** Store the device in a cool, dry place away from direct sunlight and extreme temperatures when not in use for long periods.
- **Water Resistance:** The device is not waterproof. Avoid exposure to water or high humidity.

5. TROUBLESHOOTING

- **Receiver not responding to calls:** Ensure the receiver is charged and powered on. Verify that the call buttons are functioning and properly paired with the receiver. Check for any obstructions or excessive distance between the button and receiver.
- **Display issues:** If the display is dim or unresponsive, recharge the battery. If the issue persists, contact customer support.
- **Information loss:** The receiver has an encryption function to prevent accidental data loss. If you experience unexpected data loss, ensure proper operation and consider re-registering call numbers if necessary.
- **Short battery life:** Ensure the device is fully charged. Battery life can vary based on usage and environmental factors. If battery life significantly degrades, it may indicate the need for battery replacement.

6. SPECIFICATIONS

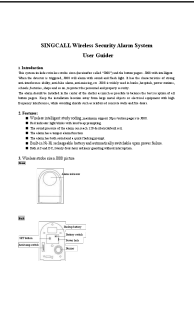


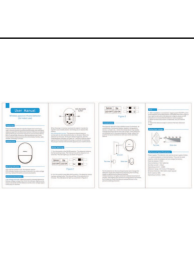

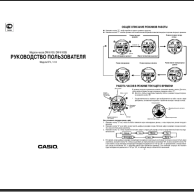
Attribute	Detail
Brand	SINGCALL
Model Number	6100
Manufacturer	EASTWAY INTERNATIONAL (HK) LIMITED
Item Weight	3.52 ounces
Product Dimensions	2.36 x 2.36 x 0.79 inches
Material Type	ABS
Date First Available	January 10, 2021

7. WARRANTY AND SUPPORT

For warranty information, technical support, or any inquiries regarding your SINGCALL Wireless Caregiver Calling System Receiver, please contact SINGCALL customer service directly. You can find more information and contact details on the official SINGCALL store page.

Visit the [SINGCALL Store](#) for additional products and support resources.

Related Documents - 6100

	<p>SINGCALL J008 Wireless Security Alarm System User Guide</p> <p>User guide for the SINGCALL J008 Wireless Security Alarm System, detailing its features, operation instructions for pairing pagers, volume and alarm time settings, and technical specifications.</p>
	<p>Philips Sonicare 6100: Gebruiksaanwijzing en Specificaties voor Patiënten</p> <p>Gedetailleerde informatie over de Philips Sonicare 6100 elektrische tandenborstel, inclusief functies zoals de A3 All-in-One en S2 Sensitive opzetborstels, 6 poetsstanden, druksensor, BrushSync-technologie en een batterijduur tot 3 weken. Ideaal voor tandheelkundige professionals die poetsinstructies geven.</p>
	<p>COLAMY HL-6100 Massage Chair: Assembly, Usage, and Safety Guide</p> <p>This comprehensive manual provides detailed instructions for assembling, operating, and maintaining the COLAMY HL-6100 Massage Chair. Includes safety precautions, warranty information, and technical specifications.</p>
	<p>SINGCALL Wireless Passive Infrared Detector User Manual</p> <p>Comprehensive user manual for the SINGCALL wireless passive infrared detector (indoor use), covering features, appearance, anti-dismantle function, mode settings, installation, testing, detection range, working indicators, functions, and technical specifications. Learn how to set up and use your PIR motion sensor.</p>
	<p>KENWOOD DVR-605/DVR-6100/DVR-6100K DVD Receiver Instruction Manual</p> <p>This instruction manual provides detailed information on the setup, operation, and features of the KENWOOD DVR-605, DVR-6100, and DVR-6100K DVD Receiver models.</p>
	<p>Casio DW-6100, DW-6100G: Руководство пользователя</p> <p>Подробное руководство пользователя для часов Casio DW-6100 и DW-6100G (модули 974, 1414), охватывающее настройку времени, термометр, секундомер, будильники, восход/заход солнца и часовые пояса.</p>

