

[manuals.plus](#) /

› [mornwell](#) /

› [Mornwell Water Flosser User Manual](#)

mornwell RST5012-BK

Mornwell Water Flosser User Manual

Model: RST5012-BK

INTRODUCTION

Thank you for choosing the Mornwell Water Flosser. This device is designed to provide effective oral hygiene by removing food debris and plaque from between teeth and below the gumline, promoting healthier gums and a cleaner mouth. Please read this manual thoroughly before use to ensure proper operation and maintenance.



Figure 1: Mornwell Water Flosser and included accessories.

IMPORTANT SAFETY INFORMATION

For your safety, please observe the following precautions:

- Keep the device away from water sources when charging.
- Be mindful of electrical safety when using the device in the shower.
- Do not add mouthwash or other solutions to the water tank that are not explicitly approved by the manufacturer, as this may damage the rubber seals.
- This device is intended for oral hygiene only. Do not use for any other purpose.
- Keep out of reach of children.

PRODUCT COMPONENTS

The Mornwell Water Flosser package includes:

- Mornwell Water Flosser Unit
- USB Charging Cable
- Multiple Jet Tips (Standard and Orthodontic/Tongue Cleaner depending on configuration)



Figure 2: Different nozzle configurations that may be included with your unit.

SETUP

Charging the Device

Before first use, fully charge the water flosser. The built-in 1400mAh lithium battery takes approximately 4 hours to fully charge and can last up to 30 days on a single charge, depending on usage.

1. Locate the charging port on the device, usually covered by a protective plug.
2. Insert the USB charging cable into the port.
3. Connect the other end of the USB cable to a compatible power source (e.g., USB adapter, computer USB port).
4. The charging indicator lights will illuminate during charging and turn off or change color when fully charged.



Figure 3: Charging the water flosser.

Filling the Water Tank

The 330ml detachable water tank provides sufficient capacity for a complete cleaning session without frequent refills.

1. Detach the water tank from the main unit by sliding it down.
2. Open the leak-proof drain lid.
3. Fill the tank with lukewarm water. For sensitive teeth, warm water is recommended as the device does not have a water heater.
4. Close the drain lid securely and reattach the tank to the main unit until it clicks into place.



Figure 4: Detachable water tank for easy filling and cleaning.

OPERATING INSTRUCTIONS

Selecting a Mode

The Mornwell Water Flosser offers three distinct cleaning modes to suit your oral care needs:

- **Normal Mode:** Provides a relatively strong water flow, suitable for daily cleaning for most users.
- **Soft Mode:** Offers a gentle water flow, recommended for first-time users or individuals with sensitive gums.
- **Pulse Mode:** Delivers an intermittent water flow, ideal for gum massage and enhanced cleaning.

Press the 'Mode' button to cycle through the available settings. The indicator lights will show the currently selected mode.

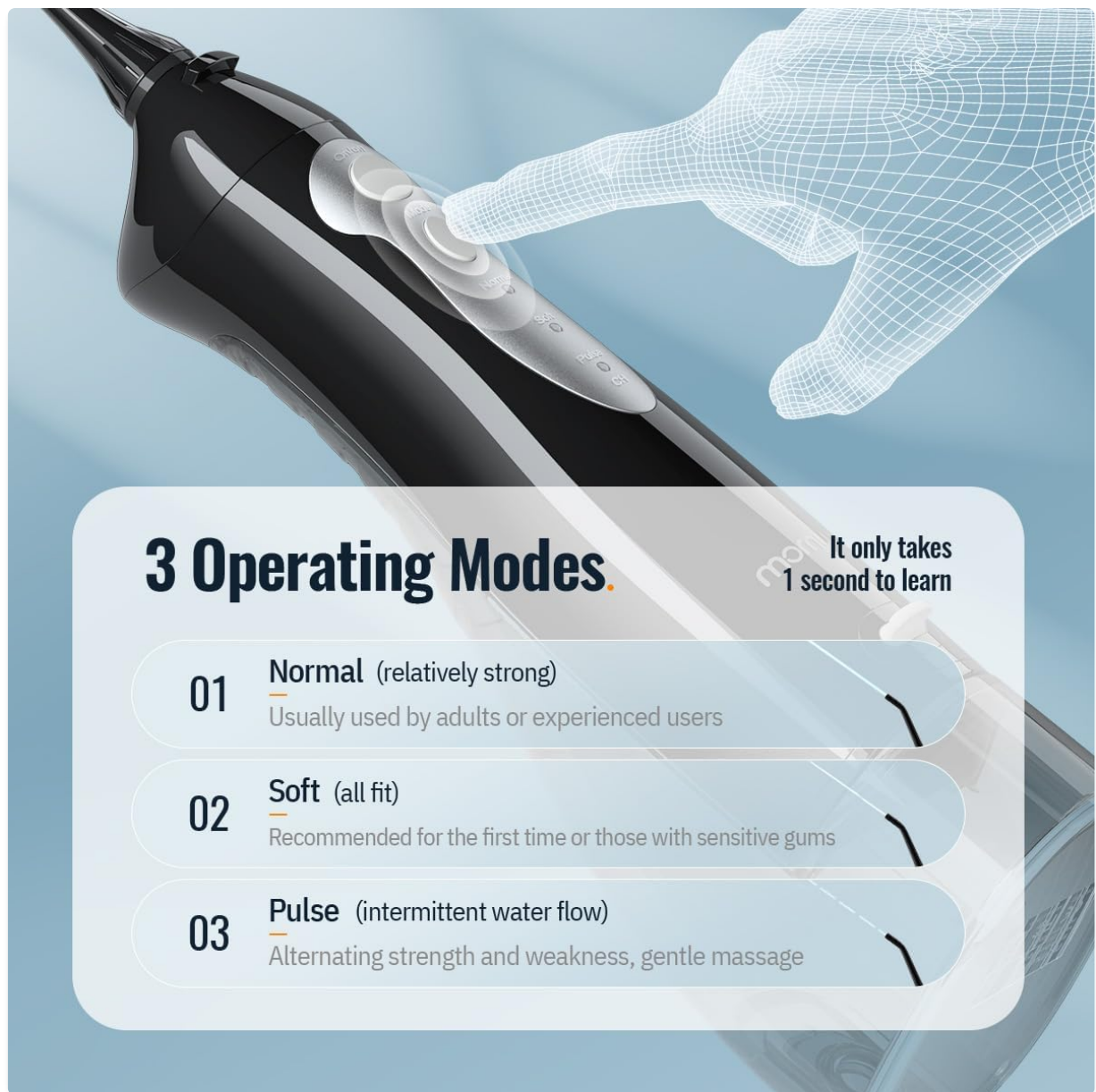


Figure 5: Three operating modes for customized cleaning.

Flossing Technique

Follow these steps for effective water flossing:

1. Select your desired mode (Normal, Soft, or Pulse).
2. Lean over a sink to allow water to drain.
3. Place the flosser tip into your mouth, aiming it at the gumline at a 90-degree angle. The 360° rotatable nozzle allows easy access to all areas.
4. Press the 'On/Off' button to start the water flow.
5. Move the tip along the gumline, pausing briefly between teeth. Ensure the water stream reaches all surfaces.
6. Keep your mouth slightly open to allow water to flow into the sink.
7. The device has an auto-timer design for optimal cleaning duration.

Effectively Remove up to **99.99%** Plaque



Figure 6: How the water flosser effectively removes plaque.

MAINTENANCE

Cleaning the Water Tank and Tips

Regular cleaning ensures hygiene and extends the life of your device.

- After each use, empty any remaining water from the tank.
- Detach the tank and rinse it thoroughly with clean water.
- Rinse the jet tips under running water.
- The device features internal and external double IPX7 waterproof protection, allowing for safe use and rinsing in the shower.
- Allow all components to air dry before reassembling or storing.



FULL BODY WATERPROOF

Use anytime, anywhere

Figure 7: The water flosser is fully waterproof for easy cleaning.

TROUBLESHOOTING

If you encounter issues with your Mornwell Water Flosser, consider the following:

- **No water flow:** Ensure the water tank is properly filled and securely attached. Check if the nozzle is correctly inserted.
- **Weak water pressure:** Verify the selected mode. If on 'Soft' mode, try 'Normal' or 'Pulse'. Ensure the device is fully charged.
- **Device not turning on:** Confirm the device is charged. If not, connect it to the charger.
- **Leaking:** Ensure the water tank is properly seated and the drain lid is tightly closed.

For persistent issues, please refer to the warranty and support section for contact information.

SPECIFICATIONS

Feature	Detail
---------	--------



Product Dimensions	9.49 x 4.33 x 2.72 inches
Item Model Number	RST5012-BK
Weight	15.17 ounces
Batteries	1 Lithium Ion battery (included)
Power Source	Battery Powered
Water Tank Capacity	330ML
Waterproof Rating	IPX7
Operating Modes	Normal, Soft, Pulse
Special Features	Portable, Timer, Rotatable Tip, Rechargeable, Multiple Tips
Wattage	3 Watt-hours




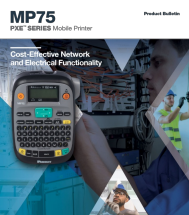
WARRANTY & SUPPORT

Mornwell Cordless Water Flosser is Certificated by CE/FCC. Mornwell electric water flossers always offer a **12-month 'replace, don't repair' warranty**.

If you have any questions or require support, please feel free to contact us. We aim to respond within 24 hours. For additional information, you may refer to the official user manual PDF:[Mornwell Water Flosser User Manual \(PDF\)](#)

Related Documents - RST5012-BK

	<p>SmarTooth STL-01-BK Dental Monitoring Device User Guide</p> <p>Comprehensive guide for the SmarTooth STL-01-BK Dental Monitoring Device, covering package contents, specifications, installation, usage, calibration, measurement, interpretation of results, and safety precautions.</p>
	<p>Oral-B iO Series IO9 M9.4B2.2A BK 9 Connected Rechargeable Electric Toothbrush Instruction Manual</p> <p>Comprehensive instruction manual for the Oral-B iO Series IO9 M9.4B2.2A BK 9 Connected Rechargeable Electric Toothbrush, covering setup, usage, features, maintenance, troubleshooting, and warranty information.</p>

	<p>MST RMX EX S PRO RC Drift Car Upgrade Guide</p> <p>Official upgrade guide for the MST RMX EX S PRO 1/10 scale RC RWD high-performance drift car, detailing recommended parts and comparisons for enhanced performance.</p>
	<p>MST RMX 2.5 RTR 1/10 Scale RC Drift Car Exploded Manual and Upgrade Guide</p> <p>Comprehensive exploded view manual and upgrade guide for the MST RMX 2.5 RTR 1/10 scale RC Drift Car. Includes part numbers, specifications, and assembly instructions.</p>
	<p>Panduit MP75 PXE Series Mobile Printer - Product Bulletin</p> <p>Detailed product bulletin for the Panduit MP75 PXE Series Mobile Printer, covering its features, specifications, and ordering information for various label cassette types, including comparisons with other models. Learn about its network and electrical functionality.</p>
	<p>Panduit MP75 PXE Series Mobile Printer: Features, Specifications, and Ordering Information</p> <p>Comprehensive product bulletin for the Panduit MP75 PXE Series mobile label printer, detailing its features, dimensions, technical specifications, and a wide range of compatible label cassette options for industrial and network identification.</p>