

## AIMENGXI AIMENGXI002

# AIMENGXI 2-in-1 Diamond Microdermabrasion and Blackhead Removal Machine

Model: AIMENGXI002

## 1. INTRODUCTION

This manual provides detailed instructions for the safe and effective use of your AIMENGXI 2-in-1 Diamond Microdermabrasion and Blackhead Removal Machine. This device is designed for facial skin care, offering both diamond microdermabrasion for exfoliation and a vacuum function for blackhead removal. Proper use can help improve skin texture, reduce the appearance of wrinkles, and promote clearer skin.

## 2. SAFETY INFORMATION AND PRECAUTIONS

Read all instructions carefully before using the device. Failure to follow safety guidelines may result in injury or damage to the device.

### General Safety

- Ensure your skin is clean before using the device.
- Do not point the diamond head directly at the eyes.
- Keep the device away from water or wet environments.
- Unplug the device when not in use or before cleaning.
- Store in a cool, dry place, out of reach of children.

### Prohibited Use Areas and Conditions

Do not use the device on the following areas or under these conditions:



- Sensitive skin.
- Skin rash.
- Damaged, inflamed, or infected skin.
- Surgical scars less than 12 months old.

- Swollen skin.
- Directly on eyeballs.
- Genital areas.


### 3. PACKAGE CONTENTS

Verify that all items are included in your package:


What you will get?


**1 \*Professional Dermabrasion Machine**




**3 \*Diamond Pens 9 \*Diamond Heads 2 \*Fuse**




**5 Blackhead Remover Tips**




**1 \*PVC Hose**

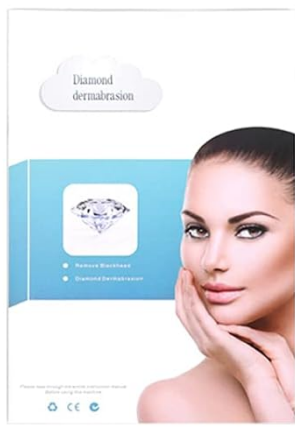

**1 \*English Manual**



**3 Pack Large Cotton Filter**



**3 Pack Small Cotton Filter**

Image: Overview of included components.

- 1 x Professional Dermabrasion Machine
- 3 x Diamond Pens
- 9 x Diamond Heads (various sizes and grits)
- 5 x Blackhead Remover Tips
- 1 x PVC Hose
- 3 x Pack Large Cotton Filters

- 3 x Pack Small Cotton Filters
- 2 x Fuses
- 1 x English User Manual

## 4. PRODUCT COMPONENTS

Familiarize yourself with the main components of the microdermabrasion machine:

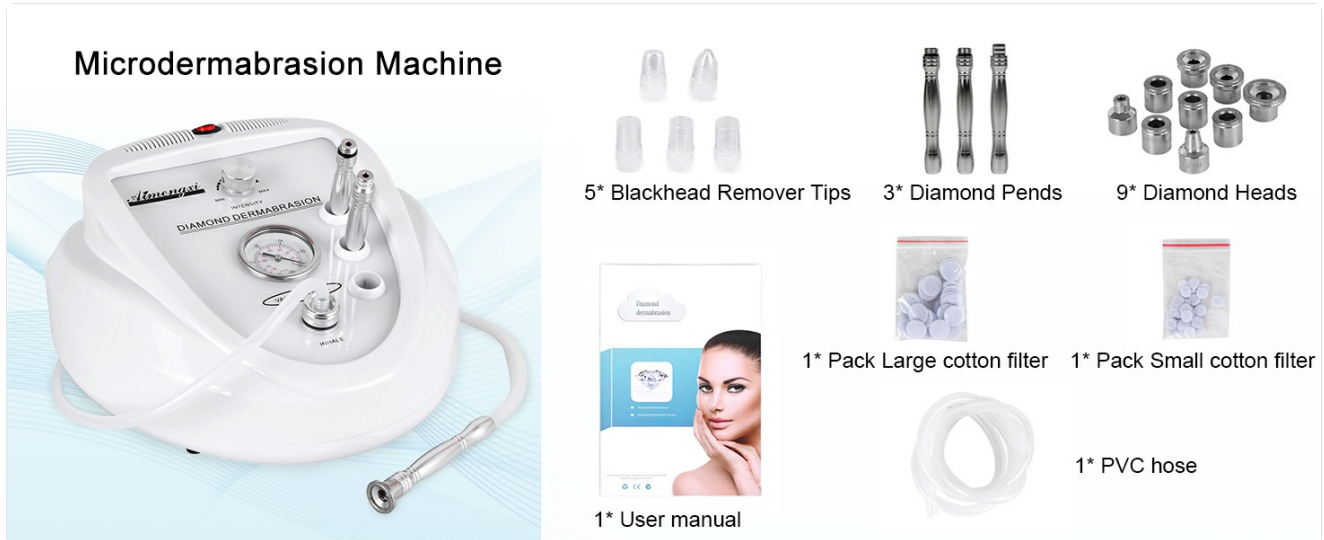


Image: Main unit with labeled controls.

### Diamond Tips Overview

The device includes various diamond tips designed for different skin areas and treatment intensities:

# Diamond Microdermabrasion Tips



D-100, D-150, D-200: Ideal for treating rough skin areas such as elbows, knees, and the back.

D-200, D-200: Suitable for use around the eyes and sides of the nose.

D-100, D-120, D-150, D-200: Perfect for facial use to improve skin texture and appearance.

Image: Assortment of diamond microdermabrasion tips.

- **D-100, D-150, D-200 (Larger Tips):** Ideal for treating rough skin areas such as elbows, knees, and the back.
- **D-200 (Smaller, Pointed Tips):** Suitable for use around the eyes and sides of the nose.
- **D-100, D-120, D-150, D-200 (Medium Tips):** Perfect for general facial use to improve skin texture and appearance.

## 5. SETUP INSTRUCTIONS

1. Unpack all components and ensure they are intact.
2. Connect the PVC hose to the designated port on the machine.
3. Attach the desired diamond pen and tip to the end of the PVC hose.
4. Ensure a cotton filter is properly inserted into the handpiece to prevent debris from entering the machine.
5. Plug the power cord into a suitable electrical outlet.

## 6. OPERATING INSTRUCTIONS

## General Operation

1. After connecting all accessories, press the power button to turn on the machine.
2. Turn the vacuum regulator knob to adjust the suction pressure and intensity to your desired level. Start with a lower setting and gradually increase as needed.
3. Move the wand gently and slowly across the skin surface.
4. After use, press the power button to stop the machine and unplug it.

## Microdermabrasion Treatment (Diamond Head)

The diamond microdermabrasion function exfoliates the skin, removing dead skin cells and promoting a smoother complexion.

### Microdermabrasion Treatment(Diamond Head)



Image: Correct application method for microdermabrasion.

- Select an appropriate diamond tip for the treatment area (e.g., finer tips for sensitive areas, coarser tips for rougher skin).
- Ensure the skin is clean and dry.
- Hold the skin taut with one hand.

- With the other hand, glide the diamond wand slowly and evenly across the skin in short, smooth strokes. Do not pull the skin or hold the wand stationary in one spot, as this can cause redness or bruising.
- Repeat strokes over the desired area, ensuring even coverage.
- After treatment, apply a soothing moisturizer or serum.

## Blackhead Removal

The vacuum function helps to extract blackheads and impurities from pores.

- Attach one of the blackhead remover tips to the handpiece.
- Ensure the skin is clean. Steaming the face beforehand can help open pores.
- Place the tip gently on the skin over a blackhead or congested area.
- Move the tip continuously across the skin; do not leave it in one spot for too long to avoid bruising.
- Adjust suction intensity as needed.
- After treatment, cleanse the skin and apply a toner or serum to help close pores.

## 7. MAINTENANCE

Regular cleaning and maintenance ensure the longevity and hygienic operation of your device.

- **Cleaning Tips:** After each use, remove the diamond tips and blackhead remover tips. Clean them thoroughly with a brush and warm soapy water, then rinse and allow to air dry completely.
- **Filter Replacement:** Regularly check and replace the cotton filters in the handpiece. This prevents debris from entering the machine and maintains optimal suction.
- **Machine Exterior:** Wipe the main unit with a soft, damp cloth. Do not use abrasive cleaners or immerse the unit in water.
- **Storage:** Store the machine and all accessories in a clean, dry place when not in use.

## 8. TROUBLESHOOTING

If you encounter issues with your device, refer to the following common solutions:

Problem	Possible Cause	Solution
No power	Power cord not connected, power switch off, or blown fuse.	Ensure power cord is securely plugged in. Turn on the power switch. Check and replace fuse if necessary.
Weak suction	Hose not properly connected, filter clogged, or vacuum regulator set too low.	Check hose connections. Replace cotton filter. Adjust vacuum regulator to a higher setting.
Device making unusual noise	Internal component issue or obstruction.	Turn off the device immediately. Contact customer support for assistance.

If problems persist after attempting these solutions, please contact customer support.

## 9. SPECIFICATIONS

- **Model Number:** AIMENGXI002
- **Brand:** AIMENGXI
- **Item Form:** Pen (for applicators)

- **Product Benefits:** Exfoliating
- **Skin Type Recommendation:** Sensitive (use with caution and lower settings)
- **Package Dimensions:** 13.5 x 13.3 x 7.5 inches
- **Item Weight:** 6.1 Pounds

## 10. WARRANTY AND SUPPORT

---

Your AIMENGXI 2-in-1 Diamond Microdermabrasion and Blackhead Removal Machine comes with a comprehensive warranty and dedicated customer support.

- **Warranty Policy:** We provide a 1-year warranty from the date of purchase.
- **Return Service:** A 30-day return service is available if you are not satisfied with the equipment.
- **Customer Support:** For any dissatisfaction, technical issues, or warranty claims, please contact AIMENGXI customer service. Refer to your purchase documentation for contact details.