



[Manuals.plus](#) /

> [AGARO](#) /

> AGARO MAGMA Air Compression Leg Massager (Model 33511) Instruction Manual

AGARO 33511

AGARO MAGMA Air Compression Leg Massager Instruction Manual

Model: 33511

[Overview](#) [Setup](#) [Introduction](#) [Safety Information](#) [Package Contents](#) [Product](#)
[Operation](#) [Maintenance](#) [Troubleshooting](#) [Specifications](#) [Warranty & Support](#)

1. INTRODUCTION

This manual provides instructions for the safe and effective use of your AGARO MAGMA Air Compression Leg Massager, Model 33511. This device is designed to provide air compression massage to the legs, feet, calves, and thighs, aiming to assist in muscle relaxation and blood circulation.

Sequential Air Compression Massage

Relieves from **fatigue**,
improve **blood circulation**

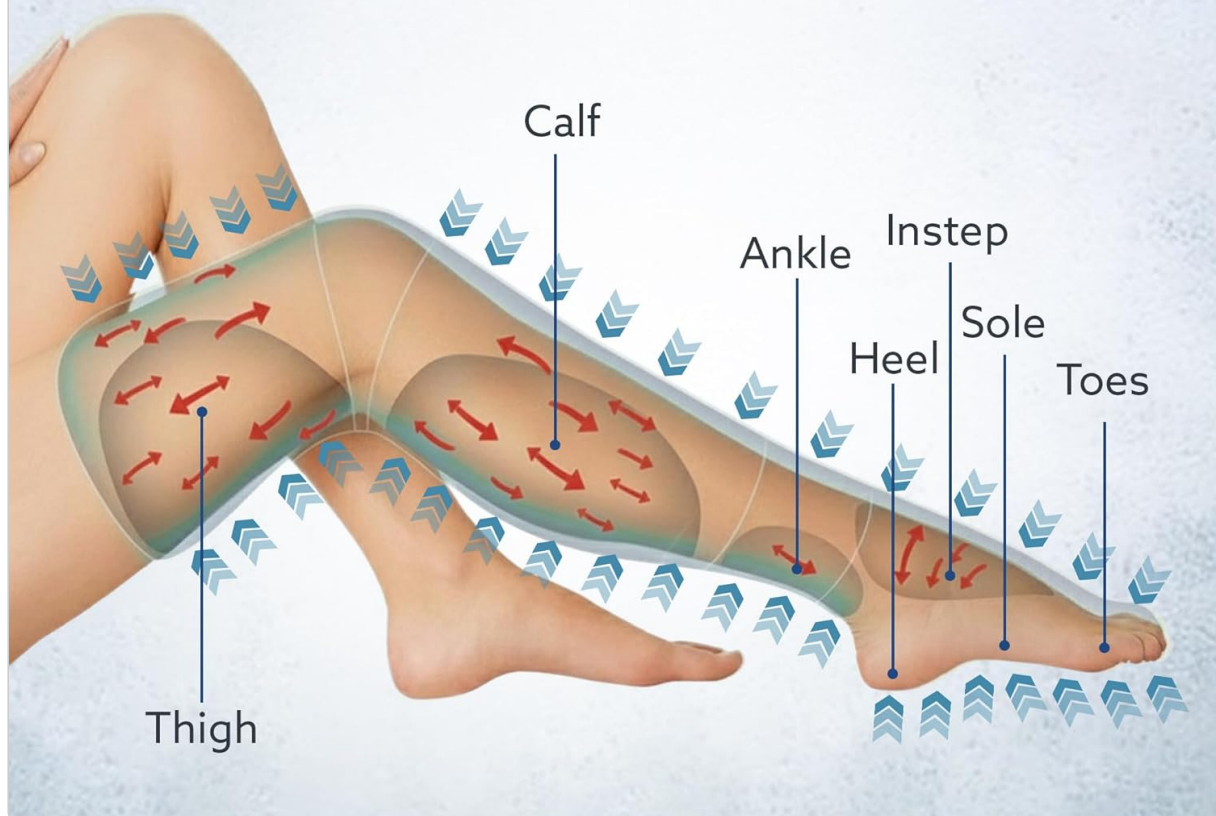


Image: The AGARO MAGMA Air Compression Leg Massager, Model Magma, showing the leg wraps and handheld controller.

2. SAFETY INFORMATION

Please read all safety warnings and instructions before using this product. Failure to follow these instructions may result in electric shock, fire, or serious injury.

- Do not use this product if you have any medical conditions that may be aggravated by massage, such as severe osteoporosis, malignant tumors, heart disease, or if you are pregnant. Consult a physician before use.
- Keep the device away from water and high humidity.
- Do not use the device if the power cord or plug is damaged.
- Unplug the device from the power outlet after use and before cleaning.
- Do not attempt to repair or modify the device yourself. Contact customer support for service.
- Keep out of reach of children.
- The device is equipped with a 20-minute automatic shut-off feature to prevent overuse. Do not override this safety feature.

3. PACKAGE CONTENTS

Verify that all components are present in the package:

- 2 x Massage Sleeves (Leg Wraps)
- 1 x Handheld Controller
- 1 x Power Adapter
- 1 x User Manual (this document)

4. PRODUCT OVERVIEW

The AGARO MAGMA Air Compression Leg Massager consists of adjustable leg wraps and a handheld controller. The leg wraps contain multiple airbags that inflate and deflate to provide compression massage.

4.1 Adjustable Straps and Dimensions



Image: The leg massager features adjustable straps to fit various body types. Maximum circumference for thigh is 30 inches, calf is 24 inches, and foot is 18 inches. Thigh length is 9.06 inches, calf length is 15.9 inches, and foot length is 9.84 inches.

The leg wraps are designed with adjustable Velcro belts to accommodate different leg sizes. The internal structure includes 11 x 2 airbags for comprehensive compression.

Adjustable straps for all body types



Max circumference:



Image: Illustration showing the full leg compression design with 11 x 2 internal airbags for effective massage.

4.2 Handheld Controller

Full leg Compression



11 x 2 airbags



Image: Detailed view of the handheld controller with labeled buttons and indicators.

The handheld controller allows you to manage the device's functions. It features:

- **Power Button:** Press and hold for 2 seconds to turn ON/OFF.
- **Power Indicator:** Green when ON, Red when OFF.
- **Mode Button (M):** Press to select massage modes (Sequence, Circulation, Whole).
- **Mode Indicators:** Lights for Sequence (M1), Circulation (M2), Whole (M3).
- **Air Pressure Button:** Press to change air pressure intensity (Max, Mid, Min).
- **Air Pressure Indicators:** Lights for Maximum, Medium, Minimum intensity.

5. SETUP

1. **Unpack:** Carefully remove all components from the packaging.
2. **Connect Leg Wraps:** Connect the air tubes from the leg wraps to the corresponding ports on the handheld controller. Ensure a secure connection.
3. **Connect Power Adapter:** Plug the power adapter into the handheld controller, then plug the adapter into a suitable power outlet.

6. OPERATING INSTRUCTIONS

6.1 Wearing the Massager

Follow these steps to correctly wear the leg massager:



Image: A visual guide demonstrating the four steps to correctly put on the leg massager wraps.

1. Sit comfortably. Identify the left and right leg wraps.
2. Wrap the foot section around your foot, securing with Velcro.
3. Wrap the calf section around your calf, securing with Velcro.
4. Wrap the thigh section around your thigh, securing with Velcro. Ensure the wraps are snug but not overly tight.

6.2 Power On/Off

- To turn **ON**: Press and hold the Power Button on the handheld controller for 2 seconds. The Power Indicator will turn green.
- To turn **OFF**: Press and hold the Power Button for 2 seconds. The Power Indicator will turn red.
- **Automatic Shut-off**: The device will automatically turn off after 20 minutes of continuous operation to prevent overuse.

Image: Text indicating the leg massager's 20-minute automatic shut-off feature for safety.

6.3 Selecting Massage Modes and Intensity

The device offers 3 massage modes and 3 intensity levels:



Image: Visual representation of the three massage modes (Sequence, Circulation, Whole) and three intensity levels (Low, Medium, High) available on the controller.

- **Massage Modes:**

- **Sequence (M1):** Airbags inflate and deflate in a sequential pattern from foot to thigh.
- **Circulation (M2):** Airbags inflate and deflate in a circulating pattern.
- **Whole (M3):** All airbags inflate and deflate simultaneously.

To change mode, press the **Mode Button (M)**. The corresponding indicator light (M1, M2, or M3) will illuminate.

- **Intensity Levels:** Low, Medium, High. To adjust intensity, press the **Air Pressure Button**. The corresponding indicator light (Min, Mid, or Max) will illuminate.

How to wear



Image: Diagram illustrating the flow of sequential air compression massage across the thigh, calf, ankle, instep, sole, and toes.

7. MAINTENANCE AND CARE

- **Cleaning:** Disconnect the device from power. Wipe the surface of the leg wraps and controller with a soft, damp cloth. Do not use abrasive cleaners or immerse in water.
- **Storage:** Store the device in a cool, dry place, away from direct sunlight and extreme temperatures.
- **Handling:** Avoid sharp objects that could puncture the airbags.

8. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Power adapter not connected properly or no power from outlet.	Ensure power adapter is securely plugged into the controller and a working power outlet.

Problem	Possible Cause	Solution
No compression or weak compression.	Air tubes are loose or disconnected; intensity level is too low.	Check that all air tubes are securely connected to the controller. Increase the air pressure intensity.
Device stops unexpectedly.	Automatic 20-minute shut-off activated.	This is a normal safety feature. You can restart the device if desired.

If you encounter issues not listed here, please contact customer support.

9. SPECIFICATIONS

Model Name	Magma
Item Model Number	33511
Product Dimensions (L x W x H)	26.5L x 56W x 73H Centimeters
Item Weight	1.55 kg (1550 Grams)
Power Source	Corded Electric
Material	Cotton
Massage Therapy Technique	Compression
Automatic Shut-off	20 Minutes
Manufacturer	AGARO, Universal Corporation Limited

10. WARRANTY AND SUPPORT

The AGARO MAGMA Air Compression Leg Massager comes with a **1 Year Warranty**. This warranty covers manufacturing defects of the product during the warranty period.

For any queries or concerns, please contact AGARO customer support:

- **Email:** info@uclindia.com
- **Customer Support:** 9821486487
- **Timings:** 10:00 AM to 6:00 PM (IST), Monday - Saturday (Excluding Holidays)

Manufacturer: Universal Corporation Limited Longgang District, Shenzhen, China.

Importer: Universal Corporation Ltd. 4/1 Middleton Street, Sikkim Commerce House, Kolkata - 700071.