

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

Manuals.plus /

> [Wutsar](#) /

> Wutsar Massager Gun Q5 Instruction Manual

Wutsar Q5

Wutsar Massager Gun Q5 Instruction Manual

Model: Q5 | Brand: Wutsar

INTRODUCTION

Thank you for choosing the Wutsar Massager Gun Q5. This percussion massage gun is designed to provide deep tissue muscle relief, enhance blood circulation, and aid in muscle recovery. Please read this manual thoroughly before use to ensure safe and effective operation.

SAFETY INFORMATION

WARNING: To reduce the risk of electric shock, fire, and personal injury, this product must be used in accordance with the following instructions:

- Do not use on injured or inflamed areas.
- Do not use for more than 15 minutes at a time.
- Keep away from water.
- Always unplug the device when not in use.
- Do not use while sleeping or in a state of reduced consciousness.
- Consult a healthcare professional before use if you have any health conditions or are pregnant.
- Always follow the instructions provided with the device.

PACKAGE CONTENTS

- 1 x Wutsar Massager Gun Q5
- 6 x Interchangeable Massage Heads
- 1 x USB-C Charging Cable
- 1 x Instruction Manual

PRODUCT OVERVIEW

6 Massage Heads Suitable For Different Body Parts

1  **Spherical head**
suitable for chest

2  **Shovel Head**
suitable for bulk muscle groups

3  **Fork Head**
suitable for achilles

4  **Flat Head**
suitable for all parts

5  **Bullet Head**
suitable for the deep tissue

6  **Taper Head**
suitable for tigger points

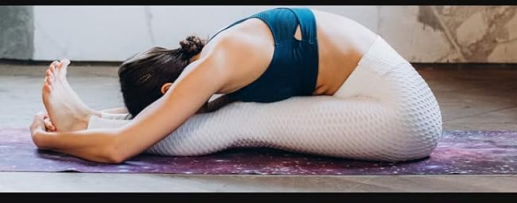
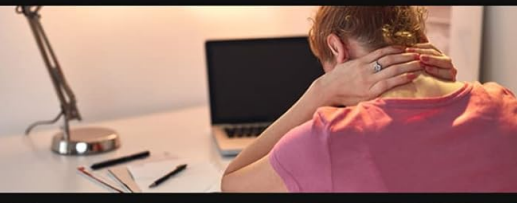


Image: Wutsar Massager Gun Q5 displaying key product information including model, dimensions, rated voltage, current, power, charging input, and color.

The Wutsar Massager Gun Q5 features a high-quality motor with quiet glide technology, an amplitude of 12mm for deep tissue penetration, and a touch LCD screen for intuitive control. It offers 8 adjustable speeds and comes with 6 specialized massage heads to target various muscle groups.

Touch Control



Image: A comparison highlighting the Wutsar Q5 Massager's features such as Liquid Crystal Display, Max 3200RPM, Noise \leq 35dB, Strong Massage Strength, and 6 Hours Battery Life.

SETUP

Charging the Device

- Connect the provided USB-C cable to the charging port on the massager and to a compatible USB power adapter (not included).
- The battery indicator on the LCD screen will show charging status.
- A full charge typically provides 4-6 hours of use.

Attaching Massage Heads

- Ensure the massager is turned off before attaching or removing any head.
- Gently push the desired massage head into the opening at the front of the gun until it is securely in place.
- To remove, simply pull the head straight out.

OPERATING INSTRUCTIONS

High-Torque Motor

deep tissue & fast cooling



Image: Close-up of the Wutsar Massager Gun Q5's touch control panel, illustrating the 8 gears for pressure monitoring and speed adjustment, with corresponding massage intensity levels.

Power On/Off

- To turn on: Press and hold the power button on the touch LCD screen for a few seconds until the device activates.
- To turn off: Press and hold the power button again for a few seconds.

Adjusting Speed Levels

- Once powered on, tap the power button to cycle through the 8 adjustable speed levels.
- The LCD screen will display the current speed setting (I-VIII).
- The device features an auto shut-off function after 10 minutes of continuous use to prevent overheating and potential muscle injury.

USAGE GUIDELINES: MASSAGE HEADS

WUTSAR Q5 Massager




Liquid Crystal Display



Max 3200RPM



Noise \leq 35dB

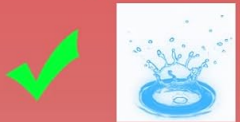

Strong Massage Strength



6 Hours Battery Life

Vs

Vs

Vs

Vs

Vs

Other Massager



No Liquid Crystal Display



Max 2400RPM



Noise $>$ 45dB



Weak Massage Strength



3 Hours Battery Life

Image: Illustration of the six different massage heads included with the Wutsar Massager Gun Q5, detailing their shapes and recommended body parts for use.

Select the appropriate massage head for your needs:

- **Spherical Head:** Suitable for large muscle groups like the chest, back, and glutes.
- **Shovel Head:** Ideal for bulk muscle groups, providing a broader massage area.
- **Fork Head:** Designed for the Achilles tendon, neck, and spine, targeting specific areas around bones.
- **Flat Head:** Versatile for all parts of the body, offering a general, comfortable massage.
- **Bullet Head:** Perfect for deep tissue massage and targeting specific knots or trigger points.
- **Taper Head:** Excellent for precise targeting of trigger points and small muscle areas.

Apply the massager to the desired muscle area, moving it slowly and gently. Avoid pressing too hard or staying in one spot for too long.

MAINTENANCE

- **Cleaning:** Wipe the device with a slightly damp cloth after each use. Do not use abrasive cleaners or immerse the device in water.

- **Storage:** Store the massager and its attachments in a cool, dry place, away from direct sunlight and extreme temperatures.
- **Battery Care:** To prolong battery life, avoid fully discharging the battery frequently. Charge the device regularly, even if not in constant use.

TROUBLESHOOTING

- **Device not turning on:** Ensure the device is fully charged. If the battery is low, connect it to the charger.
- **Weak percussion:** Check the battery level. If low, charge the device. Ensure the massage head is securely attached.
- **Unusual noise:** Turn off the device immediately. Check if any foreign objects are obstructing the massage head or internal mechanism. If the issue persists, contact customer support.

SPECIFICATIONS

Feature	Detail
Product Dimensions	7.09 x 6.5 x 2.36 inches
Item Weight	2 Pounds
Power Source	Battery Powered
Charging Input	USB-C
Amplitude	12mm
Speed Levels	8 Adjustable Speeds
Battery Life	4-6 hours (per charge)
Auto Shut-off	10 minutes

WARRANTY & SUPPORT

The Wutsar Massager Gun Q5 comes with a 3-year warranty. For any concerns regarding your order or product, please reach out via your Amazon account. Our team is committed to resolving issues promptly and ensuring your satisfaction.