

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [ZOETOUCH](#) /

› [ZOETOUCH Digital Bathroom Scale User Manual](#)

ZOETOUCH 560 lb Black (None-Bluetooth)

ZOETOUCH Digital Bathroom Scale User Manual

Model: 560 lb Black (Non-Bluetooth)

1. INTRODUCTION

Thank you for choosing the ZOETOUCH Digital Bathroom Scale. This high-capacity scale is designed to provide accurate and consistent body weight measurements for individuals up to 560 pounds (255 kg). Featuring an extra-wide platform, a large backlit LCD display, and advanced sensor technology, it offers a reliable and user-friendly weighing experience. Please read this manual thoroughly before use to ensure proper operation and maintenance.



The ZOETOUGH Digital Bathroom Scale with its clear LCD display.

2. SAFETY INFORMATION

- Do not use on carpeted or soft surfaces. Always place the scale on a hard, flat surface for accurate readings.
- Keep the scale dry. Avoid contact with water or other liquids.
- Do not jump on the scale or use it for purposes other than weighing.
- The scale's surface is made of tempered glass. Handle with care to prevent breakage.
- Keep batteries out of reach of children.
- If the scale is not used for a long period, remove the batteries to prevent leakage.

3. PACKAGE CONTENTS

Your ZOETOUGH Digital Bathroom Scale package includes the following items:

- 1 x ZOETOUGH Digital Bathroom Scale
- 4 x AAA Batteries
- 1 x Measuring Tape



Contents of the ZOETOUGH Digital Bathroom Scale package.

4. SETUP

4.1 Battery Installation

1. Locate the battery compartment on the underside of the scale.

2. Open the battery compartment cover.
3. Insert the 4 AAA batteries, ensuring correct polarity (+/-) as indicated inside the compartment.
4. Close the battery compartment cover securely.

4.2 Placement

For the most accurate readings, it is crucial to place the scale on a hard, flat, and stable surface. Avoid placing it on carpets, rugs, or uneven flooring, as this can affect accuracy.

4.3 Initial Calibration

Before your first use, or if the scale has been moved, it is recommended to calibrate it:

1. Place the scale on a hard, flat surface.
2. Lightly step on the scale to turn it on, then immediately step off.
3. The display will show '0.0' or 'CAL' and then turn off. This indicates successful calibration.
4. The scale is now ready for use.

5. OPERATING INSTRUCTIONS

5.1 Taking a Measurement

1. Ensure the scale is on a hard, flat surface.
2. Gently step onto the center of the scale with bare, dry feet.
3. Stand still until your weight is displayed and locked on the LCD screen.
4. The scale will automatically turn off after approximately 15 seconds of inactivity to conserve battery life.



The extra-wide platform provides comfort during measurement.

5.2 Unit Conversion

The scale supports three units of measurement: pounds (lb), kilograms (kg), and stones (st). To change the unit:



- Locate the 'UNIT' button on the back of the scale.
- Press the button to cycle through the available units (lb, kg, st). The selected unit will be displayed on the screen.

5.3 Weight Locking Feature

Once a stable weight reading is obtained, the scale will lock the displayed weight for approximately 15 seconds. This provides ample time to read and record your measurement before the scale automatically turns off.

5.4 Comparison Function

The scale can automatically calculate and display the weight difference compared to your last measurement. After a new measurement, the display will show '+' or '-' icons along with the numerical difference, indicating whether your current weight is higher or lower than the previous one.

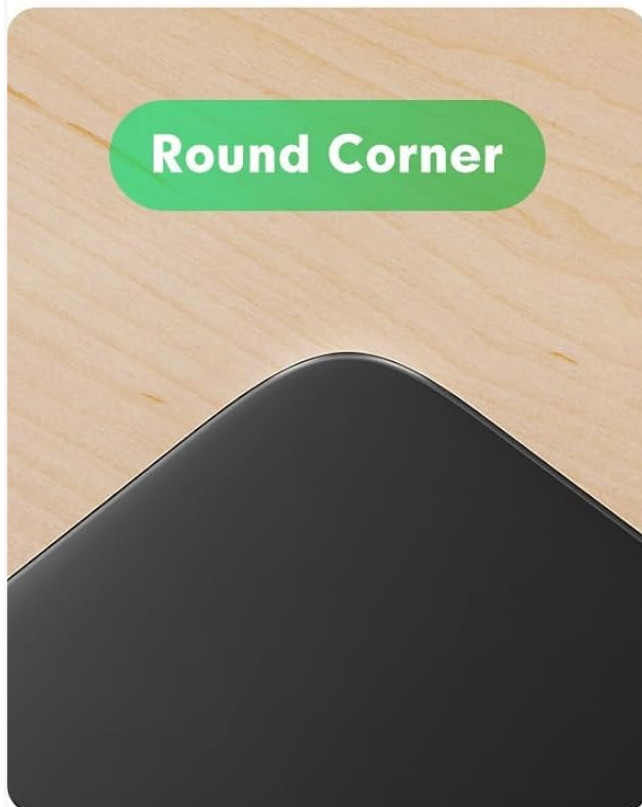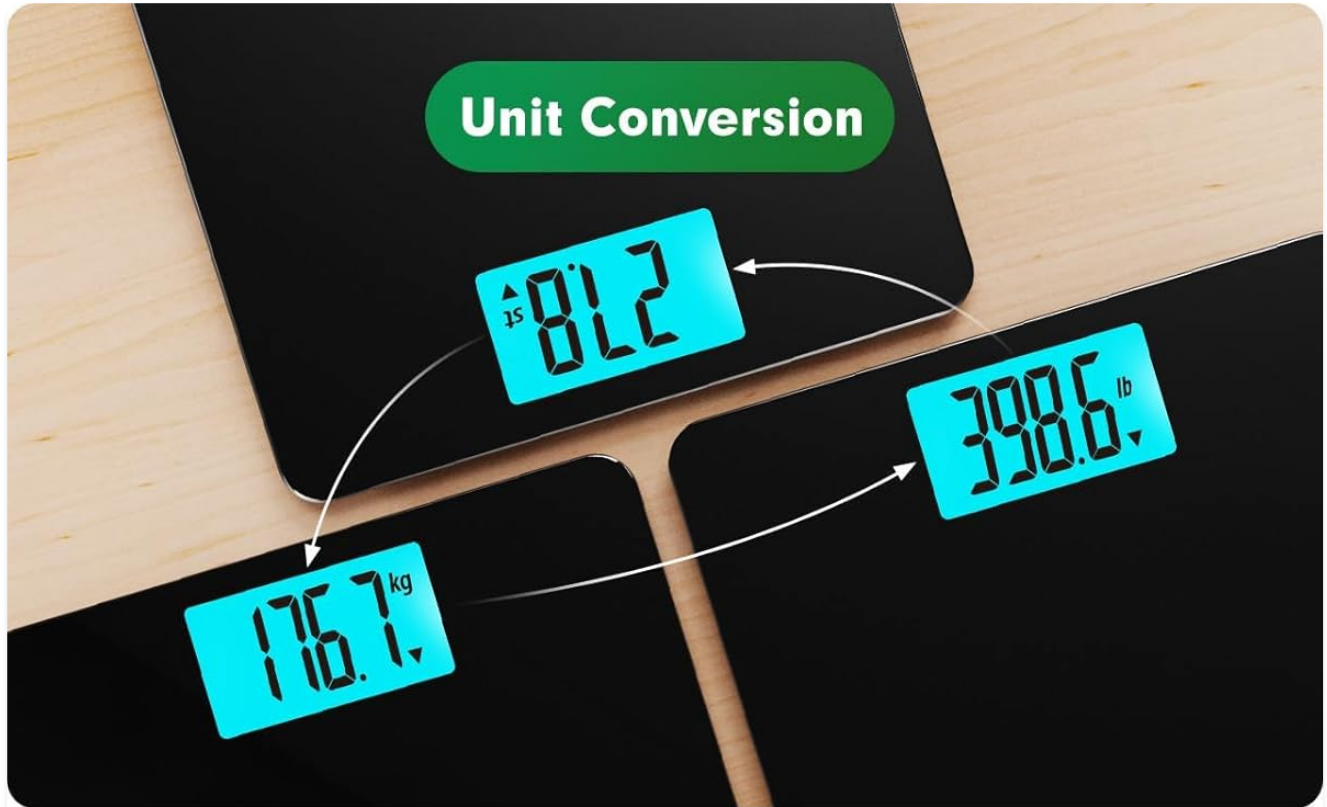

The scale features easy unit conversion and a large, clear backlit display.

6. MAINTENANCE

6.1 Cleaning

- Wipe the scale surface with a soft, damp cloth.
- Do not use abrasive cleaners or immerse the scale in water.
- Ensure the scale is completely dry before storage or next use.

6.2 Storage

- Store the scale in a cool, dry place.
- Avoid storing it vertically or with heavy objects on top, as this may affect sensor accuracy over time.

7. TROUBLESHOOTING

If you encounter any issues with your scale, please refer to the following table:

Problem	Possible Cause	Solution
Inaccurate readings	Scale on soft/uneven surface; Not calibrated; Batteries low.	Place on hard, flat surface. Perform initial calibration (Section 4.3). Replace batteries.
Display shows 'Lo'	Low battery indicator.	Replace all 4 AAA batteries.
No display/Scale won't turn on	Batteries incorrectly installed or dead.	Check battery polarity. Replace batteries.
Display shows 'Err'	Overload or internal error.	Ensure weight is within capacity (560 lbs). Step off and re-calibrate. If issue persists, contact support.

8. SPECIFICATIONS

Feature	Detail
Brand	ZOETOUCH
Model	560 lb Black (Non-Bluetooth)
Capacity	Up to 560 lbs / 255 kg

Feature	Detail
Accuracy	0.2 lb / 0.1 kg increments
Display Type	Large Backlit LCD
Platform Material	8mm Tempered Glass
Platform Dimensions	12.2 x 12.9 inches (30.98 x 32.8 cm)
Power Source	4 x AAA Batteries (Included)
Special Features	Step-on Technology, Auto Shut-off, Unit Conversion (lb/kg/st), Weight Locking, Comparison Function, High Precision Sensors, Round Corner Design

High Precision Sensors

LCD

LCD Display




Step-on

Lo


Low Battery
Indicator



The scale utilizes high precision sensors for accurate measurements.

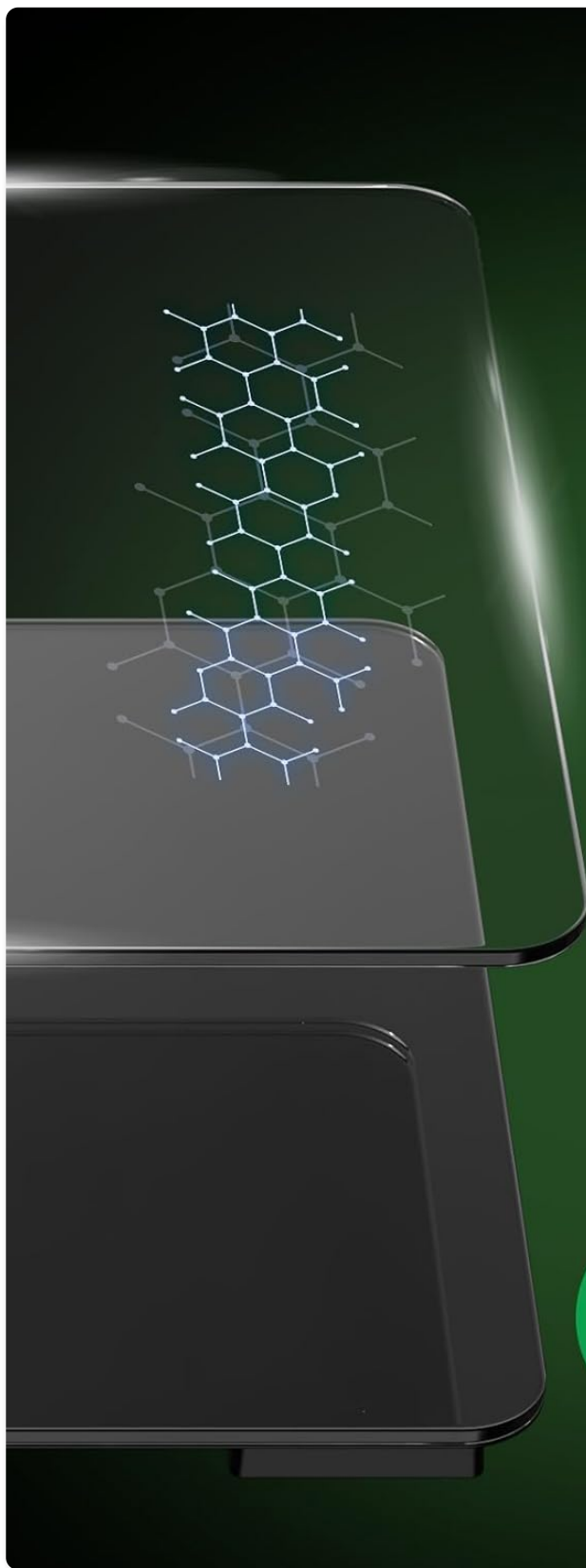

60% Much Thicker



Brand	Thickness (MM)
ZOETOUCH	8
Others	6

**Durable and Sturdy
Weighing Scale**

Constructed with durable 8mm tempered glass for enhanced sturdiness.

9. WARRANTY AND SUPPORT

ZOETOUCH products are manufactured with high-quality materials and undergo strict quality control. For any questions, concerns, or technical support regarding your scale, please refer to the contact information provided with your purchase or visit the official ZOETOUCH website. Please retain your proof of purchase for warranty claims.

