

ANLAN ANLAN

ANLAN Multifunctional Facial Massager User Manual

Model: ANLAN

Brand: ANLAN

INTRODUCTION

This manual provides comprehensive instructions for the safe and effective use of your ANLAN Multifunctional Facial Massager. This advanced beauty instrument combines hot and cold therapy, EMS microcurrent, ion introduction, vibration, and red/blue light phototherapy to address various skin concerns, including anti-aging, wrinkle reduction, deep cleansing, and overall facial care. Regular use can contribute to healthier, more radiant skin and help delay the signs of aging.

IMPORTANT SAFETY INFORMATION

Please read all instructions carefully before using the device. Failure to follow these guidelines may result in injury or damage to the device.

- Do not use on broken, irritated, or infected skin.
- Consult a physician before use if you have any pre-existing medical conditions, especially heart conditions, pacemakers, epilepsy, or are pregnant.
- Keep out of reach of children.
- Do not immerse the device in water or use it in excessively humid environments.
- Use only the provided charging cable and ensure it is properly connected.
- Discontinue use immediately if you experience any discomfort, pain, or adverse skin reactions.
- Do not disassemble or attempt to repair the device yourself.

PRODUCT OVERVIEW

Familiarize yourself with the components and features of your ANLAN Multifunctional Facial Massager.



Image: Detailed diagram of the ANLAN facial massager, highlighting its various parts such as the massage head, red/blue LED light, battery level indicator, 5-in-1 mode display, mode time, power/mode button, level button, LED button, air hole for heat dissipation, metal sensing area, and charging port.

- ✓ 5 Skincare Features
- ✓ 3 Stepped Temperatures
- ✓ 3 Level EMS Strength
- ✓ 2 Kinds Of Optical Cosmetology
- ✓ 10000(+10%)rpm Vibration Frequency



Image: The device's display screen, illustrating the various modes and settings available, including temperature levels, EMS strength, and vibration frequency.

Key Components:

- **Massage Head:** The primary contact surface with the skin, integrating hot/cold therapy and light functions.
- **Red/Blue LED Light:** Indicators for phototherapy modes, visible on the massage head.
- **Battery Level Indicator:** Displays remaining battery life on the screen.
- **Display Screen:** Shows current mode, intensity level, and remaining time.
- **Mode Button:** Cycles through the 5 operational modes.
- **Level Button:** Adjusts intensity levels for EMS and vibration.
- **LED Button:** Toggles red/blue light therapy on or off.
- **Metal Sensing Area:** Located on the sides of the device. Must be held during use for proper function of ion and EMS modes.
- **Charging Port:** Located at the bottom of the device for power connection.
- **Charging Base:** Provides a stable stand and connection for charging.

PRODUCT SPECIFICATIONS

Feature	Detail
Power Source	Battery (1 Lithium Ion, included)
Dimensions (L x W x H)	17.5 x 5.2 x 5 cm
Weight	380 g
Model Name	ANLAN
Skin Type Compatibility	Mixed
Functions	Hot/Cold Therapy, EMS Microcurrent, Ion Introduction, Vibration, Red/Blue LED Light

SETUP

Charging the Device

Before first use, ensure the device is fully charged. A full charge is recommended for optimal performance.

1. Connect the charging cable to the charging port located at the bottom of the device.
2. Place the device securely on its charging base.
3. Plug the USB end of the cable into a suitable USB power adapter (not included) or a computer USB port.
4. The battery indicator on the device's display will illuminate and show the charging status. Charging is complete when the indicator shows a full battery.



Image: The ANLAN facial massager placed on its charging base, illustrating the charging process.

OPERATING INSTRUCTIONS

The ANLAN Multifunctional Facial Massager features 5 distinct modes, each designed for specific skincare benefits. For optimal performance, ensure your hands are dry and maintain contact with the metal sensing area on the sides of

the device during use, especially for EMS and ion modes.

5 in 1 beauty device



Image: An overview of the device's five primary functions: Clean, Cool, EMS Lifting, Moisturizing, and Col Synthesis.

General Usage Tips:

- Cleanse your face thoroughly before each use to remove makeup and impurities.
- Apply your preferred skincare product (serum, essence, cream, or mask) before using the device, as recommended for each specific mode.
- Move the device gently across your skin in upward and outward motions. Avoid pressing too hard.
- Each mode has a recommended duration, typically displayed on the screen. The device may automatically turn off or prompt you to switch modes upon completion.
- Use the device for 10-20 minutes per day for effective results.

Mode Descriptions:

1. Clean Mode (5 minutes)

This mode utilizes ion export and vibration to deeply cleanse pores, effectively removing dirt, oil, and makeup residues. It prepares your skin for better absorption of subsequent skincare products.

Recommended Use: Use with a cleansing lotion or makeup remover. For enhanced cleaning, you may attach a cotton pad to the massage head.



Image: Graphic illustrating the deep cleansing action of ion export, pulling impurities from the skin.

2. Cool Mode (2 minutes)

This mode combines cryotherapy (cooling, 6-12°C) with vibration and red/blue LED light. It helps to calm and soothe the skin, reduce puffiness, shrink pores, and lock in moisture and nutrients after product application. It also aids in reducing acne and inflammation.

Recommended Use: After applying serums or moisturizers, or to soothe irritated skin. Ideal for use with sheet masks.

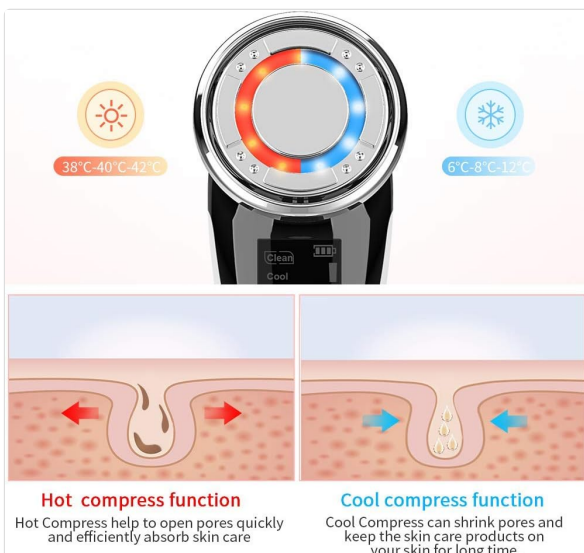


Image: Diagram illustrating the effect of cool compress on skin, helping to shrink pores and retain skincare products.

3. EMS Lifting Mode (5 minutes)

Utilizes low-frequency microcurrent technology (EMS) to stimulate facial muscles, promoting skin

firmness, improving elasticity, and reducing the appearance of fine lines and wrinkles. It helps to lift and tone the facial contours.

Recommended Use: With an anti-wrinkle cream or a conductive gel. Adjust intensity using the "Level" button to your comfort.



Image: Graphic illustrating the microcurrent stimulation on skin, showing its penetration and effect on muscle tissue.

4. Moisturizing Mode (5 minutes)

Combines ion introduction, vibration, and LED light to significantly enhance the absorption of skincare products. This mode helps active ingredients penetrate deeper into the skin, maximizing their effectiveness for hydration and rejuvenation.

Recommended Use: With serums, essences, or sheet masks. This mode can increase absorption by up to 300%.



Image: Graphic illustrating enhanced nutrient absorption through ion introduction, showing deeper penetration of skincare products.

5. CoISynthesis Mode (2 minutes)

Integrates hot compress (38-42°C), red/blue light, and vibration to stimulate collagen synthesis, firm the skin, and reduce wrinkles. The warmth helps to open pores, allowing for better product penetration and promoting skin regeneration.

Recommended Use: With anti-aging creams or collagen-boosting serums.

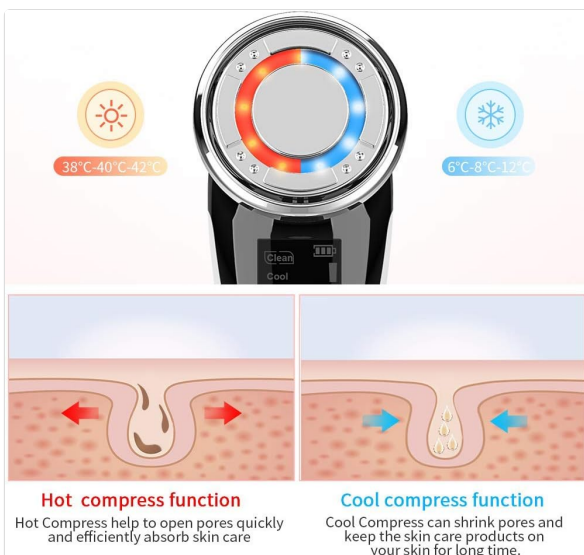


Image: Diagram illustrating the effect of hot compress on skin, helping to open pores quickly for efficient absorption of skincare products.

LED Light Therapy:

The device features integrated red and blue LED light therapy, which can be activated in various modes or independently.

- **Red Light:** Promotes collagen production, improves skin elasticity, and helps reduce the appearance of fine lines and wrinkles.
- **Blue Light:** Inhibits inflammation, reduces acne-causing bacteria, and helps balance oil production, leading to clearer skin.

Press the "LED" button to toggle between red light, blue light, or turn the light off.



Image: Close-up images of the device head displaying red and blue light therapy in action.

MAINTENANCE

Cleaning

To ensure hygiene and prolong the lifespan of your device:

1. After each use, gently wipe the massage head with a soft, damp cloth.
2. For deeper cleaning, you may use a mild, non-abrasive cleanser.
3. Do not use alcohol, harsh chemicals, abrasive materials, or immerse the device in water.
4. Ensure the device is completely dry before storing it.

Storage

Store the device in a cool, dry place, away from direct sunlight and extreme temperatures. Keep it in its original

packaging or on its charging base when not in use to protect it from dust and damage.

TROUBLESHOOTING

If you encounter any issues with your ANLAN Multifunctional Facial Massager, please refer to the table below for common problems and their solutions.

Problem	Possible Cause	Solution
Device does not turn on.	Low battery; Device not charged.	Charge the device fully using the provided cable and charging base.
No sensation from EMS/Ion modes.	Fingers not holding metal sensing area; Skin too dry; Low intensity level.	Ensure fingers are in full contact with the metal sensing area; Apply more skincare product to ensure conductivity; Increase intensity level using the "Level" button.
Hot/Cool function not working or not reaching desired temperature.	Not enough time for temperature change; Device malfunction.	Allow a few moments (approx. 30-60 seconds) for the temperature to adjust after selecting the mode; If the problem persists, contact customer support.
Device turns off unexpectedly.	Low battery; Mode cycle completed.	Recharge the device if the battery is low; This is normal behavior when a mode's timed cycle finishes.


WARRANTY AND CUSTOMER SUPPORT





For detailed warranty information, product registration, or to obtain customer support, please refer to the documentation included with your product packaging. You may also visit the official ANLAN website for further assistance.

Manufacturer: ANLAN

ASIN: B093C2S2FK

Related Documents - ANLAN

	<p>ANLAN Neck Beauty Device User Manual</p> <p>User manual for the ANLAN Neck Beauty Device, detailing its features, functions, usage instructions, and warranty information. The device utilizes positive and negative ions, vibration, warm massage, and red, blue, and green light therapy for skin cleansing, EMS care, and moisturizing.</p>
---	---

<div><div>ANLAN</div><div>NECK BEAUTY DEVICE</div><div>USER MANUAL</div><div></div><div><p>⚠ This product is not a medical device used for diagnosis or treatment, but a general household item. Please read this user manual carefully before use. Please keep this user manual after receiving the future services. Please attach this user manual when transferring this product to others.</p></div></div>	<p>ANLAN Neck Beauty Device User Manual: Features, Functions, and How to Use</p> <p>Comprehensive user manual for the ANLAN Neck Beauty Device (Model ES-081). Learn about its features, functions (EMS, warm, red/blue/green light), safety precautions, and how to use it for skin rejuvenation and care.</p>
<div><div>ANLAN</div><div>Ultrasonic Skin Scrubber</div><div>USER MANUAL</div><div></div><div><p>⚠ This product is not a medical device for diagnosis and treatment, but an auxiliary product for general household use. Thank you for purchasing ANLAN Ultrasonic Skin Scrubber. Before using this product, please read this instruction carefully. In addition, this user manual also served as a warranty card, so please keep it in a safe place for future use.</p></div></div>	<p>ANLAN Ultrasonic Skin Scrubber User Manual</p> <p>User manual for the ANLAN Ultrasonic Skin Scrubber (Model: JD-CP002), detailing its features, safety precautions, usage instructions, maintenance, and FAQs.</p>
<div><div>ANLAN</div><div>ANLAN Body Beauty Device</div><div>USER MANUAL</div><div></div><div><p>⚠ This product is not a medical device used for diagnosis or treatment, but a general household item. Please read this user manual carefully before use. Please keep this user manual after receiving the future services. Please attach this user manual when transferring this product to others.</p></div></div>	<p>ANLAN Body Beauty Device User Manual - Operation, Safety, and Specifications</p> <p>Comprehensive user manual for the ANLAN Body Beauty Device (Model OFY-2102), covering safety precautions, product details, usage instructions, troubleshooting, and warranty information.</p>
<div><div>ANLAN</div><div>ANLAN Multifunctional Skincare Device</div><div>USER MANUAL</div><div></div><div><p>⚠ This product is not a medical device used for diagnosis or treatment, but a general household item. Please read this user manual carefully before use. Please keep this user manual after receiving the future services. Please attach this user manual when transferring this product to others.</p></div></div>	<p>ANLAN Multifunctional Skincare Device User Manual</p> <p>User manual for the ANLAN Multifunctional Skincare Device (Model 01-ADRY14-10A), providing comprehensive instructions on usage, features, safety, and maintenance for effective home skincare.</p>