

## COMFIER CF-5913

# COMFIER Shiatsu Foot and Back Massager with Heat (Model CFUS-5913-CA) Instruction Manual

## INTRODUCTION

Thank you for choosing the COMFIER Shiatsu Foot and Back Massager with Heat. This device is designed to provide a soothing massage experience for your feet and back, incorporating shiatsu kneading and gentle warmth. Please read this manual thoroughly before use to ensure safe and effective operation.

## SAFETY INFORMATION

**WARNING: Always unplug the massager before cleaning or storing. Do not use the massager if the cord is damaged. Keep the massager dry and do not submerge it in water. Do not use the massager if you have any medical conditions, such as high blood pressure, heart disease, or diabetes, without consulting a doctor first. Keep the massager out of reach of children and pets. Always follow the instructions provided with the massager for safe and effective use.**

## PRODUCT FEATURES

- **2-in-1 Foot and Back Massager:** Versatile design allows for use on both feet and the lower back.
- **Flexible Nodes:** Features flexible massage nodes that adapt to different foot arches, providing a comfortable massage.
- **Gel-Node Technology:** Incorporates gel massage nodes for a gentler, professional massage experience.
- **Soothing Heat Function:** Provides optional gentle warmth to enhance relaxation during massage.
- **Toe-Touch Control:** Easily adjust settings with your foot, eliminating the need to bend over.
- **Detachable Washable Cover:** Includes a removable cover for easy cleaning and hygiene.




Figure 1: COMFIER Shiatsu Foot and Back Massager (Model CFUS-5913-CA)

## GEL AND FLEXIBLE MASSAGE NODES

Flexible massage nodes fit different foot archs while gel massage nodes deliver a gentler massage



**OTHERS**

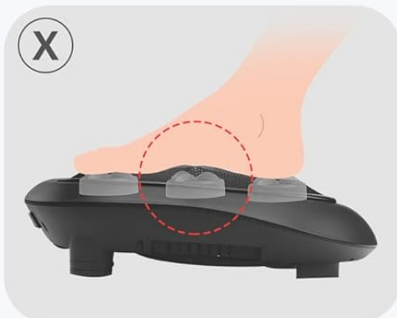

Without flexible nodes

**OTHERS**



Without gel nodes

**COMFIER**



With flexible & gel nodes

Figure 2: Illustration of flexible and gel massage nodes for varied foot arches.



## VERSATILE DESIGN

For one more option of back massage cushion



Figure 3: The massager's versatile design allows for use on feet, back, thighs, and calves.

## WHY YOU NEED A DETACHABLE DUST COVER

COMFIER



VS


OTHERS



Enhance comfort



Easy to clean

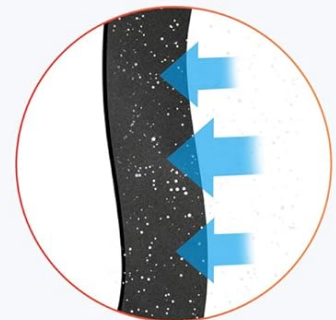

Avoid dust winding

Figure 4: The massager includes a detachable cover for enhanced comfort and easy cleaning.





Figure 5: Control buttons designed for easy toe-press operation without bending over.

Video 1: Overview of the COMFIER Shiatsu Foot Massager Machine, demonstrating its features and ease of use.

Video 2: Demonstration of the foot massager's functionality, including heat and massage modes.

Video 3: Detailed look at the massager's three massage modes, heating function, and suitability for various foot sizes.

## SETUP

1. **Unpack:** Carefully remove the massager and all accessories from the packaging.
2. **Placement:** Place the massager on a flat, stable surface. Ensure there is enough space around the unit for comfortable use.
3. **Power Connection:** Plug the power adapter into the massager's power port, then plug the adapter into a suitable electrical outlet.
4. **Optional Cover:** If desired, place the washable cover over the massage surface for added comfort or hygiene.

5. **Height Adjustment:** The massager includes extra anti-sliding pads. These can be used to adjust the height of the unit for optimal comfort during use.



Figure 6: The massager's height can be adjusted using the included anti-sliding pads.



Figure 7: Package contents include the massager, power adapter, and washable cover. Dimensions are approximately 13.4" x 14.5" x 3".

## OPERATING INSTRUCTIONS

1. **Power On/Off:** Press the Power Button ( ) to turn the massager on or off. The unit will typically start in a default massage mode.
2. **Select Massage Mode:** Press the Mode Button ( ) to cycle through the available massage modes. The massager offers 3 distinct modes.
3. **Adjust Intensity:** Press the Intensity Button ( ) to adjust the massage intensity. There are 3 intensity levels to choose from.
4. **Activate Heat:** Press the Heat Button ( ) to turn the soothing heat function on or off.
5. **Usage:** Place your feet on the massage nodes or position the massager against your lower back, thighs, or calves as desired.
6. **Automatic Shut-off:** For safety, the massager is equipped with a 15-minute automatic shut-off function.



## MAINTENANCE

---

- **Cleaning:** Always unplug the massager before cleaning. The detachable cover can be removed and machine washed. Wipe the main unit with a soft, damp cloth. Do not use abrasive cleaners or submerge the unit in water.
- **Storage:** Store the massager in a cool, dry place away from direct sunlight when not in use. Avoid placing heavy objects on the unit.

## TROUBLESHOOTING

---

- **Unit does not power on:**
  - Ensure the power adapter is securely plugged into both the massager and a working electrical outlet.
  - Check if the power outlet is functional by plugging in another device.
- **Heat function is not working:**
  - Verify that the Heat Button (🔥) has been pressed and the indicator light is on.
  - Allow a few minutes for the heat to warm up.
- **Massage feels too intense/weak:**
  - Adjust the intensity level using the Intensity Button (▶▶).
  - Ensure the washable cover is properly placed if you desire a gentler massage.

## SPECIFICATIONS

---

- **Product Dimensions:** 13.58 x 4.02 x 14.96 inches
- **Item Weight:** 5.4 Pounds
- **Model Number:** CF-5913
- **Power Source:** Corded Electric
- **Material:** Polyester
- **Brand:** COMFIER
- **Manufacturer:** Comfier

## WARRANTY AND SUPPORT

---

COMFIER products are designed for durability and performance. For information regarding warranty coverage or customer support, please refer to the documentation included with your product or contact COMFIER customer service directly.

Protection plans are available for purchase separately through authorized retailers. These plans may offer extended coverage beyond the manufacturer's warranty.

### Related Documents - CF-5913

|  |   |
|--|---|
| <div>  <div> <div>COMFIER</div> <div>Shiatsu Neck and Back Massager</div> </div> <div> <div>User Manual</div> <div>CF-2309A</div> </div> </div>   | <p><a href="#">Comfier Shiatsu Neck and Back Massager CF-2309A User Manual</a></p> <p>User manual for the Comfier Shiatsu Neck and Back Massager, model CF-2309A. Includes safety instructions, features, technical data, setup and operation, controller instructions, care and cleaning, troubleshooting, and warranty information.</p>                       |
| <div>  <div> <div>COMFIER</div> <div>Leg and Foot Massager</div> </div> <div> <div>User Manual</div> <div>CF-F25H</div> </div> </div>   | <p><a href="#">Comfier CF-F25H Leg and Foot Massager User Manual</a></p> <p>Comprehensive user manual for the Comfier CF-F25H Leg and Foot Massager, detailing safety instructions, product features, setup and operation, maintenance, troubleshooting, warranty information, and FCC compliance.</p>  |
| <div>  <div> <div>COMFIER</div> <div>Shiatsu massage pillow with heat</div> </div> <div> <div>CE</div> <div>User Manual</div> <div>CF-1228</div> </div> </div>   | <p><a href="#">Comfier CF-1228 Shiatsu Massage Pillow with Heat User Manual</a></p> <p>User manual for the Comfier CF-1228 Shiatsu massage pillow with heat, detailing safety instructions, operating procedures, features, maintenance, disposal, technical specifications, and warranty information.</p>  |
| <div>  <div> <div>COMFIER</div> <div>Shiatsu Neck and Back Massager with Heat</div> </div> <div> <div>User Manual</div> <div>CF-2303-KR</div> </div> </div>   | <p><a href="#">COMFIER CF-2303-KR Shiatsu Neck and Back Massager with Heat User Manual</a></p> <p>Comprehensive user manual for the COMFIER CF-2303-KR Shiatsu Neck and Back Massager with Heat, covering safety instructions, features, technical data, setup, operation, controller functions, care, cleaning, troubleshooting, and warranty information.</p> |
| <div>  <div> <div>COMFIER</div> <div>Shiatsu massagercushion mit wärmefunktion DE<br/>Shiatsu massager cushion with heat GB<br/>Coussin de massage Shiatsu avec chaleur FR<br/>Cushion per massaggio Shiatsu riscaldato IT<br/>Cajón masajeador Shiatsu con calefacción ES</div> </div> <div> <div>CE</div> <div>User Manual</div> <div>CF-2606MC</div> </div> </div> | <p><a href="#">COMFIER CF-2606MC Shiatsu Massage Cushion with Heat User Manual</a></p> <p>Comprehensive user manual for the COMFIER CF-2606MC Shiatsu massage cushion with heat. Includes safety instructions, operating guide, features, maintenance, and warranty information.</p>  |

COMFIER

Shiatsu neck and back massager CE  
Shiatsu neck and back massager CE  
Massleur de cou et de dos Shiatsu FR  
Massaggiatore Shiatsu per collo e schiena IT  
Masajeador Shiatsu de pes e cuello y espalda ES  
Shiatsu neck - en rughe massageur NL



CE  
User Manual  
CF-2309A

## [Comfier CF-2309A Shiatsu Neck and Back Massager User Manual](#)

Comprehensive user manual for the Comfier CF-2309A Shiatsu Neck and Back Massager, covering features, operation, safety instructions, maintenance, technical specifications, and warranty information.