



[Manuals.plus](#) /

> [Sharper Image](#) /

> SHARPER IMAGE Shiatsu Full Body Multifunction Cordless Massager User Manual

## Sharper Image 1014508

# SHARPER IMAGE Shiatsu Full Body Multifunction Cordless Massager

Model: 1014508

## INTRODUCTION

---

This user manual provides detailed instructions for the safe and effective operation, setup, and maintenance of your Sharper Image Shiatsu Full Body Multifunction Cordless Massager. Please read this manual thoroughly before using the device and retain it for future reference.



Figure 1: Sharper Image Shiatsu Full Body Multifunction Cordless Massager in use and its packaging.

## IMPORTANT SAFETY INSTRUCTIONS

To reduce the risk of electric shock, burns, fire, or injury, always follow basic safety precautions when using electrical appliances.

- Do not immerse the massager or adapter in water or other liquids.
- Do not use the massager while sleeping or drowsy.
- Keep out of reach of children.
- Do not use on sensitive skin areas, inflamed areas, or open wounds.
- Consult a physician before use if you have any medical conditions, are pregnant, or have a pacemaker.
- Unplug the adapter from the power outlet when not in use and before cleaning.
- Do not operate if the cord or plug is damaged, or if the product malfunctions or has been dropped or damaged.

## PACKAGE CONTENTS


Please verify that all items are included in your package:

- 1 x Full Body Massager
- 1 x Set of Arm Straps
- 1 x DC Adapter


## SHARPER IMAGE®




### WHAT YOU GET

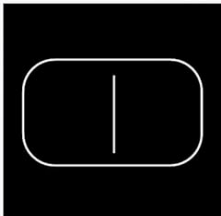

Massager

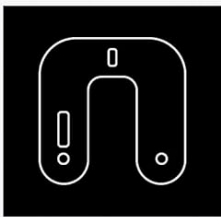

Arm strap

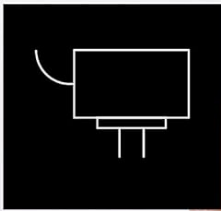

Adapter



Figure 2: Included components: Massager, Arm Strap, and DC Adapter.

## SETUP

Your Sharper Image Shiatsu Full Body Massager is designed for cordless operation once charged. Before first use, or if the battery is low, charge the device using the provided DC adapter.

1. **Charging the Massager:** Connect the DC adapter to the charging port on the massager. Plug the adapter into a standard electrical outlet. The charging indicator light will illuminate. Allow the massager to fully charge before its first use.
2. **Attaching Arm Straps (Optional):** The massager comes with removable arm straps for versatile use. Align the strap connectors with the corresponding slots on the massager and secure them firmly.

## OPERATING INSTRUCTIONS

The massager features intuitive controls for a customized massage experience.

1. **Power On/Off:** Press and hold the power button located on the control panel to turn the massager on or off.
2. **Adjusting Massage Direction:** Briefly press the direction button to reverse the rotation of the shiatsu massage nodes.
3. **Activating Heat Function:** Press the heat button to activate the warming heat function. Press it again to turn the heat off. The heat function provides soothing warmth to enhance muscle relaxation.
4. **Adjusting Speed/Intensity:** Use the speed control buttons (if available on your model) to increase or decrease the massage intensity.



Figure 3: Control panel for adjusting speed, rotation, and heat settings.

## USAGE MODES AND APPLICATION

The multifunctional design allows for versatile application across various body parts.

- **Neck and Shoulder Massage:** Place the massager around your neck and shoulders, using the arm straps to adjust pressure and position.
- **Back Massage:** Position the massager on your upper, middle, or lower back. The arm straps can be used to hold it in place, or it can be secured to a chair using the integrated chair strap (if applicable).
- **Leg and Foot Massage:** Place the massager under your calves, thighs, or feet for targeted relief. The arm straps can be used to pull the massager closer for deeper foot massage.
- **Hands-Free Use:** The massager can be used without the arm straps by simply placing it on the desired body area while relaxing.

## SHARPER IMAGE®



### 01 SHIATSU

Massage multiple muscle groups with deep kneading action

### 02 MULTIFUNCTION STRAP

Choose your style – strapless, arm straps, chair strap or hands-free strap

### 03 FULL BODY

Target any area – neck, upper back, lower back, thighs, feet and more

### 04 RECHARGEABLE

Plug in to charge to enjoy a full hour of shiatsu massage completely cord-free



Figure 4: The massager can be used on various body parts including the neck, back, and legs.



Figure 5: Example of using the multifunctional strap for lower back massage.



**MULTIFUNCTIONAL**  
Target any area



Figure 6: The massager can also be used for foot massage.

## MAINTENANCE AND CARE

Proper maintenance ensures the longevity and optimal performance of your massager.

- **Cleaning:** Always unplug the massager before cleaning. Wipe the surface with a soft, slightly damp cloth. Do not use abrasive cleaners or immerse the unit in water. The removable strap can be cleaned separately according to its care instructions.
- **Storage:** Store the massager in a cool, dry place, away from direct sunlight and extreme temperatures. Avoid placing heavy objects on top of the massager.

## TROUBLESHOOTING

Problem	Possible Cause	Solution
---------	----------------	----------

Problem	Possible Cause	Solution
Massager does not turn on.	Battery is depleted or not charged.	Connect the DC adapter and charge the massager fully. Ensure the adapter is properly plugged into both the massager and the power outlet.
Heat function is not working.	Heat function is not activated.	Press the heat button to activate the warming function. Allow a few minutes for the heat to become noticeable.
Massage feels too intense or too weak.	Incorrect pressure applied or speed setting.	Adjust the pressure by using the arm straps or repositioning the massager. If applicable, adjust the speed setting using the control buttons.

## SPECIFICATIONS

Feature	Detail
Model Number	1014508
Product Dimensions	40.13 x 18.54 x 19.81 cm
Item Weight	1.81 kg
Power Source	Battery Powered (Rechargeable)
Battery Type	1 Lithium Ion battery required
Manufacturer	MerchSource
Color	Gray/Cordless



Figure 7: Product dimensions for reference.

## WARRANTY AND SUPPORT

For warranty information or technical support, please refer to the warranty card included with your product or contact Sharper Image customer service. Contact details are typically available on the product packaging or the official Sharper Image website. Please have your model number (1014508) and purchase date ready when contacting support.