

Snefe SNEFESMOJF

User Manual

SNEFE ULTRASONIC ELECTRIC TOOTHBRUSH

Model: SNEFESMOJF

Introduction

This manual provides essential information for the safe and effective use of your Snefe Ultrasonic Electric Toothbrush. Please read these instructions thoroughly before first use and retain for future reference.

Product Overview

The Snefe Ultrasonic Electric Toothbrush is designed to provide a comprehensive oral cleaning experience. It features a unique U-shaped brush head and high-frequency sonic technology for efficient cleaning and gum protection.



Image: The Snefe U-shaped ultrasonic electric toothbrush, white in color, resting on its wireless charging base. The U-shaped brush head is illuminated with blue light, indicating its operational state or charging status.

Key Features:

- **Unique U-Shaped Design:** Engineered to fit the gum line more effectively than traditional toothbrushes, enhancing cleaning and oral protection through high-frequency sonic technology.
- **Upgraded Brush Silicone:** Constructed from food-grade silicone with soft bristles to help prevent gum recession,

ensuring a safe and comfortable brushing experience.

- **Wireless Charging & IPX7 Waterproof:** The device is IPX7 waterproof, allowing for use in wet environments like baths and showers. It fully charges in approximately 2 hours via its charging dock, providing up to one week of use.
- **Efficient Teeth Cleaning:** Achieves 360° cleaning in 45 seconds with an automatic timer. The hands-free design allows for multitasking during brushing.
- **Easy to Use:** Features a simple one-button operation with 4 distinct cleaning modes: Gentle (massage), Powerful (comfortable), Cleansing (brush), and Whitening (recommended with whitening gel).

Safety Information

Important Safety Instructions:

- Do not immerse the main unit in water for extended periods.
- Always clean and dry the toothbrush after each use.
- Avoid using excessive force while brushing to prevent damage to teeth or gums.
- Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

Setup

1. **Initial Charge:** Before first use, place the toothbrush onto its wireless charging dock. Ensure the toothbrush is properly seated. A full charge typically takes 2 hours.
2. **Attach Brush Head:** Align the U-shaped brush head with the power shaft on the main body and gently push it into place until secure.



Image: The Snefe U-shaped toothbrush is shown on its wireless charging base, connected via a USB cable. Text indicates "1 Hrs Charge, 7 Days Working," highlighting its efficient charging and long battery life.

Operating Instructions

Follow these steps for optimal cleaning:

1. **Apply Toothpaste:** This product can be used with liquid toothpaste or ordinary toothpaste. If using ordinary toothpaste, wet the U-shaped brush head first, then apply toothpaste evenly.
2. **Position in Mouth:** Place the U-shaped brush head into your mouth, allowing your upper and lower teeth to gently fit onto the alveolar (the part of the jawbone that holds the teeth).
3. **Select Mode & Start:** Press the power button to turn on the toothbrush. Cycle through the modes by pressing the button again to select your desired cleaning mode:
 - **Gentle Mode:** For massage.

- **Powerful Mode:** For comfortable cleaning.
 - **Cleansing Mode:** Standard brushing.
 - **Whitening Mode:** Recommended for use with whitening gel.
4. **Brush:** While the toothbrush is operating, gently shake your hands to move the device, ensuring it touches all left and right teeth. The recommended frequency is approximately 63 times per minute. A comfortable and fast tooth-cleaning process can be completed in 45 seconds, after which the device will automatically stop.
5. **Rinse:** After brushing, remove the toothbrush from your mouth.

STEPS TO USE



STEP 1 Fill the toothpaste into the brush head



STEP 2 Choose cleaning mode



STEP 3 Move tooth brush left and right



STEP 4 Rinse the tooth brush with water

Use Foam Toothpaste is recommended.

Image: A visual guide illustrating the four steps of using the toothbrush: 1. Filling toothpaste into the brush head, 2. Choosing a cleaning mode, 3. Moving the toothbrush left and right in the mouth, and 4. Rinsing the brush with water. Foam toothpaste is recommended.

3X Cleaning Effect than Traditional Toothbrush

High frequency vibration massages & cleans gum line

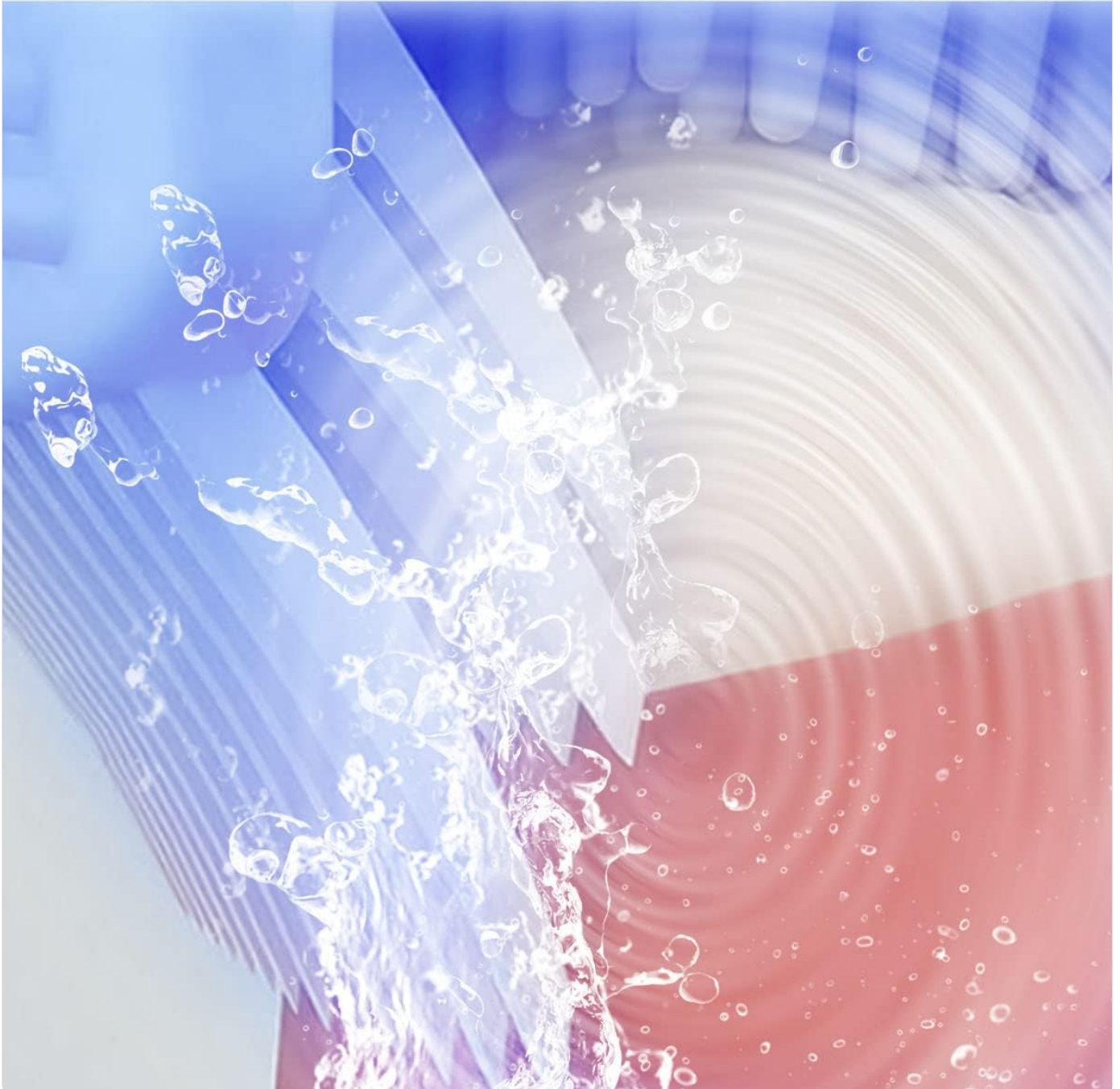


Image: A woman demonstrates the U-shaped toothbrush in her mouth. An overlay lists the four available modes: Strong, Comfortable, Massage, and Whiten, indicating the versatility of the device.

Note: If the brush cannot reach your back teeth effectively, gently move the device from side to side to ensure complete coverage.

Cleaning and Maintenance

Proper cleaning and maintenance will extend the life of your toothbrush.

- **Brush Head Cleaning:** After each use, rinse the U-shaped brush head thoroughly with water.
- **Main Unit Cleaning:** Wipe the main body of the toothbrush with a dry cloth. Do not immerse the main unit in water for prolonged periods.
- **Storage:** Store the toothbrush in a clean, dry place. Ensure the charging base is also kept dry, especially during charging.



Image: The Snefe U-shaped toothbrush is shown partially submerged in water, with bubbles around it, illustrating its IPX7 waterproof capability for easy cleaning and use in wet conditions.

Troubleshooting

If you encounter issues with your toothbrush, please refer to the following common solutions:

Problem	Possible Cause	Solution
Toothbrush does not turn on.	Low battery or not charged.	Place the toothbrush on the charging base and ensure it is fully charged (approx. 2 hours).
Toothbrush not charging.	Improper placement on charging base; charging cable/adaptor issue.	Ensure the toothbrush is correctly seated on the charging base. Check if the USB cable and adaptor are functioning correctly. Try a different USB port or adaptor.
Cleaning effectiveness seems low.	Incorrect usage technique; worn brush head; insufficient toothpaste.	Review operating instructions, ensuring proper movement and contact with all teeth. Ensure adequate toothpaste is applied. Consider replacing the brush head if it appears worn.
Device stops prematurely.	Automatic timer activated.	The toothbrush has an automatic 45-second timer. This is normal operation. If it stops before 45 seconds, ensure it is fully charged.

Specifications

Feature	Detail
Product Dimensions	2.95 x 2.17 x 4.53 inches
Item Weight	8.78 ounces
Model Number	SNEFESMOJF
Batteries	1 Lithium Ion battery required
Brand	Snefe
Age Range (Description)	Adult
Special Feature	Brushing Timer
Power Source	Battery Powered
Item Firmness Description	Soft

Upgraded Wireless charging base

Keep the toothbrush and base dry while charging



Image: A diagram illustrating the dimensions of the Snefe U-shaped toothbrush, showing its height as 11.5cm (4.52in) and width as 7.5cm (2.95in).

Warranty and Support

This product comes with a standard manufacturer's warranty. For specific warranty details, please refer to the product packaging or contact Snefe customer support.

For technical support, troubleshooting assistance, or to inquire about replacement parts, please visit the official Snefe website or contact their customer service department.

Contact information for support is typically found on the product packaging or the manufacturer's official website.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question