



Manuals.plus /

› **CallToU /**

› CallToU Wireless Caregiver Pager Call Button for Elderly, Medical Alert Systems for Seniors/Elderly Monitoring with 2 Call Bell, 500+Feet, 1 Plug-in Receiver+1 USB Receiver for Home Patients Emergency CC05-06WH 2-2

CallToU CC05-06WH 2-2

CallToU Wireless Caregiver Pager System User Manual

Model: CC05-06WH 2-2

PRODUCT OVERVIEW

The CallToU Wireless Caregiver Pager System is designed to provide reliable and convenient communication for caregivers and individuals needing assistance. This system includes a plug-in receiver, a portable USB receiver, and two wearable call buttons, ensuring alerts can be received throughout the home. It is an ideal solution for elderly care, home patients, and emergency situations, offering peace of mind through immediate notification.



Figure 1: CallToU Wireless Caregiver Pager System components.

COMPONENTS INCLUDED

- 1 x Plug-in Receiver
- 1 x USB Receiver
- 1 x Wrist-worn Call Button (with pre-installed battery)
- 1 x Pendant-style Call Button (with pre-installed battery)
- 1 x USB Charging Cable (for USB Receiver)
- 1 x Magnetic Sticker (for USB Receiver)
- User Manual (this document)

SETUP

1. Plug-in Receiver Setup

1. **Placement:** Select an electrical outlet in a central location or where the caregiver is most likely to hear the alert.
2. **Power On:** Plug the receiver directly into a standard electrical outlet. The indicator light will illuminate.
3. **Ringtone and Volume Selection:**
 - Press the "Music" button on the side of the plug-in receiver to cycle through 55 available ringtones.
 - Press the "Volume" button to adjust the volume level (5 levels, from 0dB to 110dB).



Figure 2: Plug-in Receiver with ringtone and volume controls.

2. USB Receiver Setup

1. **Charging:** Before first use, fully charge the USB receiver using the provided USB charging cable. A red light indicates charging, and it turns off when fully charged. Charging time is approximately 2 hours.
2. **Placement:** The USB receiver is portable and can be carried in a pocket or attached to metal surfaces like a refrigerator using the magnetic sticker.
3. **Ringtone and Volume Selection:**
 - The USB receiver offers 4 ringtones. Use the side buttons to select.
 - Adjust volume using the side buttons (5 levels, from 0dB to 110dB).



Figure 3: USB Receiver charging and battery indicators.



Figure 4: USB Receiver with magnetic attachment for easy portability.

3. Call Button Pairing (Pre-paired from factory)

The call buttons are pre-paired with the receivers from the factory. No manual pairing is typically required for initial setup. If a button loses connection or you wish to add more buttons, refer to the "Troubleshooting" section for re-pairing instructions.

OPERATING INSTRUCTIONS

1. Using the Call Buttons

Both the wrist-worn and pendant-style call buttons are designed for easy activation. To send an alert, simply press the large red "SOS" button on the wrist button or the blue "CALL" button on the pendant button.

- **Wrist Button:** Wear comfortably on the wrist. Ideal for quick access.
- **Pendant Button:** Can be worn around the neck using the included string, or placed near the bed, in the bathroom, or other strategic locations.

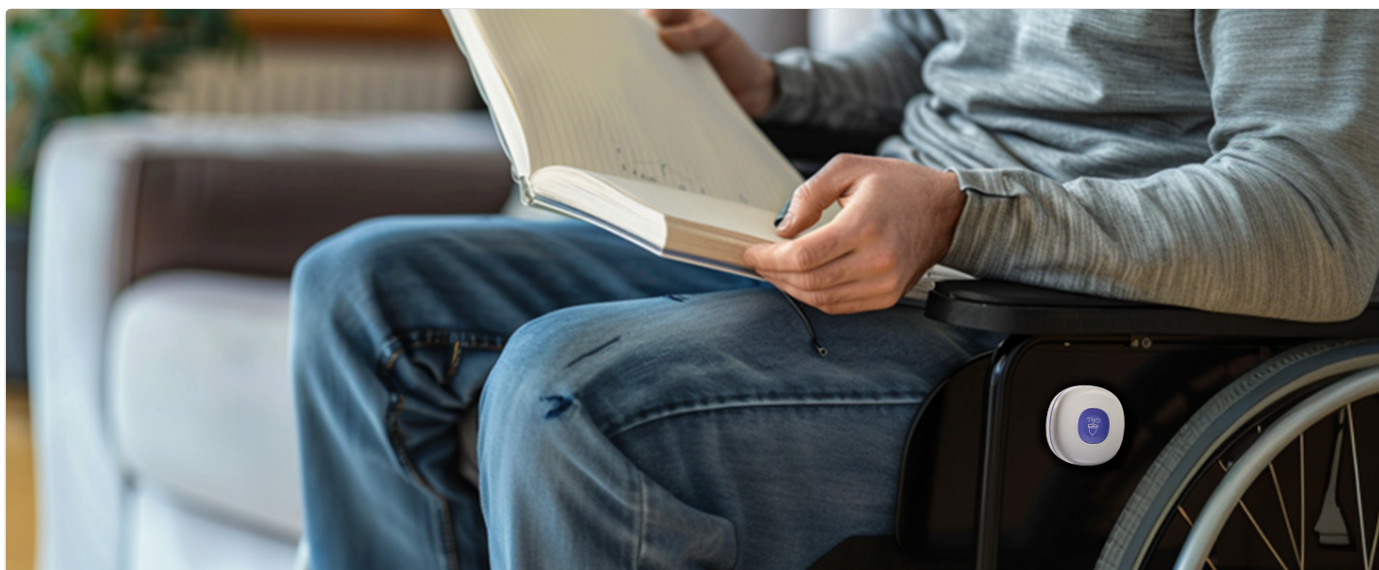


Figure 5: Various ways to use the CallToU SOS buttons.

2. Receiving Alerts

When a call button is pressed, both the plug-in receiver and the USB receiver will sound the selected ringtone and flash

their indicator lights. The system has an open area working range of 500+ feet, allowing for alerts to be received across different rooms or even outdoors within range.

Wireless 500FT

The signal of call buttons carries across 500FT easily. Never miss a call



Figure 6: Wireless range of the CallToU Caregiver Pager System.

Your browser does not support the video tag.

Video 1: Demonstration of the CallToU Wireless Caregiver Pager with Portable USB Receiver, showing setup and alert functionality.

1. Battery Management

- **USB Receiver:** Recharge the USB receiver when the low battery reminder activates. A full charge takes approximately 2 hours and provides 3-7 days of normal use or 14 days of standby battery life.
- **Call Buttons:** The call buttons use AAA batteries (included). While the product description mentions "no need to replace batteries" for the USB receiver, the specifications indicate AAA batteries for the system. Monitor button performance; if alerts become inconsistent, consider replacing the button's battery.

2. Water Resistance

The call buttons are IP55 waterproof, meaning they are protected against dust and low-pressure water jets. They are suitable for use in bathrooms or during light splashes, but should not be submerged in water.

IP55 Waterproof



Figure 7: IP55 Waterproof rating of the Call Buttons.

3. Cleaning

Wipe the devices with a soft, dry cloth. Do not use abrasive cleaners or immerse any part of the receivers in water.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Receiver does not sound when button is pressed.	<ul style="list-style-type: none"> ◦ Receiver is unplugged or uncharged. ◦ Volume is set to 0dB. ◦ Call button battery is low or depleted. ◦ Out of range. ◦ Button is unpaired. 	<ul style="list-style-type: none"> ◦ Ensure plug-in receiver is securely plugged in. Charge USB receiver. ◦ Increase receiver volume. ◦ Replace the AAA battery in the call button. ◦ Move closer to the receiver. ◦ Re-pairing: Press and hold the volume button on the receiver until it beeps. Then, press the call button. The receiver should sound, indicating successful re-pairing.
Short battery life on USB receiver.	<ul style="list-style-type: none"> ◦ Not fully charged. ◦ Frequent alerts. 	<ul style="list-style-type: none"> ◦ Ensure the USB receiver is charged for the recommended 2 hours. ◦ Battery life varies with usage frequency.
Call button is hard to press or unresponsive.	<ul style="list-style-type: none"> ◦ Battery issue. ◦ Physical damage. 	<ul style="list-style-type: none"> ◦ Replace the AAA battery in the call button. ◦ If physically damaged, contact support for replacement.

SPECIFICATIONS

Feature	Detail
Model Number	CC05-06WH 2-2
Brand	CallToU
Power Source	Battery Powered (Call Buttons), Plug-in (Plug-in Receiver), USB Rechargeable (USB Receiver)
Battery Type (Call Buttons)	1 AAA battery required (included)
Voltage	260 Volts (likely refers to receiver input)
Control Method	Touch
Noise Level	Up to 110 dB
Mounting Type	Tabletop Mount (USB Receiver), Wall Plug (Plug-in Receiver)
Maximum Range	500+ Feet (open area)
Sensor Technology	Immediate
Item Weight	9.9 ounces (total package)
Product Dimensions	2.8 x 2.4 x 3 inches (likely for a receiver)
Style	2 Receiver + 2 Call Button
Waterproof Rating (Buttons)	IP55

WARRANTY AND SUPPORT

CallToU provides comprehensive and caring service support for this product. The product comes with a **Lifetime Warranty**. For any questions or assistance regarding your CallToU caregiver pager, please feel free to contact us via Amazon's messaging system or through the contact information provided on the official CallToU website. For additional resources and product information, you may visit the [CallToU Store on Amazon](#).