

## My Weigh XL-700

# My Weigh XL-700 Talking Bathroom Scale User Manual

## 1. INTRODUCTION

---

The My Weigh XL-700 is a high-capacity talking body scale designed for ease of use and accurate weight measurement. It features advanced strain gauge technology and a comfortable, non-slip platform. This manual provides essential information for setting up, operating, and maintaining your scale to ensure optimal performance.



*Image 1: My Weigh XL-700 Talking Bathroom Scale, top view.*

## 2. KEY FEATURES

---

- **Talking Function:** Announces weight in 4 languages (English, French, German, Spanish). Features an on/off switch for the talking function.
- **High Capacity:** Measures up to 700 lbs (320 kg or 50 stones) with an accuracy of 0.2 lb (0.1 kg).
- **Multiple Weighing Modes:** Easily switch between pounds (LB), kilograms (KG), and stones (ST).
- **Extra-Large Platform:** A spacious 19" x 11.75" non-slip textured platform for stability and comfort.
- **Clear Display:** Large, reverse-lit 2" high, 4-digit LCD display for easy reading.
- **Automatic Operation:** Turns on automatically when stepped on and shuts off after 10 seconds of non-use.


# Features



Easily track your weight with My Weigh Digital Scale's convenient features.




Bright  
LCD Display

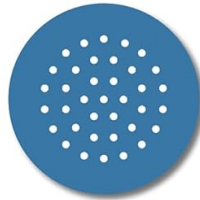

Non-slip  
Surface



Speaks in 4  
Languages

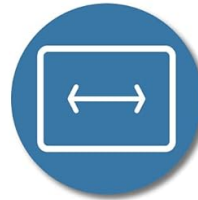

Extra Large  
Platform



Auto Shut off



Image 2: Overview of the My Weigh XL-700 scale's key features.

## 3. SETUP

### 3.1 Battery Installation

1. Locate the battery compartment on the underside of the scale.
2. Open the battery cover.
3. Insert 4 AA batteries, ensuring correct polarity (+/-). The scale includes batteries.
4. Close the battery cover securely.

### 3.2 Placement

For accurate readings, place the scale on a hard, flat, and stable surface. Avoid carpets or uneven flooring, as this can affect accuracy.

# Non-slip Surface

The sleek black comfortable non-slip textured platform measures 19" long by 11.75" wide and the reinforced bracing minimizes mechanical deviation



Image 3: Underside of the scale showing battery compartment and unit/voice controls.

## 4. OPERATING INSTRUCTIONS

### 4.1 Taking a Measurement

1. Gently tap the scale with your foot or hand to activate it. The display will show "Hello" (if voice is enabled) and then zeros.
2. Wait for the scale to display "0.0" and say "I'm Ready" (if voice is enabled).
3. Step onto the center of the scale platform with bare feet, standing still.
4. Your weight will be displayed on the LCD screen and announced verbally (if voice is enabled).
5. The scale will automatically turn off after approximately 10 seconds of non-use, announcing "Good-bye" (if voice is enabled).

### 4.2 Changing Weighing Units

The XL-700 supports three weighing modes: pounds (LB), kilograms (KG), and stones (ST).

1. With the scale off, locate the "UNIT" button on the underside of the scale (refer to Image 3).
2. Press the "UNIT" button repeatedly to cycle through the available units: LB, KG, ST.

3. The selected unit will be indicated on the LCD display.

### 4.3 Voice Function Control

The talking feature can be turned on or off, and the language can be changed.

1. **Voice On/Off:** Locate the voice switch on the underside of the scale (refer to Image 3). Slide the switch to the desired position.
2. **Language Selection:** With the voice function enabled, press and hold the "UNIT" button on the underside of the scale to cycle through English, German, French, and Spanish. Release the button when your desired language is announced.

## Friendly Talking Feature

The multi-lingual talking feature can speak in Spanish, French, English, and German



Image 4: Illustration of the multi-lingual talking feature.



# Crisp LCD Display

The reverse lit LCD display  
is clear and easy to read

Dimensions:  
LCD: 2.125 in  
Digits: 1.825 in



*Image 5: Close-up of the large, reverse-lit LCD display.*

## 5. MAINTENANCE

---

### 5.1 Cleaning

Wipe the scale with a damp cloth and mild detergent. Do not use abrasive cleaners or immerse the scale in water.

### 5.2 Battery Replacement

If the display becomes dim or the voice function is unclear, it may be time to replace the batteries. Follow the battery installation steps in Section 3.1.

## 6. TROUBLESHOOTING

---

- **No Display/No Power:** Check if batteries are installed correctly and are not depleted. Replace batteries if necessary.
- **Inaccurate Readings:** Ensure the scale is on a hard, flat surface. Avoid moving the scale between measurements. Recalibrate by stepping off, waiting for it to turn off, then stepping on again.

- **Voice Not Working:** Check the voice on/off switch on the underside of the scale. Ensure batteries are not low.
- **"Lo" or "Err" Message:** Indicates low battery or an error. Replace batteries or re-attempt measurement.

## 7. SPECIFICATIONS

Feature	Detail
Model Number	XL-700
Capacity	700 lbs / 320 kg / 50 st
Accuracy	0.2 lb / 0.1 kg
Platform Dimensions	19" L x 11.75" W
Power Source	4 AA batteries (included)
Display Type	Large, reverse-lit LCD
Special Features	Talking function (4 languages), Non-slip surface, Auto shut-off
Product Weight	8.2 Pounds

## 8. WARRANTY & SUPPORT

The My Weigh XL-700 Talking Bathroom Scale comes with a**30-year manufacturer's warranty**. For support, please refer to the contact information provided with your product packaging or visit the official My Weigh website.

## 9. PRODUCT VIDEOS

### My Weigh XL-700 Scale Overview

Your browser does not support the video tag.

*Video 1: A detailed overview of the My Weigh XL-700 Talking Bathroom Scale, highlighting its features and benefits for users.*

### Demonstration of My Weigh XL-700 Features

Your browser does not support the video tag.

*Video 2: A demonstration showcasing the functionality and user experience of the My Weigh XL-700 scale, including its talking feature and wide platform.*

### Review of My Weigh Talking Bathroom Scale

Your browser does not support the video tag.

*Video 3: An honest review of the My Weigh Talking Bathroom Scale, detailing its performance and*

*user-friendly aspects.*