



Manuals.plus /

› waxkiss /

› Waxkiss PWH-5000B Paraffin Wax Machine User Manual

waxkiss PWH-5000B

Waxkiss PWH-5000B Paraffin Wax Machine User Manual

Model: PWH-5000B

1. PRODUCT OVERVIEW

The Waxkiss PWH-5000B Paraffin Wax Machine is designed for paraffin wax treatments for hands, feet, and elbows. It features a large 5000ml capacity and precise temperature control for a comfortable and effective experience.



Image 1.1: The Waxkiss PWH-5000B Paraffin Wax Machine shown with included accessories.

Package Contents:




Image 1.2: Detailed view of all components included in the Waxkiss PWH-5000B kit.

- 1 x Waxkiss PWH-5000B Paraffin Wax Machine
- 4 x Packages of Paraffin Wax (Lavender, 0.44 lbs each)
- 1 x Heat Insulating Grid
- 1 x Pair of Cotton Mitts for Hands
- 1 x Pair of Cotton Booties for Feet
- Disposable PE Mitts for Hands
- Disposable PE Booties for Feet
- 1 x User Manual

BENEFITS IN PARAFFIN WAX SPA



Smooth and moisturize your dry skin

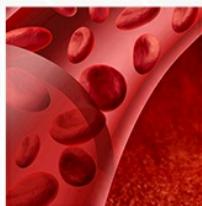

Promote blood circulation



Remove ageing corneous cells



Relieve fatigue



Heat therapy and relief the pain

2. SAFETY INFORMATION

Please read all safety instructions carefully before using the paraffin wax machine to prevent injury or damage.

- **Electrical Safety:** Ensure the machine is plugged into a grounded outlet. Do not operate with wet hands or near water. Unplug when not in use or before cleaning.
- **Hot Wax Safety:** Paraffin wax can reach high temperatures. Always test the wax temperature on a small area of your skin (e.g., wrist) before full immersion to avoid burns. Do not immerse body parts with open wounds, infections, or severe skin conditions.
- **Avoid Contact with Eyes:** If wax comes into contact with eyes, rinse immediately with cold water and seek medical attention.
- **Children and Pets:** Keep the machine out of reach of children and pets, especially when in operation or cooling.
- **Ventilation:** Use in a well-ventilated area.
- **Supervision:** Never leave the machine unattended while heating wax.
- **Material Safety:** The machine is constructed from safe materials.

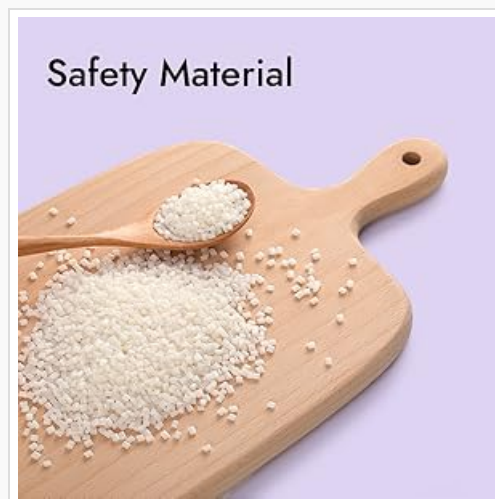


Image 2.1: Illustration of the safe materials used in the product's construction.

3. SETUP INSTRUCTIONS

1. **Unpack:** Carefully remove all components from the packaging.
2. **Placement:** Place the paraffin wax machine on a stable, heat-resistant, and level surface. Ensure there is adequate space around the unit for ventilation.
3. **Add Wax:** Place the paraffin wax blocks into the machine's basin. The machine has a large capacity suitable for hands, feet, and elbows.
4. **Power On:** Plug the machine into a suitable electrical outlet.

LARGE CAPACITY 5000ML

plenty big enough for hands and feet



Image 3.1: The machine's large capacity (5000ml) is suitable for hands and feet, with dimensions indicated.

4. OPERATING INSTRUCTIONS

Control Panel Overview:



Image 4.1: The control panel features a digital display, temperature/time adjustment buttons, a set button, and an ON/OFF button.

Melting the Wax:

1. Press the ON/OFF button to turn on the machine.
2. The machine offers two working modes: Temperature Heating Mode and Timed Heating Mode.
3. **Temperature Heating Mode:** Press the 'SET' button to toggle to this mode. Use the '+' and '-' buttons to adjust the temperature from 86°F to 185°F (30°C to 85°C). The machine will heat to and maintain the set temperature.
4. **Timed Heating Mode:** Press the 'SET' button again to toggle to this mode. Use the '+' and '-' buttons to set the heating time from 1 to 99 minutes. The machine will automatically heat the wax during this period.
5. It takes approximately 50 minutes to completely melt 2 lbs of wax at 185°F. Wait until the wax is fully melted and the digital display reaches 50°C (122°F), indicated by three 'Di, Di, Di' sounds, before proceeding.

Two Working Mode



Image 4.2: Explanation of the two operational modes: Temperature Heating Mode and Timed Heating Mode.

Applying the Wax Treatment:

1. Once the wax is melted and at a safe, comfortable temperature, place the heat insulating grid into the basin.
2. Slowly dip your hand, foot, or elbow into the melted wax. Remove it for a few seconds, allowing the wax to solidify slightly, then re-dip. Repeat this process 3-5 times to build up a thick layer of wax.
3. Immediately after the final dip, cover the treated area with a disposable PE mitt or bootie.
4. For enhanced heat retention, place a cotton mitt or bootie over the disposable layer.
5. Relax and allow the wax to work for approximately 10-15 minutes.
6. After the recommended time, remove the cotton mitt/bootie, then peel off the wax and disposable liner. The wax should come off easily, leaving your skin feeling soft and smooth.




1

Put the wax into the pot and start heating.



2

After melting place the heat insulating grid and dip your hands into the warmer.

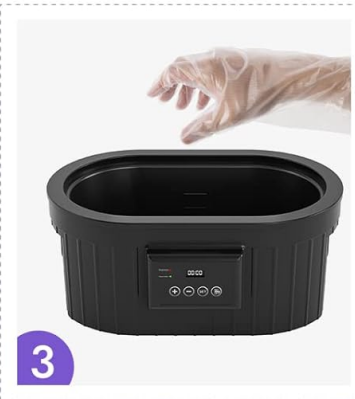

3

Warp your hand with glove liner.




4

Put on heat trapping cotton mitt.



5

Remove paraffin wax after 15mins.

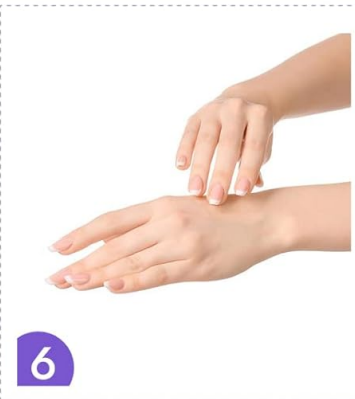

6

Your hands will be soft and fine like silk.

Image 4.3: Step-by-step guide for applying paraffin wax treatment.



PACKAGE INCLUDED

- 1 x 5000ml paraffin wax warmer
- 4 x Lavender paraffin wax, 0.44 each pack
- 1 x The heat insulating grid
- 2 x Hand thermal mitts
- 2 x Feet thermal mitts
- 100 x Plastic Gloves
- 100 x Plastic Feet Gloves
- 1 x User manual

The wax is reusable for one person using for 2-4 times.

5. BENEFITS OF PARAFFIN TREATMENT

Regular paraffin wax treatments can offer several benefits:

- **Skin Softening:** Smooths and softens dry, rough skin.
- **Joint Pain Relief:** Provides soothing heat therapy that can help relieve joint pain, particularly beneficial for arthritis in hands and feet.
- **Improved Circulation:** Promotes blood circulation.
- **Exfoliation:** Helps remove aging corneous cells, revealing smoother skin.
- **Relaxation:** Contributes to a relaxed and rejuvenated feeling.



Image 5.1: Visual representation of the benefits of paraffin wax treatments.

6. MAINTENANCE

Cleaning the Machine:

1. Ensure the machine is unplugged and completely cooled before cleaning.
2. Remove any remaining solid wax. You can reheat the wax slightly to soften it for easier removal, then turn off and unplug the machine before scraping.
3. Wipe the interior and exterior surfaces with a soft, damp cloth. Do not use abrasive cleaners or immerse the unit in water.
4. The spacious design of the basin facilitates easy cleaning.



Image 6.1: Illustration of cleaning the paraffin wax machine.

Wax Reuse:

For personal use, the paraffin wax can typically be reused 2-4 times. It is recommended to replace the wax regularly for hygiene and optimal performance.

7. TROUBLESHOOTING

- **Wax not melting or melting slowly:** Ensure the machine is plugged in and turned on. Verify the temperature setting is appropriate (e.g., 185°F for initial melting). Allow sufficient time for the wax to melt; 2 lbs of wax can take up to 50 minutes at maximum temperature.
- **Machine not turning on:** Check the power connection and ensure the outlet is functional.
- **Wax too hot/cold:** Adjust the temperature using the '+' and '-' buttons. Always test the wax temperature before application. The digital display shows the set temperature, not necessarily the actual temperature during the reservation option.
- **Machine slides during operation:** Ensure the machine is placed on a stable, non-slip surface.

8. SPECIFICATIONS

Feature	Specification
Model Number	PWH-5000B
Capacity	5000ml
Temperature Range	86°F to 185°F (30°C to 85°C)
Time Setting Range	1 to 99 minutes
Package Dimensions	16.93 x 11.46 x 7.99 inches
Item Weight	6.94 Pounds
UPC	734257215466
Manufacturer	waxkiss

9. WARRANTY AND SUPPORT

For warranty information, technical support, or any product-related inquiries, please contact Waxkiss customer service through the retailer's platform or the official Waxkiss website. Please retain your proof of purchase for warranty claims.
