

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [foreverLily](#) /

› [foreverLily Face Skin Care Machine ODM-CLC200608 Instruction Manual](#)

foreverLily ODM-CLC200608

foreverLily Face Skin Care Machine ODM-CLC200608 Instruction Manual

1. INTRODUCTION

Thank you for choosing the foreverLily Face Skin Care Machine. This device is designed to address various skin concerns through multi-functional light therapy and heat treatment. Please read this manual thoroughly before use to ensure safe and effective operation.



Image 1: The foreverLily Face Skin Care Machine in use, demonstrating its design and remote control.

Package Contents:

- 1 x Main Machine
- 1 x Adapter
- 1 x Protection Goggles
- 1 x User Manual
- 1 x Box

Package includes:

- 1*Main Machine
- 1*Adapter
- 1*Protection Goggles
- 1*User Manual
- 1*Box



Brand: Forverlily

Item Type: LED Photon Machine

Color: White

Material: ABS

Power: 5.76-63.2W

Current: 3A

Number of lamp beads: 96

Temperature range: 30°C~60°C

Control time: 1~30 minutes

Weight: 1265g

Image 2: Visual representation of the items included in the product package.

2. SAFETY INFORMATION

Please observe the following safety precautions to prevent injury or damage to the device:

- Do not use if you have any skin conditions such as open wounds, severe acne, or skin infections without consulting a medical professional.
- Always wear the provided protection goggles during light therapy sessions to protect your eyes.
- Keep the device away from water and moisture. Do not use in wet environments.
- Do not disassemble or attempt to repair the device yourself. Contact customer support for assistance.
- Keep out of reach of children.
- Ensure the power adapter is compatible with your local voltage supply.
- Discontinue use immediately if you experience any discomfort or adverse reactions.

3. SETUP

1. **Unpack the Device:** Carefully remove all components from the packaging.
2. **Connect Power:** Connect the adapter to the device and then plug it into a power outlet.
3. **Positioning:** Place the device on a stable, flat surface. The U-shape design allows for comfortable use while lying down, ensuring the light covers your face and neck area.

HOME&BEAUTY SALON USE SKIN SPA

MINI | SAFE | EASY USE

NEW-born Skin Care **20 minutes/time, 2-3 tims/week**



Image 3: The device positioned for use, highlighting its ergonomic design for home or salon settings.

4. OPERATING INSTRUCTIONS

Follow these steps for optimal use of your foreverLily Face Skin Care Machine:

1. **Clean Your Skin:** Ensure your face is clean and free of makeup or other products before starting a session.
2. **Apply Skincare Products (Optional):** For enhanced results, you may apply a mask, essence, or lotion to your skin before use.
3. **Connect Power and Turn On:** Connect the power supply and short press the power button on the control panel to turn the device on.

4. **Select Time:** Use the 'TIME+' or 'TIME-' buttons to set your desired session duration (5, 10, 15, 20, 25, or 30 minutes).
5. **Select Light Mode:** Choose your preferred light therapy mode using the dedicated buttons:
 - **R/Y (Red/Yellow Cycle):** Cycles between Red and Yellow light.
 - **B/P (Blue/Purple Cycle):** Cycles between Blue and Purple light.
 - **G/C (Green/Cyan Cycle):** Cycles between Green and Cyan light.
 - **CYCLE (White Light/7-Color Cycle):** Cycles through all available colors, including White light.
6. **Adjust Temperature (Optional):** Use the 'TEMP+' or 'TEMP-' buttons to activate and adjust the heat treatment function, ranging from 30°C to 60°C (86-140°F).
7. **Wear Goggles:** Always wear the provided protection goggles before the light therapy begins.
8. **Enjoy the Session:** Lie down comfortably and allow the device to work for the selected duration. A typical session is 15-20 minutes.

1 Time+ **2** Time-

Time selection: 5/10/15/20/25/30 min

3 Switch
Short press: turn on/turn off

4 Red/yellow/red-yellow cycle

5 Blue/purple/blue-purple cycle

6 Green/cyan/green-cyan cycle

7 White light/7-color cycle



HOW TO USE

-  1. Clean your skin
-  2. Apply a mask or apply skin care products
-  3. Connect the power supply, turn on and select the light and time
-  4. Apply a mask or apply skin care products
-  5. Enjoy a 15-20 minute spa

Image 4: A visual guide detailing the usage steps, from skin preparation to enjoying the treatment.

Light Therapy Modes and Benefits:



Red(640nm)
Deep repair
Increase skin elasticity



Blue(415nm)
Inhibit bacteria
Remove acne



Yellow(590nm)
Improve wrinkles
Increase immunity

LIARTY



Green(532nm)
Lighten spots
Calms the skin



Purple
Improves fine lines



Cyan
Deep metabolism
Oil control

Image 5: An overview of the various light colors and their specific benefits for skin treatment.

- **Red Light (640nm):** Promotes deep repair and increases skin elasticity.
- **Blue Light (415nm):** Inhibits bacteria and helps remove acne.
- **Yellow Light (590nm):** Improves the appearance of wrinkles and supports skin immunity.
- **Green Light (532nm):** Helps lighten spots and calms the skin.
- **Purple Light:** Improves fine lines.
- **Cyan Light:** Supports deep metabolism and oil control.

Color Spectrum Cycling IMPROVE YOUR SKIN

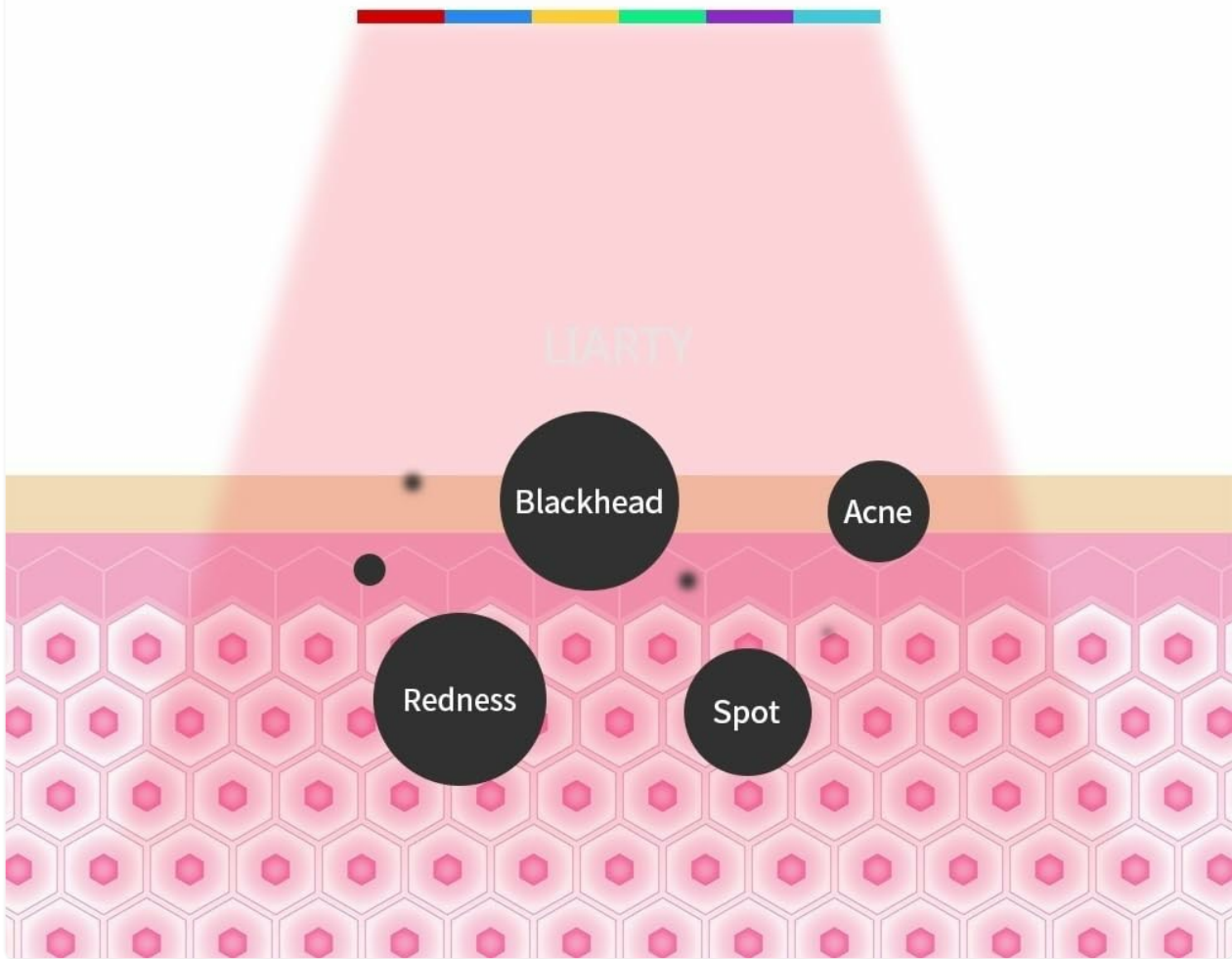


Image 6: Illustration of how the device's color spectrum cycling technology targets various skin concerns.

Heat Treatment:

The device features a heat treatment function that can be adjusted between 30°C and 60°C. Heat helps to dilate pores, promote the absorption of skincare products, and improve blood circulation for enhanced results.

HEAT TREATMENT

Dilates pores / Promotes absorption / Promotes blood circulation



30°C~60°C
(86-140°F)



Image 7: Detail of the heat treatment feature, showing temperature control via the remote.

5. MAINTENANCE

- **Cleaning:** After each use, gently wipe the surface of the device with a soft, dry cloth. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the device in a cool, dry place, away from direct sunlight and extreme temperatures. When not in use, you can fold the device for compact storage.
- **Adapter Care:** Ensure the power adapter is not bent or damaged. Unplug it from the wall when not in use for extended periods.

6. TROUBLESHOOTING

- **Device not turning on:**
 - Ensure the power adapter is securely connected to both the device and a working power outlet.
 - Check if the power outlet is functional by plugging in another appliance.

- **Lights not working:**

- Verify that the device is powered on and a light mode has been selected.
- If only some lights are off, contact customer support.

- **No heat or insufficient heat:**

- Ensure the heat function is activated and the temperature is set to your desired level using the 'TEMP+' button.

- **Remote control not responding:**

- Ensure there are no obstructions between the remote and the device.
- Check if the remote's batteries need replacement (if applicable).

If you encounter any issues not listed here or if troubleshooting steps do not resolve the problem, please contact customer support.

7. SPECIFICATIONS

Feature	Detail
Brand	foreverLily
Model Number	ODM-CLC200608
Item Type	LED Photon Machine
Material	ABS
Color	White
Number of Lamp Beads	96 PCS
Power	5.76-63.2W
Current	3A
Temperature Range	30°C ~ 60°C
Control Time	1-30 minutes
Item Weight	1.93 kg (1930 g)
Package Dimensions	38.28 x 27.89 x 21.69 cm

Note: Specifications are subject to change without prior notice.

8. WARRANTY AND SUPPORT

foreverLily is committed to providing a satisfactory shopping experience. While specific warranty details are not provided in this manual, please retain your proof of purchase for any potential warranty claims.

For any questions, concerns, or technical support regarding your foreverLily Face Skin Care Machine, please do not hesitate to contact our customer service team. Refer to the product packaging or the retailer's website for contact information.

